



Artisan Technology Group is your source for quality new and certified-used/pre-owned equipment

- FAST SHIPPING AND DELIVERY
- TENS OF THOUSANDS OF IN-STOCK ITEMS
- EQUIPMENT DEMOS
- HUNDREDS OF MANUFACTURERS SUPPORTED
- LEASING/MONTHLY RENTALS
- ITAR CERTIFIED SECURE ASSET SOLUTIONS

SERVICE CENTER REPAIRS

Experienced engineers and technicians on staff at our full-service, in-house repair center

*InstraView*SM REMOTE INSPECTION

Remotely inspect equipment before purchasing with our interactive website at www.instraview.com ↗

WE BUY USED EQUIPMENT

Sell your excess, underutilized, and idle used equipment. We also offer credit for buy-backs and trade-ins. www.artisanng.com/WeBuyEquipment ↗

LOOKING FOR MORE INFORMATION?

Visit us on the web at www.artisanng.com ↗ for more information on price quotations, drivers, technical specifications, manuals, and documentation

Contact us: (888) 88-SOURCE | sales@artisanng.com | www.artisanng.com

1700°C Controller, Programmable, With Communications

Lindberg/Blue M 1700°C Programmable Controllers provide multiple programs and multiple segments for ramp (up and down) and dwell (timed hold) temperature control. The controller visually displays ramp rate, dwell time, program segment and percent power output. A holdback feature allows the operator to set a “process vs setpoint” temperature value which, when exceeded, holds the program to allow the process to catch up. Please see page 35 for additional information.


The controller includes a selectable self-tuning feature which sets the best PID settings for the thermal process. Two digital displays simultaneously indicate actual temperature and

setpoint temperature. High limit overtemperature protection is standard. The control console includes a circuit breaker, power module, transformer and cooling fans.

Controllers include RS485 data port (communications card and port) for connection to remote computer, allowing modification, interrogation and data transfer of all instrument control and configuration parameter. Up to 30 units can be connected to one PC. Software is not included, but is available as an option. Please see page 35 for additional options and information.

Option B Overtemperature Control (OTC)

Adjustable digital overtemperature control, factory installed on selected control consoles with “B” suffix designation; see chart. Protects furnace and load in the event of primary control circuit failure. Overrides main controller and shuts off power to furnace if high limit is reached. Manual re-set required for safety. Operates via magnetic contacts through signal from independent thermocouple.

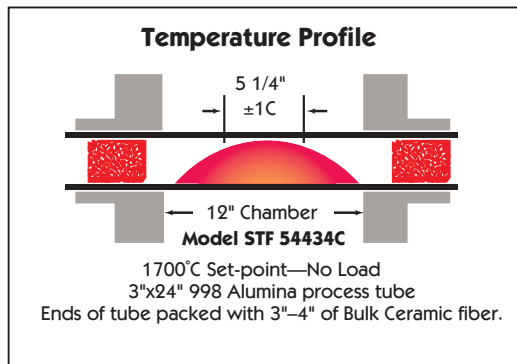

Uniformity Profile, Model STF54434C,
1700°C, No Load.



Programmable Controller, 1700°C
Model CC59256PCOMC.

1700° Controller, Programmable, with Communications

Controller Model No.	Digital	With Programmer	With Overtemp Control	Electrical Volts, Hz	Exterior Dimensions H x F-B x W in" (mm)	Ship Weight lbs (kg)
CC59256PCOMC	■	■		208/240V, 50/60Hz	10" (254) x 15" (381) x 21" (533.4)	130 (59)
CC59256PBCOMC	■	■	■	208/240V, 50/60Hz	10" (254) x 15" (381) x 21" (533.4)	130 (59)
CC59256PCOM2CKTC	■	■		208/240V, 50/60Hz	30" (762) x 15" (381) x 22" (558.8)	230 (105)
CC59256PBCOM2CKTC	■	■	■	208/240V, 50/60Hz	30" (762) x 15" (381) x 22" (558.8)	230 (105)

Note: Required power cord, hardwiring and interconnecting wiring are not included.

Vestibules

Optional Moldatherm® vestibules permit operation with 1", 2" and 3" O.D. process tubes for increased versatility. Two vestibules are required for each furnace. One 3" vestibule set is included and installed with each furnace.

Moldatherm Sleeves

Tube sleeves may be placed over customer supplied process tubes to reduce thermal shock to the process tube. Sleeves must be ordered separately. All tube sleeves are 3" long.

Furnace Model Number	Vestibule Catalog No. (Ordered Separately)		
	1" dia.	2" dia.	3" dia.
STF54434C	7219-2147-001	7219-2147-002	7219-2147-003*
STF54454C	7219-2147-013	7219-2147-012	7219-2147-011*

* One set of (2) included with furnace

Furnace Model Number	Sleeve Catalog No. (Ordered Separately)		
	1" dia.	2" dia.	3" dia.
STF54434C	7219-2134-001	7219-2134-002	7219-2134-003*
STF54454C	7219-2134-013	7219-2134-012	7219-2134-011*

* One set of (2) included with furnace

Actual performance may vary depending on load, chamber size, sample placement, ambient temperature and environmental conditions.



Artisan Technology Group is your source for quality new and certified-used/pre-owned equipment

- FAST SHIPPING AND DELIVERY
- TENS OF THOUSANDS OF IN-STOCK ITEMS
- EQUIPMENT DEMOS
- HUNDREDS OF MANUFACTURERS SUPPORTED
- LEASING/MONTHLY RENTALS
- ITAR CERTIFIED SECURE ASSET SOLUTIONS

SERVICE CENTER REPAIRS

Experienced engineers and technicians on staff at our full-service, in-house repair center

*InstraView*SM REMOTE INSPECTION

Remotely inspect equipment before purchasing with our interactive website at www.instraview.com ↗

WE BUY USED EQUIPMENT

Sell your excess, underutilized, and idle used equipment. We also offer credit for buy-backs and trade-ins. www.artisanng.com/WeBuyEquipment ↗

LOOKING FOR MORE INFORMATION?

Visit us on the web at www.artisanng.com ↗ for more information on price quotations, drivers, technical specifications, manuals, and documentation

Contact us: (888) 88-SOURCE | sales@artisanng.com | www.artisanng.com