ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *BrandTech / Brinkmann Transferpette-8* at our website: *Click HERE*

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Operating Manual

Please read carefully before use, and follow all operating and safety instructions!

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Contents

	Page
Safety Instructions	4
Purpose	4
Operating Limitations	5
Operating Elements	6
Handling	6
Pipetting	8
Maintenance	11
Adjustment	12
Checking the volume	13
Components	16
Cleaning	17
Repairs	20
Autoclaving	28
Trouble Shooting	30
Specifications	31
Ordering Data	33
Spare Parts and Accessories	25
Return for Repair	37
Disposal	35
Warranty	38

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Read carefully before use!

Safety Instructions

1. Observe general hazard prevention instructions and safety regulations; e.g., wear protective clothing, eye protection and gloves.
2. Observe the reagent manufacturer‘s information.
3. Read this Operating Manual carefully.
4. Always use the instrument in such a way that neither the user nor any other persons are endangered.
5. Only use manufacturer‘s original parts.
6. Operate only with tips attached.
7. Use original PLASTIBRAND® pipette tips, or make sure that the tips to be used are compatible.

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Purpose

The multichannel pipette is an air-displacement system for pipetting aqueous solutions with medium density and viscosity into microtiter plates.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Operating Limitations

The user has to ensure the compatibility of the instrument with the intended application.

This instrument cannot be used:

- for liquids incompatible with polypropylene
- for liquids of a very high vapour pressure
- in Reverse Mode (ISO 8655/2)

Viscous and wetting liquids may impair the accuracy of the volume. The same applies to liquids of a temperature deviating from the ambient temperature by more than ± 5 °C/41 °F (permissible operating range: 15-40 °C/59-104 °F).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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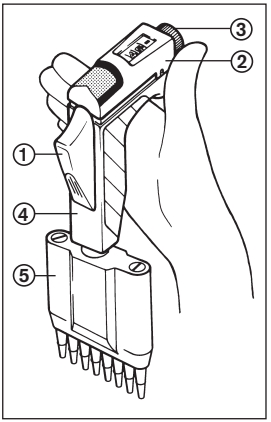
Handling

Use proper pipette tips

Pipette performance can only be guaranteed if high-quality pipette tips are used. We consequently recommend the use of PLASTIBRAND® pipette tips. Pipette tips are disposables.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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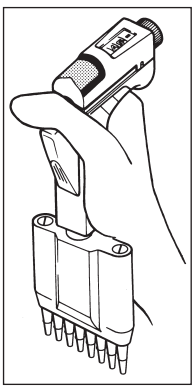
Adjusting the manifold

The **manifold** can be rotated counterclockwise to the most convenient pipetting position.

Counterclockwise rotation
for ease of handling

Note:
The manifold must be completely mounted to the handle for accurate pipetting. During the counterclockwise rotation, a clicking noise must be heard.
Turning in the opposite direction will loosen the manifold.

Hand posture while pipetting

Securely attach tips to nose cone assembly. Put your thumb **across** the pipetting key – unlike with conventional pipettes.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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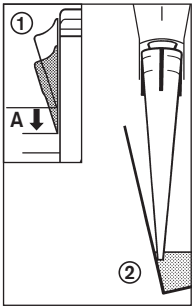
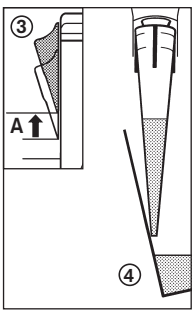
Pipetting

Aspirating the sample

Rinse the tips once with the sample liquid (at normal pipetting speed).

1. Press the lateral pipetting key to **the first stop (A).**
2. Immerse the pipette tips 2-3 millimeters into the sample.

Note:
Leave the tips immersed in the liquid for about another second to avoid the intake of air.

3. Let the pipetting key slide back **slowly**.
4. Lightly wipe the tips against the wall of the vessel.

Note:
Don’t lay the instrument horizontally when the tips are filled, or liquid may flow inside the instrument.


Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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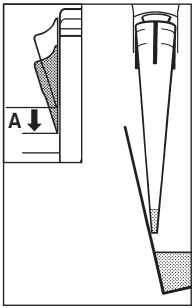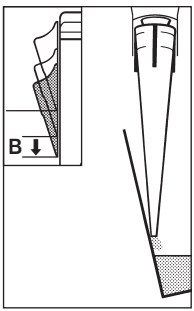
Discharging the sample

Lightly place the pipette tips against the wall of the vessel.

– Press the pipetting key down to the **stop (A)** and hold it down.

Note:
For serums and liquids of high viscosity or low surface tension, observe adequate waiting time to improve accuracy.

– The blow-out stroke empties the tips completely: Press down to the **stop (B)**.

Note:
To improve accuracy when working at low volumes: rinse the tips with the reagent contained in the vessel

– Lightly wipe the pipette tips against the wall of the vessel.
– Let the pipetting key slide back.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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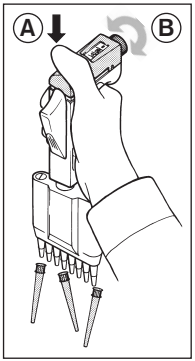
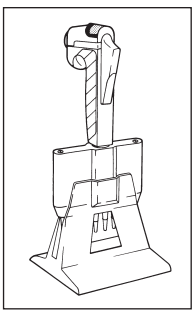
Ejecting the tips (A)

Press the ejector cap down firmly in the region of the color code (see arrow).

Tips that are jammed on too tightly may need to be pulled off separately by hand.

Volume setting (B)

Set the desired volume by turning the volume setting knob (➤ “Volume setting knob”, page 6).

After use

Always store the instrument upright on the **stand** supplied.

Note:
To reduce strain on the sealing rings (➤ “V-rings”, page 24), always store the instrument without the tips attached.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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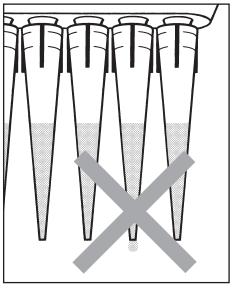
Maintenance

Check for tightness

At regular intervals, check the pipetting system – consisting of Transferpette®-8/ -12 and pipette tips – for tightness.

After drawing in the sample:
Hold down the instrument vertically for about 10 seconds. If a drop forms at one of the tips,

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Lubricating the plunger

After autoclaving the manifold several times, it is necessary to lubricate the pistons to maintain smooth operation:
Unscrew nose cones, push pistons downward, and lubricate (➤ “Repairs”, page 20, 21).


To lubricate the pistons, use only the silicone grease supplied with the instrument. The silicone grease is also available separately as accessory (➤ “Ordering Data”, page 33).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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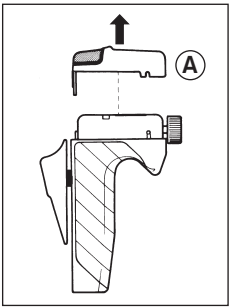
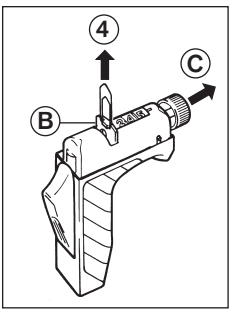
Adjustment

The instrument is permanently adjusted for aqueous solutions. It can be adjusted for solutions of other density and viscosity.

1. Check the volume, determine actual value (➤ page 15).
2. Release lateral closures by turning into position ◍.
3. Remove ejector cap (A).
4. Use a paper clip to extract the red distance plate (B) from the top.
5. Pull out volume setting knob (C) to the stop in the direction of its axis (approx. 2 mm).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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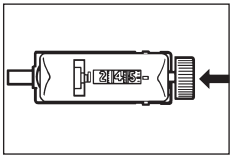6. Set to the previously measured actual value. Push back volume setting knob to the stop (approx. 2 mm), without turning it.
7. Replace the red distance plate across the axis.
8. Mount ejector cap (A).
9. Secure ejector cap by turning the lateral closures into position ⊜.

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Checking the volume

The instrument should be calibrated at least once a year. A more frequent interval may be dictated by individual lab rules, applications and conditions.
The gravimetric volume test of the pipette requires the following steps (page 14):
(Testing procedures are described e.g., in ISO 8655, Part 6, ➤ “Testing Instruction”, page 36.).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Adjust volume to nominal capacity.

Conditioning the pipette:

To condition the pipette before testing, mount new pipette tips and pipette five times with testing liquid (H_2O dist.). Discard the pipette tips.

Carrying out the test (observe instructions by the balance manufacturer):

- Mount new pipette tips and pre-rinse once with testing liquid.
- Take up testing liquid and pipette into the weighing vessel.
- Weigh the pipetted quantity with an analytical balance.
- Calculate the volume, taking the temperature into account.
- The number of tests per channel should be selected according to the required accuracy. We recommend 5-10 tests per channel in 3 volume ranges (100 %, 50 %, 10 %).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Variationskoeffizient

VK% $= \frac{100\ s}{\bar{V}}$

Calculation of accuracy (A%) and coefficient of variation (CV%):

Accuracy (A%) and coefficient (CV%) are calculated according to the equations used in statistical quality control.

Calculations (for the nominal volume)

Mean value $\bar{x} = \frac{\sum x_i}{n}$ x_i = results of weighings, n = number of weighings

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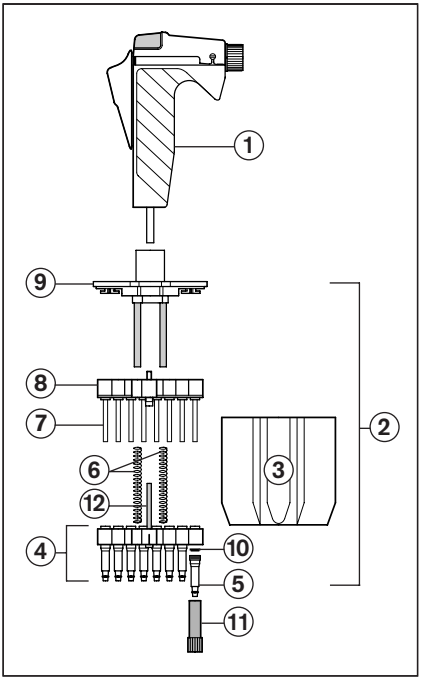
Components

1. Handle
2. Manifold
3. Manifold housing
4. Nose cone assembly
5. Individual nose cone
6. Piston spring
7. Piston
8. Piston support bar
9. Manifold housing cover
10. Seal
11. Mounting tool
12. Central guide rod

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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
Cleaning

Note: The Transferpette®-8/-12 requires regular maintenance, depending on frequency of use and pipetted media.

Dismantling the instrument

Note: Work over a tray to collect any part that may fall out during dismantling.

- Eject the tips.
- Set the instrument to nominal volume.
- Carefully remove the manifold (②) from the handle (①).
- Turn both closures of the manifold housing cover (⑨) with a coin (rotation through 90°) and slide off the housing.

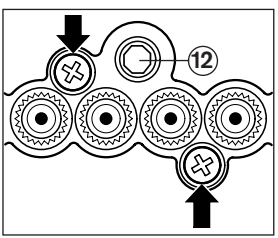- **Remove both outer screws (Phillips). Do not loosen the central guide rod** (⑫)**!**
- Remove the nose cone assembly (④), piston springs (⑥), and piston support bar (⑧).
- Unscrew single nose cones with seals from the nose cone support bar (➤ "Repairs", page 20).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Cleaning the components

- Clean single nose cones (⑤), pistons (⑦), piston support bar (⑧) and piston springs (⑥) **(these components only)** with an appropriate cleaning agent (e.g., water, alcohol. Do not use organic solvents!).
- Let these parts dry completely. Residual moisture in the nose cone assembly may result in a loss of accuracy.
- Lubricate the pistons lightly (only with the supplied silicone grease).

Important:
To lubricate the central guide rod, use only the recommended fluorstatic grease (➤ page 33)!

Note:
To protect against contamination, filter tips can also be used.

Assembling the instrument

- Reassemble the nose cone assembly without cones (④), piston springs (⑥), piston support bar (⑧) and manifold housing cover (⑨). Tighten all screws.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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- Set volume to minimum.
- Push the seals on the pistons with the flat side facing upward (direction of the handle) (➤ Fig. 4, page 21); use the mounting tool to mount and tighten the single nose cones clockwise.
- Slide on the manifold housing (③). Turn the closures into place by rotation through 90°.
- Push the handle and the manifold together and turn until they click into place.

Important:
Do not press on the ejector cap during this procedure.

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After cleaning

- Check if pipetting and ejector functions are correct. If this is not the case, dismantle the instrument and carefully reassemble it. If necessary, lubricate the pistons (➤ “Maintenance”, page 11).
- Check for proper fit of tips.
- Check for correct volume delivery (➤ page 13).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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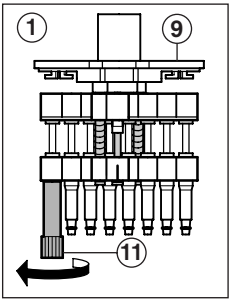
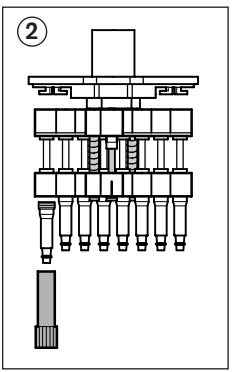
Repairs

➤ **“Accessories and spare parts”, pages 33-35.**

Replacing of single nose cones and seals

Removing a single nose cone


- Set the instrument to nominal volume.
- Turn both closures of the manifold housing cover (⑨) with a coin (rotation through 90°) and slide off the housing.
- Push the mounting tool (⑪) on a single nose cone (Fig. 1).
- Turn counterclockwise to loosen and remove the nose cone (Fig. 2).


Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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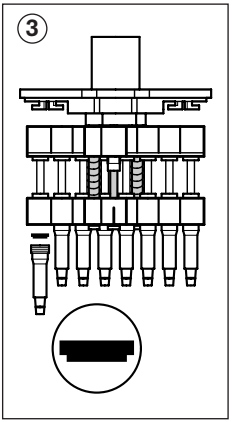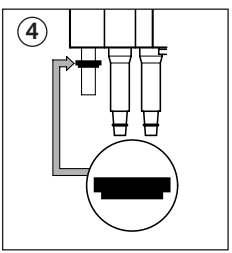Note:
The seal will either remain inside the nose cone, or will stay on the piston after the nose cone is removed.
By pressing the pipetting key it can easily be removed and cleaned (Fig. 3).

- Set volume to minimum.
- Push down the pipetting key completely.
- If required, lubricate the piston lightly with the **supplied** silicone grease.
- Push the seal on the piston with its flat side facing upward (direction of the handle) (Fig. 4).

 Important: The 300 µl instrument has a pressure ring above the seal.
- Use the mounting tool to mount and tighten the new nose cone clockwise.
- Reassemble the manifold (➤ page 18).
- Check instrument for tightness (➤ “Check for tightness”, page 11).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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

Replacing the pistons

- Unscrew the manifold from the handle.
- Turn both closures of the manifold housing cover (⑨) with a coin (rotation through 90°) and slide off the housing.
- **Remove both outer screws (Phillips)** which hold the nose cone assembly (④) to the metal support rods of the manifold housing cover (⑨).
 Do not loosen the central guide rod (➤ page 16)!

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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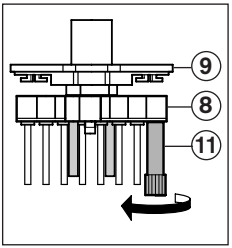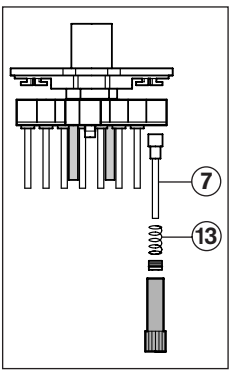- Push the mounting tool (⑪ on a piston (⑦).
- Turn counterclockwise to loosen and remove the piston (⑦) along with the spring (⑬).
- Insert a new piston.
- Use the mounting tool to mount and tighten the new piston.
- Reassemble the instrument (➤ page 18) and check for tightness (➤ page 11).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Replacement of nose cone V-rings (O-rings)

Required materials:

1 V-Ring removal tool (square)

1 V-Ring mounting tool (cone shaped)*


Replacement V-rings (O-rings)

* The mounting tool is only necessary for replacing the V-rings on instruments of the volume ranges 2.5-25 µl, 5-50 µl, 10-100 µl, 20-200 µl and 30-300 µl.


Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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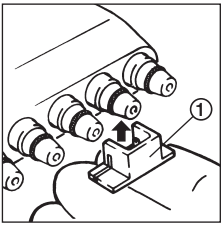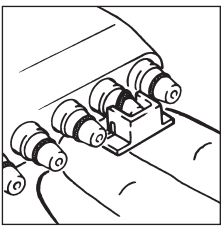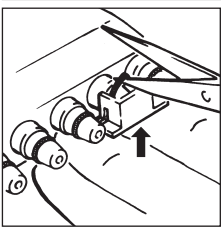
To remove an old V-ring (O-ring):

1. Position the square removal tool ① under one of the nose cones.


2. Move the tool up against the nose as far as possible to push the V-ring (O-ring) up, away from the nose cone.

3. If necessary, cut the dislodged portion of the V-ring (O-ring) with a scissor and discard. Do not damage the nose cone.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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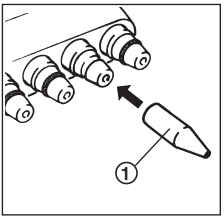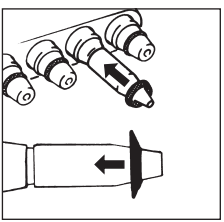
Mounting a new V-ring (O-ring):


4. Only for instruments of size 2.5-25 µl, 5-50 µl, 10-100 µl, 20-200 µl and 30-300 µl: Slide the cone-shaped mounting tool ① over the nose cone.
 For instruments of size 0.5-10 µl and 2-20 µl, no mounting tool is required; the O-ring can be placed directly onto the nose cone.

Note:
Moisten mounting tool, nose cone and V-ring (O-ring) with dist. H_2O in order to facilitate mounting.

5. Place the new V-ring (O-ring) onto the mounting tool. The widest side of the V-ring should face towards nose cone manifold. Push the V-ring (O-ring) over the mounting tool into the groove on the nose cone.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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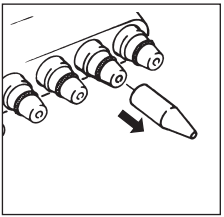6. Remove the mounting tool. Check correct seat of the V-ring (O-ring). If twisted, turn it several times around the nose cone to restore the V-ring (O-ring) to its correct shape.

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After any repair:

- Check if pipetting and ejector functions are correct.
- Check for proper fit of the tips.
- Check for correct volume delivery (➤ page 13).

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After every change of the complete manifold:

Check for correct volume delivery (➤ page 13). Readjust instrument if necessary.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Autoclaving

Attention:
Only the manifold may be autoclaved. Eject the tips before autoclaving! For disinfection the handle you may use commercially available disinfecting solutions as long as they do not attack PP.

- Remove the manifold from the handle.
- For sterile work or for decontamination, the complete manifold is autoclavable at 121 °C/250 °F, 2 bar absolute, for 20 minutes exposure time, according to DIN. It is the user's responsibility to ensure effective sterilization and correct autoclave temperature.

After autoclaving

- Allow the manifold to cool down completely.
- Push the manifold and the handle together and turn until they click into place. **Important:** Do not press the ejector cap during this procedure.
- If necessary, tighten the nose cones with a sterile pipette tip of proper size.
- Check for smooth function of pipetting and ejector keys.
- Check for proper fit of pipette tips.
- Check for correct volume delivery (➤ page 13).

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Trouble Shooting

Problem	Possible Cause	Corrective Action
Pipette tip leaks or volume too low	- Tip does not sit properly	- Press tip on firmly. Use only suitable high-quality tips, page 6.
	- Damaged V-ring (O-ring) on single nose cone	- Replace V-ring (O-ring) on nose cones. “Repairs”, pages 24-26.
	- Piston damaged	- Replace piston. “Repairs”, pages 22, 23.
	- Seal dirty, damaged	- “Repairs”, pages 20, 21.
	- Instrument is dirty	- Clean instrument. “Cleaning”, page 17.
	- Non-aqueous solutions pipetted	- “Adjustment”, page 12.
Poor sample aspiration	- Single nose cone clogged	- Clean nose cone. Remove cones only; do not disassemble instrument! “Cleaning” page 17.
Volume too large	- Pipetting key pressed too far down into the blow-out stroke before aspirating liquid.	- Operate pipette properly. “Operating the pipette”, page 8.
Pipetting key moving stiffly	- Key movement reduced	- Lubricate piston and centrale guide rod with the appropriate grease, pages 11, 15.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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* Test values determined at the same temperature of instrument, liquid (H_2O dist.) and ambience (20 °C) and smooth, jerk-free operation.

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When pipetting liquids with different viscosity and/or flow properties than H_2O dist., deviations from the values stated can occur.

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Accuracy and coefficient of variation refer to each single channel.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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PLASTIBRAND® Tip-Box, PC, empty Cat. No. 7023 00

Reservoirs for reagent aspiration (set of 10), autoclavable at 121 °C Cat. No. 7034 59

Stand for one instrument Cat. No. 7034 40

Silicone grease for piston (tube, approx. 1 g) Cat. No. 7036 77

Fluorstatic grease for central guide rod (tube, approx. 1 g) Cat. No. 7036 78

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Individual nose cone (2 pieces, incl. seals and mounting tool)

For capacity, µl	Cat. No.
0.5 - 10	7033 31
2 - 20	7033 32
2.5 - 25	7033 33
5 - 50	7033 34
10 - 100	7033 35
20 - 200	7033 36
30 - 300	7033 37

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V-Rings (Set of 24), incl. removing and mounting device (except 10 µl and 20 µl capacity: no mounting tool needed)

For capacity, µl	Cat. No.
10, 20 (O-ring)	7034 13
25, 50, 100, 200, 300	7034 91

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Note:

ISO 9001 and GLP-guidelines require regular examinations of your volumetric instruments.We recommend checking the volume ever 3-12 months. The interval depends on the specific requirements on the instrument. For instruments frequently used or in use with aggressive media, the interval should be shorter.

The detailed testing instruction can be downloaded from www.brandtech.com.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Return for Repair

Attention:
Transporting of hazardous materials without a permit is a violation of federal law.

BrandTech Scientific, Inc. will not accept instruments that are not appropriately cleaned and decontaminated.

➤ Therefore, contact BrandTech Scientific, Inc. and obtain return authorization **before** sending your instrument for service.

Return the instrument, with the Return Authorization-Number prominently displayed on the outside of the package to the address provided with the Return Authorization Number. Include an exact description of the type of malfunction and the media used.

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Disposal

For the disposal of instruments and tips, please observe the disposal regulations valid in your area.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

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Warranty

We shall not be liable for the consequences of improper handling, use, servicing or unauthorized repairs of the instrument or the consequences of normal wear and tear especially of wearing parts such as pistons, seals, valves and the breakage of glass as well as the failure to follow the instructions of the operating manual. We are not liable for damage resulting from any actions not described in the operating manual or if non-original parts have been used.

Subject to technical modification without notice. Errors excepted.

Artisan Scientific - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisan-scientific.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP