ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Measurement Computing / Data Translation DT2802-4* at our website: *Click HERE*

DATA TRANSLAT

UM-1

DT777
User Manual

Screw Terminal Panel
for Chromatography Boards

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

First Edition
June, 1991

Copyright © 1991 by Data Translation, Inc.

All rights reserved. No part of this publication may b
reproduced, stored in a retrieval system, or transmit
form by any means, electronic, mechanical, by phot
recording, or otherwise, without the prior written per
Data Translation, Inc.

Information furnished by Data Translation, Inc. is be
be accurate and reliable; however, no responsibility
assumed by Data Translation, Inc. for its use; nor fo
infringements of patents or other rights of third parti
may result from its use. No license is granted by im
or otherwise under any patent rights of Data Transla

Data Translation, Inc.
100 Locke Drive
Marlboro, MA 01752-1192
(508) 481-3700
Easylink MBX 62825999
Telex 951646
(DATATRANS MARO)
Fax (508) 481-8620

Data Translation® is a registered trademark of Data
Translation, Inc.

IBM® and Personal Computer AT® are registered tra
of International Business Machines Corporation.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Getting Starte

Using the DT7

Attenuation Resisto

Connector Pin Assignmen

Specificatio

Troubleshooti

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Table of Contents

How to Use This Manual v
Where to Get Help v

Chapter 1: Getting Started 1
About the DT777 1
Required Equipment 1
Accessory Products 1
Unpacking 2
Connecting to a Chromatography Board 3

Chapter 2: Using the DT777 5
Screw Terminal Assignments 5
Detector Connections 6
Detector Input Attenuation 6
External Event Device Connections 7
Input Connections 7
Output Switch Connections 8
Other Connections 9

Reference 1: Attenuation Resistors 11
About the Resistors 11

Reference 2: Connector Pin Assignments 13
J1 Connector 13

Reference 3: Specifications 15
Electrical Specifications 15
Physical Specifications 16

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reference 4: Troubleshooting . . . . . .
General Checklist . . . . . .
Service and Support . . . . . .
If Your Board Needs Factory Service . . . . . .

Index

Warranty Policy

Reader's Comment Form

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Figures

Figure 1: Connecting to a Chromatography Board . . . . . .
Figure 2: DT777 Screw Terminal Assignments 1
Figure 3: Preferred Detector Connections 1
Figure 4: Alternate Detector Connections 1
Figure 5: Resistor Network Location 1
Figure 6: Sensing Vref 1
Figure 7: Supplying Vref 1

Tables

Table 1: Attenuation Resistors (-0.500V to 5.000V) 2
Table 2: Attenuation Resistors (-1.000V to 10.00V) 2
Table 3: J1 Pin Assignments 2

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

How to Use This Manual

This manual describes Data Translation's DT777 screw termin panel. This panel facilitates user connections to Data Translation's chromatography boards.

Listed below is an overview of the information presented in thi manual.

Chapter 1, "**Getting Started**," introduces you to the major features of the DT777 screw terminal panel and explains how connect it to chromatography boards.

Chapter 2, "**Using the DT777**," describes how to connect the DT777 to the various subsystems on chromatography boards.

Reference 1, "**Attenuation Resistors**," lists the resistor netwo for attenuating detector signals.

Reference 2, "**Connector Pin Assignments**," lists the pin assignments for the DT777's J1 connector.

Reference 3, "**Specifications**," lists the electrical and physical specifications for the DT777.

Reference 4, "**Troubleshooting**," explains what to do if you encounter difficulty using your screw terminal panel.

Where to Get Help

If you have difficulty installing or using your screw terminal par see Reference 4, "Troubleshooting."

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Chapter 1

Getting Started

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Chapter 1, "**Getting Started**," describes the DT777 scr terminal panel and explains how to connect it to a Data Translation chromatography board. This chapter is org follows:*

- *About the DT777*
- *Unpacking*
- *Connecting to a Chromatography Board*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

About the DT777

The DT777 screw terminal panel is an accessory product whic permits the user connections to a chromatography board to be made on screw terminals. The function of each screw terminal on the panel is indicated by a silk-screened label.

The DT777 can accommodate signal attenuation on the detect inputs, effectively extending the detector input range of the chromatography board.

The DT777 provides transient overvoltage protection on the inp switch detectors and output switches, protecting both the chromatography board and the user's interface from transient overvoltages.

The DT777 screw terminal panel has an optional rugged fully-shielded plastic enclosure. Data Translation recommends you use this enclosure to shield the screw terminal panel from flow and electrical noise.

Required Equipment

The following equipment is required to install and use the DT77 screw terminal panel:

- Data Translation chromatography board
- EP777 cable assembly (available from Data Translation)
- Flat blade screwdriver

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Accessory Products

The EP777 cable assembly (available from Data Transla consists of individually shielded twisted pairs housed in flexible cable. This cable assembly has a D-shell conne each end. One end plugs into connector J1 on the chromatography board. The other end plugs into the J1 connector on the DT777.

Data Translation also has a resistor pack which provides additional resistors required to attenuate the detector inp signals. Consult the factory for ordering information for resistor pack.

The EP246 shielded enclosure (available from Data Tra helps preserve the integrity of analog signals by shieldin DT777 from electrical noise (RFI and EMI) and from am airflow. Data Translation recommends using the EP246 enclosure for best performance.

Note *If you do not use the EP246, protect your connections to DT2802 following one of these two methods:*

(1) Enclose the DT777 in a metal box and connect the box to SGND on the DT777. The metal box will pro shielding from air flow and electrical noise.

(2) Enclose the DT777 in a cardboard box or wrap it in pack to shield it from air flow. This method does no the connections from electrical noise.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Unpacking

The screw terminal panel package includes the following items

- DT777 screw terminal panel
- *DT777 User Manual*

Note *Save the original packing material in the unlikely event that you panel requires servicing in the future.*

If an item is missing or damaged, call Data Translation's **Customer Service Department** at **(508) 481-3700**. They will guide you through the appropriate steps for replacing missing o damaged items.

Note *If you are located outside the USA or Canada, call your local distributor. The name and telephone number of your nearest distributor are listed in your Data Translation catalog.*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Connecting to a Chromatogra Board

The following steps explain how to connect the DT777 s terminal panel to a chromatography board.

Note *The chromatography board must be installed in your co before you begin. Refer to the user manual supplied wi board for installation instructions.*

1. Turn off the computer and all accessories.

2. Position the computer so that you can easily access chromatography board's rear panel connector (J1).

3. Attach one end of the cable to the J1 connector on t terminal panel.

4. Attach the other end of the cable to the rear panel c (J1) of the chromatography board.

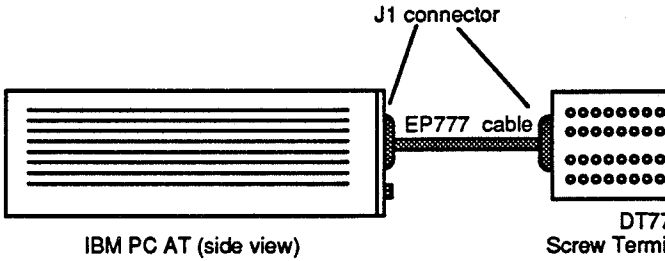


Figure 1: Connecting to a Chromatography Bo

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Chapter 2

Using the DT777

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Chapter 2, "**Using the DT777**," explains how to connect*
to the various subsystems on the chromatography board
screw terminal panel. This chapter is organized as follo

- *Screw Terminal Assignments*
- *Detector Connections*
- *External Event Device Connections*
- *Other Connections*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Screw Terminal Assignments

Connections for the chromatography board's data acquisition subsystem and external event switches are made by connecting the board's J1 connector pins to the screw terminals. Each J1 signal name is indicated directly on the panel next to the corresponding screw terminal.

The DT777 screw terminal panel accommodates the differential analog input signal from a detector. Blank circuit pads provided onboard the panel allow you to install a resistor network to attenuate the detector signal. Tables 1 and 2 on page 21 list the resistors necessary for attenuating a 5V and 10V detector input signal.

Figure 2 illustrates the screw terminal and resistor locations for the DT777. Refer to this figure while performing the connection procedures listed in the sections that follow.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

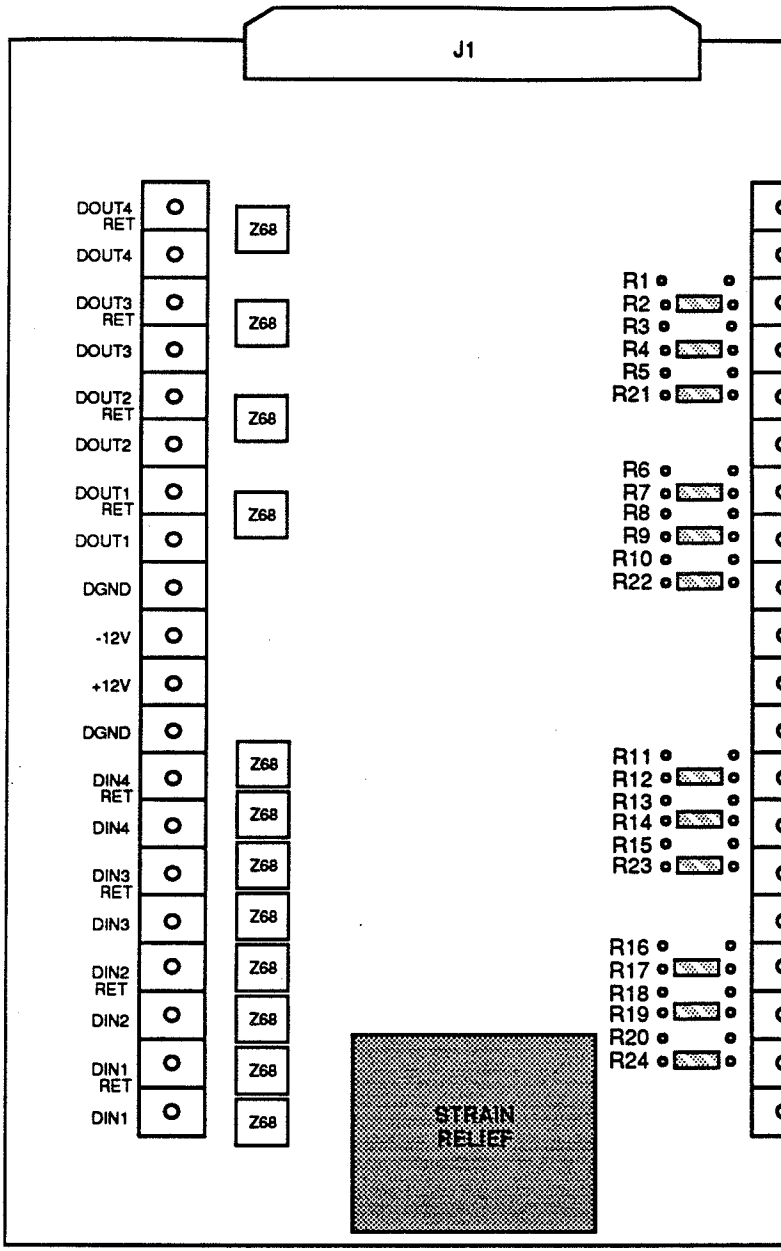


Figure 2: DT777 Screw Terminal Assignments

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Detector Connections

The DT777 screw terminal panel accommodates connections
up to four detectors to the chromatography board. Each integrator on the board supports two detector input channels. Integrator A supports channel 1 and channel 2; integrator B supports channel 3 and channel 4.

The detector connections are differential. Each channel has a separate return. Differential inputs are relatively insensitive to electrical noise since the integrator sees only the difference in voltage between the high and return sides of each channel. TI DT777 provides a 100kΩ resistor from each channel return to AGND. This resistor provides a return path for the bias curren input.

The chromatography board is capable of converting analog inp voltages ranging from –0.207V to +2.000V. The detector you connecting need not have a voltage output that ranges to +2.000V. If your detector's output voltage exceeds +2V, however, you must reduce the voltage using a resistor divider circuit. The DT777 screw terminal panel can accommodate 5V and 10V detector outputs. Refer to "Detector Input Attenuation on page 13 for more information.

In order to reduce the noise pickup in the cabling, use a shield twisted pair cable for each detector input connection. The two wires that make up the pair are used for the signal HI and sign LOW connections. The shield protects the signal wires from outside electrical influences. For further information about making detector connections, see Chapter 4, "Detector Connections," in the *DT2802 User Manual.*

To connect the detector to the DT777 (no signal attenuation):

1. Connect the screw terminal panel to the chromatography board.

2. Connect the voltage output from the detector to the desired channel terminal on the screw terminal panel.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

3. Connect the low side of the detector to the correspo terminal for that channel.

4. Connect the shield to the SGND terminal on the scr panel. *Do NOT connect the other end of the shield detector! (See Figure 3.)*

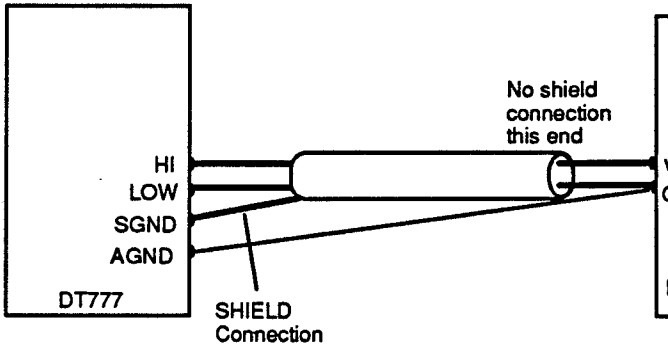


Figure 3: Preferred Detector Connection

Note *If you cannot connect the shield at the DT777, leave tha unconnected and connect the other end of the shield to detector as shown in Figure 4.*

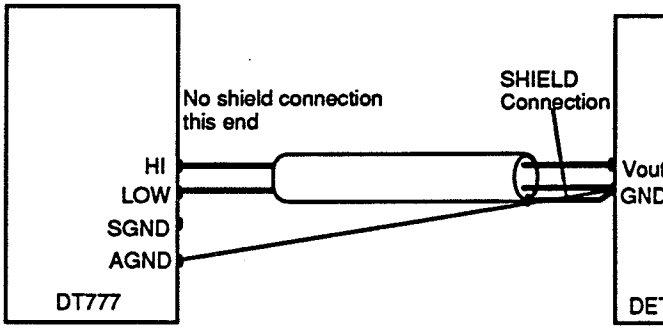


Figure 4: Alternate Detector Connection

When the shield is connected, it carries no electrica This means that the entire shield is at the same volt making it an equipotential or electrostatic shield.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

5. Connect a wire between the AGND terminal on the screw terminal panel and the detector GND terminal.

6. Connect each detector input to the DT777 in this manner. Shielding the input channels from each other minimizes crosstalk, particularly when one input channel is near full-scale and the other is near zero.

Detector Input Attenuation

The DT777 provides circuit pads for installing resistors for attenuating the detector input signal. Each detector input channel can be configured with resistors to increase the –0.20 to +2.000V input range of the chromatography board to –0.500 to +5.000V or –1.000V to +10.00V.

Each detector input channel requires a network of three individual precision, low temperature coefficient resistors. Zero ohm resistors are factory-installed at the appropriate locations provide a factory configuration of no attenuation. When you attenuate the detector signal, all four detector input channels must be configured for the same level of attenuation.

To install the resistor network:

1. Cut the 0Ω resistors from the DT777. See Tables 1 and 2 page 21 for the resistors to remove.

2. Install the attenuation resistors for each input channel at th appropriate locations. See Tables 1 and 2 on page 21 for the resistors to use and the resistor locations for each inpu channel.

3. Make each detector input connection as described in the previous section.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

External Event Device Connections

The DT777 provides screw terminal connections for the switch detectors and output switches of the chromatogr board. The chromatography board protects the input sv detectors to ±50V. The DT777 provides transient (surg protection with a metal oxide varistor (MOV) connected parallel at each input switch detector terminal. The MO protection devices enhance the input protection already by the chromatography board.

The chromatography board's output switches are optica isolated switch closures having a voltage rating of 24VD There is an MOV connected in parallel at each output sv terminal for transient (surge) voltage protection.

Input Connections

Many external devices in a chromatographic system out signal when an event occurs at that device. You can co this signal output to one of the input switch detectors on chromatography board.

To connect an input signal to the DT777:

1. Connect the high side of the input signal to the desi switch terminal on the DT777.

2. Connect the low side of the input signal to the corre input switch low terminal on the DT777.

Note *Use a shielded cable for minimum noise on the input sig the detector.*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Output Switch Connections

Many external devices in a chromatographic system accept a switch closure signal from another device. These devices use this signal as a trigger to perform a particular activity during the separation.

To connect an output switch to an external device:

1. Connect one lead of the cable between the desired switch low terminal on the DT777 and a GND terminal on the device.

2. Connect the other leads on the cable to the corresponding switch high terminal on the DT777 and the input terminal o the device.

Note *Use a shielded cable for minimum noise on the input signal to the detector.*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Other Connections

±12V Supplies The chromatography board provides the user access to supplies of the PC/AT bus. These lines are brought out terminal connections on the DT777. Returns for the sup provided through DGND. The amount of power available from each supply.

Vref I/O The chromatography board provides access to the inter voltage reference of integrator A. You may use the Vref connection to input a stable, precision voltage reference the integrator's own reference voltage as an output.

To sense integrator A's internal reference, you must buf Vref I/O and Vref RETURN signals with low drift, low bia op amps. (Data Translation recommends OP-27 op am must make a differential amplifier with a gain of 1.00 (se 5) where the positive input is the buffered Vref I/O and t negative input is the Vref RETURN. The output of the d amplifier is then an "image" of the +5.000V internal refe with respect to your ground system.

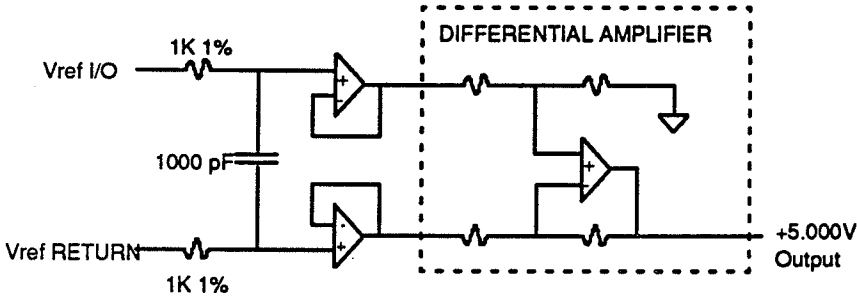


Figure 5: Sensing Vref

To supply a voltage reference to the chromatography bo must provide a voltage source of +5.000V ±1%. This vo source must be an op amp or other low impedance outp respect to the Vref RETURN signal. See Figure 6 on th following page.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

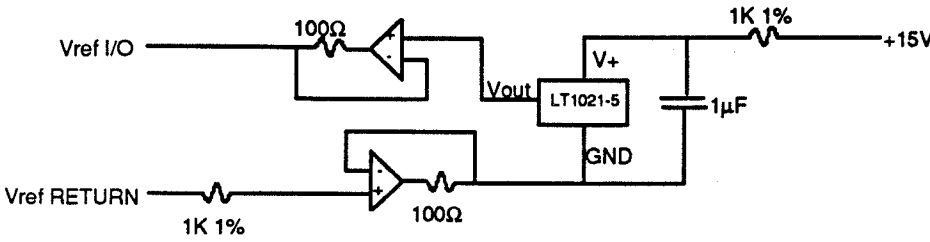


Figure 6: Supplying Vref

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reference 1

Attenuation Resistors

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Reference 1, "**Attenuation Resistors,**" lists the resistor for attenuating the detector input signal.*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

About the Resistors

To ensure the accuracy of the chromatography board's converter, the attenuation resistors should be precision (±0.01 or better tolerance), high stability (tempco of 5 ppm/°C matche MAR5 or wirewound resistors of at least 1/10 watt power dissipation. The resistors provided in the resistor pack access meet these requirements.

Each detector input channel requires a network of three resistors.

Tables 1 and 2 list the attenuation resistors to install and remo on the DT777 for each detector input channel.

Table 1: Attenuation Resistors (-0.500V to 5.000V)

	Install			Remove	
Channel	Resistor 1	Resistor 2	Resistor 3	Resistor 1	Resistor 2
1	R1: 300kΩ	R3: 400kΩ	R5: 300kΩ	R2: 0Ω	R4: 0Ω
2	R6: 300kΩ	R8: 400kΩ	R10: 300kΩ	R7: 0Ω	R9: 0Ω
3	R11: 300kΩ	R13: 400kΩ	R15: 300kΩ	R12: 0Ω	R14: 0Ω
4	R16: 300kΩ	R18: 400kΩ	R20: 300kΩ	R17: 0Ω	R19: 0Ω

Table 2: Attenuation Resistors (-1.000V to 10.00V)

	Install			Remove	
Channel	Resistor 1	Resistor 2	Resistor 3	Resistor 1	Resistor 2
1	R1: 400kΩ	R3: 200kΩ	R5: 400kΩ	R2: 0Ω	R4: 0Ω
2	R6: 400kΩ	R8: 200kΩ	R10: 400kΩ	R7: 0Ω	R9: 0Ω
3	R11: 400kΩ	R13: 200kΩ	R15: 400kΩ	R12: 0Ω	R14: 0Ω
4	R16: 400kΩ	R18: 200kΩ	R20: 400kΩ	R17: 0Ω	R19: 0Ω

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reference 2

Connector Pin Assignments

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Reference 2, "**Connector Pin Assignments**," lists the pi*
assignments for the DT777's J1 connector.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

J1 Connector

Table 3 lists the pin assignments for the DT777's 37-pin J1 connector.

Table 3: J1 Pin Assignments

Description	Pin		Description
Channel 1 High	19	37	Channel 1 Low
Analog Ground	18	36	Analog Ground
Channel 2 High	17	35	Channel 2 Low
Voltage Reference High	16	34	Voltage Reference Low
Channel 3 High	15	33	Channel 3 Low
Analog Ground	14	32	Analog Ground
Channel 4 High	13	31	Channel 4 Low
Analog Ground	12	30	Analog Ground
Digital Ground	11	29	Digital Ground
+12V[1]	10	28	−12V[1]
Input Switch 1 High	9	27	Input Switch 1 Low
Input Switch 2 High	8	26	Input Switch 2 Low
Input Switch 3 High	7	25	Input Switch 3 Low
Input Switch 4 High	6	24	Input Switch 4 Low
Switch 1 High	5	23	Switch 1 Low
Switch 2 High	4	22	Switch 2 Low
Switch 3 High	3	21	Switch 3 Low
Switch 4 High	2	20	Switch 4 Low
Digital Ground	1		

[1] *No more than 10mA should be drawn from the ±12V terminals (pins 10 ar 28). The return for the ±12V is Digital Ground.*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reference 3

Specifications

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Reference 3, "**Specifications**," lists the electrical and ph*
specifications for the DT777. All specifications are typic
+25°C and rated voltage unless otherwise specified. Th
section is organized as follows:

- *Electrical Specifications*
- *Physical Specifications*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Electrical Specifications

Analog Input Attenuation

Note *Resistors must be supplied by the user. Resistor kits are available from Data Translation.*

Resistor Value: -0.500V to 5.000V

Channel 1	R1 and R5 = 300k, R3 = 400k
Channel 2	R6 and R10 = 300k, R8 = 400k
Channel 3	R11 and R15 = 300k, R13 = 400k
Channel 4	R16 and R20 = 300k, R18 = 400k

Resistor Value: -1.000V to 10.00V

Channel 1	R1 and R5 = 400k, R3 = 200k
Channel 2	R6 and R10 = 400k, R8 = 200k
Channel 3	R11 and R15 = 400k, R13 = 200k
Channel 4	R16 and R20 = 400k, R18 = 200k

Resistor tolerance ±0.01% or better

Resistor power dissipation
1/10 watt

Resistor temperature coefficient
±5 ppm/°C TCR (matched)

Resistor type MAR5 or precision wirewound

External Event I/O Protection

Protection device Metal oxide varistor

Rated voltage 45V DC

Varistor voltage 56V ±10%

Rated peak single pulse transient current
500A

Rated single pulse transient energy
5.6 joules

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Rated transient average power dissipation
0.4 watts

Maximum clamping voltage
70V @ 10A

Typical capacitance 1800pF

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Physical Specifications

PCB dimensions	0.875"H x 6.00"W x 4.59.00"L (2.22 cm x 15.24 cm x 11.43 cm)
Weight	16 ounces
Screw terminal	40 barrier strip connections: 4 sections with 10 connections each
Barrier strip screws	
Size	6-32 screws, cadmium plated steel
Terminal material	Bright acid, tin plated brass
Barrier strip material	Thermoplastic polyamide, UL grade, recognized flammability rating 94V-0
Acceptable wire size	14 to 22 AWG

EP246 Enclosure

Enclosure dimensions	2.95"H x 6.50"W x 5.00"L (7.49 cm x 16.51 cm x 12.71 cm)
Enclosure material	Aluminum

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reference 4

Troubleshooting

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

*Reference 4, "**Troubleshooting**," describes what to do if have difficulty using your screw terminal panel. This sect organized as follows:*

- *General Checklist*
- *Service and Support*

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

General Checklist

Should you experience difficulty using your screw terminal pa
the following sections describe some possible problems and
suggest steps you can take to correct them.

If you are experiencing a problem not described here, or have
followed the steps suggested and are still having problems,
please follow the instructions provided in "Service and Suppo

Go through the following checklist when troubleshooting to he
you isolate your problem.

1. Make sure that all connections are made properly.

   *Refer to Chapter 2, "Using the Screw Terminal Panel," fo
   the appropriate connection scheme.*

2. Verify that all cables associated with the board are firmly
   connected.

3. Make sure that all unused analog input channels are
   referenced to analog ground.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Service and Support

Please be prepared to answer the following questions wh
call for technical support. This information helps to ident
specific system and configuration-related problems.

- What is the serial number of the product?
- What other boards or applications are installed in yo
system?
- Have you verified that all your connections are made
properly and are secure?
- Have you been able to isolate the source of your pro
input or output device, board, software?
- Are all unused input channels grounded?

In addition, if your problem appears to be software-relate
specify whether or not it has occurred more than once a
possible, document the sequence of steps that caused th
problem. This information helps us to reproduce the prob
in-house.

If your problem is not solved by the steps outlined in the
sections, the Technical Support staff is available to provi
prompt, reliable assistance. They can be reached by cal

(508) 481-3700

If you are located outside the USA, call your local distribu
The name and telephone number of your nearest distribu
listed in your Data Translation catalog.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

If Your Board Needs Factory Service

If your screw terminal panel must be returned to Data Translation:

1. Call the Customer Service Department and obtain a Retu Material Authorization (RMA) number.

   All return shipments to Data Translation must be marked the correct RMA numbers to ensure proper processing.

2. Package the board as follows:

   - Wrap the board in an electrically conductive plastic material, handling it with ground protection. *A static discharge can destroy components on the board.*
   - Pack the board in its original shipping container, if available and in good condition. If the original contain is not available, use any secure shipping container.

3. Return the board to the following address, making sure th RMA number is visible on the outside of the box.

   Customer Service Dept.
   Data Translation, Inc.
   100 Locke Drive
   Marlboro, MA 01752-1192

Outside the USA

If you are located outside the USA, call your local distributor fo authorization and shipping instructions. The name and telephone number of your nearest distributor are listed in your Data Translation catalog.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Index

+/-12V supply connections, 16

A

Accessories, 4
  EP246 enclosure, 4
  EP777 cable, 4
  resistor network, 4

Attenuation resistors, 21

C

Connecting digital inputs, 14

Connecting output switches, 15–18

Customer service, 5, 36

D

Detector connections, 11–13
  alternate, 12
  input attenuation, 13
  preferred, 12

DT777
  accessories, 4
  connecting to the DT2802, 6
  detector connections, 11–13
  diagram, 10
  electrical specifications, 29
  external event device connections, 14–15
  general description, 3
  J1 connector pin assignments, 25
  other connections, 16–18
  physical specifications, 31–32
  screw terminal assignments, 9-
  unpacking, 5

E

Electrical specifications, 29–30

EP246 enclosure, 4
  physical specfications, 31

EP777 cable, 4

External event device connection 14–15
  digital inputs, 14
  output switches, 15

External event devices, I/O prote 29

F

Floating channels, 35

H

Help, ix

I

Input connections, troubleshootin

Installation procedure, 6

J

J1 connector, 25

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

M

Minimizing noise, 11

N

Noise protection, 4

O

Other connections, 16–18
  +/-12V supplies, 16
  Vref I/O, 16

P

Physical specifications, 31

R

Required equipment, 3

Resistor network, 4, 21
  installing, 13

Resistors, list of functions, 21

Returning the board for serv

S

Screw terminal assignments

Service and support, 36

Signal attenuation, 3, 9, 13
  specifications, 29

T

Troubleshooting, 35

U

Unpacking, 5

V

Vref
  sensing, 16
  supplying, 17

Vref I/O connections, 16

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

WARRANTY AND SERVICE POLICY

WARRANTY

Data Translation, for the effective period of the warranty set out below, warrants that its standard hardware products and software media will be free from defects in materials and workmanship under normal use and service. Data Translation's obligation under this warranty shall not arise until the Buyer returns the defective product, freight prepaid, to Data Translation's facility or another specified location. The only responsibility of Data Translation under this warranty is, at its option, to replace or repair, free of charge, any defective component part of such products.

EFFECTIVE PERIOD OF WARRANTY

All boards, modules, systems, and software media are warranted for one year from date of invoice. Products repaired or replaced within this warranty period maintain their original warranty. Products repaired or replaced outside the warranty period are warranted for 30 days from invoiced date of return.

RETURN OUTSIDE USA AND CANADA

Contact the local sales representative or factory for authorization and shipping instructions.

LIMITATIONS ON WARRANTY

The warranty set forth above does not extend to and shall not apply to:

1. Products which have been repaired or altered by other than Data Translation personnel, unless Buyer has properly altered or repaired the products in accordance with procedures previously approved in writing by Data Translation.

2. Products which have been subject to misuse, neglect, accident, or improper installation.

THE WARRANTY AND REMEDIES SET FORTH ABOVE ARE IN LIEU OF ALL OTHER WARRANTIES EXPRESSED OR IMPLIED, ORAL OR WRITTEN, EITHER IN FACT OR BY OPERATION OF LAW, STATUTORY OR OTHERWISE, INCLUDING WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE WHI TRANSLATION SPECIFICALLY DISC TRANSLATION, INC. NEITHER ASS AUTHORIZES ANY OTHER LI CONNECTION WITH THE SALES, INS OR USE OF ITS PRODUC TRANSLATION, INC. SHALL HAVE N FOR INCIDENTAL OR CONS DAMAGES OF ANY KIND ARISING SALE, INSTALLATION, OR US PRODUCTS.

CUSTOMER SERVICE POLICY

a. Data Translation, Inc. will repair or option, any faulty item within ten day of said part regardless of its warran

b. A Quick Turnaround Plan is availab servicing. Contact the Data Translat Service Department for details.

c. If a product should fail during the wa it will be repaired free of out-of-warranty repairs, the custo invoiced for repair charges.

RETURN/REPAIR PROCEDURES

Upon determining that repair services or to return a product to Data Translatio Data Translation Customer Service D (508) 481-3700, extension 394 to obt Material Authorization (RMA) numbe following information ready when you c product model number, product serial nu address, and telephone number of per product; and any special repair instruct

Carefully package the product i packaging, making sure the RMA numb the outside of the package. Ship prepa

Customer Service Departme
Data Translation, Inc.
100 Locke Drive
Marlboro, MA 01752-1192 U

On completion of required service, a issued stating charges, when applica work that has been completed. A pa enclosed in the return shipment.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Reader's Comment Form

Please use this form to let us know about your experience using this manual. Your suggestions will help us continue to improve the accuracy, organization a overall quality of our user publications.

Manual title: _____

Document UM– _____

What is your general reaction to this manual?

Accuracy	☐ excellent	☐ very good	☐ good	☐ fair	☐ po
Completeness	☐ excellent	☐ very good	☐ good	☐ fair	☐ po
Usability	☐ excellent	☐ very good	☐ good	☐ fair	☐ po

What parts of the manual do you use most?

1. _____ 2. _____

3. _____ 4. _____

What errors have you found?

_____

_____

_____

Additional comments or suggestions:

_____

_____

_____

If we have questions about your comments, may we call you? ____Yes ____N

Name: _____ Title: _____

Company: _____

Address: _____

_____

Phone: _____

Please mail this form to:

Technical Publications Department
Data Translation, Inc.
100 Locke Drive
Marlboro, MA 01752–1192

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP