

ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.


- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Visit our website - *Click HERE*

CPS Disc Centrifuge

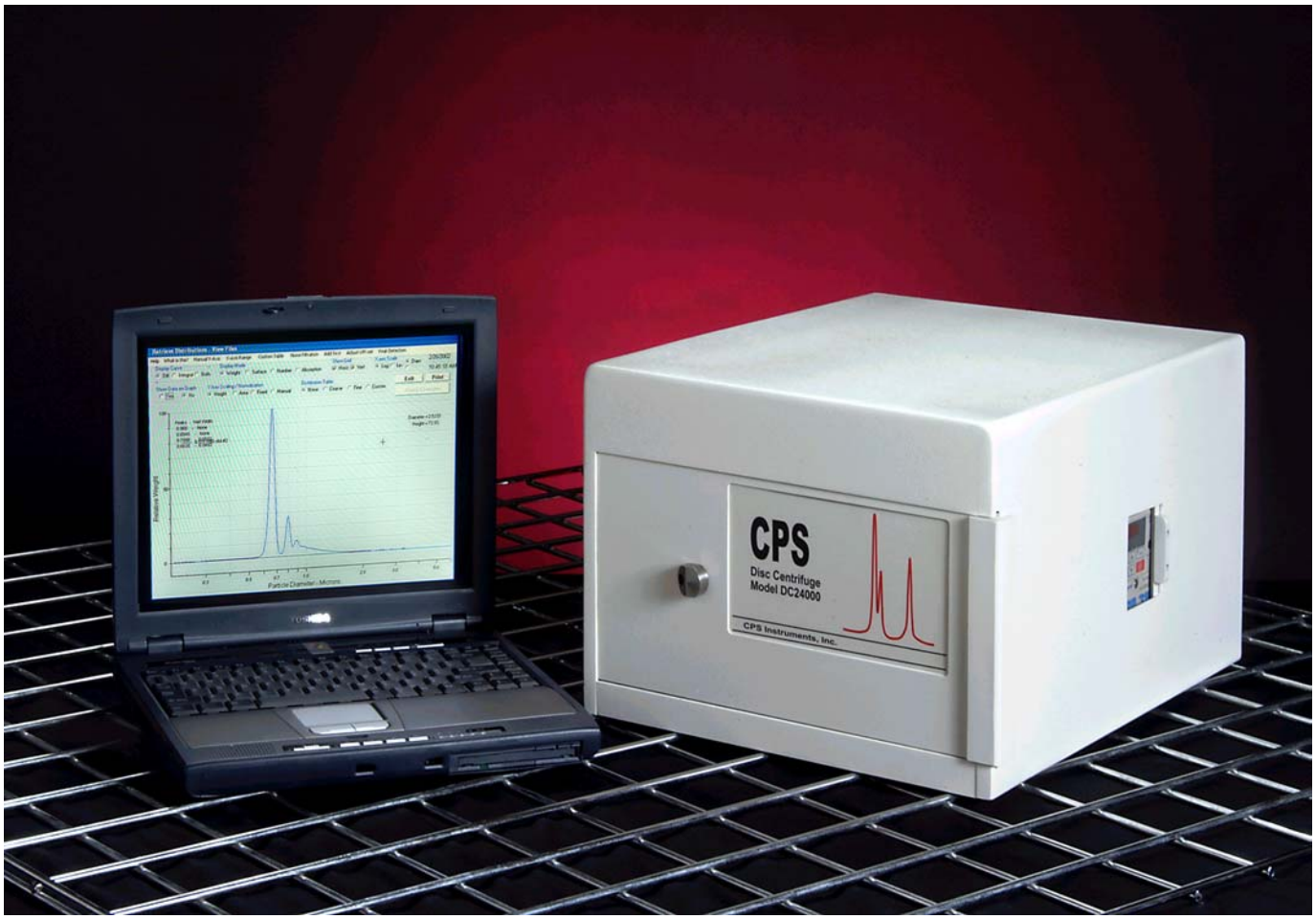

ARE YOUR PARTICLE SIZE RESULTS REAL?

Get realistic results with CPS Particle Size Analyzers.

CPS Instruments, Inc.
7349 SE Seagate Lane
Stuart, Florida 34997 (USA)
Tel: 772 221 7892 Fax: 772 221 7893
http://www.cpsinstruments.com
Email: stevefitzpatrick@cpsinstruments.com
ronellis@cpsinstruments.com

The CPS Disc Centrifuge can measure virtually any material between <0.02 and 40 microns.

Some example applications:

Polymer latexes and emulsions

Microabrasives

Printer/copier toner powders

Microspheres

Agglomeration patterns

Impact modifier particles

Protein clusters

Fillers ($CaCO_3$, clay, barites, etc.)

SiO_2 dispersions

Oil emulsions

Carbon black

Magnetic iron oxide

Liposomes

Pigments - water and oil based

Inkjet inks

Particles in diagnostic tests

Microencapsulated drugs

Virus particles/virus-like particles

Cells (culture) and cell fragments

Abrasives (of all types)

CMP compounds for integrated circuits

Micro-fiber paint viscosity modifiers

Starch/flour particles

Disclaimer
The descriptions, tables, and illustrations in this document are supplied for general information only. The content is presented in good faith and is believed to be completely accurate. However no warrantee or guarantee as to its accuracy is made; the content of this document does not commit CPS or its representatives to specific instrument performance guarantees.


Differential Sedimentation Method

The CPS Disc Centrifuge measures particle size distributions using sedimentation, a well known and reliable method of particle size analysis. Particles settle in a fluid under a gravitational field according to Stokes' Law. Sedimentation velocity increases as the square of the particle diameter, so particles that differ in size by only a few percent settle at significantly different rates.


In differential sedimentation, all the particles in a sample begin sedimentation as a thin band. If all the particles are of the same size, they settle at the same speed and arrive at a detector beam as a thin band. The time needed to reach the detector is used to calculate the size of the particles. A broad distribution of sizes separates during sedimentation into a broad band, while a mixture of narrow sizes separates into separate and distinct narrow bands, each of which arrives at the detector at a different time.

Particles between <0.02μ and >40μ can be measured. The actual size range that can be measured depends on the difference in density between the particles being measured and the density of the fluid in which the analysis is run.


The CPS Disc Centrifuge can routinely separate particles that differ in size by less than 5%, and separation of particles that differ by as little as 2.5% is possible. This resolution is much higher than available from other particle sizing techniques.

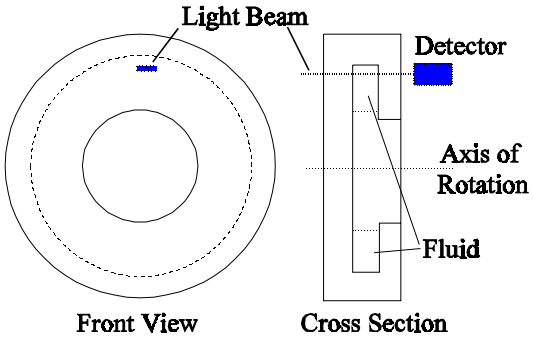

Hollow Disc Centrifuge Design

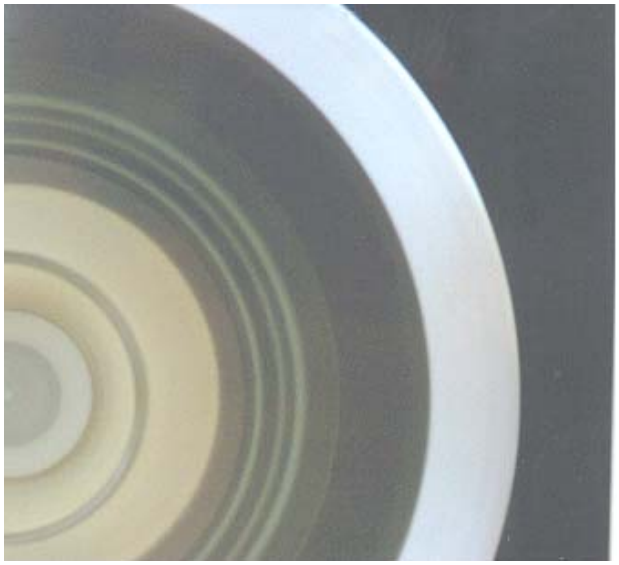


Three Bands of Particles During Sedimentation

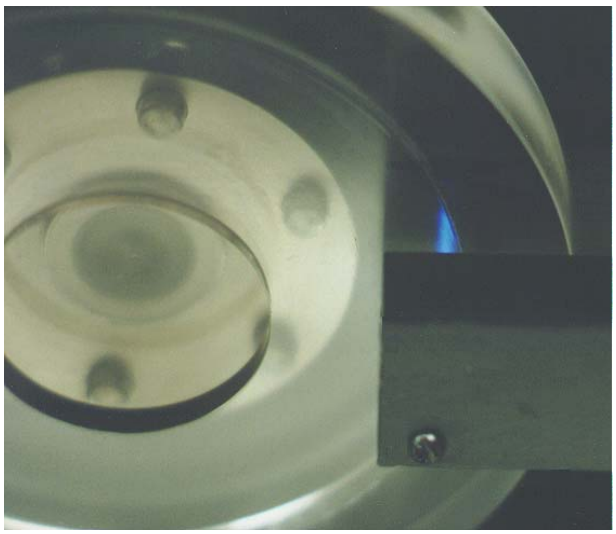

Closeup of Disc and Light Source

CPS Disc Centrifuge particle size analyzers provide high resolution, accurate results, even with very non-ideal samples that completely mislead other particle sizing methods.

Resolution
The size distribution graph below shows typical resolving power of the CPS Disc Centrifuge. This distribution is from a mixture of narrow polystyrene calibration standards. The 0.366μ and 0.343μ peaks differ in size by only 6.5%, yet they are almost completely resolved. Most of the width of these peaks is real, not from the instrument. Very narrow peaks that differ by as little as 3% can be completely resolved, while narrow peaks that differ by as little as 2% can be partly resolved. No other particle sizing instrument can match this resolving power.

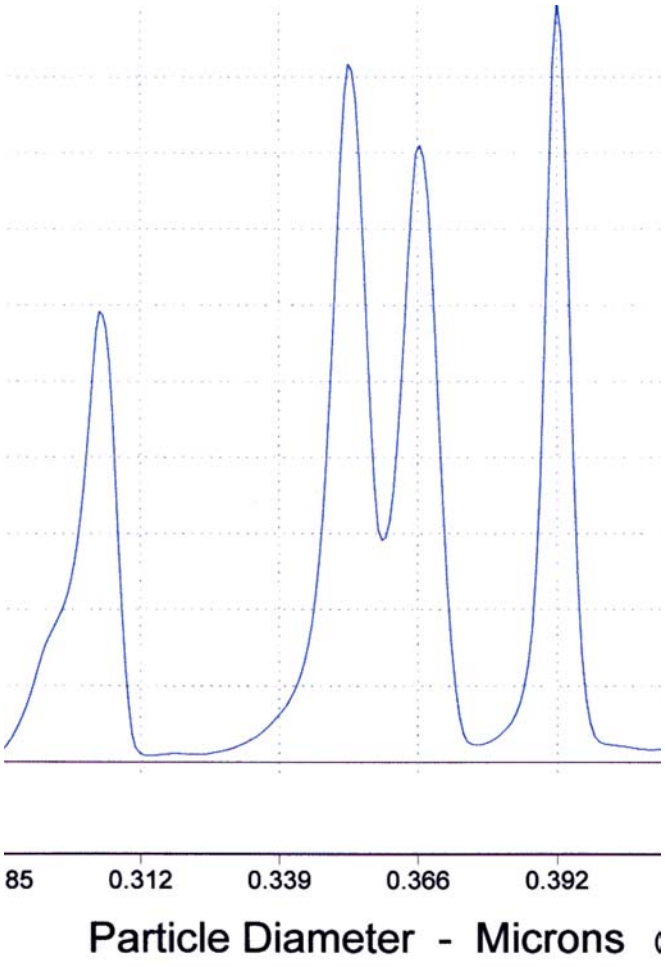

Accuracy and Precision
All analyses are run against a known calibration standard, so high accuracy is assured. Calibration can be either external (calibration standard injected before the unknown), or internal (calibration standard mixed with the unknown). Typical precision of reported sizes with an external standard is about ± 0.5% (95% confidence), and better than ± 0.25% with an internal standard. Replicate runs of the same sample produce virtually duplicate results in all cases.

Sensitivity
Even at 10^{-6} gram active sample weight, the CPS Disc Centrifuge reports an accurate particle size distribution. The lower detection limit for narrow samples is **well below 10^{-8} gram**, so even trace quantities of many kinds of particles can be detected. This high sensitivity allows accurate analysis of microgram samples on a routine basis.

Dynamic Range
Dynamic size range is the ratio of largest to smallest sizes that can be measured in one analysis. The CPS Disc Centrifuge has a practical dynamic range of up to 60 using a fixed centrifuge speed, and up to ~1000 using centrifuge speed ramping. The broad dynamic size range with speed ramping allows ultra high resolution analysis of very broad size range samples, which in the past could only be measured at lower resolution using light scattering methods.

Want to Learn More About Performance?
Visit the CPS Instruments web page:

http://www.cpsinstruments.com

Compare Particle Sizing Methods

The table below shows how the CPS Disc Centrifuge compares in some important characteristics with other particle sizing instruments. In cases where more than one competitive instrument of the same basic type is available, typical values are used.

The CPS Disc Centrifuge offers an exceptional combination of resolution, operating range, and measurement accuracy, at a far lower price than competitive instruments.

Sizing Method	Resolution[1]	Medium	Total Range[2]	Dynamic Range[3]	Analysis Time[4]	Precision	Price[5]
Electrozone Particle Counter	15-20%	Water	500µ - 0.5µ	~30	~3 minutes	±2.5%	1.4
Low Angle Laser Light Scattering	15-20%	Most Fluids	2,000µ - 0.1µ	~20,000	3 - 10 minutes	±0.5%	1.7
Photon Correlation Spectroscopy	~70%	Most Liquids	~4µ - 0.004µ	~1,000	~3 minutes	±1%	1.3
Capillary Hydrodynamic Chromatography	~ 15%[6]	Water[7]	1.1µ - 0.015µ	73	~10 minutes	±1%	1.3
Optical Particle Counter	~20%	Most Liquids	2,500µ- 0.5µ	~200	3 - 5 minutes	±2%	1.3
Competitive Disc Centrifuge	~4%	Many Liquids	40µ - 0.03µ	~40	9 - 300 minutes[8]	±1%	1.0
CPS Disc Centrifuge	**~2%**	**Most Liquids**	**40µ - 0.01µ**	**~1000[9]**	**3 - 40 minutes[10]**	**±0.25%[11]**	**1.0**

Notes

1. Resolution is the minimum difference in size between two perfectly narrow peaks that allows 95% separation.
2. Total range is from the largest practical size to the smallest practical size.
3. Dynamic range is the ratio of largest to smallest sizes in a single analysis.
4. Analysis time is total time, sample to sample.
5. Relative purchase price.
6. 15% resolution only after applying mathematical deconvolution; physical resolution is ~ 40%.
7. Only operates using a proprietary mixture of surfactants in water; solution must be purchased.
8. Maximum dynamic range for a 10 minute analysis is ~15.
9. With disc speed ramping; dynamic range is ~ 60 for fixed speed.
10. Maximum dynamic range for a 10 minute analysis is ~150.
11. ±0.5% precision with external calibration.

Designed for Reliability

- The CPS Disc Centrifuge is designed for long life without maintenance. The instrument has only a few moving parts: the drive motor, the disc, the cooling fan, and the door safety lock.

- The system is mechanically robust and durable. The permanently lubricated, precision bearings have a projected lifetime of >30,000 hours; that's 10 years of operation 8 hours per day at maximum speed. If the centrifuge is used at less than maximum speed, the bearing lifetime will be much longer.

- All critical electronics are encased in epoxy potting compound, virtually eliminating the possibility of damage from mechanical fatigue.

- All electronic components operate at 10% to 80% of their rated capacity, to insure extra long life.

- There is absolutely no scheduled maintenance, except to clean the disc chamber after use (normally with soap and water).

- The CPS Disc Centrifuge comes with a 2 year free replacement warrantee on **any part** that fails in normal use.



18,000 RPM, Solvent Resistant Disc - Reinforced with High Strength Alloy


Software


Operating software for the CPS Disc Centrifuge is completely graphical (Windows native), and is easy to learn and to use. The software handles all data collection and analysis, and has the following capabilities:


- Background data collection
- Publication quality graphs; full page or half page, on any available printer
- Weight, surface, number, and absorption distributions
- Log or linear scaling
- Graphical overlays of up to 20 distributions
- Averaging of multiple distributions
- Overlay of two averages
- Overlay average and individual distributions
- Notations can be added to graphs before printing, 8 point to 24 point letters
- Custom distribution summary tables
- Password protected file manager utilities
- On-screen centrifuge speed control
- Import and export of distributions
- Adjustable data functions, including noise filtration.


When there is a software function you need that is not available, CPS can usually customize the operating software to meet your needs.

The software also controls the CPS Autosampler, which allows completely unattended instrument operation, with automatic injection of up to 20 samples and calibration standards. The screen images below show a few examples of the operating software. The software functions under MS Windows 95, 98, NT, 2000, and XP.

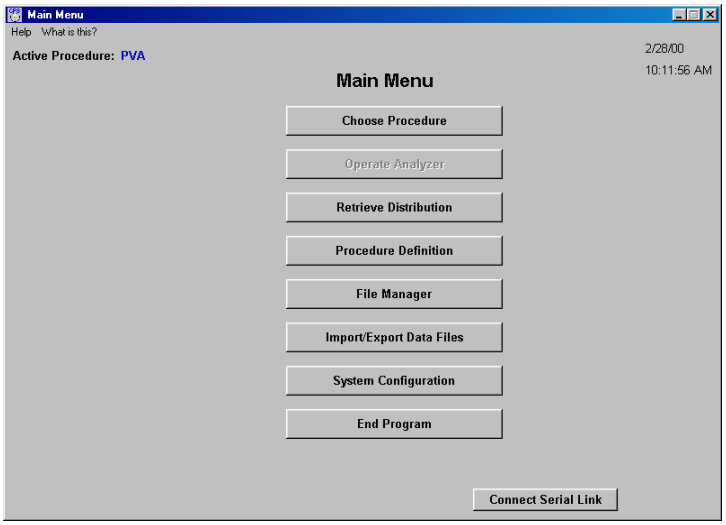

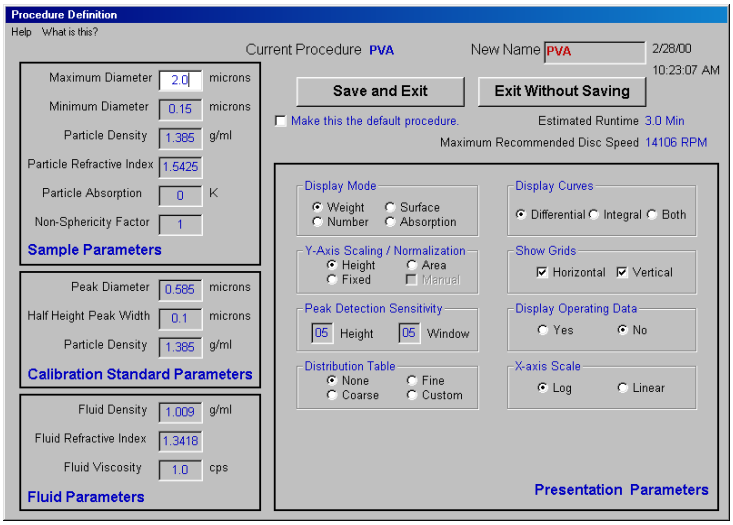

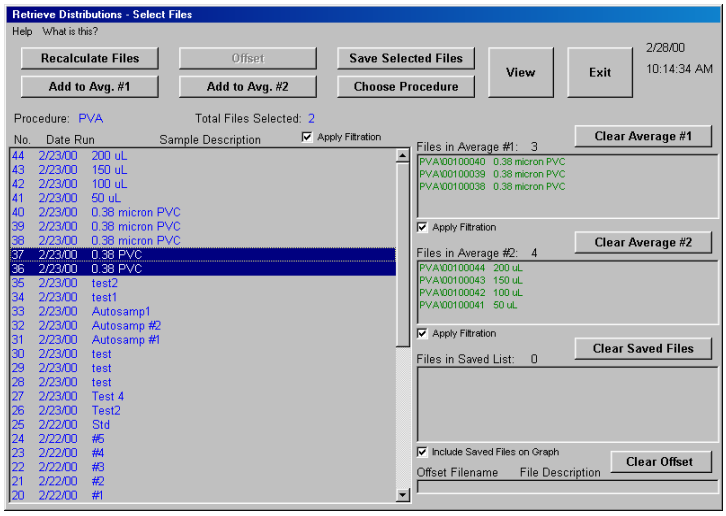

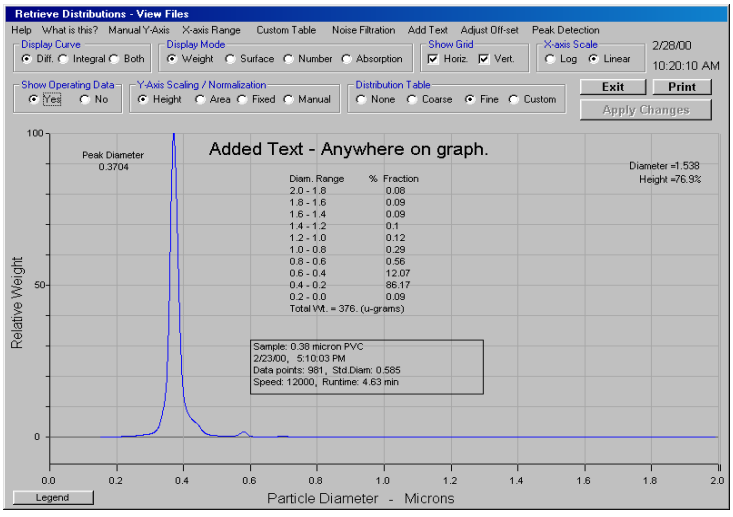

The Leading Differential Sedimentation Technology

Centrifuge Speed

The CPS Disc Centrifuge operates at up to 24,000 RPM. Analyses of very small particles are completed three times faster than competitive disc centrifuges. Analyses of samples with very wide size ranges can be also be completed rapidly, through the use of centrifuge speed ramping. Many analyses that are impractically long (many hours) on competitive sedimentation instruments are completed rapidly on the CPS Disc Centrifuge.

High Accuracy Calibration

The CPS Disc Centrifuge is always calibrated against known particle size standards. Internationally recognized size standards insure that reported size distributions are both consistent and accurate. You can add a trace of calibration standard to an unknown so that analyses are *internally calibrated.* Internal calibration yields typical accuracies of better than ±0.25% for reported peak sizes.

Measure Low Density Particles

Differential sedimentation has historically been limited to particles more dense than the fluid in which they are suspended. CPS has developed a new differential technique (US Patent 5,786,898), where the sample is delivered to the **bottom** of the centrifuge chamber, rather than to the fluid surface. This method allows particles that are lower in density than the fluid to float toward the fluid surface. Even neutral buoyancy particles in water (density ~1.0 g/ml) can be measured by conducting the analysis in deuterium oxide (density 1.108 g/ml) instead of water. Many materials that were difficult or impossible to measure by differential sedimentation can now be measured easily and with high resolution.

Speed Ramping

CPS has developed a special disc design that allows the centrifuge speed to be changed during an analysis, without creating turbulence in the sedimentation fluid. The centrifuge speed can be ramped up and down over a 20 fold range of speeds without disruption of the fluid. A 20 fold range of speeds in a single analysis expands the dynamic range for the analysis by almost a factor of 20. Other disc centrifuges are limited to a dynamic range of ~40, but the CPS Disc Centrifuge has a dynamic range of **up to 1000.** Many kinds of samples that in the past could only be measured by light scattering methods can now be measured by differential sedimentation, but at much higher resolution.

Low Noise Light Source/Detector

The light source and electronics in the CPS Disc Centrifuge are carefully designed to minimize noise in the signal. Typical signal to noise ratio is in the range of +/- 1 part in 85,000. Low noise in the signal gives very high sensitivity: **the detectible weight for a peak is oftem less than 10^{-8} gram.** High resolution analysis of a few micrograms (active) weight is routine.

Fast, High Resolution A/D Interface

The usable signal resolution (what the software works with) is only as good as the analog to digital (A/D) conversion process. The A/D conversion in the CPS Disc Centrifuge is better than **20 bits at 31 readings per second;** a conversion error of less than ±1 part in 10^6. Most competitive instruments work with 12 bits (±1 part in 4,096) or 16 bits (±1 part in 65,536).


Rapid Detector Response

The detector and electronics in the CPS Disc Centrifuge operate with an effective response time (to 95% accuracy) of ~0.1 second, while competitive sedimentation instruments have a response time of ~1 to 4 seconds. Rapidly moving peaks, which pass the detector beam during a period of only 1 to 2 seconds, are accurately recorded by the CPS Disc Centrifuge. The ability to accurately record rapidly moving peaks expands dynamic range, and allows much faster analyses of most samples. For example, analysis of polyvinyl chloride latex between 3μ and 0.07μ requires ~ 20 minutes with the CPS Disc Centrifuge (~8,000 RPM), but over an hour and 20 minutes with competitive sedimentation instruments (~4,000 RPM).

Resolution

No other particle size analyzer (using any analysis method) matches the resolving power of the CPS Disc Centrifuge. Clean separation of peaks that differ by 5% is routine. The particle size distribution is shown more clearly than by any other instrument. When the instrument is set up for maximum resolution, the measured peak widths for NIST tracable standards are nearly identical to those reported for electron microscopy.

The resolution is so good that for most analyses the instrument's contribution to the measured distribution is essentially zero... you see the true sample distribution rather than a mixture of the sample and the instrument.

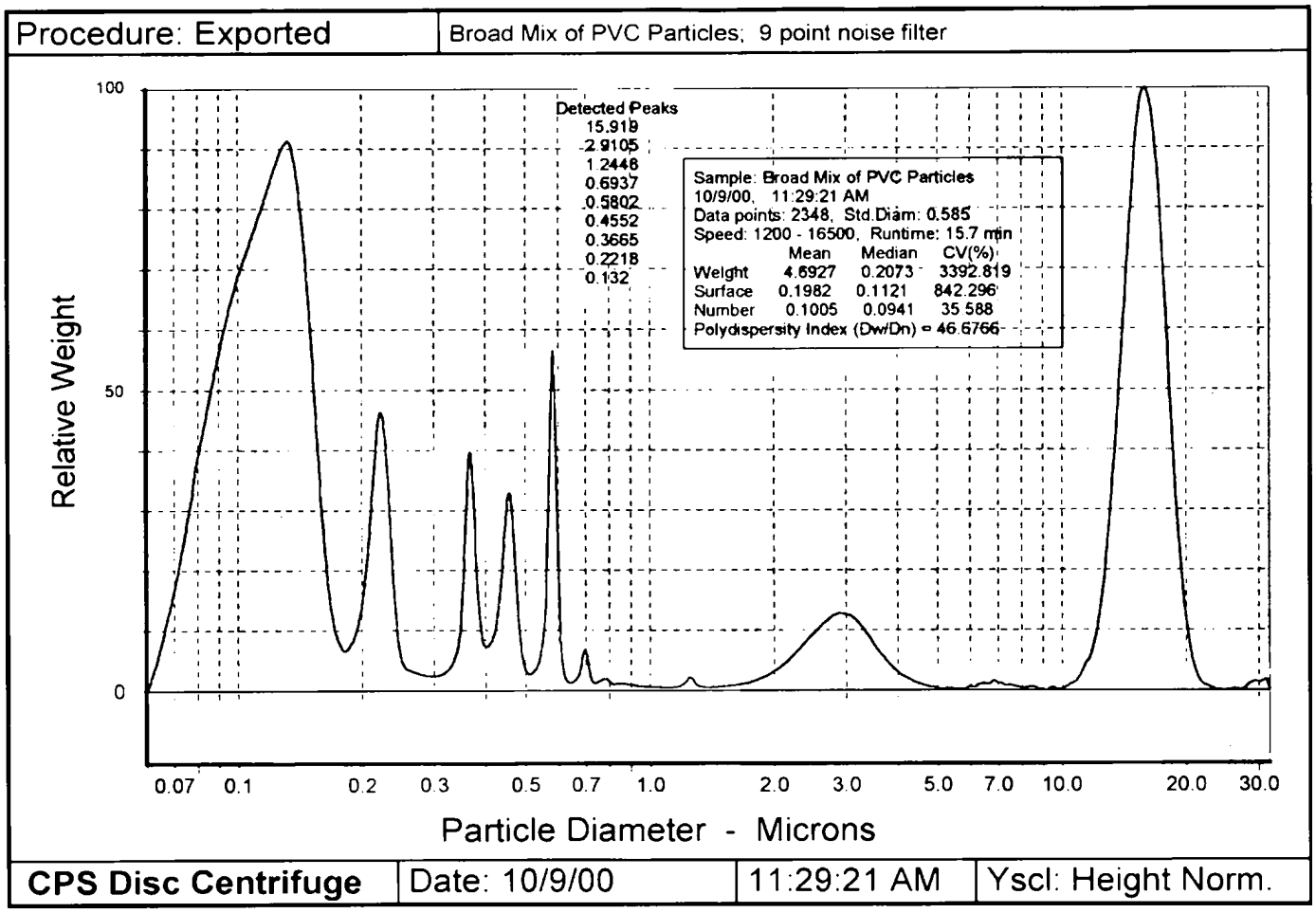

System Configuration

Basic Systems

Model DC12000 - 12,000 RPM Maximum Speed

Model DC18000 - 18,000 RPM Maximum Speed

Model DC24000 - 24,000 RPM Maximum Speed

Operating Software[1]
Manual
Two NIST Traceable Calibration Standards
Starter set of sedimentation fluids and sample preparation fluids

Disc Options[2]

Solvent Resistant Disc - allows analyses in nearly all organic solvents

Disc for Low Density Particles

Disc for Speed Ramping

Other Options

Automatic Density Gradient Builder

Automatic Sample Injector

Notes:

1. **A minimum computer is 233 MHz Pentium with several gigabytes available hard disc space. A color printer is recommended.**

2. **All disc options may be combined in a single disc if desired.**

For More Information

In the United Kingdom:
Analytik Ltd.
Unit 35 Byards Green
Potton
Sandy, Bedfordshire SG19 2SB
Tel: 44(0) 870 991 4044
Email: ian.laidlaw@analytik.co.uk

Germany, Switzerland, Austria:
L.O.T. Oriel GmbH & Co. KG
Im Tiefen See 58
D-64293 Darmstadt
Tel: (0)61 51 8806 63
Fax: (0)61 51 8966 67
Email: Wittmer@LOT-Oriel.de

The rest of Europe:
Inventech Benelux BV
Beneluxweg 37
Postbus 161, 4900 AD Oosterhout
The Netherlands
Tel: 31 (0)16 2472478
Fax: 31 (0)16 2472460
Email: steinmetz@inventech.nl

In Japan:
Nihon Rufuto Co. Ltd.
5-1-8 Higashiueno
Taitoh-ku, Tokyo 110- 0015
Japan
Tel: 81 3 3874 6880
Fax: 81 3 3847 6890

USA:
Michael G. Brown & Associates, Inc.
14 Locust Lane
Newtown, PA 18940
Tel: 215 497 1060
Fax: 215 497 1059
Email: mgba@voicenet.com

The rest of the World:
CPS Instruments, Inc.
7349 S.E. Seagate Lane
Stuart, Florida 34997
USA
Tel: 001 772 221 7892
Fax: 001 772 221 7893
stevefitzpatrick@cpsinstruments.com
ronellis@cpsinstrumetns.com



Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP