



Artisan Technology Group is your source for quality new and certified-used/pre-owned equipment

- FAST SHIPPING AND DELIVERY
- TENS OF THOUSANDS OF IN-STOCK ITEMS
- EQUIPMENT DEMOS
- HUNDREDS OF MANUFACTURERS SUPPORTED
- LEASING/MONTHLY RENTALS
- ITAR CERTIFIED SECURE ASSET SOLUTIONS

SERVICE CENTER REPAIRS

Experienced engineers and technicians on staff at our full-service, in-house repair center

*InstraView*SM REMOTE INSPECTION

Remotely inspect equipment before purchasing with our interactive website at www.instraview.com ↗

WE BUY USED EQUIPMENT

Sell your excess, underutilized, and idle used equipment. We also offer credit for buy-backs and trade-ins. www.artisanng.com/WeBuyEquipment ↗

LOOKING FOR MORE INFORMATION?

Visit us on the web at www.artisanng.com ↗ for more information on price quotations, drivers, technical specifications, manuals, and documentation

Contact us: (888) 88-SOURCE | sales@artisanng.com | www.artisanng.com



horizon

Model 755

ELITE

24-Place Horizontal • Variable Speed

A laboratory centrifuge in a class by itself.

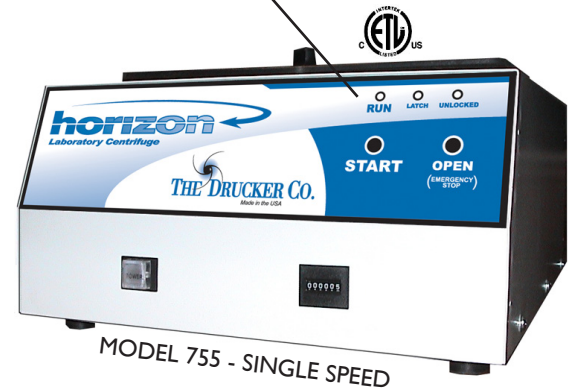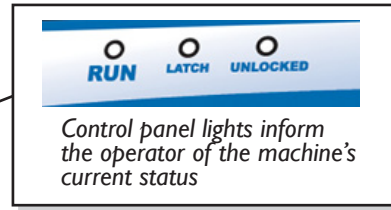

- Large Capacity- up to twenty-four (24) 10mL Test Tubes, twelve (12) 15mL Test Tubes, or six (6) 50mL tubes
- Brushless-DC Motor
- Full Horizontal Separation for Best Results
- Automatic Braking Reduces Cycle Time
- Precise Variable Speed Control (model 755V)
- Audible End-of-Run Indicator
- Push-Button Operation/Safety Lid-Lock
- ETL Compliant to ANSI/UL Standards
- Cycle Counter to Monitor Machine Usage

The Horizon Elite, model 755, truly delivers. Its brushless DC motor provides for extremely smooth operation and precise speed control. Select from several rotor options to suit your needs, (see reverse) and choose either single-speed or variable.


SPECIFICATIONS:

General Specifications for the Horizon Elite, Model 755 Centrifuge

Overall Dimensions:	Height:	9.0 in. (23 cm)
	Width:	14.5 in. (37 cm)
	Depth:	17 in. (43 cm)
Max Acceleration Time:		15 seconds
Max Deceleration Time:		15 seconds
Weight:		39 lbs. (17.7 kg)
Centrifuge Motor:		1/30 H.P. AC
Protection Breaker:		4 Amp. Re-settable Electronic,
Timer:		1 to 30 minutes +/- 1%
Max. Power Requirements:		200 Watts
Voltage:		115/230 (+/-10%) Volts
Frequency:		50/60 Hz

Why Horizontal Separation?

A horizontal rotor produces better separation with a tighter gel line, eliminating remixing and re-spins.

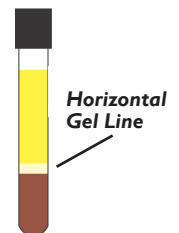

10mL test tubes up to 17mm x 100mm, 15mL test tubes up to 17mm x 125mm, with the Standard Rotor UL Std. 61010A-1 and UL Std. 61010A-2-20, CAN/CSA C22.2 No. 1010.1 and CAN/CSA C22.2 No. 1010.2.20



OPERATING FEATURES

- **ELECTRONIC TIMER:** Time is selectable from 1 to 30 minutes.
- **BRUSHLESS DC MOTOR:** For smooth and quiet operation.
- **ELECTRONIC BRAKING:** Braking reduces the deceleration time to approximately 15 seconds.
- **END-OF-RUN AUDIBLE INDICATOR:** Several beeps are sounded once the rotor has stopped to inform the operator that a cycle is completed.
- **VARIABLE SPEED:** The precision variable speed control gives the operator the ability to select from 500 RPM to the max. variable rotor speed (see table).
- **QUIET OPERATION WITH LOW VIBRATION:** The Double-encased motor and steel construction provide quiet performance and low out-of-balance vibration.
- **COOL RUNNING:** Unique design of the motor housing and rotation chamber allows for cool operation even under the heaviest workloads

SAFETY FEATURES

- **AUTOMATIC LID LOCK:** The centrifuge will not run with the lid open and the lid cannot be opened while samples are spinning.
- **POSITIVE LID LATCHING:** Firm mechanical latching of the lid prevents the escape of debris in case of breakage or an accident.
- **WELDED STEEL CONSTRUCTION:** For strength and durability.
- **OUT-OF-BALANCE PROTECTION:** If an attempt is made to run the centrifuge while unbalanced, the run is terminated and the rotor is safely brought to a stop.

ROTOR / CARRIER OPTIONS


The model 755 is available with many different rotor and carrier packages to choose from. Select the package which best suits your needs. Consult your authorized dealer or The Drucker Company for additional information.



	PART NO.	# PLACES	MAX TUBE SIZE	SPEED RANGE ±10%	FORCE RANGE
6-Carrier Horizontal (Standard)	7786022				
100mm 4-pl carrier	7713023	24	17mm x 100mm	500 to 3,150 RPM	40 to 1,760 xg
100mm 4-pl carrier	7713023	12	17mm x 125mm	500 to 3,150 RPM	40 to 1,760 xg
50mL horizontal carrier (optional)	7713037	6	34mm x 115mm	500 to 3,150 RPM	40 to 1,760 xg



12-Place Horizontal	7786024				
75mm tube holder	7713033	12	17mm x 75mm	500 to 3,200 RPM	40 to 1,280 xg
100mm tube holder	7713031	12	17mm x 100mm	500 to 3,200 RPM	40 to 1,540 xg
125mm tube holder	7713032	12	17mm x 125mm	500 to 3,200 RPM	40 to 1,850 xg

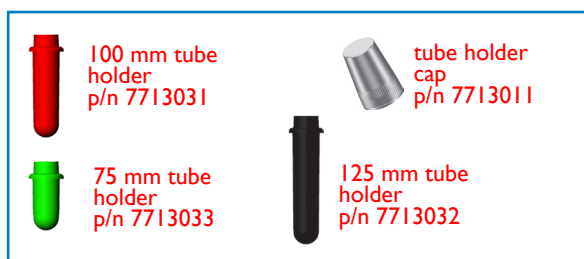

Standard Accessories




THE DRUCKER CO.
Providing Quality Centrifuges since 1932

Drucker, The Drucker Company, and Horizon are trademarks of The Drucker Company, Inc.
 © The Drucker Company 2005-2011. All Rights Reserved

MADE IN THE USA

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantech.com | www.druckercompany.com

168 Bradford Drive, Port Matilda, PA 16870 U.S.A.
 +1-814-692-7661 • Fax: +1-814-692-7662



Artisan Technology Group is your source for quality new and certified-used/pre-owned equipment

- FAST SHIPPING AND DELIVERY
- TENS OF THOUSANDS OF IN-STOCK ITEMS
- EQUIPMENT DEMOS
- HUNDREDS OF MANUFACTURERS SUPPORTED
- LEASING/MONTHLY RENTALS
- ITAR CERTIFIED SECURE ASSET SOLUTIONS

SERVICE CENTER REPAIRS

Experienced engineers and technicians on staff at our full-service, in-house repair center

*InstraView*SM REMOTE INSPECTION

Remotely inspect equipment before purchasing with our interactive website at www.instraview.com ↗

WE BUY USED EQUIPMENT

Sell your excess, underutilized, and idle used equipment. We also offer credit for buy-backs and trade-ins. www.artisanng.com/WeBuyEquipment ↗

LOOKING FOR MORE INFORMATION?

Visit us on the web at www.artisanng.com ↗ for more information on price quotations, drivers, technical specifications, manuals, and documentation

Contact us: (888) 88-SOURCE | sales@artisanng.com | www.artisanng.com