ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

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Contents

Safety Symbols	3
Safety precautions	4
Warranty terms	5
Product	6
Short description	6
Motor	7
Pump-head	7
Spring-loaded tube-bed	7
Mains voltage	8
Inserting the tubing	9
Starting your pump	10
Running the tubing in	11
Service life of tubing	11
Tubing quality	11
Type of tubing	12
Pumping against back-pressure	12
Pump tubing	12
When the pump is not in use	13
Tube connectors	13
Adhesives for tubing connections	13
Tubing life	13
Flow-rate chart	14
Maintenance	15
Changing a defective fuse	15
Spare parts	16
Disposal	16
Standards and Specifications	16
Technical Specifications	17

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Safety Symbols

Below you find the following symbols used throughout this manual:

This "STOP Sign" denotes a hazard. The user is instructed to refrain from a certain practice or procedure. Failure to comply with such warnings can result in personal injury and/or damage to the product.

This "Caution Sign" denotes a hazard. Certain practices or procedures may only be carried out by taking the necessary precaution.

Caution: rotating pumphead creates a pinch and crush hazard.

Caution: risk of electrical shock.

Protective earth: (Ground.)

General information about something important.

➡ **Wherever you see these signs, you are requested to follow the specific warnings and recommendations given.**

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Caution

Rotating pumphead creates a pinch and crush hazard.

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Safety precautions

ISMATEC® Pumps are designed for pumping applications in laboratories and industry. As such it is assumed that Good Laboratory Practice (GLP) and our following recommendations will be observed:

- The circuit between mains supply and pump has to connect to earth ground.
- The pump must not be operated outside the designated operating and environmental conditions.
- Ensure that the pump is switched OFF and the power cable disconnected before changing pump tubing.
- The pump must not be used:
 - for medical applications on human beings
 - in explosion proof chambers or in the presence of flammable gases or fumes.
- The pump must be switched OFF when pump-heads, cassettes or tubing are inserted or changed.
- The permeability of tubing depends on the material used and pressure conditions. Tubing can also become electrostatically charged. Please be aware of possible hazards when laying tubing in explosion-proof chambers.
- Pump-heads consist of rotating parts. Therefore, the pump must not be operated before the cassettes are fully snapped-in or the tube-bed completely shut.
- Do not manipulate the pump-head before the pump is switched off and disconnected from the mains supply.
- Be particularly cautious that no parts of your body such as fingers, long hair, etc. or jewellery, or loose objects such as cables or tubing, etc. can be trapped by the rotating roller head.
- Tubing can tear and burst during operation. If this could cause damage, the necessary safety measures based on the specific situation must be taken. When pumping corrosive media we recommend placing the pump in a collecting basin.
- The permeability of pump-tubing depends on the material used and pressure conditions. Moreover, tubing can also become electrostatically charged.

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Safety precautions

- Only new fuses, according to the specifications stated on Page 15 in this manual, may be used.
- The fuse-holder must not be short-circuited.
- Do not open or remove the housing while the pump is operating.
- Repairs may only be carried out by a skilled person who is aware of the hazard involved.
- For service and repairs carried out by the customer or by third-party companies, ISMATEC® denies any responsibility.

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Warranty Terms

For all ISMATEC® panel-mounted tubing pumps excluding consumables:

from date of delivery: 2 years

We warrant the perfect functioning of our products, provided they have been installed and operated correctly according to our operating instructions.

If production or material faults can be proved, the defective parts will be repaired or replaced free of charge at our discretion.

The duration of the warranty is not affected by making a claim for warranty service.

Further claims are excluded. Shipping costs for instruments under warranty are charged to the customer.

Our warranty becomes invalid in case of:

- improper operation by the user, or if the pump is diverted from its proper use,
- unauthorized modification or misuse by the user or by a third-party,
- improper site preparation and maintenance,
- operation outside environmental and electrical specifications for the product,
- use of third-party software, hardware, interfacing or consumables purchased by the user and which do not comply with our specifications.

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Product

➡ Please check the package and its contents for transport damage. If you find any signs of damage, please contact your local ISMATEC® representative immediately.

Complaints can only be accepted within 8 days from receipt of the goods.

Check if the following items are included in the delivery:

- Mini-S Pump
- Power Cord
- Operating Manual

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Short description

The **Mini-S** is a compact, stand-mounted peristaltic pump designed to easily meet the requirements of numerous pumping applications at low flow-rates.

This **3-channel fixed-speed pump system** allows you to deliver liquids or gases in the clockwise and counter-clockwise direction.

The following **6 models** are available:

Modell/Model	Bestell-Nr./Order #	Rollen/Rollers	Kanäle/Channels	Fließrate/Flow-rate
Mini-S 620	ISM857	6	3	0.1 - 9.0 ml/min
Mini-S 640	ISM858	6	3	0.1 - 20.0 ml/min
Mini-S 660	ISM859	6	3	0.2 - 30.0 ml/min
Mini-S 820	ISM860	8	3	0.06 - 8.0 ml/min
Mini-S 840	ISM861	8	3	0.1 - 14.5 ml/min
Mini-S 860	ISM862	8	3	0.2 - 21.0 ml/min

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All six pump models have the following features:

Fixed-speed drive
Spring-loaded tube-bed

Intended area of application:

The Mini-S is a standard tubing pump designed for suctioning, feeding and circulating liquids and gases in chemical and pharmaceutical laboratories and similar areas of applications. Depending on the tubing material chosen, it is suitable for almost all liquids (including liquids which contain particulates), corrosives and abrasive substances and gases. When using sterilizable tubing, the pump can also be used for pumping under sterile conditions.

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This instrument must not be used for medical applications on human beings.

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Motor

The Mini-S pump has a synchronous motor with an excellent start/stop behavior. The power consumption is 11 Watts.

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Pump-head

This pump is equipped with a 3-channel pump-head with 6 or 8 self-lubricating pump rollers.

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Spring-loaded tube-bed

The Mini-S pump series features a 3-channel tube-bed with a spring-loaded pressure lever. The pressure on the tubing can easily be adjusted with the knurled-head screw (see figure below).

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➡ The Mini-S pumps are equipped with a 12 mm rod and must be mounted to a stand (see picture above).

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Before operating your pump:

Please check that the voltage indicated on the identification plate (on the rear panel) agrees with the local mains voltage available.

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The mains socket used has to be connected to earth ground (protective conductor contact);

Otherwise a safe operation is not ensured and the instrument may be damaged. (For the power connection, it is recommended that the original power cord be used.)

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
230 V (220-240 VAC) / 50 Hz
oder
115 V (110-120 VAC) / 60 Hz

Mains voltage

The Mini-S pump models are available in two voltage versions. Depending on your order it is supplied for a voltage of

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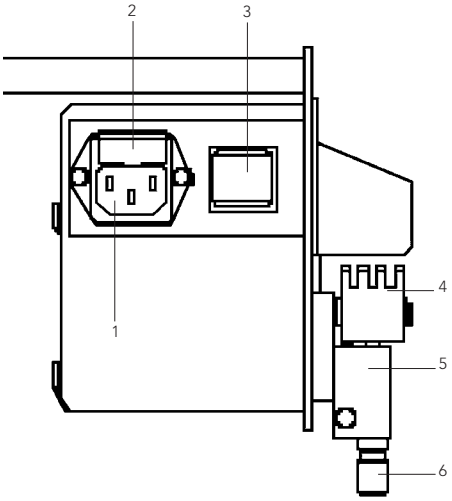
Offenes Schlauchbett
Open tube-bed

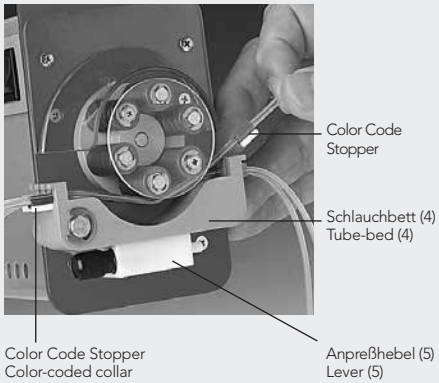Geschlossenes Schlauchbett
Closed tube-bed

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Inserting the tubing

➡ The Mini-S pumps take 3-stop collared tubing, allowing the same tubing to be used twice, and hence extending service life by 100%.

For loading or changing tubing you proceed as follows:

1. ⚠ Switch the pump OFF and unplug it.
2. Open the tube-bed by turning the white lever (5) to the right or the left. The spring-loaded 3-channel tube-bed (4) is lowered and provides free access.
3. Insert one collar of the tubing into the tube-holding slit of the tube-bed. Feed the tubing over to the other end of the tube-bed and insert the second collar into that slit (see figure below).


➡ – Insert the first tubing in the back of the tube-bed.
– The tube must lie correctly in the tube-bed (collars pointed upwards). It must not be twisted.

4. Do the same with the other 2 pump tubes.
5. Close the pump-head by turning the white lever (5) into the vertical position.
6. Start the pump as described on Page 10.
7. Adjust the correct pressure to the tubing by turning the knurled-head screw (6) so that a smooth flow is achieved.

 ➡ It is most important to adjust the tube-pressure under operating conditions, i.e. when all the tubes are connected and the pump system is completely installed. Apart from the tube size and the material, the correct pressure also depends on other parameters such as back-pressure, viscosity, etc.

 ⚠ The knurled-head screw (6) must not be tightened more than necessary. Otherwise, tubing and pump are exposed to unnecessary strain, which reduces the service life considerably.

 ➡ Depending on the application and tubing material, re-adjustment of the tube pressure may be necessary from time to time.

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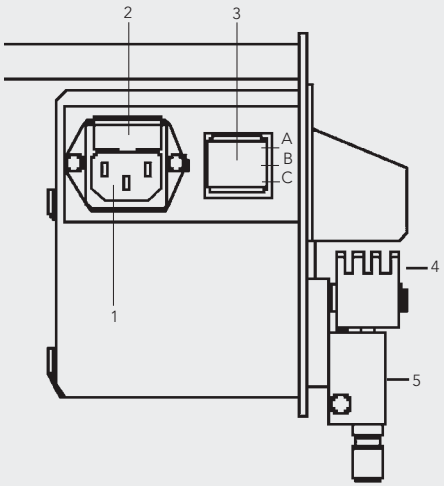Starting your pump

When the pump is running, the pump-head must be closed. A rotating, open pump-head is dangerous and is capable of trapping fingers, hair and clothing, causing personal injury or damage to the instrument.

1. Ensure that the pressure lever (5) is locked and the tubes lie correctly in the tube-beds. Check that all tubing connections are tight.
2. ⚠ Check that the location of the pump, as well as the connection and safety measures, agree with the "Warranty Terms" and "Safety Precautions" stated on Pages 4 and 5.
3. ⚠ Connect the pump socket (1) to the mains supply by using a protective earth cable. (It is recommended that the original power cord be used).
4. Press on the top or bottom of switch (3), depending on the flow-direction required (see figure at left).

 If the pump is connected to the mains, the pump-head will turn in the selected direction.
5. To stop the pump, press switch (3) again so that it is back in the middle position (see figure at left).

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Running the tubing in

➡ Every new tubing needs a certain run-in time. If **constant and reproducible flow-rates** are required, it is recommended that **new tubing be run-in for at least 1 hour using water or the medium to be delivered.**

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Service life of tubing

The service life of tubing depends mainly on the application (chemicals being used), tubing material, pressure, temperature, viscosity, etc. (For more information about the service life, please refer to the current ISMATEC® catalog.

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Tubing quality


➡ The pump-heads of ISMATEC® peristaltic pumps are designed for ISMATEC® tubing. The use of ISMATEC® tubing guarantees reproducible results and long service life.

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Pumping against back-pressure

The Mini-S pumps **must not be used for applications with a differential pressure (back-pressure) exceeding 1.5 bar.** (When using any combination of three tubes with the largest inner diameters (**2.54, 2.79 or 3.17 mm**), the back-pressure **must not exceed 0.5 bar**).

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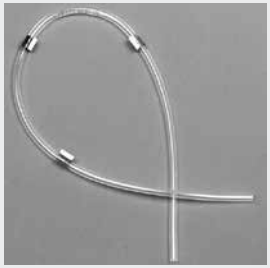Type of tubing

These cassette pumps take tubing with 3 color-coded collars. Thanks to these collars the tubes are always inserted under the same conditions, which results in a high repeatability of the flow rate. What is more, the different color codes help to easily recognize the various tubing diameters.

Tubing with 3 collars can be used twice. Thanks to this advantage the service life of the tubes can be doubled.

The selection of the right tubing material, especially with regard to chemical compatibility, is entirely the customer's responsibility. For chemical compatibility, please refer to the guidance chart in our current pump tubing selection guide.

The pump-heads are designed for ISMATEC® tubing. The use of ISMATEC® tubing guarantees a long tubing life and reproducible results. For further information about pump tubing, please refer to the current pump tubing selection guide.

Pump tubing

These cassette pumps take tubing with 3 color-coded collars. Thanks to these collars the tubes are always inserted under the same conditions, which results in a high repeatability of the flow rate. In addition, the different color codes help to easily recognize the various tubing diameters.

Tubing with 3 collars can be used twice. Thanks to this advantage the service life of the tubes can be doubled.

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Syphoning effect

When the tubing is un-occluded, the fluid can flow back to the reservoir.

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When the pump is not in use

➡ When the pump is idle, it is recommended that pressure on the tubing be released. This helps to protect the tubing from unnecessary strain and prolongs its service life.

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Tube Connectors

A number of different tube connectors are available. Please refer to the latest ISMATEC® catalog.

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Tubing life

The service life of the tubing depends on the application and the tubing material used.

Example: chemicals, rotation speed, differential pressure, temperature, viscosity, pressure on tubing, etc.

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Please note

The above mentioned flow rates are approximate values, determined with water, at ambient temperature, without back-pressure and with Tygon® tubing.

Exact flow rate values depend on various paramters such as quality and age of the tubing, back pressure, temperature and viscosity of the liquid, etc.

We will be pleased to help the user select the tubing material. However, the user has the final responsibility for the selection of the correct tubing material.

For reproducible results, it is recommended that ISMATEC® tubing be used.

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Please refer to our detailed tubing documentation or our website

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Flow rate charts

The Mini-S pump models can be used with 26 different tube sizes which provide the following flow rates per pump-channel:

Modell / Model	Mini-S					
Typ / Type	620	640	660	820	840	860
Pumpenrolle /Pump Rollers	6	6	6	8	8	8
Kanäle / Channels	1 - 3	1 - 3	1 - 3	1 - 3	1 - 3	1 - 3
Schlauchtyp / Tubing type	3 Stopper			3 Stopper		
Schlauch IØ / Tube I.D.	ml/min			ml/min		
0.13	auf Anfrage/on request			auf Angrage/on request		
0.19	0.1	0.1	0.2	0.06	0.10	0.20
0.25	0.2	0.3	0.4	0.14	0.23	0.37
0.38	0.2	0.5	0.7	0.20	0.40	0.60
0.44	0.3	0.7	1.0	0.30	0.60	0.80
0.51	0.5	1.0	1.2	0.40	0.78	0.97
0.57	0.6	1.1	1.8	0.48	0.95	1.50
0.64	0.7	1.3	2.1	0.56	1.10	1.80
0.76	0.9	1.7	2.6	0.82	1.40	2.30
0.89	1.3	2.4	3.5	1.10	2.00	3.00
0.95	1.5	2.7	4.1	1.25	2.30	3.50
1.02	1.5	2.9	4.4	1.40	2.60	3.94
1.09	1.6	3.1	4.8	1.50	2.80	4.20
1.14	1.7	3.5	5.4	1.60	3.00	4.81
1.22	1.9	3.8	6.2	1.75	3.30	5.30
1.30	2.4	4.2	7.1	1.90	3.60	6.00
1.42	2.8	5.7	8.2	2.30	4.00	6.40
1.52	3.1	6.0	9.3	2.50	4.50	7.50
1.65	3.6	6.4	10.1	2.98	5.30	8.65
1.75	4.1	7.8	11.8	3.30	6.20	9.70
1.85	4.7	9.0	13.5	3.80	7.20	11.20
2.06	5.3	9.5	15.3	4.40	8.60	12.80
2.29	5.8	12.5	18.0	5.00	9.50	14.00
2.54	7.0	14.0	21.0	5.40	11.00	15.50
2.79	7.8	16.2	24.3	6.00	12.00	19.28
3.17	9.0	20.0	30.0	8.00	14.50	21.00

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
The selection of the correct tubing material, especially with regard to chemical compatibility, is entirely the customer's responsibility.

For chemical compatibility, please refer to the guidance chart in our current pump tubing selection guide.

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Ensure that the new fuses are also slow blow and the rating is correct as described on this page.

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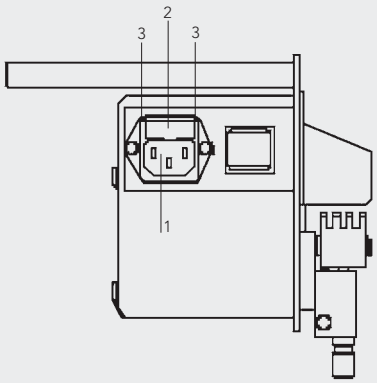Maintenance

Your ISMATEC® peristaltic pump does not require any special maintenance. Provided that you operate it properly and in accordance with the guidelines stated in this instruction manual, the tubing is the only part that needs replacing from time to time depending on the application and the operating conditions.

Changing a defective fuse

The fuse-holder is included in the mains supply socket (see drawing below).

1. ⚠ **Make sure that the pump is disconnected from the mains.**
2. Press the two clips (3) on both sides against the fuse-holder (2) and, at the same time, pull the fuse-holder out.
3. Two slow blow fuses will be found in this fuse-holder.
 – T100mA/230 V
 or
 – T200mA/115 V

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Disposal

Please retain packaging materials until the product warranty ends. Afterwards please discard packaging materials in an environmentally-friendly manner according to local regulations.

Once the useful life of the product has ended, please ensure proper disposal according to local laws. Plastic and electronic components should be disposed of at a recycling facility. Please refer to local regulations regarding proper disposal.

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Spare parts

After the warranty period, your ISMATEC® representative will be pleased to send you on request:

- spare parts
- parts lists
- wiring diagrams

Please give information on defect, date of purchase, serial-no., and model.

➡ **In case of any queries or if anything is unclear regarding the handling of your ISMATEC® pump, please contact your local ISMATEC® representative.**

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Standards and Specifications

The ISMATEC® pump supplied with this instruction manual complies with the CE-specifications and has the **CE-mark**. Additionally, this instrument meets the following standards:

Electromagnetic emission:	**EN 50081-1** **EN 55011, Class B**
Electromagnetic immunity:	**EN 50082-1**
Safety standards:	**EN 61010-1** **IEC 61010-1**

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Schlauch IØ
0.13 – 3.17 mm

Technical Specifications

Motor type
Synchronous motor with fixed-speed

Operating conditions

Temperature	+5 up to +40°C
Rel. humidity	max. 80 %

- not condensing, at normal environmental conditions

Differential pressure
1.5 bar for tubing with an I.D. up to 2.29 mm
0.5 bar for tubing with an I.D. 2.54, 2.79 or 3.17 mm, when all three tube-channels are used simultaneously

Dimensions/Weight

DxWxH	125 x 80 x 135 mm
Weight	1.1 kg

Mains connection
230 V (50 Hz)
115 V (60 Hz)

Fuse rating
2 x T100mA/230V
or
2 x T200mA/115V

Speed

20 rpm	(Mini-S 620 and 820)
40 rpm	(Mini-S 640 and 840)
60 rpm	(Mini-S 660 and 860)

Flow rates
See Page 14

Power consumption
11 Watts

Pump channels
3

Pump rollers
6 (Mini-S 620, 640 and 660)
8 (Mini-S 820, 840 and 860)

Type of tubing
3 Color-coded collars

Tubing I.D.
0.13 – 3.17 mm

ISMATEC®

Learn more about the Ismatec® product line by visiting: www.ismatec.com

For ordering and technical support, please contact:

North America sales@ismatec.com | 1-800-323-4340 | 1-847-549-7600

Europe sales.europe@ismatec.com | +49 (0) 9377 9203-0

Cole-Parmer®

©2015 COLE-PARMER INSTRUMENT COMPANY, LLC. 14-022, REV. C

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP