ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *IKA CC3 308B* at our website: *Click HERE*

huber

Compatible Control

The World of Thermoregulation

Simply efficient thermal regulation.

Operating instructions

for units with

Compatible Control

CC1, CC2, CC3

Peter Huber Kältemaschinenbau GmbH
Werner-von-Siemens-Strasse 1 · D-77656 Offenburg / Germany
Tel. +49-781-96030 · Fax +49-781-57211
E-Mail: info@huber-online.com · Internet: www.huber-online.com

Version 1.0/05

Symbols

Safety
Warning! A potentially hazardous situation. Identifies hazards sufficient to cause death or severe injuries if the safety instructions are disregarded.
Caution! A potentially hazardous situation. Identifies hazards sufficient to cause light injuries if the safety instructions are disregarded.

Definitions from our Huber Glossary and fundamental technical knowledge.

Entries at the Compatible Control controller

Service
This is where you get help: the Huber Hotline

EXTRA
Additional information

Preface

Dear Customer,

Congratulations! Units and devices manufactured by Peter Huber Kältemaschinenbau GmbH are always a good choice. Thank you very much for your trust.
To meet your demands as a user, we have revolutionized the user interface of our thermostats and implemented a uniform method of operating almost the entire HUBER product range. Many devices, ranging from small immersion thermostats to large Unichillers are operated via a single controller generation: **Compatible Control**.
The units operated with the Compatible Control are listed in the following table with a remark for those which also have a supplementary instruction manual.

On the type label on the rear of your device you will find important information such as:

CC2-505wl vpc	-50°... +200°C	SNr.: 59845/05
[Device name]	[Temp.range]	[Serial number]

Free choice in operation:
All the listed Units can be operated with three different controller versions:
Compatible Control CC1, the simple one
Compatible Control CC2, the convenient one
Compatible Control CC3, the one with dialog capabilities

You controller version is identified by the label above the display.
For details on your type of controller, please refer to Chapter 2.2

Please consider only those instructions in the present documentation that apply to your device type and controller version.

To avoid confusions:

- **Compatible Control CC1, Compatible Control CC2** or **Compatible Control CC3** means the removable controller only!

- **CC1-E,CC2-E** or **CC3-E** means the complete immersion thermostat including the Compatible Control.

All units with *no* separate instruction manual are listed in the following table. They are included in the present manual.

Combinations of immersion thermostats with baths and Compatible Control Thermostats with Compatible Control CC1, CC2, CC3

Controller CC1, CC2 and CC3	**Please note: In the present manual the 3 controller options (CC1, CC2, CC3) are designated as "CCx".**
CCx-E	
CCx-E-BX	
CCx-105A - CCx-118A	Compatible Control immersion thermostats with Polycarbonate bath
CCx-208B - CCx-225B	Heating thermostats with stainless steel bath
CCx-205B	Heating circulation thermostat
CCx-202C	Heating circulation thermostat
CCx-304B vpc	Compatible Control heating thermostats
CCx-308B vpc	Compatible Control heating thermostats
CCx-315B vpc	Compatible Control heating thermostats
CCx-300B2X	Compatible Control bridge heating thermostats
CCx-300B4X	Compatible Control bridge heating thermostats
CCx-K6	Cooling and heating bath thermostats
CCx-K12 - CCx-K25	Cooling and heating bath thermostats
Minichiller 1-3 /-H1	chiller

All units listed in the following table have a suplementary instruction manual. Please read the instructions in the other manual as well.

Compatible Control Thermostats with Controller CC1, CC2 and CC3

Variostat CC	Cooling and heating bath thermostat for irregularly shaped objects
CCx-305 vpc	Cooling and heating thermostats
CCx-405 vpc / CCx-405w vpc	Cooling and heating thermostats
CCx-415 vpc / CCx-415wl vpc	Cooling and heating thermostats
CCx-410wl vpc	Cooling and heating thermostats
CCx-505 vpc / CCx-505wl vpc	Cooling and heating thermostats
CCx-510 vpc / CCx-510w vpc	Cooling and heating thermostats
CCx-515vpc / CCx-515w vpc	Cooling and heating thermostats
CCx-520w vpc	Cooling and heating thermostats
CCx-525w vpc	Cooling and heating thermostats
CCx-705 vpc	Cooling and heating thermostats
CCx-805 vpc	Cooling and heating thermostats
CCx-815 vpc	Cooling and heating thermostats
CCx-820 vpc / CCx-820w vpc	Cooling and heating thermostats
CCx-905 vpc / CCx-905w vpc	Cooling and heating thermostats

Unichiller with Compatible Control CC1, CC2, CC3	
UC006A / UC006B	Circulation coolers
UC015 / UC015w	Circulation coolers
UC022 / UC022w	Circulation coolers
UC020 / UC020w	Circulation coolers
UC025 / UC025w	Circulation coolers
UC006T	Unichiller air cooled bench units
UC009T	Unichiller air cooled bench units
UC012T	Unichiller air cooled bench units
UC016T	Unichiller air cooled tower units
UC020T	Unichiller air cooled tower units
UC025T	Unichiller air cooled tower units
UC040T	Unichiller air cooled tower units
UC045T	Unichiller air cooled tower units
UC055T	Unichiller air cooled tower units
UC060T	Unichiller air cooled tower units
UC080T	Unichiller air cooled tower units
UC006Tw	Unichiller water cooled bench units
UC009Tw	Unichiller water cooled bench units
UC012Tw	Unichiller water cooled bench units
UC015Tw	Unichiller water cooled bench units
UC023Tw	Unichiller water cooled bench units
UC016Tw	Unichiller water cooled tower units
UC020Tw	Unichiller water cooled tower units
UC025Tw	Unichiller water cooled tower units
UC030Tw	Unichiller water cooled tower units
UC040Tw	Unichiller water cooled tower units
UC055Tw	Unichiller water cooled tower units
UC060Tw	Unichiller water cooled tower units
UC080Tw	Unichiller water cooled tower units
Unichiller ETH WL	Circulation cooler to operate under the workbench to be connected to an air-extract duct.
Hotbox HB1	Compact recirculators to maintain the temperature of externally open systems.
Hotbox HB2	Compact recirculators to maintain the temperature of externally open systems.
Hotbox HB3	Compact recirculators to maintain the temperature of externally open systems.
Hotbox HB4	Compact recirculators to maintain the temperature of externally open systems.
Hotbox HB5	Compact recirculators to maintain the temperature of externally open systems.
BFT 1-2	Beer "force ageing" thermostats (with Compatible Control CC2 only)
BFT 1w-2	Beer "force ageing" thermostats (with Compatible Control CC2 only)
BFT 2-2	Beer "force ageing" thermostats (with Compatible Control CC2 only)
BFT 2w-2	Beer "force ageing" thermostats (with Compatible Control CC2 only)

When contacting our service hotline please remember that the machine name and the serial number have to be indicated. This information is to be found on the type label on the rear of your device.

Important note:

Grouping our products into "families" to best present them for an easy overview has always been a priority. Now, with the expansion of the product range and the introduction of new technologies such as "VPC" (Variable Pump Control") we have made some changes to simplify the range to make it easier for you to navigate your way to the best unit for your application.

The "Polystat" and "Compatible Control" ranges have been merged into one family called the "Compatible Control Thermostats". The following table shows the new names for each product together with the old name.

(„x" stands for the controller variants CC**1**, CC**2** und CC**3**)

New names	**Old names**
Compatible Control CCx	(Polystat Control ccx) (**controller**)
CCx-E	(Polystat ccx)(**immersion thermostat**)
CCx-105A - CCx-118A	(Polystat A5 – A18)
CCx-208B - CCx-225B	(Polystat B8 – B25)
CCx-205B	(Polystat 202)
CCx-202C	(Polystat 201)
CCx-304B vpc	(CC 301-x)
CCx-308B vpc	(CC 302-x)
CCx-315B vpc	(CC 303-x)
CCx-300B2X	(CC 302BX-x)
CCx-300B4X	(CC 303BX-x)
CCx-305 vpc	(CC 130-x)
CCx-405 vpc / CCx-405w vpc	(CC241-x / Cc 241w-x)
CCx-415 vpc / CCx-415wl vpc	(CC 240-x / CC 240wl-x)
CCx-410wl vpc	(CC 245wl-x)
CCx-505 vpc / CCx-505wl vpc	(CC 250-x / CC 250wl-x)
CCx-510 vpc / CCx-510w vpc	(CC 150-x / CC 150w-x)
CCx-515 vpc / CCx-515w vpc	(CC 155-x / CC 155w-x)
CCx-520w vpc	(CC 156w-x)
CCx-525w vpc	(CC 157w-x)
CCx-705 vpc	(CC 175-x)
CCx-805 vpc	(CC 180-x)
CCx-815 vpc	(CC 185-x)
CCx-820 vpc / CCx-820w vpc	(CC 181-x / CC 181w-x)
CCx-905 vpc / CCx-905w vpc	(CC 190-x / CC 190w-x)

Contents

Preface

1. Safety Instruction

1.1. Intended Use 11-1-2
General Safety Instructions

2. Device Description

2.1. Structure 21-1-5
2.2. Controller, Remote Control 22-1

3. Commissioning

3.1. Safety Instructions and Principles 31-1-4
3.2. Connecting Main Power and Media 32-1-3
3.3. Filling Thermofluid 33-1-5
3.4. Major Pre-settings 34-1

4. Thermoregulation via Controller

4.1. Safety Instructions and Principles 41-1-2
4.2. Main Menu 42-1
4.3. Editing Default Settings 43-1-10
4.4. Limiting the Thermoregulation Range 44-1
4.5. Calibration 45-1-2
4.6. Utilities 46-1
4.7. Enter a Setpoint-Start 47-1
4.8. Convenient Thermoregulation - Programs 48-1-6
4.9. Terminating the Thermoregulation Process 49-1

5. Shut Down

5.1. Safety Instructions and Principles 51-1
5.2. Draining, Deactivating and Dismantling 52-1
5.3. Maintenance and Service 53-1-2

6. Appendix

6.1. Presettings 61-1
6.2. Interface specification, Data 62-1
6.3. Device Messages 63-1-4
Huber Glossary L-1-8

Attached

Technical data
Spare Parts

1. Safety

1.1. Intended Use

General Safety Instructions

The thermostat is designed for industrial applications.
The thermostat is used for direct and indirect thermoregulation, i.e. for heating or cooling external substances through suitable thermal fluids.
It must be operated strictly in compliance with the operating instructions.
The thermostat must not be modified by the plant operator or any operating personnel.
The thermostat must not be used for purposes other than thermoregulation in compliance with the operating instructions.Unintended use or use not in compliance with the operating instructions may lead to severe personal injury or property damage.

Your device has been designed and constructed according to the state of the art and in compliance with the generally accepted safety rules. Nonetheless, your device may constitute an imminent or unexpected hazard. For this reason, your device has been equipped with safety devices. Deactivating these safety devices bears high risks and invalidates the warranty.

Use the device only if it is in good order and condition.
Shut down the device immediately in the case of malfunctions or failures.
Only qualified personnel is permitted to perform repairs.
Do not bypass, bridge, dismount or deactivate any safety devices.

The manufacturer assumes no liability for damage due to technical modifications, improper handling or use of the device disregarding the operating instructions.
The manufacturer assumes no liability for damage due to technical modifications, improper handling or use of the device disregarding the operating instructions.

Warning! Risk of injuries!

While operating at high temperatures, the bath lid and the housing could become very hot.
Only touch the housing and the lid by the grips otherwise there is a *risk of burns*!

Never lift the bath's lid during operation at high temperature:
- Risk of *scald /burn* through thermal fluid overflow.
Depending on the type of thermal fluid used:
- Risk of *caustic vapours* causing injuries to the respiratory tract and/or skin!
(For further information about the chosen thermal fluid please read the material safety data sheet delivered with it.)

This warning is only applicable for units with this  **warning sign.**

Important: transport damage!
When unpacking the device, inspect it for transport damage.
Please revert to the haulage contractor or shipping agent for settlement of claims.
Commission a damaged device only after the damage has been repaired or you have ascertained the full effects of damage and the insurance agent/haulage contractor/shipping agent has given their permission.

1. Safety

1.2. Intended Use

General Safety Instructions

Duties of the Plant Operator

The operating instructions must be kept readily available in the immediate vicinity of the thermostat.
Only sufficiently qualified operating personnel are permitted to use the thermostat.
The operating personnel must be trained in handling and using the thermostat.
Verify that operating personnel have read and understood the operating instructions.
Precisely define the fields of responsibility of the operating personnel.
Provide protective clothing for the operating personnel.

Requirements to be Met by the Operating Personnel

Only personnel assigned and trained by the plant operator may handle and operate the thermostat.
The minimum age for operating personnel is 16 years.
Within the workspace, the device operator is responsible for third parties.
The device operator must be sufficiently qualified.

Duties of the Operatorating Personnel

The operating personnel must read the operating instructions thoroughly before handling or using the thermostat.
The operating personnel must heed all the safety instructions.
The operating personnel must wear protective clothing when handling or using the thermostat.

Workspace

The workspace is defined to be at the control panel in front of the thermostat.
The workspace is further defined by the peripheral equipment connected by the customer. The customer is responsible for taking suitable safety measures.

Safety Devices

Over-temperature protection device
Low liquid level protection

Protection -Devices

Mains failure automatic
Alarm functions

Emergency Plan – Switch off the Power Supply!

Hazardous emission of fluid/vapor from the thermostat or connected pipes/hoses (very hot, very cold, hazardous chemical compositions) and or fire/explosion/implosion:
Strictly heed the safety instructions of the plant operator relating to the risk of injury and danger to life as well as to the limitation of damage.
Observe the instructions included in the safety data sheet of the respective thermal fluid!

2

Device Description

2. Device Description

2.1. Structure

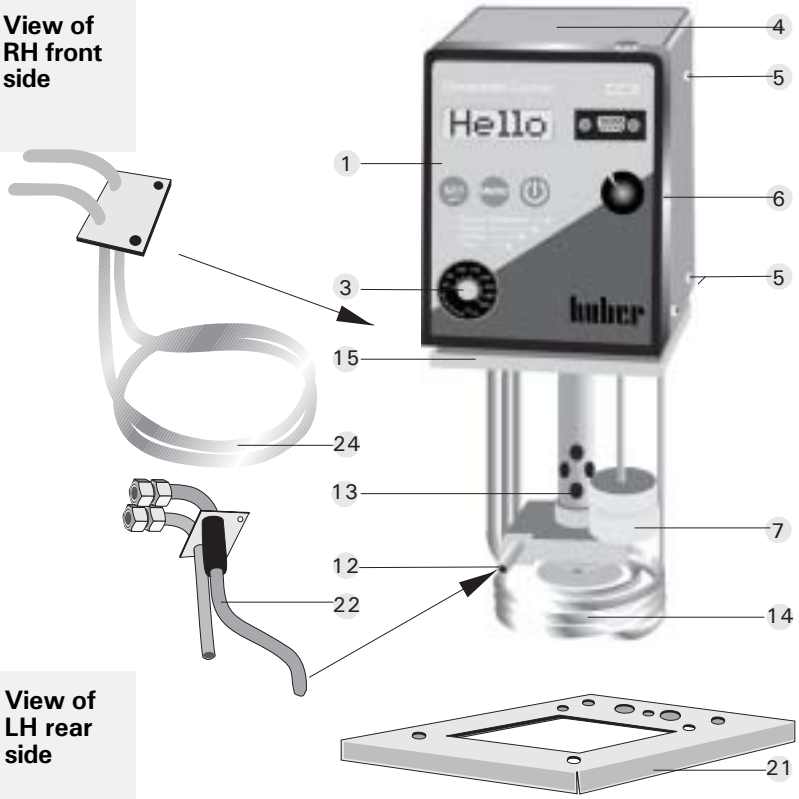
Bath thermostats and circulators used as heating thermostats
Models: compatible Control CC1-3, CC1-3-300B2/4X

1 Removable, exchangeable Compatible Control
2 Power switch
3 Over-temperature protection device
4 Hood
5 Hood fixing screws
6 Controller fixing screws
7 Float (level protection)
8 Internal control sensor
9 Temperature sensor – over-temperature protection device
10 Pressure pump with:
11 Stirrer
12 Hose pressure connection for circulating thermal fluid from the thermostat to the external device
13 Inlet for the pressure pump
14 Heater
15 Base plate
16 Power supply
17 Socket for cooling bath control
20 ID plate
21 Bath bridge #6309*
22 M16 x 1 Pump adapter for lateral attachment (hose connecting piece)*, for additional external thermoregulation
23 M16 x 1 Pump adapter for rear attachment (hose connecting piece)*, for additional external thermoregulation
24 Cooling coil for mounting with lateral cooling water inlet and outlet connections*
25 Cooling coil for mounting with rear cooling water inlet and outlet connectors*
26 Screw clamp for fastening the Polystat cc1-3 to a bath tank.

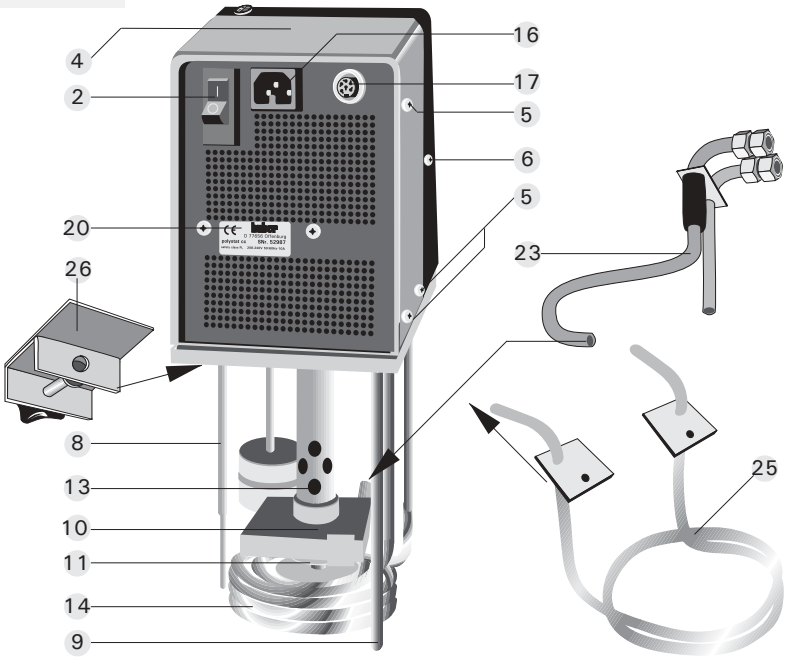
Within the scope of technological progress and practicality, we reserve the right to change the location of individual components.
*Can be ordered optionally for Compatible Control CC1-3-E; for other accessories, please refer to the product catalog.

The devices heat up to 200 °C.
They are suitable for cooling with cooling water and through an externally connected cooler to -30 °C.

View of RH front side

View of LH rear side

2. Device Description

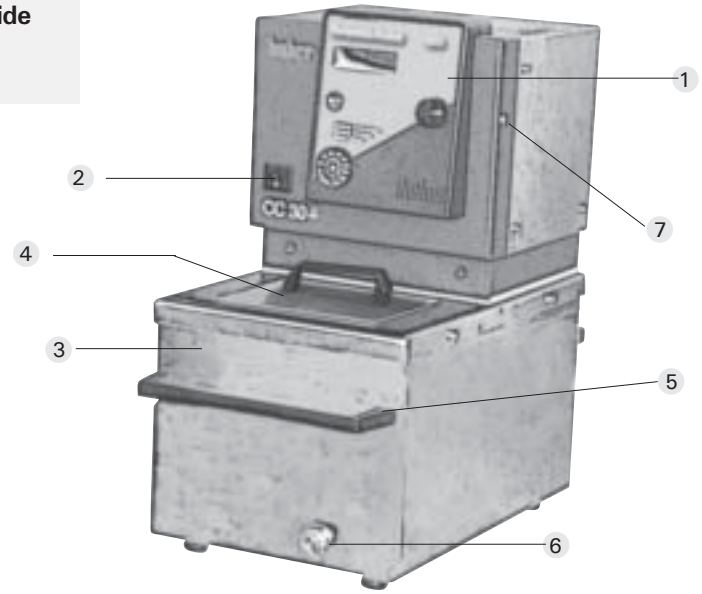
2.1. Set-Up

Compatible Control Thermostat CCx-304B vpc
Bath and Circulation Thermostat

Working Temperature Range (-20) 28...300°C.
With Refrigerator -20...300°C.

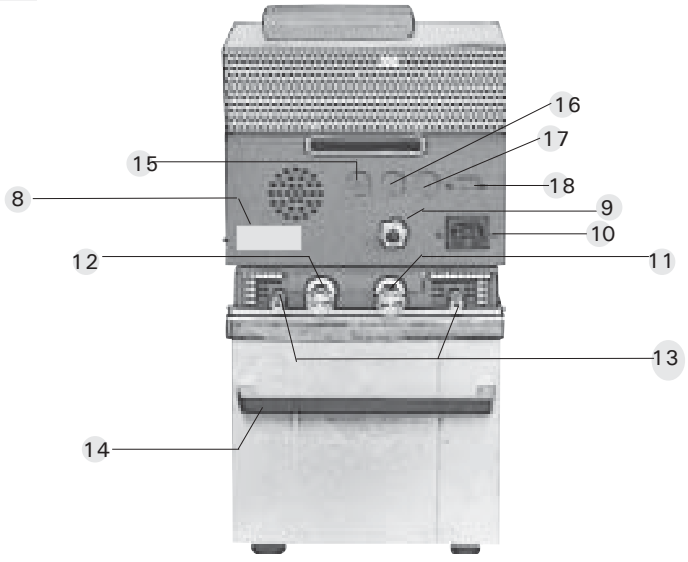
1 Removable, exchangeable Compatible Control
2 Main switch
3 Bath
4 Bath lid
5 Front grip
6 Bath outflow
7 Fixing screw of the controller

View of RH front side

8 ID plate
9 CAN-Interface
10 Power supply
11 Flow (circulation)
12 Return (circulation)
13 Connection for cooling coil
14 Back grip
15 POKO (Potential-Free Contact)
16 AIF (Analog Interface)
17 Enable
18 RS232/485

View of rear side

2. Device Description

2.1. Structure

Unichiller (chiller)

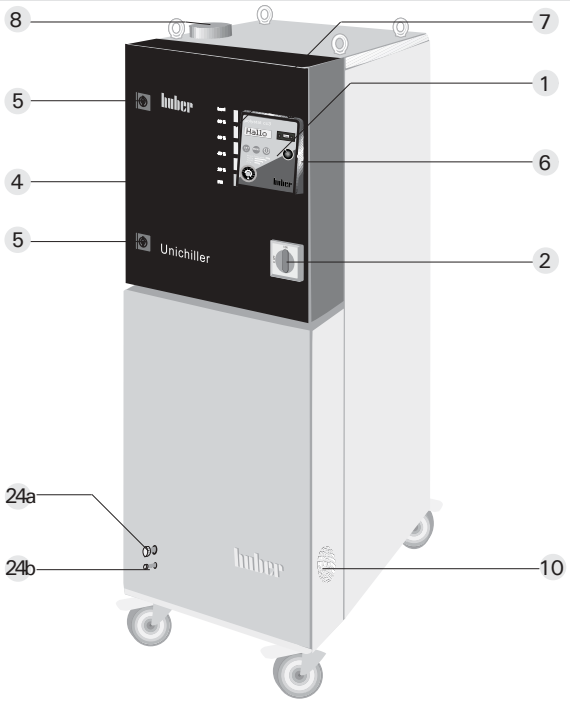
Models: UC016T-3(-H), UC020T-3(-H), UC025T-3(-H), UC040T-3(-H)

The devices cool to -10 °C/-20°C with air-cooled chiller H-devices are additionally equipped with a heater. H-devices heat up to 100 °C.

1 Removable, exchangeable Compatible Control
2 Main switch
3 Overheating protection (not visible)
4 Electric system box with swing door
5 Lock of the electric system box door
Keep the key at a safe place!
6 Controller attachment screws (on either side)
7 Level indicator
8 Expansion vessel with filler connection

10 Air supply for pump motor

View of RH front side

20 ID plate
21 Power supply without plug

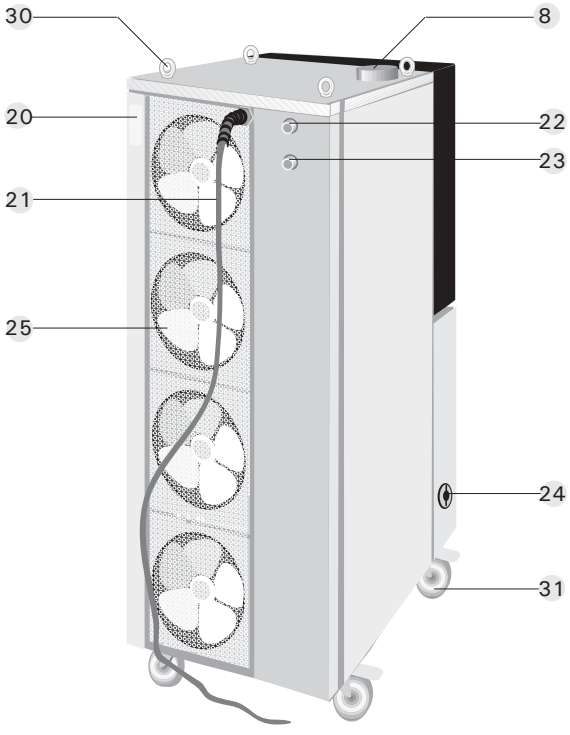
22 Hose connection for circulating thermal fluid from the thermostat to the external application

23 Hose connection for circulating thermal fluid from the external application to the thermostat
24 Thermal fluid drain valve:
24a fluid circulation pump draining connection
24b for the evaporator draining connection
25 Air grill with forced ventilation behind it

30 Ring bolt (detachable eyebolt)

31 Lockable wheels

View of LH rear side

Within the scope of technological progress and practicallity, we reserve the right to change the location of individual components and to modify the appearance of individual components.

2. Device Description

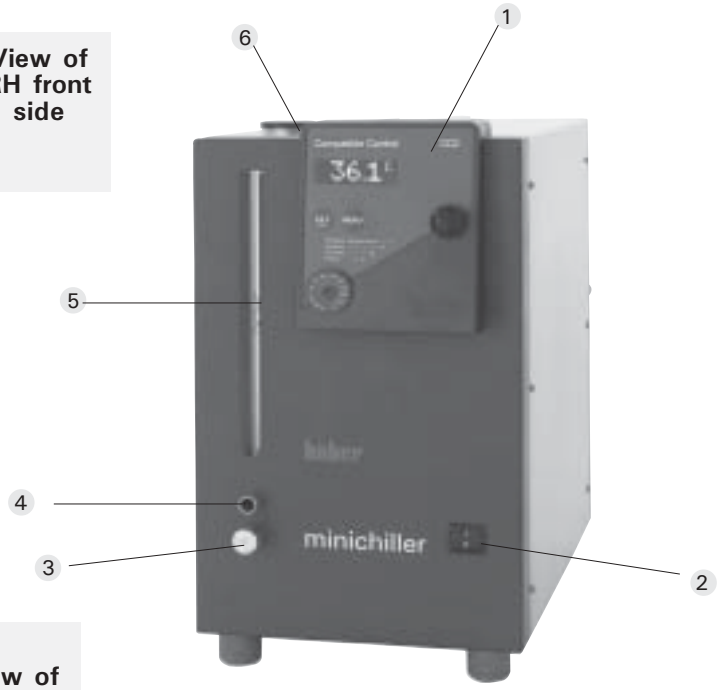
2.1. Set-Up

Circulation Thermostat
minichiller 1-3 /-H1
minichiller w-1-3

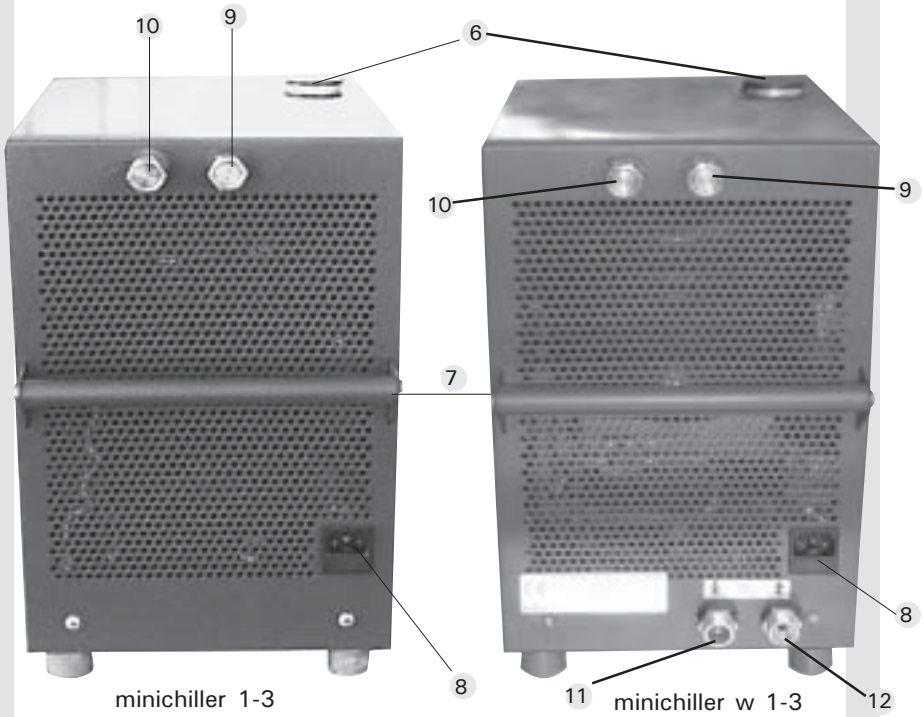
Working Temperature Range minichiller 1-3: **-20...40°C.**
Working Temperature Range minichiller 1-3 with H1: **-20...100°C**
Working Temperature Range minichiller w-1-3: **-20...40°C**

1. Removable, exchangeable Compatible Control
2. Main switch
3. Draining connection
4. Overflow connection
5. level indicator
6. Filler connection

View of RH front side

View of rear side

7. Back grip
8. Power supply
9. Return (circulation output)
10. Flow (circulation input)
11. Cooling water inlet
12. Cooling water outlet

2. Device Description

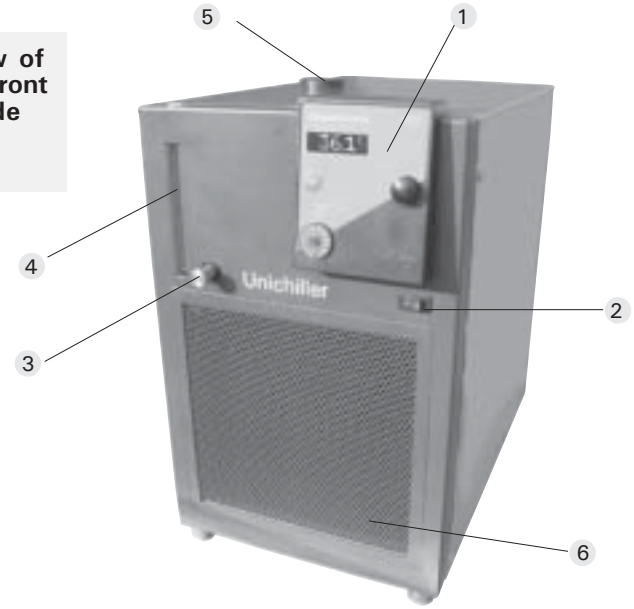
2.1. Set-Up

Circulation Thermostat
UC007-1-3 /-H

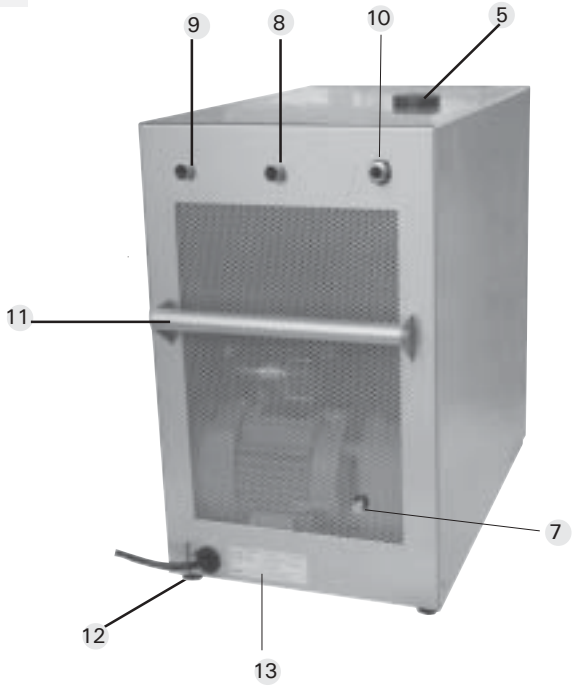
Working Temperature Range UC007-1-3 -20...40°C
UC007-1-3 mit H: -20...80°C

1 Removable, exchangeable Compatible Control
2 Main switch
3 Draining connection
4 Level indicator
5 Filler connection
6 Guard

View of RH front side

View of rear side

7 Draining pump
8 Overflow connection
9 Flow (circulation input)
10 Return (circulation output)
11 Back grip
12 Power supply
13 ID plate

2. Device Description

2.2. Controller

Compatible Control CC1
Compatible Control CC2
Compatible Control CC3

Displays, control panel

R1 Digital status display

R2 Prompt for a temperature set-point, input confirmation (data transfer)

R3 Call of the user menu for convenient handling

R4 On/Off key to start/stop thermoregulation

R5 Activity indicator LEDs

R6 Over-temperature protection device

R7 Encoder/ rotate: Entry of program numbers, step numbers, parameters (e.g. temperature set-point)

Encoder/ press: Input confirmation (data transfer)

Connections

R8 Pt100 sensor socket

R9 RS 232/485 interface

R10 Fixing screw (on either side) for connecting the Compatible Control to the thermostat

Activity indicator LEDs

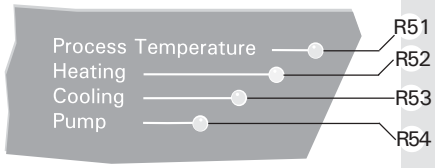
Activity indicator LEDs

R51 Process temperature control active (green LED), only with CC2, CC3:
Temperature is measured by a process sensor located at the point of control e.g. in a reactor.

R52 Heating active (yellow LED)

R53 Cooling/compressor active (blue LED)

R54 Pump active (green LED)
Thermal fluid is being pumped through the connected application, e.g. around the reactor jacket.

3

Commissioning

3. Commissioning

3.1. Safety Instructions and Principles

Preparatory Measures for Commissioning

Plan the thermoregulation target and procedure.
Determine the device configuration and system structure.
Select an appropriate thermal fluid.
Selection criteria for thermal fluids: Temperature range of the thermostat, application restrictions building safety regulations, your projected working temperature, viscosity, flash point.Position the thermostat and external devices on a stable and even surface. Ensure that the surface can safely hold the weight of the thermostat and external equipment.

Positioning

Caution! Potential risk of injury and material damage:
Keep the device upright during transport.
Place the device in an upright and stable position and make sure that it cannot tilt over.
Keep the vicinity of the device clean: Prevent slip and tilting hazards.
Lock the wheels of floor-mounted devices once you have positioned them as desired!
Useful note: Place appropriately large collecting trays under the thermostat and the application.

Location

The use of the device is only allowed in normal environment conditions according to DIN EN 61010-1:2001:
- use only in the inside;
- Hight up to 2000m;
- Place on a dense, level, non-slip and non-flammable surface.
- Place the device at a distance to walls and the ceiling that permits sufficient air circulation (heat dissipation, supply of fresh air for the thermostat and the workspace). A water-cooled unit requires a minimum of 10 cm and an air-cooled unit needs a minimum of 20 cm.
- Ambient temperature min- 5°C to max. 32 °C;
- Short distances to the supplies;
- The device souldn't be placed such that the access to the mains could be obstructed;
- Mains voltage fluctuation do not exceed $\pm$ 10% of the stated input voltage;
- Transient overvoltage like they usualy appear in the mains supply;
- Applicable degree of pollution: 2;
- Class of overvoltage: II.

The workspace of the thermostat must comply with local safety regulations.

Please Note:
All the safety instructions are vital and must be considered during the operation of the unit in compliance with the operating instructions.

3. Commissioning

3.1. Safety Instructions and Principles

The operating instructions contain additional safety instructions. These are identified through a triangle with an exclamation mark. Thoroughly read and heed the instructions Non-observance may involve considerable consequences such as device damage, physical damage or personal injury with fatal consequences.

Workspace
The workspace is defined to be at the control panel in front of the thermostat.
The workspace is further defined by the peripheral equipment connected by the customer.
The customer is responsible for taking suitable safety measures.

Safety devices

^ Overtemperature protection
^ Low liquid level protection
^ Mains failure automatic
^ Alarm functions.

Hazardous emission of fluid/vapour from the thermostat or connected pipes/hoses (very hot, very cold, hazardous chemical compositions) and or fire/explosion/implosion:
Strictly heed the safety instructions of the plant operator relating to the risk of injury and danger to life as well as to the limitation of damage.
Observe the instructions included in the safety data sheet of the respective Thermofluid!
Classification according to DIN12876:

Class designation	Fluid for thermoregulation	Technical requirements	Designation[d]
I	Non-combustible [a]	Overheating protection [c]	NFL
II	Combustible [b]	Adjustable overheating protection	FL
III	Combustible [b]	Adjustable overheating protection and additional low level protection	FL

a Generally water, other fluids only if they are not combustible in the temperature range of an individual fault.
b The fluid for thermoregulation must have a burning point of $\geq$65°C; this means that ethanol can only be used under supervision.
c The overheating protection can be achieved e.g. through an appropriate fill level sensor or appropriate temperature control devices.
d Optional, depending on the manufacturer's choice.

Your thermostat complies with the class designation III / FL

We cannot guarantee the suitability of our components (individually or combined) for use in a potentially explosive environment.
Not suitable for use as a medical appliance (e.g. in vitro diagnosis methods).

3. Commissioning

3.1. Safety Instructions and Principles

Thermal fluid

Only non-flammable thermal fluids are permissible for unattended operation!
Requirements for thermal fluids classified as III / FL:
EN 61010-1: Maximum permissible working temperature 25°C below the flash point!
Maximum viscosity at the lowest working temperature: 50 mm^2/s.
Maximum density of the thermal fluid: 1 kg/dm^3.
Possible thermoregulation range within the range of the planned minimum and maximum working temperature.

Do not use thermal fluids with any of the additives ether, ester or amine. Do not use demineralized water, mineral water, sea water or $CaCl_2$ brines!
Compatibility with the materials of the thermostat (stainless steel 1.4301 (V2A) and Viton) and all materials of the connected application.

Hazards during thermo-regulation

Caution! Potential risk of injury and material damage during thermoregulation:
In the course of operation of the thermostat, extreme changes in temperature and pressure and the specific characteristics of the thermal fluids used may constitute hazards.

3. Commissioning

3.1. Safety Instructions and Principles

Hazards through emission of fluids

Caution! Potential risk of injury and material damage:
The floor will be slippery when fluids have been spilt!
Thermal fluids with a low flash point constitute a fire hazard!
Hazard of scalding/burning when touching exposed or defective connections that are hot.

Prevent overflow of the sight glass. (A site glass is not fitted to open-bath units. This refers to the UCO range and other Chillers)
Prevent leaking fluid-conveying pipes/tubes and connections.
Always remove any liquids spilt on the floor immediately.
Always clean contaminated devices immediately.
Place an appropriately large collecting tray under your external application.

Hazards through emission of gases

Caution! Potential risk of injury and material damage:

Risk of causticization of your respiratory tracts and skin through vapors!
Prevent leaks on closed external devices.
Ensure good aeration and ventilation in the vicinity of the thermostat.
Choose thermal fluids for thermoregulation that are not detrimental to health.

3. Commissioning

3.2. Connecting Mains Power and Media

Release the Transport securing srcrews

Table machines do not have transport securing screws.

Connecting Power

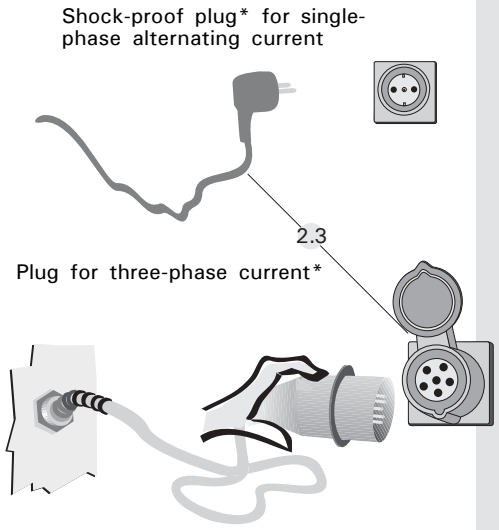
2. Check the fuse, power and voltage ratings according to the Data Sheet (Appendix) and the type plate (on the rear of the device).

2.2 Thermostats designed for use with three-phase current are supplied without a power plug.
In the country of use, an appropriately trained and qualified electrician must fit the power cable with a suitable power plug.

Caution!
Please check the rotating field of your mains supply before connecting the plug.
Your thermostat must be wired in accordance with a dextrorotatory field and the cable markings. Incorrect connection will lead to false results when working with the thermostat!
Interchange the two phases in your power outlet
to correct the dextrorotatory field. (Only an appropriately trained and qualified electrician may perform this work!)

2.3 Connect the power plug to the power outlet.

* Figure applies in the country of manufacture (Germany) only.

Warning! Hazard through electric current!
Exclusively appropriately trained electricians are permitted to replace plugs, change the voltage in compliance with the Data Sheet and perform permanent installations.
Otherwise there is a risk of injury and danger to life.

3. Commissioning

3.2. Preparing the thermostat for use

Applies to bath thermostats with external pump connection

For internal thermoregulation, please proceed as follows:
close the pump connections with plugs and sleeve nuts (M16x1 / 19 spanner/wrench), at the same time hold against it at the pump connection with 17 spanner/wrench.

For external thermoregulation, please proceed as follows:
remove the plugs and sleeve nuts on the pump connections. Replace them with appropriate hose connection to your external application.

For more information please see chapter 4 „thermoregulation via controller".

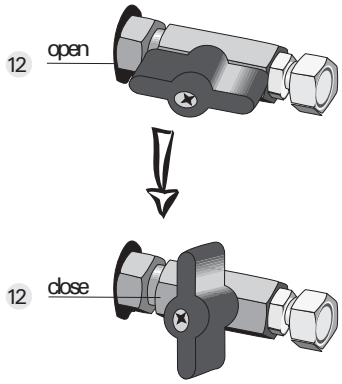
Thermal Fluid
Preparatory measures for thermal fluid circulation

Applies to chillers (thermostats with closed external circulation)

2.4 Make sure that the draining valve is closed. (upright position or „right stop").

2.5 Connect hoses at connecting pieces 10 and 11 for circulating thermal fluid to and from the external device (if fitted).

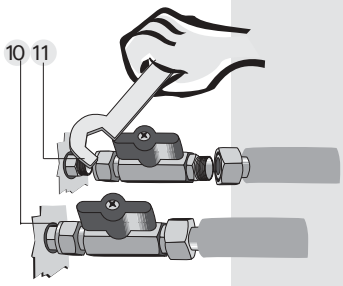
Note:
Recommended tightening torques for hose couplings when using silicone oil as thermal fluid:
M16 x 1 30 Nm
M30 x 1.5 85 Nm

Alternative

Hoses can be connected via shut-off valves. (Please refer to your Huber catalogue under „accessories")

Hazard!
If fitted, close valves only when performing work on the reactor. Remember that thermal fluid expands and contracts with changes in temperature. Sealing the external application will expose the application to these forces!

Control

Verify the following:
Make sure that all connections are correct and that there are no leaks.
For water-cooled thermostats:
Unobstructed inflow and outflow of the cooling-water exists
Pressure difference between the cooling-water inflow and outflow is 3 to 5 bar. Cooling-water supply temperature is 15 °C to 25 °C.

3. Commissioning

3.2. Preparing the thermostat for use

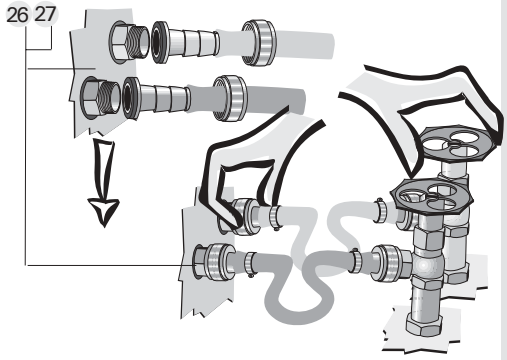
Media

Preparatory measures for water cooling

Applies to water-cooled thermostats only

4.1 Using suitable hosing, connect the Huber water-cooled unit to the cooling water supply via connections 26 and 27

4.2 Open all the water supply valve slowly and check for leaks.

Attention!
When working **above 40°C** the **cooling water** supply should be **turned off** to avoid wastful counter cooling!

Preparatory measures for cooling baths **CCx-K12 bis CCx-K25**

1 Mains supply *
(from Compatible Control CC-E to cooling bath)

2 Control connection *
(from Compatible Control CC-E to cooling bath)

3 Main switch
I = On (continuous duty of the cooling machine)

0 = Off

II = „Auto." the compressor is automatically controlled

4 Mains supply
(fuse protected)

Rear side of Compatible Control CC-E (immersion thermostat)

Rear side of CC-K12

* These connectors have to be plugged in even if the cooling bath is not being automatically controlled, otherwise the error message „SENSOR COND.!!" will be displayed (ref. to error code -11, page 63-3).

3. Commissioning

3.3. Filling Thermofluid

Overtemperature protection

Require ments

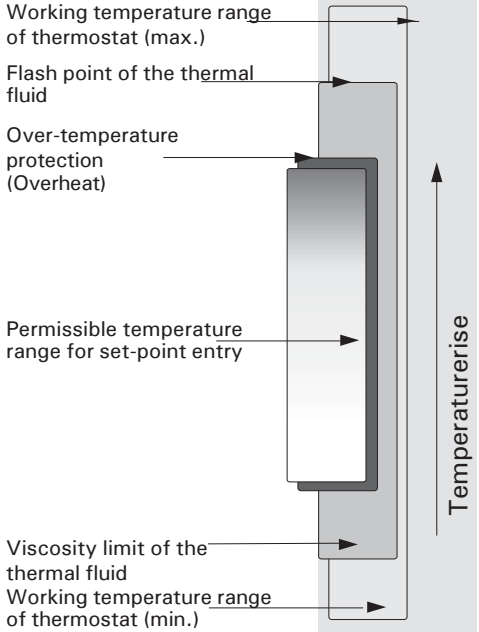
Prepare the thermostat for thermoregulation and take safety measures as described below.

Setting the overtemperature protection device.

Requirements:

A suitable thermal fluid has been selected for the process requirements.
The flash point of the thermal fluid is known.

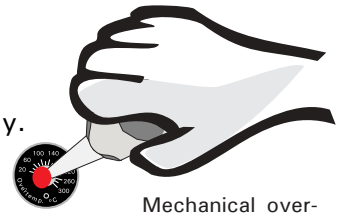
Procedure:
The over-temperature protection is set to at least 25 °C below the flash point of the thermofluid.

Caution!
The overtemperature protection is an especially important safety device of your thermostat. It should always be operable and be tested periodically!

Mechanical over-temperature protection
CC1-E, CC2-E, CC3-E, 105A to 118A 208B to 225B, 202C,205B, K6 to K25

2 For Compatible Control CC1-E, CC2-E, CC3-E (immersion thermostat), CCx-E with polycarbonate bath 105A to 118A, heating bath thermostats 208B to 225B, heating thermostat 202C and 205B, cooling heating bath K6 to K25,the over-temperature protection is set mechanically. It is independent of the controller.
Use a suitable tool (screw driver or the like) to *turn* the red button in the center of the over-temperature protection device (section 2.2) to the desired temperature value.

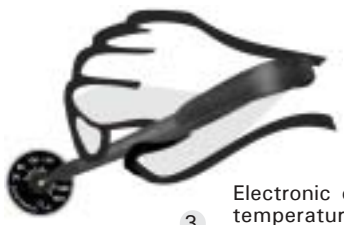
2 Mechanical over-temperature protection

Electronic over-temperature protection
for all units not listed above

3 For all other units not listed above, the over-temperature protection is set electronically. It is independent of the controller.
Use a suitable tool (screw driver or the like) to *press* the button in the center of the over-temperature protection device.
The over temperature protection menu will be displayed.

```
Over temp. protection
   OT Setpoint
   OT temp . display
   Exit
```

OT setpoint : over temp. value (adjustable on the controller)
OT temperature display: display of the over temperature actual value

3 Electronic over-temperature protection

3. Commissioning

3.3. Filling Thermofluid

Bath and heating thermostat
CCx-304B vpc

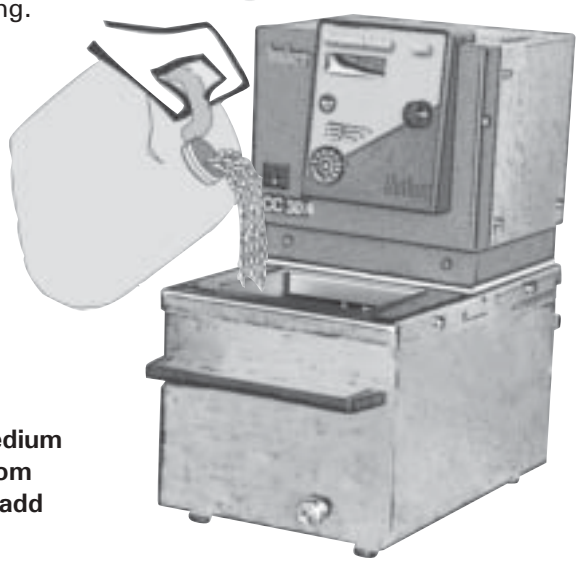
Filling the bath

L2 Lift the bath lid.
Fill in the thermofluid.

Optimum liquid level:
3 cm below bath opening.

Close the bath lid.

L2

Please note that the medium should cool down to room temperature bevor you add thermofluid!

Warning! Risk of injuries!

While operating at high temperatures, the bath lid and the housing could become very hot.
Only touch the housing and the lid by the grips otherwise there is a *risk of burns*!

Never lift the bath's lid during operation at high temperature:
- Risk of *scald /burn* through thermal fluid overflow.

- Risk of *caustic vapours* causing injuries to the respiratory tract and/ or skin!

(For further information about the chosen thermal fluid please read the material safety data sheet delivered with it.)

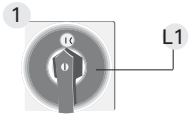
3. Commissioning

3.3. Filling Thermofluid

Unichiller (Chiller)
Models: UC016T-3(-H), UC020T-3(-H), UC025T-3(-H), UC040T-3(-H) etc.

Switch on

Condition:
make a hose connection to your external application.

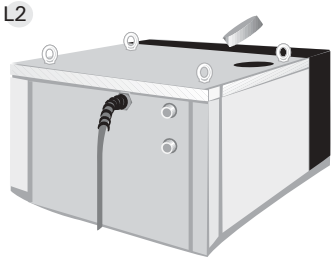
L1 Switch on the main switch.

Open, Filling

UC016T-3, UC020T-3, UC025T-3, UC040T-3 etc. also all UCs with heating(-H)

L2 Lift the bath lid.

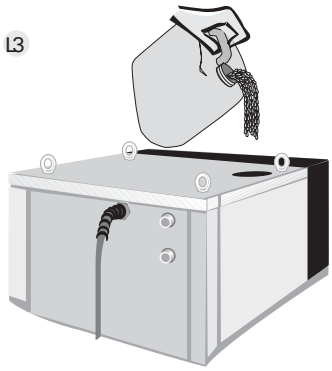
L3 Fill in the Thermofluid.
The thermofluid flows from the expansion vessel into the thermostat than through the hose connections to the external application.

L4 Support the distribution process with the program „Venting“. (46-1)

Continue filling procedure and „Venting“, until Gases are fully displaced off the System and the level on the luminescent display stays constant at 20-40%.

Optimal level:
20-40%.

Terminate filling, close

L5 The liquid level on the luminescent display stays constant on reaching of the optimal level.

L6 Stop „Venting“ Program.

L7 Close the bath lid.

Caution! Potencial risk of injury and material damage:
A proper venting is the condition for a troublefree operation!
Commissioning often means a change of liquid.
The manufacturer mostly uses other thermofluids as the user during the test run. Despite thorough draining, smallest residual could remain in the circuit and degasify after filling a new fluid and the beginning of the thermoregulation.
The concequence of the above described and an incomplete venting could be an increasing presure *that leads to sudden discharge of fluid and vapours*!
This danger can be prevented to a large extend with the „Venting“ program
Change of liquid refer to 46-1., program „venting“!
Absolutely pay attention to a thorough degasification at the beginning of the thermoregulation!

Please note!
For detailed information on the range of thermofluids incl. technical data and price quotations, please refer to the Huber catalog or contact your Huber agent.

3. Commissioning

3.3. Filling Thermofluid

Circulation thermostat
minichiller 1-3 /-H1
minichiller w-1-3

Filling the thermostat

Please note:
The thermostat is fitted out with an overfill protection.
when the fluid level exceeds the max mark in the sight-glass, the fluid begins to outflow from the overflow connection (4).
An overflow of the internal tank (wich could flood the machine internals) can thus be avoided.

Warning:
The overflow connection (4) should never be closed, the fluid should be able to flow out without hindrance. However it is recommended with frequent change of the liquid or the application and/or with the first commissioning to put a hose on the overflow connection and to place the other end into a suitable receiver.

Make a hose connection to your external application and make sure that everything is leakproof.
Take off the lid of the filler connection (6) and slowly fill the thermal fluid.
The fluid level should allways be between MIN and MAX on the level indicator.
Switch the unit on with the main switch (2).
The thermostat starts circulation and the fluid is pumped into your application. If there is no fluid in the application yet, the level in the Minichiller will lower continualy, please refill with fluid untill the level stays constant.

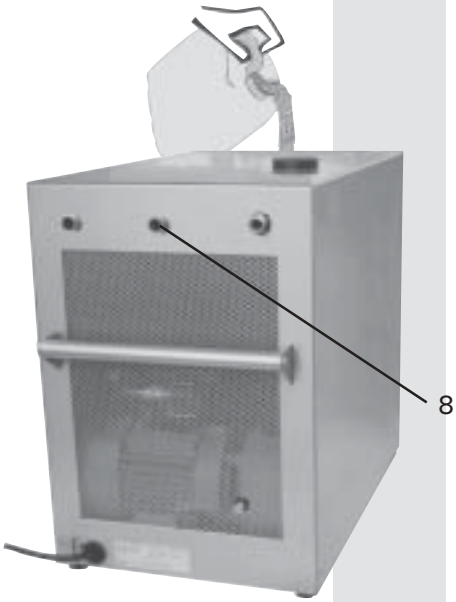
3. Commissioning

3.3. Filling Thermofluid

Circulation thermostat
Unichiller CC007 1-3

Filling the thermostat

Please note:
The thermostat is fitted out with an overfill protection.
when the fluid level exceeds the max mark in the sight-glass, the fluid begins to outflow from the overflow connection (8).
An overflow of the internal tank (wich could flood the machine internals) can thus be avoided.

Warning:
The overflow connection should never be closed, the fluid should be able to flow out without hindrance. However it is recommended with frequent change of the liquid or the application and/or with the first commissioning to put a hose on the overflow connection and to place the other end into a suitable receiver.

Make a hose connection to your external application and make sure that everything is leakproof.
Take off the lid of the filler connection (5) and slowly fill the thermal fluid.
The fluid level should allways be between MIN and MAX on the level indicator.
Switch the unit on with the main switch (2).
The thermostat starts circulation and the fluid is pumped into your application. If there is no fluid in the application yet, the level in the Unichiller will lower continualy, please refill with fluid untill the level stays constant.

3. Commissioning

3.4. Major Pre-settings

Language

Deutsch

When delivered the controller displays will be in German.

Other options can be selected in the „Language“ menu (refer to 4.3.2): English, Francais

Setpoint

The thermostat controls the temperature to the predefined setpoint of 20 °C. Use the SET key and the encoder to select other options (refer to 4.1 and 4.7): Any temperature value within the performance parameters of the thermostat (refer to the Data Sheet 2.4 Appendix, type plate!) and the safety limits (refer to 3.1!).

Minimum set-point

Maximum set-point

The permissible setpoints are limited to the range of the minimum and maximum setpoints.

The option „setpoint limit“ allows this range to be modified.

Notes to temperature range limit:
- Do not select the limits greater than needed.
- Observe the characteristics of the thermal fluid (viscosity, flashpoint).
- Consider the defaults of your application’s manufacturer.

The programs for operating the thermostat are described in Chapter 4 of the Operating Instructions.

4

Thermoregulation via Controller

At work

4. Thermoregulation via controller

4.1. Safety Instructions and Principles

Checks to be Performed: Test Arrangement

Compare the device configuration, the system structure and the selected thermofluid to the thermoregulation target.
Verify the stability of the thermostat and external devices.
Make sure there are no leaking connections.
Make sure the shut-off valves for thermofluid and cooling water (if applicable) have been opened.
Check the connection to the power supply.

Caution! Potential risk of injury and material damage:
Slip hazard! Floor and workspace may be slippery (spilt fluids)!
Tilt hazard! check stable position of the thermostat.
Shock hazard! check connection to the power supply
Scalding and burning hazard! extreme temperatures!
Causticization hazard! Risk of causticization of your eyes, your skin and your respiratory tracts through hazardous vapors.

Safety Measures

Setting the overtemperature protection:
(Applies to Compatible Control Thermostats):
When: Immediately after filling the system with thermofluid!
Recommendation:
Periodically test the function of the overtemperature protection by entering a higher setpoint, e.g. as follows:
Set the overtemperature protection to 30 °C.
Enter the maximum setpoint of 40 °C (since the overtemperature protection is independent of the controller, you can enter this excessive value without problems).
Enter the new setpoint of 33 °C (former setpoint 20 °C); the heating will heat to 33 °C, then the thermoregulation process will be stopped automatically with the error message: „Temp“
The error message will be displayed until the error has been remedied.
Remedy the error. Set the overtemperature protection to a temperature above 40 °C.Switch the thermostat off and then on again.

Entering the minimum and maximum setpoints:
In combination with the working temperature range limits, the minimum and maximum setpoints provide additional safety for the thermoregulation process. This means, accidental entry of a setpoint that is too low or too high will be rejected.

Low liquid level protection:
Monitor the liquid level during operation.
Applies to bath thermostats (Compatible Control Thermostats):
Fill level to approx. 60 – 80 % of the bath height;
for chillers: fill level to approx. 60 – 80 % on the level indicator.
Thermofluid level too low: Risk of the thermostat pump running dry.
The controller will report an error and stop the thermoregulation process.
Thermofluid level too high: Overflow, soiling, slip and causticization hazard!

Change of fluid:
Rinsing fluid and thermofluids come into contact with stainless steel (V2A), Viton and Perbunan and must be compatible with these materials.

Room ventilation:
Sufficient aeration and venting in the vicinity of the thermostat minimizes the risk of overheating and the accumulation of harmful gases and vapors.

Please Note: All the safety instructions are vital and must thus be considered on the job in compliance with the present operating instructions.

4. Thermoregulation via Controller

4.1. Safety Instructions and Principles

Principles of displays and entries

Salutation

Flash displays on the controller

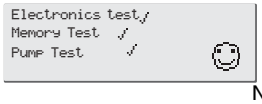
Manufacturer

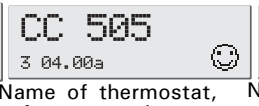
Salutation

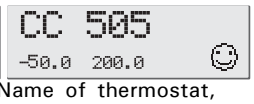
System test

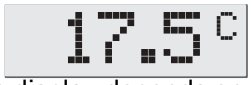
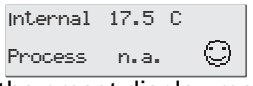
Name of thermostat, software version

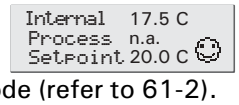
Name of thermostat, working temperature range

Status display

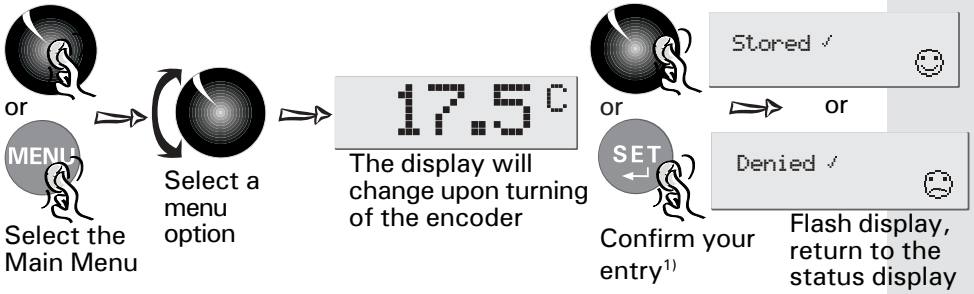
The status display depends on the preset display mode (refer to 61-2).

Operating the controller

[1)] If you fail to confirm your entry within 4 minutes, your selection will not be saved. The program returns to the status display.
To exit the menu instantaneously at any point, use the **Break function***: Press the SET and MENU keys simultaneously. Your selection will not be saved.

Messages

During operation, ad-hoc messages may be displayed on the controller. They provide information on irregularities and hazards in the thermostat. In the case of imminent danger, the controller will display a message and stop the thermoregulation process/switch off the thermostat at the same time.

On/Off

*Press the On/Off key to start/stop the thermoregulation process (thermoregulation combined with circulation).
Alternative:
MasterClear function**: Pressing the MENU and TEMP keys simultaneously stops a thermoregulation program in progress and switches off the analog interface or switches from the digital interface

Note:
To learn more about the menu and the individual menu options, please also read sections 4.2 to 4.9 of the present operating instructions.

*applies to Compatible Control CC2 and CC3 only
**applies to Compatible Control CC3 only

4. Thermoregulation via Controller

4.2. Main Menu

Contents

The main menu provides menu options and submenu options including all the settings and selections required to operate the thermostat. (Refer to 42-1!)

[1] This menu option is not available for Compatible Control CC1.
[2] This menu option is not available for Compatible Control CC2.

Selection

MENU

Page	Main Menu	Submenus
46-1	I⇨ 2nd Setpoint [1), 2)]	
43-3	Alarm Clear	
43-3	Alarm Confie.	I⇨ ALARM CONFIGURATION: I⇨ Alarm Mode Lower Alarm Limit Upper Alarm Limit Level Alarm Delay Exit
43-6,7	Analog-Interface [1), 2)]	
43-2	Display	I⇨ Display I⇨ Display modes optimise display Options Exit
43-8	Digit. Interface [1), 2)]	I⇨ Digit. Interface : I⇨ Hardware RS Baud rate Protocol slave address Exit
46-1	Venting	
43-5	Machine Options	
43-4	Max. Heat Power	
45-1	Calibration Proc.	I⇨ 2-P.CALIBR.: I⇨ Edit TCal1 Edit TCal2 Control to TCal1 Control to TCal2 Exit
43-4	Compressor Auto	
43-2	Mains Failure Auto	
45-2	Offset Calibration	I⇨ OFFSET CALIBRATION: I⇨ internal sensor process sensor Exit
43-9	PI-Parameters	
48-2_5	Edit Program [1)]	
48-6	Program Start/Stop [1)]	
48-1	Start Ramp [1)]	
43-10	Acoustic Alarm	
43-4	Software version	
44-1	Setpoint Limits	
43-1	Language	
43-1	Temperature Scale	
43-1	Control Mode [1), 2)]	
46-1	Circulation	
43-10	Select Usermenu	
43-10	Confie Usermenu	
43-1	Factory Default	I⇨ INITIALISE: I⇨ Unit Data User menus Programmer All together Exit
43-2	Time Scale	
	Exit	

4. Thermoregulation via Controller

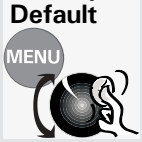
4.3. Editing Default Settings

1. Restoring the factory default
2. Selecting a language for the controller display
3. Temperature scale
4. Control mode – internal, process (cascade)

All the factory defaults can be customized in the menu.

All the factory defaults can be restored via the „Factory Default“ option (with the thermoregulation function switched off). For this purpose, proceed as described below. (Also refer to „Salutation“ 4.1).

Factory Default

```
MAIN MENU:
   Config User
I->Factory default
   Time Scale
```

„Factory Default“ option

```
INITIALISE:
I-> Unit Data
    Programmer
    User menus
    All together
    Exit
```

Submenu selection:
„Unit Data“: Important if you replaced device components or accessories.

„Programmer“
Deletion of all the thermoregulation programs incl. ramps entered by the customer.

```
  INITIALISE:
    User menus
I-> All together
    Exit
```

„User menus“
Restores the „User menu“ factory default.
„All together“
Restores all the factory defaults.
Caution! All the thermoregulation programs incl. ramps entered by the customer will be deleted!
„Exit“ option:
The factory defaults will not be restored.

Language

```
MAIN MENU:
   Set-point limits
I-> Language
   Temperature Scale
```

```
Language:
I-> Deutsch
   English
   Francais
```

For selecting the language that is to be displayed on the controller.

Temperature Scale

```
MAIN MENU:
   Language
I-> Temperature Scale
```

```
Temperature Scale:
I-> Celsius
   Fahrenheit
   Kelvin
```

For selecting the temperature unit for display on the controller.

Control Mode

```
MAIN MENU:
   Temperature Scale:
I-> Control mode
```

Not available for Compatible Control CC1!

```
Control mode:
I->Internal
   Process (Cascade)
```

For selecting the control mode, either Internal or Process (Cascade).

For definitions of internal and external thermoregulation, please refer to Appendix 5.1. Huber Glossary, keyword Control Mode – Internal, Process (Cascade)!

4. Thermoregulation via Controller

4.3. Editing Default Settings

5. Display
6. Time scale
7. Mains failure auto

Display

```
MAIN MENU:
   Digit. Interface
I→ Display
   Edit Program
```

```
DISPLAY:
I→Display mode
   optimise display
   display
   Exit
```

```
DISPLAY MODES:
  Standard
I→ Double
   Double1+Setp.
   Double1+OT
   Options
   Service1
   Service2
   Service3
   exit
```

„Display“ option

Selection in the „Display Modes“ submenu:
„Standard“:
Single-line status display, the actual value of the controller temperature (according to the control mode (internal or process temperature) is displayed in maximum font size.
„Double“:
Double-line status display, the actual values of the internal temperature and the (external) process temperature are displayed in medium font size.
„Double + Setp“:
Three-line status display, the actual values of the internal temperature and the process temperature as well as the setpoint are displayed in small font size.
„Double + OT“:
Three-line status display, the actual values of the internal temperature and the process temperature as well as the overtemperature setpoint are displayed in small font size.
Options:
Additionally in the third line of the display, depending on the unit and the mode of operation, the pump speed or the resulting value from the analogue output will be displayed.
„Service1“ , „Service2“ and „Service3“:
Service related variable status display.

```
DISPLAY:
   Display modes
→ optimise display
   Exit
```

```
Please Enter:
Optimise display
2
```

Selection in the „Display Angle“ submenu:
Selection of a value to change the display angle by turning the encoder.

Time Scale

```
MAIN MENU:
    Temperature Scale
I→ Time Scale
   Ventine
```

```
TIME SCALE:
   Seconds
I→ Minutes
```

„Time scale“: option
Select the unit on which the timing of the thermoregulation programs is to be based.

Mains Failure Auto

```
MAIN MENU:
  PI-Parameters
I→ PowerOff AutoStart
  Select Usermenu
```

```
POWEROFF AUTOSTART.:
I→Off
  On
```

„Mains Failure Auto“ option
„Off“ After mains failure, manual input is required to continue the thermoregulation process.
„On“ After mains failure, the thermoregulation process is continued automatically. The setpoint programmed last will be used for thermoregulation.

Refer to 4.1 Safety Instructions!

4. Thermoregulation via Controller

4.3. Editing Default Settings

8. Alarm Configuration
9. Alarm Clear

Alarm: read the alarm message - eliminate the cause - acknowledge the message (alarm clear).

```
MAIN MENU:
   Alarm clear
I> Alarm confie.
   Analoe Interface
```

Selection of the parameters that are to trigger an alarm monitored by the software.

```
ALARM CONFIG.:
I> Alarm mode
   Lower Alarm Limit
   Upper Alarm Limit
   Level alarm delay
   Exit
```

„Run Mode“ (ref. to „Alarm clear“)
„Stop Mode“ (ref. to „Alarm clear“)

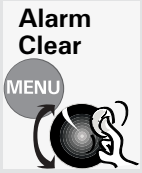
Submenu selection:
„Lower Alarm Limit“ / „Upper Alarm Limit“
The lower and upper alarm limits define the temperatures that trigger an alarm and stop the thermoregulation process, depending on the Alarm Mode settings.
Refer to 4.1 Safety Instructions!

```
ALARM MODE:
   Run Mode
-> Stop Mode
```

```
Please enter:
Upper Alarm Limit
          40.0°C
```

„Level Alarm Delay“
The level alarm delay is defined by entering the delay time in seconds.
At commissioning, a delay of 60 seconds is recommendable to provide for uninterrupted distribution.

```
Please enter:
Level alarm delay
40
```

```
MAIN MENU:
   2. setpoint
I> Alarm clear
   Alarm confie.
```

„Alarm Clear“ option
Acknowledge the alarm.

In the case of software-monitored alarm messages, the thermoregulation process is continued.
The alarm message will persist until the cause auf the alarm has been eliminated and the alarm message has been acknowledged.
„Alarm Mode“ when „stop mode“ has been previously set:
The thermoregulation process will be stopped.
Once the cause of the alarm has been eliminated and the alarm message has been acknowledged, it can be restarted manually.

```
ALARM MODE:
   Please switch unit
   off and on
```

```
ALARM CLEAR:
I> Restart
   Alarm clear
```

„Alarm Mode“ when „Run mode“ has been previously set after the alarm has been eliminated:
„Restart“: The controller is restarted.The message disappears.
„alam clear“:
The message disappears.

The following applies for each alarm mode:
Faults caused by hardware failures (overtemperature protection / overheating protection, low fluid level protection, pump winding thermostat) result in cut-out of the thermoregulation process.
Alarm messages due to hardware failures may be triggered by the overtemperature protection/overheating protection, by the low liquid level protection and by the pump winding thermostat.
Refer to 4.1 Safety Instructions!

4. Thermoregulation via Controller

4.3. Editing Default Settings

10. Compressor Automatic
11. Maximum Heating Power
12. Software Version

Compressor Auto

```
MAIN MENU:
   Calibration Pros.
I> Compressor Auto
   Mains Failure Auto
```

```
COMPRESSOR AUTO:
I> Automatic
   Always on
   Always off
```

Selection of the compressor starting mode

„Automatic“
The compressor will start automatically depending on the topical demand. The chiller will work on demand only.
Benefit: Saving of energy.
Drawback: Heating time (idle time) in the case of a sudden demand.

„Always On“
The compressor remains switched on, continuous chiller operation.

„Always Off“
The compressor remains switched of, the chiller is not in operation.

Heating Capacity

```
MAIN MENU:
   Machine options
I> Max. heat power
   Calibration Pros.
```

```
Please enter:
Max. heat power (%)
100
```

„Max. Heat Power“ option

Enter the desired maximum heating power in percent.
This is required for devices equipped with a heater and a compressor (chiller).
The heating power must be reduced to enable simultaneous operation of heater and compressor with a view to the fusing of the device.

Software version

```
MAIN MENU:
   Acoustic alarm
I> Software version
   Set-Point limits
```

```
SOFTWARE:
Series No.: 0
Version     04.00s
Date   19.02.03.  15.41
```

„Software Version“ option

For example:
04.00s vom 19.02.2003, 15.41 PM.

4. Thermoregulation via Controller

4.3. Editing Default Settings

13. Machine Options

Machine features

The possible settings of the unit are dependent on the features of the unit.
The settings can be selected in the submenu. The integrated / embeded features of the unit will be displayed.

```
MENU FEATURES:

I→ Pump speed
```

The desired pump speed can be set if the device is equiped with a pump with variable speed drive.

```
MENU FEATURES:

I→ External Signal
```

An external control signal
can be used to achieve different results. For this purpose, a submenu for selection of the resulting action on activation of the external control signal will be displayed.

```
MENU OPTIONS:

I→ External Signal
   -Off-
```

„Off“: the external control input has no function.

```
MENU OPTIONS:

I→ External Signal
   -Standby-
```

„Standby“: when the external control signal is activated the thermoregulation will be switched on and stays active untill the external control signal is deactivated.

```
MENU OPTIONS:

I→ External Signal
   -Act.2.Setpoint-
```

„Act. 2nd setpoint“ on activation of the external control signal input the unit will switch to the 2nd setpoint.

```
MENU OPTIONS:

I→ External Signal
   -Exit-
```

„Exit“: no change of setting, exit menu.

```
MENU OPTIONS:

I→ Reserved
```

has no significance

```
MENU OPTIONS:

I→ Exit
```

no change of setting, exit menu.

```
MENU OPTIONS:
   Venting
I→ Machine features
   Max. heat power
```

„Machine features“
This menu option is available in special cases only.

4. Thermoregulation via Controller

4.3. Editing Default Settings

14. Analog Interface – Parameter Input
Applies to all thermostats with the **Compatible Control CC3**

Anlog Interface

```
MAIN MENU:
   Alarm mode
I> Analog Interface
   Display
```

„Analog Interface“ (AIF) option: The thermoregulation process is controlled via an analog signal (currents from 4 to 20 mA), the strength of which represents the setpoint. The temperature range can be set by the user. The difference between *Zero* and *Span* must be at least 10 K but must not exceed 320 K.

```
ANALOG INTERFACE:
I> Temp. T1 (Zero)
   Temp. T2 (Span)
   Conf. Input
   Conf. Output
   Exit
```

Submenu selection: „Temp. T1 (Zero)“

```
Please enter:
Temp. 1 (zero)
            5.0°C
```

Default „Temp.1“: Lower limit of the temperature range

```
ANALOG INTERFACE:
   Temp. T1 (zero)
I> Temp. T2 (Span)
   Conf. Input
   Conf. Output
   Exit
```

Submenu selection: „Temp. T2 (Span)“

```
Please enter:
Temp. 2 (Span)
           35.0°C
```

Default „Temp.2“: Upper limit of the temperature range

```
ANALOG INTERFACE:
   Temp. T1 (Zero)
   Temp. T2 (Span)
I> Conf. Input
   Conf. Output
   Exit
```

Submenu selection: „Param. Input“

Default: Parameter input disabled/enabled or configuration.

```
ANALOG INPUT:
   AIF-Input Off
   AIF -> Set-point
I>  Adjust
   Exit
```

Submenu selection: „Configuration“

```
ANALOG INPUT:
I> AD-value at T1
   AD-value at T2
   On error Analog
   Exit
```

Default: Measured value of the analog-digital converter at T1/T2: „AN. CURR OK?- Yes“ if the analog device (provided by the customer) is synchronized with the Huber controller.
„AN. CURR OK?- No“ if the analog device (provided by the customer) must be resynchronized with the Huber controller. Automatic return to the menu.

```
ANALOG INPUT:
   AD-Value at T1
   AD-Value at T2
I> On error Analog
   Exit
```

Selection in the „If Error at Analog“ submenu:
Response to errors: Cut-out or enabling of the 2nd setpoint.

```
ANALOG INPUT:
   AD-Value at T1
   AD-Value at T2
I> On error Analog
   Exit
```

When the AIF is enabled, the input current determines the setpoint. If a setpoint is entered via the CC3 keyboard in this period, this setpoint will be enabled only after the AIF has been disabled. The definition of the setpoint via the AIF can be aborted with the MasterClear function. The setpoint defined prior to enabling the AIF will then be used for thermoregulation.

4. Thermoregulation via Controller

4.3. Editing Default Settings

15. Analog Interface – Parameter Output
Applies to all thermostats with the **Compatible Control CC3**

```
MAIN MENU:
   Alarm mode
I> Analog Interface
   Display
```

```
ANALOG INTERFACE:
   Temp. T1 (Zero)
   Temp. T2 (Span)
   Conf. Input
I> Conf. Output
   Exit
```

```
ANALOG OUTPUT:
I> Output value
   Current at T1
   Current at T2
   Exit
```

```
SOURCE ANALOG OUT:
   No output
   Internal temp.
   External temp.
   Programmer
I> Manual value
   Exit
```

Select „Analog Interface“ (AIF): The thermoregulation process is controlled via an analog signal (currents from 4 to 20 mA), the strength of which represents the setpoint. The temperature range can be set by the user. The difference between *Zero* and *Span* must be at least 10 K but must not exceed 320 K.

Submenu selection: „Output Value“ Definition of the measuring point for temperature output values.

„No Output“ Constant output of 4 mA as a power supply, e.g. for an external thermometer.

„Internal Temp.“: Measurement using an internal sensor in the thermostat. In the case of a bath thermometer, the sensor is inside the bath.

„External Temp.“: Measurement with an external sensor, e.g. in an external bath.

„Programmer“: In a thermoregulation program, each segment can be assigned an output current, e.g. for adapting the speed of an external pump to individual program segments.

```
ANALOG OUTPUT:
   Source for output
I> Output value
   Current at T1
   Current at T2
   Exit
```

„Manual Value“ Any percentage from 0 %...100 % analog to 4...20 mA with encoder.

```
Please enter:
DAC-Output at T1
591
```

Submenu selection: „Current at T1“ Presetting of a converter output value for T1 (Zero): The encoder is used to set a value that results in the desired current value for T1, e.g. 4 mA at the analog setpoint encoder.

```
ANALOG OUTPUT:
   Output value
   Current at T1
I> Current at T2
   Exit
```

Submenu selection: „Current at T2“

```
Please enter:
DAC-Output at T2
4104
```

Presetting of a converter output value for T2 (Span): The encoder is used to set a value that results in the desired current value for T2, e.g. 20 mA at the analog setpoint encoder.

Caution! The controller module may be destroyed if currents exceeding 20 mA are used and/or if the polarity is confused!

4. Thermoregulation via Controller

4.3. Editing Default Settings

16. Digital Interface

Applies to all thermostats with the **Compatible Control CC3**

Digital Interface

```
MAIN MENU:
   Control Mode
I> Digit.Interface
   Display
```

```
DIGIT. INTERFACE.:
I> Hardware RS
   Baud rate
   Slave address
   exit
```

```
HARDWARE RS:
I> RS 232
   RS 485
```

```
DIG. SCHNITTST.:
   Hardware RS
I> Baud rate
   Slave address
   Exit
```

```
BAUD RATE:
    1200 Baud
   2400 Baud
   4800 Baud
I> 9600 Baud
   19200 Baud
```

```
DIGIT. INTERFACE:
   Hardware RS
   Baud rate
I> Slave address
   exit
```

```
Please enter:
Slave address
1
```

„Digit. Interface“ option
The controller is equipped with a bidirectional RS232 interface and an RS485 interface. These digital interfaces enable remote control via a PC (Remote mode).

Submenu selection: „Hardware RS“
Preselection of the RS232 (for 1 PC) or RS485 (for up to 32 PCs) interface.

Submenu selection: „Baud rate“

Preselection of the data transfer rate between the thermostat and the connected PC. You can select one of five baud rates.

Factory default: 9,600 Baud

Submenu selection: „Slave address“

The Huber thermostat is assigned an „address“, i.e. an assignment across the entire device system of the user. Selection range: 0 to 99.

If the analog input is enabled as the setpoint source in the menu, then this setpoint has a higher priority than the setpoint sent to the controller via the digital interface.

Keyboard entries are not possible in remote mode. There is only one exception: the MasterClear function (press the MENU and TEMP keys simultaneously). In this case, the program exits the remote mode and the controller can be operated via the keyboard again. At the same time, the controller setpoint active prior to selecting „RS232“ or „RS485“ will be reactivated (auxiliary setpoint).

4. Thermoregulation via Controller

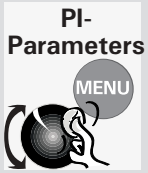
4.3. Editing Default Settings

17. PI-Parameters

All the factory defaults can be customized in the menu.

All the factory defaults can be restored via the „Factory Default“ option (with the thermoregulation function switched off). (Also refer to „Salutation“ 4.1).

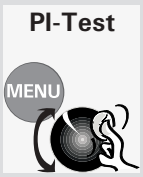
Factory Default:
The P-parameter (proportional parameter) and the I-parameter (integral parameter) influence the thermoregulation behavior of your thermostat. The factory default of the P-parameter is 5,000, that of the I-parameter is 1,000.
The factory default is well suited for a large number of applications.
New settings:
The PI-parameters can be edited as desired.
Value range of the P-parameter: 50 ... 30,000
Value range of the I-parameter: 0 ... 30,000

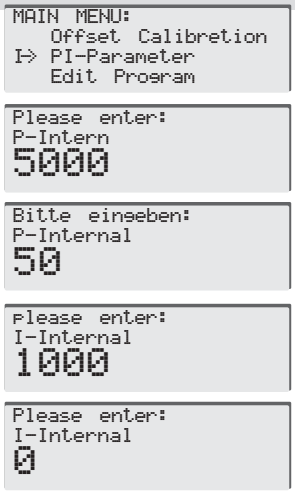
Test your thermoregulation process for optimum settings by entering new value pairs for the thermoregulation mode (internal or process (cascade)) you intend to use.

Test Part 1	External thermoregulation Temperature change by 20 °C after each setting up to U* (ballistic effect)
P - portion	50 1000 2000 3000 ... 10.000 ...up to Ü*e.g. 12.000
I - portion	0 0 0 0 0
Test Part 2	External thermoregulation Temperature change by 20 °C after each setting up to O**
P - Portion	12.000 12.000 12.000 12.000
I - Portion	1.000 2.000 5.000 ... up to O** e.g. 15.000

* U: Ballistic effect: At the thermoregulation target, the actual value oscillates about the setpoint.
**O: Optimum ratio between the accuracy of the controller and your desired speed.
Rule: Fast thermoregulation due to a high P and a low I-parameter results in a high ballistic effect.

Thermoregulation:
1. Heating by 20 °C with P/I = 50/0, then cooling by 20 °C
2. Heating by 20 °C with P/I = 1000/0, then cooling by 20 °C
3. Continue as described for Test Part 1 (table above) until all the relevant P-parameters have been tested.
4. Heating by 20 °C with P/I = 12,000/1,000, then cooling by 20 °C
5. Heating by 20 °C with P/I = 12,000/2,000, then cooling by 20 °C
6. Continue as described for Test Part 1 (table above) until all the relevant I-parameters have been tested.

4. Thermoregulation via Controller

4.3. Editing Default Settings

15. Acoustic Alarm
16. Select Usermenu
17. Configure Usermenu

Acoustic Alarm

```
MAIN MENU:
   2nd Stpoint
I> Aucoustic alarm
   Alarm Clear
```

```
AUCOUSTIC ALARM
I> OFF
  ON
```

„Acoustic Alarm“ option
„Off“ Alarm signals and error messages without acoustic alarm.
„On“ Alarm signals and error messages with acoustic alarm.

Select Usermenu

```
MAIN MENU:
   PowerOff AutoStart
I> select user menu
   Setpoint Limit
```

```
USER MENU:
I>Administrator
  User menu 1
  User menu 2
  User menu 3
  User menu 4
  User menu 5
  User menu 6
  User menu 7
  Exit
```

„Select Usermenu“ option:

Only the administrator can select this submenu after entering the correct password.

Config Usermenu

```
MAIN MENU:
   Compressor Auto
I> Confie. User menu
   Control Mode
```

```
Please enter:
User menu
1
```

The „Config Usermenu“ option is visible to the administrator only (default see above „Select Usermenu“).
For this reason, only the administrator can select this submenu after entering the correct password.
Via the „Config. Usermenu“ option you define the menu options that are to be visible in the Main Menu, i.e. the factory default is customized.

The administrator password is forwarded to the user separately (on request).

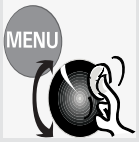
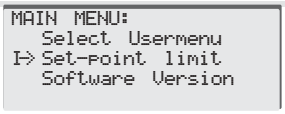
4. Thermoregulation via Controller

4.4. Limiting the Thermoregulation Range

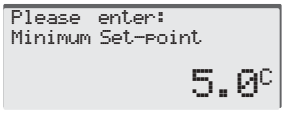
1. Setpoint Limits

Setpoint Limit

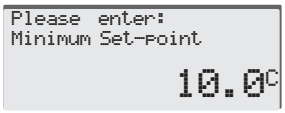
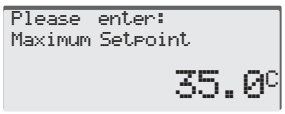
„Setpoint Limits“ option:
Enter the desired minimum setpoint and confirm your entry (e.g. 10.0 °C).

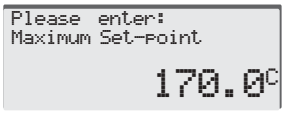
The minimum setpoint is a safety limit for thermoregulation. Where:
The lowest permissible temperature value > = minimum setpoint (3.1, 4.1).
It is not possible to enter a setpoint lower than the minimum setpoint.

Enter the desired maximum setpoint and confirm your entry (e.g. 170.0 °C).

The maximum setpoint is a safety limit for thermoregulation. Where:
The highest permissible temperature value < = maximum setpoint (3.1, 4.1).
It is not possible to enter a setpoint higher than the maximum setpoint.
The maximum setpoint should not be set to a value lower than 5-8 degrees below the value set for the overtemperature protection. In this way, the controller can tolerate a minimum ballistic effect of the attained temperature when thermoregulating to the maximum setpoint.

Warning
The overtemperature protection is an especially important safety device of your thermostat. It should always be operable and thus be tested **periodically**!
If the actual value exceeds the set overtemperature limit, an alarm is triggered and the thermostat will cut-out the thermoregulation process. This process can be restarted only after the cause of the alarm has been eliminated and the alarm message acknowledged.
Refer to 3.1 Principles and Safety Instructions!

4. Thermoregulation via Controller

4.5. Calibration

1. Calibration Program

Calibration Prog.

```
MAIN MENU:
   Analog Interface
I> Calibration proc.
   Circulation
```

„Calibration Prog.“ option
Is used exclusively to calibrate the internal sensor.
Use a calibrated reference thermometer as a second temperature sensor for calibration.

```
2-P.CALIBR. INTERNAL:
I> Edit T-Cal1/2
   Control to T-Cal1
   Control to T-Cal2
   exit
```

```
Please enter:
Set-point TCal1
              6 C
```

Submenu selection: „Edit TCal1/2“

```
Please enter:
Set-point TCal1
             10 C
```

Entry of the 1st of two calibration temperatures, e.g. change from 6°C to 10°C.

```
Please enter:
set-point TCal2
            100 C
```

Entry of the 2nd of two calibration temperatures, e.g. change from 100°C to 40 °C.

```
Please enter:
set-point TCal2
             40 C
```

Submenu selection: „Control to TCal1“

```
2-P.calibration:
   Edit TCal1/2
I> Control to T-Cal1
   Control to T-Cal2
   exit
```

Start thermoregulating until the 1st calibration temperature has been reached.
If your reference thermometer indicates the set temperature reliably, compare the value to the actual value display of the Compatible Control controller. Deviations can be corrected using the encoder.
Submenu selection: „Control to TCal2“
Proceed as described for TCal1.
Exit the menu once you have completed the calibration procedure.

```
2-P.calibration:
   Edit TCal1/2
   Control to T-Cal1
I> Control to T-Cal2
   Exit
```

The two-point calibration includes all the values between T1 and T2.
With offset calibration, on the other hand, the entire temperature level is shifted by one value on a linear level.

4. Thermoregulation via Controller

4.5. Calibration

2. Offset calibration

```
MAIN MENU:
   Mains Failure Auto
I→ Offset Calibration
   PI-Parameters
```

```
OFFSET CALIBRATION:
I→ Internal sensor
   Process sensor
   Exit
```

„Offset Calibration“ option
Suitable for calibrating all the sensors used.
Use a calibrated reference thermometer as a second temperature sensor for calibration in the area of the sensor to be calibrated.
Compare the actual controller temperature displayed to the reference thermometer reading.

```
Please enter:
Internal sensor
              0.0°C
```

Any deviation can be corrected via the options provided in the „Internal Sensor“ submenu by entering the difference with the versed sign.

Example:
Controller display (Internal Sensor) 10 °C, reference thermometer reading 8.5 °C, difference = 1.5 K.
Calibrate by entering „-1.5“. The output level of the controller display will decrease by 1.5 K. The controller display will now correspond to the reference thermometer reading.

```
Please enter:
Internal sensor
             -1.5°C
```

```
OFFSET CALIBRATION:
   Internal sensor
I→ Process sensor
   exit
```

Submenu selection: „Process Sensor“
Calibrate the process sensor (e.g. in an external bath)
Proceed as described for the „Internal Sensor“.

Offset calibration shifts the entire temperature level by one value on a linear level.
Two-point calibration, on the other hand, includes all the values between T1 and T2.

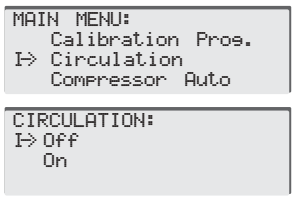
4. Thermoregulation via Controller

4.6. Utilities

1. Circulation
2. Venting
3. 2. Setpoint (for **Compatible Control CC3** only)

Circulation

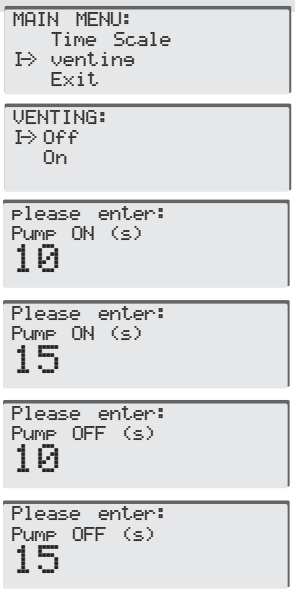
„Circulation“ option
„Off“
No pump operation (related to the thermoregulation process) or stop of the alternating pump operation.
„On“
Start of the alternating pump operation (without thermoregulation) e.g. to enhance the filling procedure.

Venting

„Venting“ option:
Can be selected only with thermoregulation switched off.
Using the venting option, the pump can be operated in intervals in alternating mode, e.g. for enhanced venting of external applications.

„On“
The default settings for the time intervals for pump operation/pause may have to be edited (take into consideration the viscosity of the thermofluid and the system dimensions) and – at the same time – start of the alternating pump operation in intervals.

„Off“ Stop the alternating pump operation in intervals.

(Refer to 3.4!)

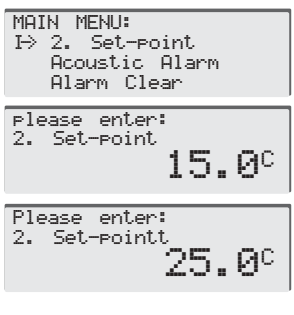
Caution! Potential risk of injury and material damage!

Please take into consideration the capacity and fill level of the thermostat and the connected systems as well as the viscosity and expansion characteristics of the thermofluid used.
Please prevent overflow of the fluid. Refer to 3.4!

2. Setpoint

„2nd Setpoint“ option *(exclusively available for Compatible Control CC3)*

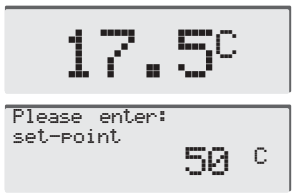
Entry of the 2nd setpoint. This setpoint is enabled only if an error occurs in the analog control. Refer to „Analog Interface“!
When entering the 2nd setpoint, the same applies as to the „standard“ setpoint: the characteristics of the thermofluid, the thermoregulation objective and the safety measures must always be taken into consideration.

4. Thermoregulation via Controller

4.7. Enter a Setpoint - Start

Press the SET key.
The setpoint will be shown on the status display.
Enter the new setpoint by turning the encoder (e.g. to 50 °C).
Confirm your entry by pressing the encoder or the SET key.

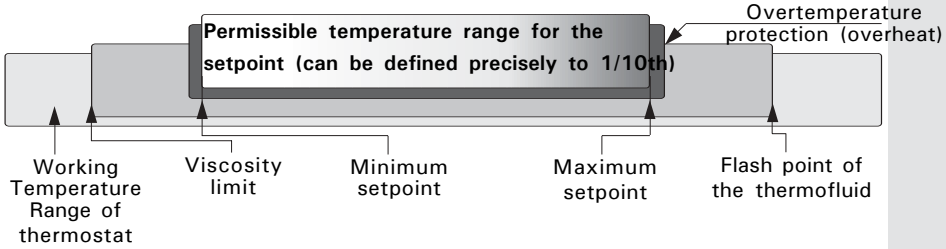
Start thermoregulation only after all the prerequisites have been met (refer to 3 Commissioning), especially:

Suitable location (3.1)
Correct connections (3.2)
Ambient temperature max. 30 °C (3.1)
Correct setting of the overtemperature protection (4.1)
Correct setting of the setpoint limits (4.4)

Compatible Control CC1 and CC2 controllers:
Start the thermoregulation process by turning the encoder.

Compatible Control CC3:
The thermoregulation process to the new setpoint is started by pressing the TEMP key.

To stop the thermoregulation process: Refer to section 4.9!

It is not possible to enter a setpoint beyond the setpoint limits.

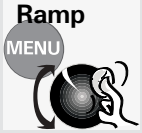
4. Thermoregulation via Controller

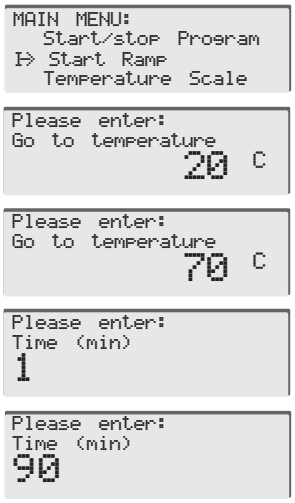
4.8. Convenient Thermoregulation – Programs

1. Start Ramp (for Compatible Control CC2/CC3 only)

Start Ramp

If you want to change the working temperature slowly and smoothly instead of suddenly, you should implement the setpoint change via a ramp.

„Start Ramp“ option:

Enter the desired final temperature of the ramp (ramp setpoint), e.g. 70 °C.

Enter the time (in minutes) the thermoregulation to the ramp setpoint should take, e.g. 90 minutes.

Starting the ramp:
The ramp will automatically be started once you have confirmed the time parameter.

Completion of the ramp:
Once the ramp setpoint has been reached (e.g. after 90 minutes), the thermostat will keep the new actual temperature constantly on the new value (e.g. 70 °C

Interrupting the ramp:
The ramp can be interrupted in the course of the process by entering a new setpoint (SET/encoder).

Caution!
The setpoint must be 25 degrees below the flash point of the thermofluid and above the temperature, at which the thermofluid attains a viscosity higher than 50 mm^2/s. (3.1.!)

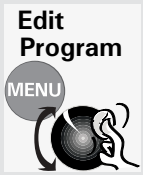
4. Thermoregulation via Controller

4.8. Convenient Thermoregulation – Programs

2. Edit Program
(Exclusively for **Compatible Control CC2/ CC3**)

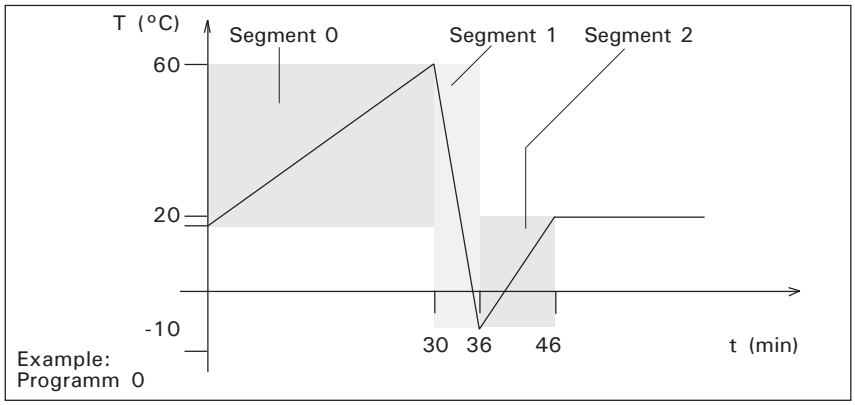
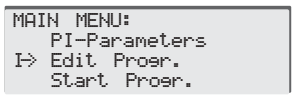
All the factory defaults can be customized in the menu. All the factory defaults can be restored via the „Factory Default“ option (with the thermoregulation function switched off). (Also refer to „Salutation“ 4.1).

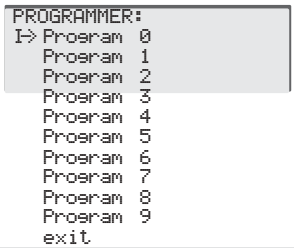
Creation of a separate thermoregulation program
Compatible Control CC2: 1 thermoregulation program for editing (can always be overwritten) with as many as 5 segments.
Compatible Control CC3: 10 thermoregulation programs for selecting and editing (can always be overwritten). A total of 50 segments are available.
One segment can maximally span a period of 54 hours.
Programming is done in steps according to the „Edit Program“ menu item.

Planning (Example):
A fluid in an external bath is to be heated and cooled in 3 steps. The heating process is to be temperature-stable and the cooling process time-stable. For the time it takes to cool the bath, an agitator is to be controlled via the potential-free contact.
At the end of the program, the thermostat is to maintain the bath temperature attained last until the process is terminated manually.

„Edit Program“ option

Options in the „Programmer“ (PR) submenu:
10 Programs
„Program 1“ option
Upon initial commissioning, all the programs will still be „empty“, i.e. they have not been assigned any segments. In the course of the application, these programs may be completely filled with segments.
For further programs, you may edit existing ones.

Caution!
The setpoint must be 25 degrees below the flash point of the thermofluid and above the temperature at which the thermofluid attains a viscosity higher than 50 mm^2/s. (3.1.!)

4. Thermoregulation via Controller

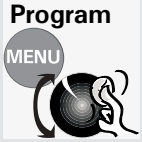
4.8. Convenient Thermoregulation – Programs

3. Edit Program
(exclusively for Compatible Control CC2/ CC3)

All the factory defaults can be customized in the menu.
All the factory defaults can be restored via the „Factory Default“ option (with the thermoregulation function switched off). (Also refer to „Salutation“ 4.1).
Select „Exit“ to return to the superordinate level.

Edit Program

```
FUNCTIONS PR:
   Attach Segment
I> insert Segment
   modify Segment
   Delete Segment
   show Segment
   Delete Program
   Exit
```

Submenu selection:
„Programmer“ (PR)/ „Program 1“/
„Functions PR“:
„Program 1“ options
„Attach Segment“ (for the 1st segment, this corresponds to „Insert Segment“). The parameters for the 1st segment are defined in the following. (Refer to the example of program 1 on page 48-2!)

```
PROGRAM NO.1:
I> Set point SegEnd
   Segment period
   Modify Segment
    Control mode
    Options
   Save & exit
   Exit
```

```
Please enter:
Set point SegEnde
          0.0 C
```

```
Please enter:
Set point SegEnde
         60.0 C
```

„Setpoint SegEnd“ option: entry of the setpoint for the 1st segment of the 1st program, e.g. 60 °C.

```
PROGRAM NO.1:
   Set point SegEnd
I> Segment period
   Modify Segment
    Control mode
    Options
    Save & exit
    Exit
```

```
Please enter:
Segment period (s)
1
```

```
Please enter:
Segment period (s)
1800
```

„Segment period“ option: entry of the time period for the 1st segment of the 1st program, e.g. 1800 seconds.

```
PROGRAM NO.1:
   Set point SegEnd
   Segment period
I> Control mode
   Options
   Save & exit
   Exit
```

```
CONTROL MODE:
I> Internal
   External
```

„Control Mode“ option: example: selection of external thermoregulation.

```
OPTIONs PR:
   Pot.free Contact
   Analoe Output
   End condition
-> Stability
   exit
```

```
STABILITY:
I> Time-stable
   Temperature-stable
```

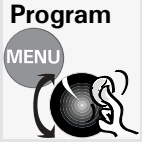
„Stability“ option, e.g. „Time-Stable“ for the 1st Segment of the 1st program.
Time-Stable The segment period entered has priority for the segment end.
Temperature-Stable: The „Setpoint SegEnd“ has priority for the end of the segment.

4. Thermoregulation via Controller

4.8. Convenient Thermoregulation – Programs

3. Edit Program
(Exclusively for Compatible Control CC2/ CC3)

Edit Program

```
OPTIONs PR:
   Pot.free Contact
   Analoe Output
   End condition
   Stability
 I→exit
```

Select „Exit“ to return to the superordinate level „Program 1“.

```
PROGRAM NO.1:
    Set point SesEnd
    Seament period
    Control mode
    Options
I→ Save & exit
    Exit
```

„Save & Exit“ option: Saving of all the data for the 1st segment in the 1st program. The 1st program has thus been created.
If you exit the „Program 1“ level without „Save & Exit“, all the data entered for this segment up to now that have not yet been saved will be discarded.
After „Save & Exit“ the program will return to the superordinate level „Functions PR“.

```
FUNCTIONS PR:
   Attach Seament
   insert Seament
   modify Seament
   Delete Seament
I→ show Seament
   Delete Proaram
   Exit
```

```
Pros:0      Ses:0
Temp:60     Contr:Cas
Time:1800   Stab:Time
Poco:0      An0:0
```

„Functions PR“ option, „Show segment“ option: this is where you can check your entries. Press the MENU key to exit the display.
The program will return to „Functions PR“, where you select „Attach Segment“ as shown in the example. The entries now correspond to those for segment 1 taking the values from our example. In this segment, you must additionally select the potential-free contact for controlling the agitator in segment 2.

```
FUNCTIONS PR:
I→ Attach Seament
   insert Seament
   modify Seament
   Delete Seament
   show Seament
   Delete Proaram
   Exit
```

```
Please enter:
Set point SesEnd
          60.0
```

```
Please enter:
Set point SesEnd
         -10.0 C
```

```
PROGRAM NO.1:
    Set point SesEnd
    Seament period
    Control mode
    Options
I→ Save & exit
    Exit
```

```
Please enter:
Seament period (s)
1800
```

```
Please enter:
Seament period (s)
360
```

```
PROGRAM NO.1:
    Set point SesEnd
    Seament period
 I→ Control mode
    Options
    Save & exit
    Exit
```

```
CONTROL MODE:
I→ Internal
   External
```

```
CONTROL MODE:
   Internal
I→ External
```

Potential-free contact

```
PROGRAM NO.1:
    Set point SesEnd
    Seament period
    Control mode
 I→ Options
    Save & exit
    Exit
```

```
OPTIONS PR:
I→ Pot.free contact
   Analoe output
   Stability
   Exit
```

```
POCO ACTIVE:
   No
I→ Yes
   Unchanaed
```

```
OPTIONS PR:
   Pot.free contact
I→ Analoe output
   Stability
   Exit
```

```
Please enter:
An. output (%)
100
```

```
OPTIONS PR:
   Pot.free contact
   Analoe output
I→ Stability
   Exit
```

```
STABILITY:
   time-stable
I→ Temperature-stable
```

```
OPTIONS PR:
   Pot.free contact
   Analoe output
   Stability
I→ Exit
```

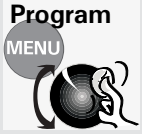
```
PROGRAM NO.1:
    Set point SesEnd
    Seament period
    Control mode
 I→ Options
    Save & exit
    Exit
```

4. Thermoregulation via Controller

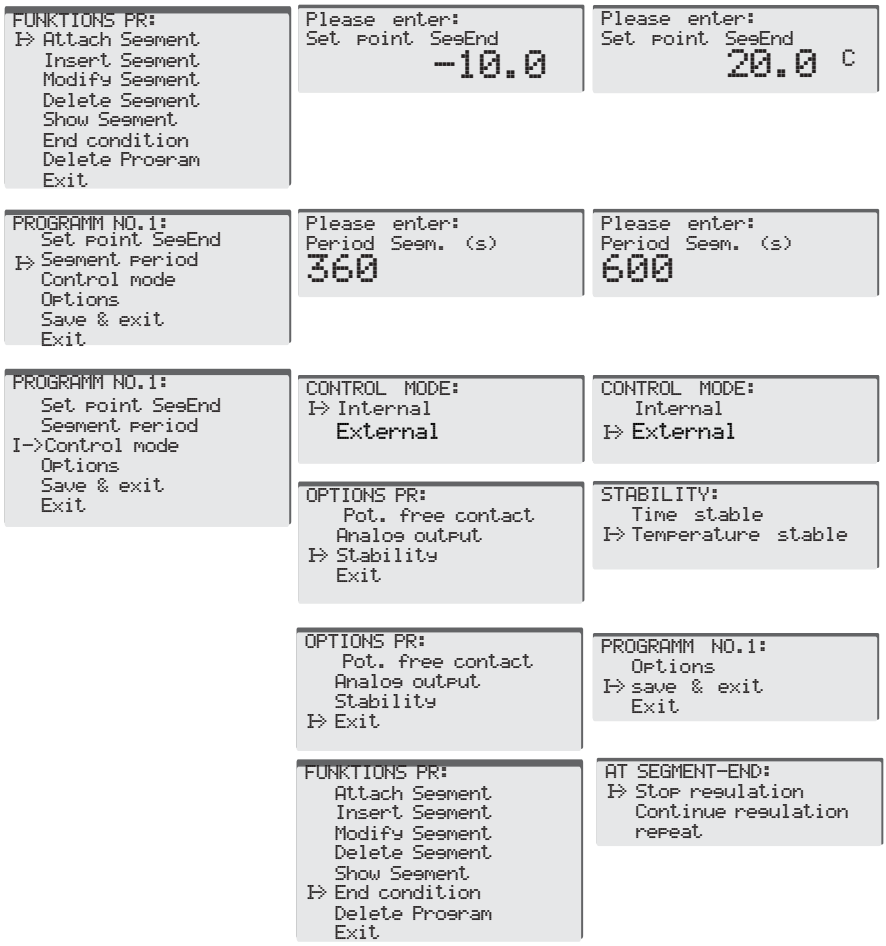
4.8. Convenient Thermoregulation – Programs

3. Edit Program
(Exclusively for Compatible Control CC2/ CC3)

Edit Program

„Functions PR“ option, „Show segment“ option: this is where you can check your entries. Press the MENU key to exit the display.
The program will return to „Functions PR“, where you select „Attach Segment“ as shown in the example. The entries now correspond to those for segment 1 taking the values from our example. In this segment, you must additionally select the potential-free contact for controlling the agitator in segment 2.

You have now created a thermoregulation program with three segments according to our example.

Caution!
Before activating the program, make sure you are using a suitable thermofluid!
The setpoint must be 25 °C below the flash point of the thermofluid and above the temperature at which the thermofluid attains a viscosity higher than 50 mm^2/s. (4.1.!)

Please Note:
You cannot create thermoregulation programs with the Compatible Control CC1.
You can create one thermoregulation program with the Compatible Control CC2.
You can create 10 thermoregulation programs with the Compatible Control CC3.

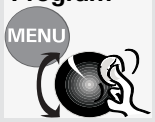
4. Thermoregulation via Controller

4.8. Convenient Thermoregulation – Programs

4. Start Program
(Exclusively for **Compatible Control CC2/ CC3**)

For alarm messages, please refer to 4.3!

Start/ Stop Program

```
MAIN MENU:
   Software Version
I> Start/Stop Program
   Start Ramp
```

„Start/Stop Program“ option

```
PROGRAMMER:
I> Program 0
   Program 1
   Program 2
   Program 3
   Program 4
   Program 5
   Program 6
   Program 7
   Program 8
   Program 9
   Exit
```

If ***no*** program has been started:
Submenu selection:
„Program 0“ example Compatible Control CC3:
10 programs / 50 segments
(However one program cannot contain more than 40 segments)

Compatible Control CC2: 1 program / 5 segments

```
SERVICE-PROGRAM:
   Program Pause
   Program Continue
   Go To Segment No.
   Program stop
   Exit
```

If a program has ***already*** been started:
Submenu selection:
„Program Pause“
„Program Continue“
„Go To Segment No.“
„Program Stop“
„Exit“

```
Progr.:1        Segm.:0
Internal          21.2
Process           17.5
Setpoint          17.6
```

Program 1 will be started.

„Internal“:
development of the internal actual value

„Process“:
development of the external actual value

„Setpoint“: Calculated topical setpoint, development until the segment end value has been attained

Aborting the thermoregulation program
With Compatible Control CC2: Power switch or new menu (using„Program Stop“).
With Compatible Control CC3:
Break function (press the MENU and On/ Off keys simultaneously).
Alternative: Power switch or new menu.

Standard end of the thermoregulation process once the thermoregulation program has been completely executed, according to:
- programmed segment end (the temperature of the last segment setpoint is maintained (Continue) or thermoregulation to the last setpoint entered outside the program) or repeating of Temp. Program.
- Stability: Time-Stable (i.e. after the programmed segment period has elapsed) or Temperature-Stable (i.e. after the segment setpoint has been attained).

```
   60.0 C
```

Once the program has been completely executed, the status display will be shown.

Caution!
Before activating the program, make sure you are using a suitable thermofluid!
The setpoint must be 25 °C below the flash point of the thermofluid and above the temperature at which the thermofluid attains a viscosity higher than 50 mm^2/s. (3.1.!)

4. Thermoregulation via Controller

4.9. Terminating the Thermoregulation Process

Entry New setpoint

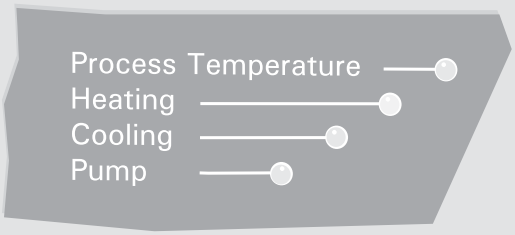
The thermoregulation process is automatically continued once the setpoint has been reached, i.e. the thermostat keeps the temperature at the measuring point at the temperature defined by the setpoint.
The thermoregulation process is continued until it is terminated manually.
This thermoregulation process can be changed at any time by entering a new setpoint (refer to 4.7).

Abort

Compatible Control CC1/CC2: The thermoregulation process can be aborted at any time by pressing the power switch. You can also turn the encoder until the „Thermoregulation“ menu is displayed, in which you can select „Off“ or „On“.

Compatible control CC3: The thermoregulation process can be aborted at any time by pressing the On/Off key on the controller.

When you press the On/Off key or the power switch to stop the thermoregulation process, all the LEDs will go out.

Actuate the power switch of the thermostat to interrupt the power supply.

Caution!
Do not stop the thermoregulation process by pulling the power plug. When the controller is switched back on, various device messages may be displayed and faults may occur.

5

Shut Down

5. Shut Down

5.1. Safety Instructions and Principles

Caution! Potential risk of injury and material damage:

Slip hazard! The floor and the workspace may be slippery when fluids have been spilt!
Tilt hazard! Make sure the thermostat and external devices are in a stable position.
Shock hazard! Make sure the connection to the power supply is undamaged and in perfect working order.
Scalding and burning hazard! Always be aware of extreme temperatures. **Causticization hazard!** Risk of causticization of your eyes, your skin and your respiratory tracts through hazardous vapors (depends on the thermofluid used).

Draining thermo-fluid

Avoid damage caused by extreme temperatures:
Before draining it, allow the thermofluid to cool or warm to an acceptable temperature that does not constitute a hazard.
Avoid contamination in the vicinity of the thermostat through spilt or splashed thermofluid and through escaping hazardous vapors.
Use appropriate drain tubes and collecting vessels.

Draining cooling water

For water-cooled thermostats:
Slip hazard! Avoid contamination in the vicinity of the thermostat through spilt or splashed cooling water.
Use appropriate drain tubes and collecting vessels.

Electrical connections

Make sure that the power connections are in good condition. Plugs, sockets and cables must be in good order and of the correct rating.
Avoid contact with fluids.

When replacing fuses, always use fuses with the correct ratings (refer to 5.3.!)

Cleaning

Clean and service the thermostat periodically (refer to 5.3!).
Keep the vicinity of the thermostat clean!
Always clean contaminated areas immediately.

Transport

Keep thermostats upright during transport.

Transport thermostats on clean, level and non-slip surfaces only.

Keep all the thermostat connectors/fittings closed.
Close the valves for fluids!
Protect plug-type connections using the caps provided for this purpose.

Please Note:

All the safety instructions are vital and must be considered during use in compliance with the present operating instructions!

5. Shut Down

5.2. Draining, Deactivating and Dismantling

Requirements

Keep at hand a collecting vessel, a suitable drain tube, a screwdriver, compressed air and a cleaning cloth.

Terminate the thermoregulation process via menu selection or by pressing the On/Off key (R4 with cc3). (The final temperature should be as close as possible to the ambient temperature. Highly viscous thermofluids may need to be above ambient temperature).

Draining process

Thermo-fluid

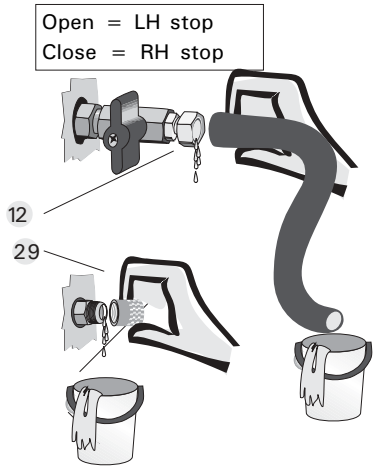
Applies to chillers (thermostats with closed fluid circuit)

Remove the plug on the draining connection 12.

Attach the drain tube onto the draining connection. Hang the free end of the drain tube into the collecting vessel.
Open the valve of draining connection 12.

Open draining connection 29 to empty any remaining thermofluid.

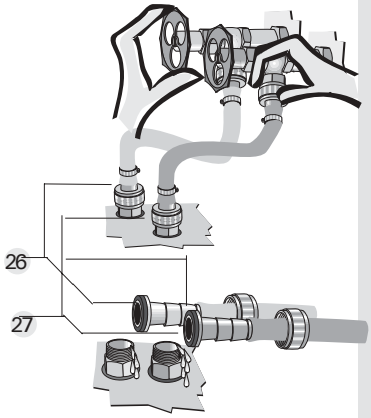
Close draining valve 12

Draining process

Cooling water

For water-cooled thermostats: Close the water shut-off valves supplying the cooling water to the machine.Unscrew hose connections from the device to the water supply at connections 26 and 27. Completely drain the cooling water via the cooling water draining connections with the aid of compressed air.
(Observe safety procedures when using compressed air)
(Take preventive measures against frost damage during transport and storage caused by cooling water residues)

Switching off

If the thermostat features a main switch, switch off the thermostat at the main switch (1, position „0“).Otherwise, switch off the unit switch (2).

Disconnect the power plug.

Note:The drained thermofluid must be disposed of in compliance with the instructions of the fluid manufacturer.The drained cooling water can be disposed of through the standard sewage.

5. Shut Down

5.3. Maintenance and Service

Require-ments

Keep at hand a collecting vessel, a suitable drain tube, a screwdriver, compressed air and a cleaning cloth.The thermoregulation process has been terminated by selecting a menu or by pressing the On/Off key (R4 on cc3).
(The final temperature should be as close as possible to the ambient temperature, with highly viscous thermofluids slightly above the ambient temperature).

Replacing fuses

In case of a fuse failure trace the source of the problem and eliminate it. Thereafter change the fuse.

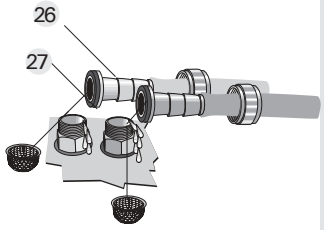
Cleaning cooling water filter

For water-cooled units:

Place the collecting vessel below cooling water draining connection 27.

Close the water supply valves to the unit. Unscrew hose connections 26 and 27

Remove the cooling water filter from the connectors (water inlet). Clean and replace it. Refit the hose connections.

Cleaning surfaces

Cover unused electrical connectors using the protective caps supplied.

Clean the stainless-steel surfaces using steel cleaning spray.

Carefully clean painted surfaces using a mild detergent.

Caution!
to prevent accidents:
Keep the working area clean!
Always clean contaminated areas immediately.
Service the equipment periodically!

5. Shut Down

5.3. Maintenance and Service

Disposal

For the avoidance of environmental damage, leave „retired“ thermostats exclusively by certified specialized enterprises despose to.

Decalcifying

If the water supply is „hard“, water-cooled thermostats must be decalcified periodically. We recommend that this procedure is carried out with De-calcifying unit „SEK 28“available from your Huber agent.

Remove the cooling water hoses.
Close the water supply connection to the unit (26) with the cap.
Connect the hoses for decalcifying fluid from the „SEK 28“ to the water draining connection and the water supply connection of the thermostat.
Follow the user instructions supplied with the SEK 28.

Service Hotline

Germany: (+49) 7 81 – 96 03 – 2 44
Hotline Service Center USA: (+1) 7 40 – 3 73 – 68 09
Hotline Service Center NL: (+31) 4 85 – 542 - 811

Caution!
to prevent accidents:
Keep the working area clean!
Always clean contaminated areas immediately.
Service the equipment periodically!

6

Appendix

6. Appendix

Programs for changing the factory default settings

6.1. Presettings

Standard Parameter Set - Factory Default

4.3.3.	Alarm Configuration – Upper Alarm Limit	305 °C
4.3.3.	Alarm Configuration – Lower Alarm Limit	-35 °C
4.3.3.	Alarm Configuration – Level Alarm Delay	4 s
4.3.3.	Alarm Mode	Stop Mode
4.3.6-7	Analog Interface – Temp. T1	5 °C
4.3.6-7	Analog Interface – Temp. T2	35 °C
4.3.6-7	Analog Interface – Parameter Input	disabled
4.3.6-7	Analog Interface – Parameter Output / Source Analog Out / Internal Temp.	
4.3.6-7	Analog Interface – Current at T1	DAC 593
4.3.6-7	Analog Interface – Current at T2	DAC 4104
4.3.2.	Display – Display Modes	Standard
4.3.2.	Display – Display Angle	100
4.3.8.	Digital Interface – Hardware RS	RS 232
4.3.8.	Digital Interface – Baud Rate	9600
4.3.8.	Digital Interface – Slave Address	1
4.6.1	Venting	Off / 10 s
4.3.4.	Max. Heat Power	100 %
4.5.1.	Calibration Prog. – TCal 1	0 °C
4.5.1.	Calibration Prog. – TCal 2	100 °C
4.3.4.	Compressor Auto	Automatic
4.3.2.	Mains Failure Auto	Off
4.5.2.	Offset Calibration – Internal Sensor	0 °C
4.3.9.	PI-Parameters – P-Internal	5000
4.3.9.	PI-Parameters – I-Internal	1000
4.3.10.	Edit Program	0
4.3.10.	Acoustic Alarm	On
4.7.	Setpoint	20 °C
4.4.1.	Setpoint Limits – Minimum Setpoint	5 °C
4.4.1.	Setpoint Limits – Maximum Setpoint	35 °C
4.3.1.	Language	Deutsch
4.3.1.	Temperature Scale	°C
4.3.1.	Control Mode	Internal
4.6.1.	Circulation	Off
4.3.10.	Usermenu	Administrator
4.3.2.	Time Scale	Minutes
4.6.1.	2nd Setpoint	15 °C

EXTRA

All the factory defaults can be restored via the „Factory Default“ option (refer to 42-1).

6. Appendix

6.2. Interface Specification, Data

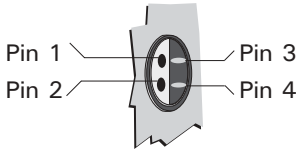
Note:
Before using the interfaces provided on the Huber device, please review the Data Communication Manual for information on the correct handling of interfaces and PC/ PLS data communication.

Please adhere to the following general principles:

- Always switch off all associated devices before connecting any cables.
- Only low voltages may be applied to the interfaces of the Huber devices.
- Always use an process sensor with a shielded cable. Otherwise the sensor may be electro-statically charged. The charge on the sensor could damage the measuring circuit when the sensor is connected.
- Before connection of plug-type connectors, ensure that they are in good and proper condition.

Pt100 external

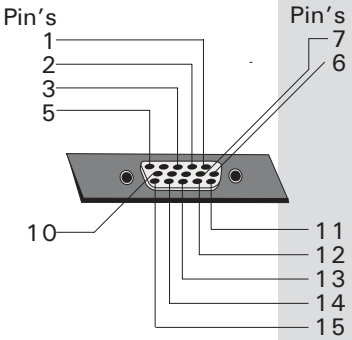
Pin 1: I+
Pin 2: U+
Pin 3: U--
Pin 4: I- (four-conductor technology)
Refer to 2.2.!

Potential-free contact Poko

R232/RS485

Pin 1: Terminating resistor 120 Ù, RS4851
Pin 2: Data input RS232 (RXD)
Pin 3: Data output RS232 (TXD)
Pin 4: Reserved, please do not use!
Pin 5: Earth (DGND)
Pin 6: Potential-free contact (Poko)
max. 1 A / 24 V DC
Pin 7: Terminating resistor 120 Ù, RS485
Pin 8: Reserved, do not use!
Pin 9: Reserved, do not use!
Pin 10: Potential-free contact (Poko)
Pin 11: A (RS485)
Pin 12: Analog output for AIF
Pin 13: Analog earth (AGND for AIF)
Pin 14: Analog input for AIF
Pin 15: B (RS485)

Please bridge pin 1 and 7 if the terminating resistor between A and B has to be activated.
POKO: the potencialfree contact signalises the state of the thermostat. The contact is closed when the unit is in working order. In case of an alarm or power failure the POKO contact is open.

Data Communication

Various protocols are available for complex parameterization to customize a Huber device to an automated application.
Examples: Point-to-Point commands for direct operation with a PC; LAI instruction set for RS485 bus systems, Modbus connection.
For more detailed information, please refer to: **Data Communication Manual.**

6. Appendix

6.3. Device Messages

1. Basic information

The alarm and warning messages of the Compatible Control controllers are divided into 3 categories:

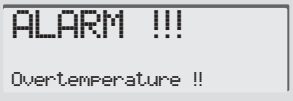
- Hard alarms, which are released by fault switches
- Soft alarms, which are released by software functions and plausibility checks
- Warnings, which are released by software functions and plausibility checks

Hard alarms have top priority, soft alarms have medium priority and warnings low priority. Warnings are comparable with soft alarms, the tempering is, however, **not** deactivated with warnings.

As long as no alarm or warning is pending, the program control can be assigned to an occurring event of all three groups. If a warning or an alarm is pending, only an alarm with higher priority can take over the program control.

Hard and soft alarms can only be deactivated if the device is switched off, the fault cause is rectified and the device is re-started. Unlike * alarms, warnings can be deactivated in the „alarm clear“ menu.
*(refer to 43-3: menu option „stop mode“).

An occurring event of all three groups is signalled both on the display and via the signal transmitter (with default setting).
The acoustic alarm can be disabled via the „Acoustic Alarm“ option.

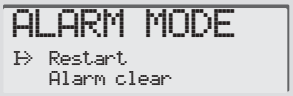
The display shows the alarm message in the form of a „flash“ display. In intervals of 2 seconds, the display alternates with the status display.

„Alarm!!!“ is displayed in the upper line of the alarm message.
The cause of the alarm is displayed i the bottom line of the alarm message.

To acknowledge the alarm message:
Select the „Alarm mode“ option.
Display in Stop Mode (refer to 43.9).
Follow the displayed prompts.

Display in Stop Mode (refer to 43.9).
Select „Restart“ or „Alarm clear“.

The following applies to all automatic messages:
Please check the external conditions of your thermostat and your controller inputs according to thbe operating intructions and your process arrangement.
If required. re-program the condition.
If the message persists after you have repeated this action, please contact our service hotline.

Hotline Deutschland:	**+49 - 781 - 9603 - 244**
Hotline Service-Center USA:	**+1- 740 - 373 - 6809**
Hotline Service-Center NL:	**+31 - 485 - 542 - 811**

6. Appendix

6.3. Device Messages

2. Hard Alarms

Hardware error

During the controller boot up before starting thermoregulation. Hardware error: Defect in the thermostate or system does not recognize the controller.

The circulation and thermoregulation process cannot be started.

Please call our customer service!

Device ID error

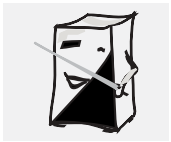
During the controller boot up before starting thermoregulation. Hardware error: Defect in the thermostate or unsuitable software (possible after a controller replacement).

The circulation and thermoregulation process cannot be started.

Please call our customer service!

The FSW?[CR][LF] command transmits a zero in fault-free conditions; in error cases error codes in accordance with the following tables are signalled.

Error	Released by	Error code
OT error „OVERTEMPERATURE !!“	Independent overtemperature protection, such as electronic overtemperature protection, mechanic overtemperature protection, pill	-1
Level Error „LEVEL!!“	Independent level monitoring, e.g. floater or pressure cell	-2
Winding thermostat error - with hook-in unit and UC006, 015, 020, 022 a. 025 „WT PUMP MOTOR !!“ - with devices with busboard „pressostat“ (UC006T-015T, CC240-250 a. BFT2) or „WT error chain“ with all other devices	Independent winding thermostats (hook-in unit) or error chain, which contains, among others, pressostat(s) or overcurrent relays (CCs and UCs with busboard)	-3

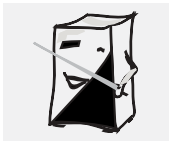
6. Appendix

6.3. Device Messages

3. Soft alarms

Error	Released by	Error code
Upper alarm temperature internal“ALARM TEMPERATURE !!“	The alarm temperature (not to be confused with the overtemperature protection) was reached in the internal bath	-4
Lower alarm temperature internal“ALARM TEMPERATURE !!“	The alarm temperature (not to be confused with the overtemperature protection) was reached in the internal bath	-5
Upper alarm temperature external“ALARM TEMPERATURE !!“	The alarm temperature (not to be confused with the overtemperature protection) was reached in the external process	-6
Lower alarm temperature external“ALARM TEMPERATURE !!“	The alarm temperature (not to be confused with the overtemperature protection) was reached in the external process	-7
Internal temperature above max. setpoint „Int.Temp >> Max. SetP.“	The internal temperature has exceeded the upper value of the setpoint limitation	-8
Line break internal sensor „INTERNAL SENSOR !!“	Detection of a fault in the internal sensor or supply line	-9
Line break external sensor „EXTERNAL SENSOR !!“	Detection of a fault in the external sensor or supply line	-10
Line break condensation sensor „SENSOR COND.!!“	Detection of a fault in the condensation sensor or supply line	-11
Cooling water error in CC40X „NO COOLING WATER“	The temperature at the expansion tank has considerably increased despite the water valve’s opening command	-12
Watchdog Alarm in RS „RS WATCHDOG !!“	The time for a watchdog operation at Labworld has expired due to a missing actuation	-13
Condensation alarm „COND. TEMPERATURE“	The temperature at the condenser and hence of the pressures in the refrigeration cycle has further increased despite a reduction when the warning temperature was reached. Switch off as otherwise a possibly installed pressostat trips.	-14
Pump error „PUMP SPEED“	10 seconds after the pump start, the actual speed is still below 250 rpm. Jammed or faulty pump?	-15
Circuit breaker error „Circuit Breaker“	A fault in the actuation of the circuit breaker has been detected	-16
Error in OT communication „OT Communication“	Despite the return of the OT actual values, the controller is not able to detect the receipt of OT-CAN command every 300msec	-17
Considerably deviating actual temperatures of OT and internal sensors „Diff. OT-actual value“	Except for the CC40X, the internal sensor and the OT sensor are positioned very close to each other. The deviation of the measured values exceeds 20K.	-18
Error in PCP Communication „PCP Communication“	CC40X: The controller cannot detect any CAN communication in the required time frame	-19
PCP-Firmware Error „PCP-Software Version“	CC40X: Firmware in controller and ControlPanel is not identical	-20
Error in Controller Communication „Controller Communication“	CC40X: Disturbed PCP – controller communication.	-21
Error during OT Test „OT Plausibility“	Controller checks OT setpoints in EEPROM with RAM setpoints. In case of deviations, an error is generated.	-22

6. Appendix

6.3. Device Messages

4. Warnings

Error	Released by	Error code
Upper water temperature internal „ALARM TEMPERATURE !!“	The warning temperature was reached in the internal bath	-258
Lower warning temperature internal „ALARM TEMPERATURE !!“	The warning temperature was reached in the internal bath	-259
Upper warning temperature external „ALARM TEMPERATURE !!“	The warning temperature was reached in the external process	-260
Lower warning temperature external „ALARM TEMPERATURE !!“	The warning temperature was reached in the external process	-261
Error in analogue input „SETPOINT DEVIATION“	The measured values at the analogue input are not plausible	-262
2nd error at analogue input „ANALOGUE INPUT!!“	The measured values at the analogue input are too high	-263
Hardware error „Controller replacement“ „detected, please execute“ „alarm elimination“ „menu point“	The configuration of a device with AC board (and code resistances) stored in the controller can no longer be found. The controller was probably plugged onto a different device with different code resistances and the respective initialisation (factory setting) has not yet been carried out	-264
Condensation warning “COND.TEMPERATURE“ “EXCESSIVE!!“ “WATER SUPPLY?“ “CLEAN CONDENSATOR?“	The temperature at the condenser and hence of the pressures in the refrigeration cycle has considerably increased; therefore, the refrigerating capacity has to be reduced.	-265
Cooling water errror in CC40X „EXP. TEMPERATURE“ „EXCESSIVE!!“ „WATER SUPPLY?“	The temperature at the expansion tank has considerably increased despite the water valve’s opening command	-266
PUMP-DAC error „No feedback“ „from D/A converter“ „on the busboard“	The module responsible for the analogue setpoint for the pump speed is not identified by the CPU	-267
Busboard-EEPROM error „No feedback“ „from EEPROM“ „on the busboard“	The memory containing the device configuration for CCs, ministats and some UCs cannot be identified	-268
Stepping motor communication error „No feedback“ „from SM module“	No trouble-free CAN communication with the stepping motor module can be identified	-269
DigPot error „No feedback“ „from digital-pot.“	The module outputting the perspective / contrast setting cannot be identified	-270
Busboardeepromdataerror „Error at R/W“ „of EEPROM“ „on the busboard“	Error when reading or writing busboard data	-271
Busboarddataerror „Busboard data“ „older than controller“ „software, please“ „update busboard“	The firmware contained in the controller requires a later version of busboard data	-272

Appendix

Huber Glossary

AIF – Analog Interface The analog interface is used to transfer constantly changing values, mostly temperatures, in the form of a signal current (generally 4 to 20 mA). Refer to „External Control"! Refer to 4.3!

Alarm LimitThe defined absolute temperature for the thermofluid or product depending on the control mode.When the temperature exceeds or drops below the defined value, the thermostat will respond as defined in the alarm mode.

Alarm Mode
The response of the thermostat when the temperatureerature has exceeded or dropped below the defined alarm limits.

Ambient Temperature Range
This is the temperature range permissible in the area surrounding the device for the device to function properly.
For all devices by HUBER, this is the temperature range from 5 to 30 °C. The cooling performance quoted is always based on an ambient temperature of 20 °C.

Automatic Mains Failure Control (Mains Failure Auto)
Determines the response of the thermostat in the case of a power failure.
The corresponding settings are made in the main menu of the controller.
Default setting:
Mains Failure Auto is disabled.Once the supply of power has been restored, manual intervention is required to continue the thermoregulation process.
Alternative setting:
Mains Failure Auto is enabled.Once the supply of power has been restored, thermoregulation is automatically continued.
Refer to 4.3!

Automatic Cooling Capacity Control
A method for saving energy and resources.The microprocessor control determines whether the required cooling capacity can be reduced according to the operating temperature.The adaptation is performed continuously and – apart form saving energy – it also protects the compressor, decreases heat dissipation into the environment and reduces the consumption of cooling water.
Refer to Compressor Auto, 4.3!

Automatic Compressor Control (Compressor Auto)
Supports an energy-saving method for cooling thermostats in specific thermoregulation situations.A control logic in the microprocessor helps decide whether the chiller is needed at all or whether it can be switched off.
Refer to 4.2.! Refer to Automatic Cooling Capacity Control!

Actual Temperature
The temperature at the measuring point at the time of display.The measuring point is either in the thermostat, in the thermofluid (factory default: internal thermoregulation) or external to the thermostat, e.g. in the core of a reactor (external thermoregulation).

Bath, Closed
Double-wall, closed „shell" with inflow and outflow connections.
Indirect thermoregulation of the customers' products in the „core".Designed in metal or glass.Refer to „Reactor", refer to 3.3!

Bath Opening
This is the active surface available for direct thermoregulation, generally across the entire usable depth.

Bath Thermostat
A bath thermostat is a thermostat equipped with a circulation pump (force pump in Polystats) and a bath tank for holding the object to be thermoregulated.The circulation pump is primarily used for mixing the bath fluid, if required, however, it can also pump the fluid through a closed external circuit, e.g. when flow water coolers for cooling heating thermostats are connected.

Bath Thermostat/Circulator
Thermostats with a bath opening large enough to accommodate objects for direct thermoregulation in the bath, including a circulation pump (force or suction pump with Compatible Control thermostats) for closed (force pump) or open external (suction pump) circuits.

Bath Volume (also referred to as fill volume)
The volume of the thermofluid inside the bath, required for the intended operation of the thermostat, excluding the thermofluid volume in external fluid circuits.
Especially in closed external applications, the expansion tank must be dimensioned with care since the circulator must additionally handle the expansion of the fluid in the external circuit.

Baud Rate
The data transfer rate for serial communication (mandatory setting).
Refer to:4.2 Digital Interface!

Bath Tank
An open tank that holds the thermofluid (fluid for thermoregulation).

Appendix

Huber Glossary

Calibration Thermostat (CAL)
A bath thermostat with an especially high degree of temperature stability and especially uniform temperature distribution across the bath fluid.

Calibration
Verification of the measured values displayed and – if required – recording of the values by which the displayed values deviate from the actual temperatures measured.Refer to Offset Calibration!

Cascade Control
Refer to Process Control!

Characteristic Temperature
This term denotes the operating temperature of a heating thermostat attained in stationary status with the heating switched off and the pump in operation.It depends on the motor capacity of the pump installed, the bath fluid used (viscosity, density) and the thermostat insulation, e.g. on whether the bath cover is used or not.

Compressor
A compressor is a machine for compressing gases and vapors.

Condenser
Apparatus in chillers for condensing refrigerant vapors (liquefier).

Cooler (Special Feature of the Unistat Tango)
This type of cooler is a heat exchanger through which water flows during the cooling process. This cooling water leads off just enough heat from the superheated refrigerant to prevent condensation.
The cooling water may flow only while the compressor is in operation, otherwise the cooling unit will not work!
It is important that only a tiny runnel of cooling water is routed through the cooler.

Cooling Capacity
Heat abstraction from the thermofluid via heat exchangers, using either ambient air or cooling water.
Refer to HT-Cooler!

Cooling/Heating Thermostat
A thermostat with a working temperature ranging from below to above the ambient temperature that is capable of abstracting heat from and transferring heat to the thermofluid.

Cooling Thermostat
A thermostat, the working temperature range of which is generally below the ambient temperature and that preferably abstracts heat from the thermofluid.HUBER cooling thermostats are actually cooling/heating thermostats since their working temperature ranges from below to above the ambient temperature and they are capable of abstracting heat from or transferring heat to the thermofluid.

Cooling Water
Water routed through the thermostat via a system of pipes for cooling the chiller.It should contain only minimized amounts of lime and corrosive substances that limit the life cycle of the system.

Cooling Circulators (IC, UC)
These are cooling thermostats designed in the form of circulators.Due to their special design (desktop, tower, without a readily accessible bath), their cooling performance and pumping capacities, they constitute a separate device group.

Circulation Pump
The circulation pump is used to circulate the thermofluid in a closed circuit.

Circulator (Unistat)
A circulator is a thermostat, in which the thermofluid is routed through an open or closed external circuit.Unistats feature a thermically decoupled, active surface (expansion tank). In this case, however, the surface temperature differs from the operating temperature. Unistats are not equipped with a bath.

Circulation (operating mode)
Circulation means that the thermofluid circuit is driven by the circulation pump.The heater and/or cooling unit are switched off.
Refer to 4.6!“Circulation“!

Controller
In a long tradition at Huber’s, „controller“ refers to an exchangeable electronics unit, including the required control and thermoregulation software.Refer to 2.3!

Control Mode – Internal, Process (Cascade)
Internal thermoregulation:The temperature measured inside the thermostat (flow temperature) is used as the controlled variable for thermoregulation.**Process (Cascade)**:With the aid of an external Pt100 sensor, a measuring point outside the thermostat in an external application, e.g. in a reactor, is used for controlling the thermoregulation process.
Refer to:4.4!

Appendix

Huber Glossary

Default Parameters
Refer to Factory Default!

Delivery Pressure
The overpressure of the circulation pump of a thermostat directly at the pressure connecting piece, measured using water.On charts, the delivery pressure is expressed as a function of the flow rate.

Delivery Suction
The suction of the circulation pump (force/ suction pump or duplex pump) directly at the pressure connecting piece, measured using water.On charts, the delivery suction is expressed as a function of the flow rate.

Digital Interface
The digital interface is used to transfer data in the form of bits and bytes.

DW-Therm
DW-Therm is a thermofluid (fluid for thermoregulation) developed especially for Unistats (Huber thermostats with a closed thermofluid circuit) with an extremely large temperature range (-90°C... + 200°C).

Effective Heating Capacity
Heat flow delivered to the thermofluid (fluid for thermoregulation) by heat sources.

Extended Working Temperature Range
This is the working temperature range extended to lower temperature levels achieved through the use of a manufacturer-defined cooling coil with cooling water.

External Control
The thermostat is not controlled through the controller but rather through external equipment. The external interface (multifunction socket) on the Polystat controller is used for this purpose.
Refer to:2.3, 4.3!

Factory Default
The manufacturer has set the thermostat parameters such that the temperature range constitutes only a minimum potential hazard and that the easiest and most probable thermoregulation process can be performed. Factory defaults ensure safe commissioning by the customer and can be changed using the respective programs according to the requirements of the customer
All the factory defaults can be restored with Program 52 with the thermoregulation function switched off.
Refer to 4.2 Restart!

Fill Volume
Refer to Bath Volume!

FL
Refer to Safety Classification!

Flash Point
The term flash point denotes the temperature to which a fluid in an open bath must be heated until the vapor/air mixture at the surface inflames when a flame is held near and keeps on burning once the flame is pulled away again.
Refer to:3.1.EN 61010-1EN 61010-1

Flash Display
Controller display shown for seconds only to provide information on a variable intermediate program status or the status of the thermostat.

Flow Cooler (DC)
Flow coolers are „uncontrolled" coolers without circulation pump that are interposed in an external circuit, thus expanding the functionality of a heating thermostat to that of a heating/ cooling thermostat.On the one hand, they substitute water cooling, on the other hand they are used to achieve lower temperatures.

Flow Rate
The volume of fluid delivered by the circulation pump per time unit, measured using water.In charts, the flow rate is expressed as a function of the delivery pressure (backpressure).

Flow Temperature
Denotes the temperature measured directly at the point the thermofluid exits the circulator.

Force Pump
The force pump is used for circulating the thermofluid in a closed external circuit and for mixing the fluid inside the bath.

Force/Suction Pump
This pump features a pressure and a suction stage driven by a common motor.The pressure stage forces the thermofluid from the thermostat into the circuit, the suction stage draws the fluid back into the thermostat.A force/suction pump can be used in a closed circuit.Benefit:In the external circuit, the pressure drops from positive values (pressure) to negative values (vacuum) and is almost zero in the consumer.This facilitates thermoregulation even of pressure-sensitive glass vessels.In addition, an open external circuit (bath tank) can also be connected using a force/suction pump.A level constant ensures a constant fluid level in the bath tank.

Appendix

Huber Glossary

Gas Venting
A special thermoregulation regime common to all Unistats with the objective of perpetuating thermoregulation despite the evaporation of part of the thermofluid and removing elements with a low boiling point from the circulating fluid.

Heat Exchanger
Refer to Thermofluid!

Heating Capacity
The maximum installed electrical power of the heating element.The heating capacity depends on the voltage of the thermostat. It is continuously controlled and reduced when the temperature approaches the programmed setpoint.

Heating Thermostat
A thermostat, the working temperature range of which is primarily above the ambient temperature and that preferably transfers heat to the thermofluid.

High-Pressure Cut-Out
Protective switch installed in Unistats and Unichillers. Shuts down the thermoregulation process if the pressure on the compressor side gets too high.
Applies to Unistats:Once the Unistat is ready for operation again, the high-pressure cut-out must be pressed to continue the thermoregulation process.
Applies to Unichillers:Once the Unichiller is ready for operation again, the high-pressure cut-out is automatically reset to continue the thermoregulation process.(The high-pressure cut-out in Unichillers is not readily accessible.)

High-Temperature Cooler (with Unistats):
A high-temperature cooler is a heat exchanger that cools a thermofluid from a high temperature to ambient temperature using air or water.It is installed on the upstream side of chillers in Unistats, thus reducing the load on and power consumption of the chiller.

High Temperature Stage (HT) (with Unistats):
The ***high temperature stage*** is the top temperature stage in multi-stage chillers in Unistats.Sie dient der It is used for decreasing the temperature (to approx. −60 °C) and thus prepares ***the low temperature stage*** (NT) in dual-stage chillers (to approx. −90 °C) and the medium temperature stage (MT) (to approx. −90 °C) in triple-stage chillers.
The medium temperature stage, in turn, prepares the low temperature stage (currently to approx. −120°C) in triple-stage chillers.

Housing Volume
The volume resulting from the outer dimensions of the thermostat.

Immersion Cooler (TC)
A cooling unit with a flexible tube and a cooling coil (evaporator) for immersing in baths.

Immersion Thermostat (Polystats cc)
This is a thermostat that is combined with a bath tank that forms an independent unit. Immersion thermostats are equipped with a screw clamp for fastening to the walls of any desired bath tank. Using a bath bridge, they can also be mounted permanently on the top of a bath or on a tripod.

Industrial Thermostats (IC-Hx)
Industrial thermostats are cooling circulators (intelligent chillers) with a heater installed ex works.High cooling, heating and pumping capacities and small fluid volumes provide for fast cooling and heating rates.They are ideal for thermoregulation in process engineering and in narrow temperature ranges (-20..120 °C).

Interface, Analog (AIF)
This interface is used for entering the temperature setpoint or for the output of the actual temperature as an analog value in the form of a signal current (4 – 20 mA).Refer to: 4.3 Main Menu, Analog Interface!

Interface, Digital
This interface is used for digital data exchange between connected devices via the data communication line.Primarily, the temperature setpoint and the actual temperature are transmitted.Refer to:RS232 Interface and RS485 Interface!
Refer to:4.3 Main Menu, Digital Interface!

Low Temperature Stage (NT)
Refer to High Temperature Stage!

Medium Temperature Stage (MT)
Refer to High Temperature Stage!

Net Refrigerating Capacity
The effective capacity the cooling thermostat or chiller provides for the application.The frictional heat generated by the circulation pump and the heat entering due to insufficient insulation has already been taken into consideration.The cooling capacity data correspond to those of the net cooling capacity.

NFL
Refer to Safety Classification!

Appendix

Huber Glossary

Offset Calibration
Single-point correction of a temperature sensor at a specific temperature.

Operating Temperature Range
This is the temperature range that is limited by the permissible minimum and maximum operating temperature.
The operating temperature is the temperature to which the thermostat is permitted to heat/cool the respective thermofluid.

Overheat Point
Refer to Overtemperature Protection!

Overheating
The positive difference between the condensation temperature and the topical gas temperature of the thermofluid in the refrigeration process.(Does not refer to overheating of the thermostat!)

Overtemperature Protection
Each thermostat with a heater is equipped with overtemperature protection.The overtemperature protection operates completely independently of the controller, i.e. it cannot be influenced through the controller.If at least one temperature sensor reports that a temperature limit has been exceeded, the thermoregulation process is stopped at once.
In the course of the thermoregulation process, the sensors are tested for short-circuit and failure.The overtemperature protection is equipped with an independent protective device.This prevents very reliably that the OK status is set while there actually is a failure.The user sets the upper temperature limit (overheating point) at the thermostat to ensure operating safety considering the respective thermofluid used.For unattended operation, the overheating point must be selected at least at 25 °C below the flash point of the respective thermofluid.
Refer to:4.1.!

Physical Mass
The physical mass describes the P/I parameters (P = proportional, I = integral) of the thermostat.
Refer to 4.3!

PLC (Programmable Logic Controller)
A device frequently used in industrial environments to control operations and processes.

POKO Potential-Free Contact
The POKO is a change-over contact led through the POKO relay in the controller.The potential-free contact is designed for an ohm resistive load up to 30 V and a maximum of 0.1 A.
The following applies due to the working current principle:The OK status is linked to the flow of current through the winding of the POKO relay.The POKO relay is actuated by the controller.

Presettings ex Works, Factory Default
Refer to Factory Default!

Process Control (often referred to as cascade control)
Thermoregulation of a connected application.A temperature sensor (Pt100) that forms part of this application is connected with the thermostat.The actual temperature of the external process is sensed and the operating temperature of the thermostat is permanently calculated and adapted.Depending on the operating temperature, losses through insulation and exothermic reactions, the operating temperature at the consumer may deviate considerably from the setpoint and the actual temperature of the application.
Caution!
Observe the safety-relevant limits of the respective thermofluid!
Refer to: 3.1!

Process Temperature (core temperature)
The process temperature is the temperature measured in the core of a connected external application (in the case of external thermoregulation), e.g. in the core of a reactor.Refer to Reactor, Control Mode!
Refer to 2.3!

Program (Thermostat Control)
The controller provides programs for the user to operate the thermostat, to monitor and control statuses and processes, i.e. to perform the thermoregulation process.

Program (Thermoregulation)

Refer to Thermoregulation Program!

Appendix

Huber Glossary

Protocols
Protocols are used in digital data exchange.The term protocol in this context refers to a series of rules drawn up to facilitate uniform coding and decoding of the bits and bytes in the data exchange between the thermostats and a **PC, PCS** (process control system) or the like. For the time being, three protocols are available, two of which have been implemented in all the product families.
PPP – point-to-point protocol
LAI protocol – for communication on the bus
The third protocol is currently only available for UNISTATS on demand:
Modbus protocol (Modbus RTU by Gould)

Ramp
A desired change in temperature within a given time period.A ramp is a simple and frequently used programming of a temperature pattern.
A ramp is defined through the increase in temperature and the thermoregulation period.Once the ramp has been entered, the thermoregulation process thus programmed will start automatically and does not require any accompanying inputs.
The ramp can be defined through different combinations and sequences of the setpoint, time period and slope parameters.
Where:Time period = slope x setpoint
Refer to 4.3!

Reactor
A cylindrical double-wall tank permitting the flow of thermofluid through the space between the two walls (jacket).The temperature of the thermofluid is transferred to the hollow space inside (core) via the inner wall.The user fills his/her reagents into this core with separate inlets and outlets.The reagents are then indirectly thermoregulated via the inside reactor wall.

Refrigerant
Refrigerant is used in the circuit of the refrigerating unit. It abstracts heat from the thermofluid when the compressed gas expands and evaporates in the evaporator.Since 1992 HUBER has been using exclusively refrigerants free of CFC and since 1994 even refrigerants free of HCFC (e.g. R22) that do not harm the ozone layer, i.e. with an ODP (Ozone Depletion Potential) = 0 and minimized GWP (Global Warming Potential, e.g. greenhouse effect).

Refrigerating Capacity
Heat flow led off from the thermofluid by a cooling thermostat.

Remote Control Panel
An external device permitting operation of the Huber Thermostat (possibly to a limited extent).

RS232 Interface
Digital controller interface for digital data exchange between connected devices via the data communication line.RS232 is a serial interface for establishing a point-to-point connection.This means, only two parties, e.g. a thermostat and a PC, can communicate via the interface at a time.
Refer to:4.3. Main Menu, Digital Interface!

RS 485 Interface
Digital controller interface for digital data exchange between connected devices via the data communication line.
As many as 32 parties can be connected to the RS 485 interface.Each party connected to this bus system is assigned a separate address.
Refer to:4.3 Main Menu, Digital Interface!

Safety Classificatio
In thermostats, either flammable (FL) or non-flammable (NFL) thermofluids (fluids for thermoregulation) may be used.The respective safety requirements have been laid down in DIN EN 61010-2-010.Accordingly, the safety classes NFL and FL have been defined:
NFL (non-flammable) thermostats with integrated overheating protection exclusively for non-flammable fluids.
Thermostats classified as FL (flammable) with adjustable overtemperature protection and low liquid level protection for flammable fluids (all the thermostats by HUBER).Refer to: 3.1. Thermofluid!

Segment
A segment is a section of a thermoregulation program that creates part of the temperature graph.
Refer to Thermoregulation Program!

Setpoint
The temperature to be attained at the point of measurement.
Maximum Setpoint
The highest permissible temperature to be attained at the point of measurement.
Minimum Setpoint
The lowest permissible temperature to be attained at the point of measurement.
Refer to Control Mode! Refer to 4.4!

Self Test
Test procedure in the thermostat after power on.It ensures the operability of safety-relevant thermostat components.
(Except for overtemperature protection!)

Appendix

Huber Glossary

Slave Address
Identification number the user assigns to a device to be able to allocate data to a device during data transfer via a bus protocol.
Refer to RS485 Interface!

Standards
Safety requirements for electrical laboratory equipment, and especially for thermostats, have been laid down in European Standards EN 61010-1 and EN 61010-2-010.These standards replace DIN 12879 (among others).The terms and characteristic data are defined in DIN 12876-1 and DIN 12876-2.

Status Display
Display on the UniCop that provides information on the topical program sequence or the topical status of the thermostat.

Temperature – internal, external
Refer to Control Mode, Flow Temperature!

Temperature Stability
This term denotes half the difference in temperature between the highest and lowest temperature measured in a thermostat within 30 minutes of attaining a stable value at a specific setpoint.
This value is determined at a temperature of 70 °C (with water) for a heating thermostat and at – 10 °C (with ethanol) for a cooling thermostat.
Also refer to DIN 12876!

Temperature-Stable
Refer to Thermoregulation, Temperature-Stable!

Tempmove
Possibility of selecting the input variable for a temperature.

Thermofluid (fluid for thermoregulation)
Heat transfer fluid. Transfers energy from the thermostat to the application or vice versa.Siehe 3.1.!Refer to 3.1!

Thermoregulation Fluid
Refer to Thermofluid!

Thermoregulation
The active manipulation of the temperature of a substance to the desired level by way of heating and/or cooling

Thermoregulation, Temperature-Stable!
Thermoregulation will proceed until the defined setpoint has been attained (if required, the desired thermoregulation period will be exceeded).**Thermoregulation, Time-Stable**
Thermoregulation will proceed according to the time targets entered (independent of the actual temperature value attained).

Thermoregulation Program
Arrangement of segments to create a reproducible temperature graph.
A thermoregulation program is made up of a program header and segments.
The program header is independent of the number of active segments and the length of the program.
The program header defines the beginning and sequence of the thermoregulation program (TP) through the following parameters:
Start temperature of segment n, start ramp slope (K/min),
thermoregulation priority: temperature-stable or time stable, program end.
A segment is defined through the following parameters:
Start temperature of segment n (= end temperature of segment n-1)
Segment period n
Refer to:4.4!

Transparent Bath
A bath thermostat with transparent walls permitting direct monitoring of the object to be thermoregulated (Polystat models A5 – A18).

Appendix

Huber Glossary

Usable Depth
The term usable depth denotes the depth of fluid available in the bath thermostat for direct thermoregulation.

Venting
A special circulation regime with the objective of driving out air trapped in closed thermofluid circuits.

Viscosity
Also referred to as „internal friction“. The characteristic of a fluid or gaseous substance that causes frictional tension at deformation in addition to the thermodynamic pressure, which counteracts the displacement of fluid or gaseous particles relative to each other.
This means, the viscosity of a substance increases with decreasing temperature, the available cooling capacity decreases with decreasing temperature.
The viscosity of thermofluids (fluids for thermoregulation) used in HUBER thermostats should not exceed 50 mm^2/s.

Watchdog
A safety device in the electronic control, similar to the principle of the dead-man switch in an engineThe functionality of the system is monitored through the input signals.

Warning
Device messages pointing to irregularities in the operation of the thermostat that do not cause a shut-down of the thermoregulation process. Programs 10, 11.

Working Temperature Range
The working temperature range is defined with reference to an ambient temperature of 20 °C. It denotes the temperature range the thermostat attains by itself without the aid of any auxiliary means*, utilizing electric power only.
Due to the heat transferred by the pump motor and the insulation, the working temperature range of heating thermostats starts above the ambient temperature and ends at the upper limit of the operating temperature.
In the case of cooling/heating thermostats, the working temperature ranges from the lower limit of the operating temperature to the upper limit of the operating temperature. This temperature is permissible in continuous operation with the use of a chiller.
The working temperature range of cooling thermostats starts at the lower limit of the operating temperature and ends at ambient temperature.*Auxiliary means are additional heating elements, coolers, heat exchangers or fans.

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Immersion Circulators a powerful pressure pump, moistened parts in stainless steel or high-resistant plastics. Plug & Play Technology – new generation of microprocessor controlled Compatible Control. Easy handling and digital display. Alarms and warnings are given both visually and acoustically. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and calibrateable. 3 year warrantee (registration necessary).
Simple user guidance using a 4-line graphical LCD display and a high performance processor. Administrator level and user levels are available.
CC1: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2, similar to **CC1**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function.
CC3, similar to **CC2**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko). This allows it to be easily integrated into a laboratory automation system.

Technische Daten	**Technical Data**		**CC1-E**	**CC2-E**	**CC3-E**
Arbeitstemperaturbereich	Operating temperature range		25...200°C		
mit Wasserkühlung	with water cooling		20...200°C		
mit Kühlgerät	with refrigeration chiller		-30..200°C		
Temperaturkonstanz bei 70°C (15l)	Temperature stability at 70°C (15 l bath)		0,02 K		
Temperatureinstellung	Temperature adjustment		digital		
Temperaturanzeige	Temperature indication		digital		
Absolutgenauigkeit	absolute accuracy		kalibrierbar / setup for calibration		
Temperaturfühler	Temperature sensor		Pt100		
Anschluss externer Fühler	connection external sensor	---	Pt100	Pt100	
Analoges Interface In/Out	Analogue interface in/out	---	---	4-20mA alternativ 1-5V	
Sicherheitsklasse	Safety classification		FL		
Heizleistung	Heating capacity		2,0 kW		
Druckpumpe Anschluss 12mm	Pressure pump adapter nom 12 mm		10 l/min max.		
Anschluss 8 mm	adapter nom 8 mm		7 l/min max.		
Förderhöhe (Druck)	Delivery pressure (head)		0,2 bar max.		
Pumpenanschluss	Pump connection		M16x1		
Pumpenanschluss für Schlauch	Pump connections for hose		NW8/ NW12		
Max. zulässige kin. Viskosität	Max. permissible kin. viscosity		$50mm^2/s$		
Abmessungen B x T x H	Overall dimension w x l x h		120 x 135 x 300 mm		
Eintauchtiefe	immersion depht		150 mm		
Gewicht, netto	Net weight		4,0 kg		
Netzanschluss	Power supply requirement		230 V 1~ 50/60 Hz		
max. Stromaufnahme	max. current		10 Amp.		
Absicherung	Fuse		16 Amp.		
min./max. Umgebungstemperatur	min./max. ambient temperature		5°C / 32°C		
Bestell-Nr.	**Order-No.**	**688.0051**	**688.0052**	**688.0053**	

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Accessories and periphery: Cooling coil for water-cooling, pump adapter for external thermoregulation, adapter nom. dia 8 mm/ 12 mm, dummy plugs and sleeve nuts thread M16x1, connection tubes. Waterbathes in polycarbonate and stainless steel and units for cells and glasses.

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Immersion Circulators a powerful pressure pump, moistened parts in stainless steel or high-resistant plastics.
Plug & Play Technology – new generation of microprocessor controlled Compatible Control. Easy handling and digital display. Alarms and warnings are given both visually and acoustically. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and calibrateable. 3 year warrantee (registration necessary).
Simple user guidance using a 4-line graphical LCD display and a high performance processor. Administrator level and user levels are available.
CC1: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2, similar to **CC1**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function.
CC3, similar to **CC2**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko). This allows it to be easily integrated into a laboratory automation system.

Technische Daten	**Technical Data**	**CC1-202C**	**CC2-202C**	**CC3-202C**
Arbeitstemperaturbereich	Operating temperature range		45...200°C	
Mit Wasserkühlung	With water cooling		20...200°C	
Mit Kühlgerät	With refrigeration chiller		-30..200°C	
Temperaturkonstanz bei 70°C (15 l)	Temperature stability at 70°C (15 l)		0,05 K	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Absolutgenauigkeit	Absolute accuracy	kalibrierbar / setup for calibration		
Temperaturfühler	Temperature sensor		Pt100	
Externer Temperaturfühler	External Temperature sensor	---	Pt100	Pt100
Analoges Interface In/Out	Analogue interface in/out	---	---	4-20mA alternativ 1-5V
Sicherheitsklasse	Safety classification		FL	
Heizleistung	Heating capacity		2,0 kW	
Druckpumpe Anschluss 12mm	Pressure pump adapter nom 12 mm		10 l/min max.	
Anschluss 8 mm	adapter nom 8 mm		7 l/min max.	
Förderhöhe (Druck)	Delivery pressure (head)		0,2 bar max.	
Pumpenanschluss	Pump connection		M16x1	
Pumpenanschluss für Schlauch	Pump connections for hose		NW8/ NW12	
Max. zulässige kin. Viskosität	Max. permissible kin. viscosity		50mm^2/s	
Füllvolumen max.	Bath capacity max.		1,8 lit.	
Badöffnung Durchmesser / Tiefe	Bath opening diameter / depth		25 mm / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		150 x 215 x 350 mm	
Arbeitshöhe Bad	Heith of bath opening		190 mm	
Gewicht, netto	Net weight		7,5 kg	
Netzanschluss	Power supply requirement		230 V 1~ 50/60 Hz	
Max. Stromaufnahme	Max. current		10 Amp.	
Absicherung	Fuse		10 Amp.	
Min./max. Umgebungstemperatur	min./max. ambient temperature		5°C / 32°C	
Bestell-Nr.	**Order-No.**	**685.0011**	**685.0012**	**685.0013**
				V2.2/05

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Accessories and periphery: Adapter nom. dia 8/12 mm*, dummy plugs* and sleeve nuts thread M16x1*, adapter nom. dia 8 mm, connection tubes, displacement rack for reducement of the bath volume.
* standard equipment

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Accessories and periphery: Adapter nom. dia 8/12 mm*, dummy plugs*, sleeve nuts thread M16x1*, connection tubes, external sensor (not with CC1), connection cable, stopcock, displacement rack for reducement of the bath volume, calibration insert. * standard equipment

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Accessories and periphery: Adapter nom. dia 8/12 mm*, dummy plugs*, sleeve nuts thread M16x1*, connection tubes, external sensor (not with CC1), connection cable, stopcock, displacement rack for reducement of the bath volume, calibration insert.

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Bridge Thermostat with telescopic arm (up to 884 mm). Housing and all moistened parts in stainless steel.
With cooling coil for water-cooling, pressure- and suction pump.
Plug & Play Technology – new generation of microprocessor controlled Compatible Control. Easy handling and digital display. All controllers (CC1, CC2, CC3) are backward compatible, interchangeable and calibrateable. 3 year warrantee (registration necessary).
Simple user guidance using a 4-line graphical LCD display and a high performance processor. Administrator level and user levels are available.
CC1: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2, similar to **CC1**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function.
CC3, similar **CC2**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko).

Technische Daten	**Technical Data**	**CC1-300B4X vpc**	**CC2-300B4X vpc**	**CC3-300B4X vpc**
Arbeitstemperaturbereich	Operating temperature range		28...300°C	
mit Kühlgerät	with refrigeration chiller		-20...300°C	
Temperaturkonstanz bei 70°C	Temperature stability at 70°C		0,02 K	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Analoges Interface In/Out	Analogue interface in/out	---	---	4-20mA alternativ 1-5V
Temperaturfühler	Temperature sensor		Pt100	
externer Temperaturfühler	external Temperature sensor	---	Pt100	Pt100
Sicherheitsklasse	Safety classification		FL	
Heizleistung umrüstbar	Heating capacity converible		3 kW oder/or 4kW	
Druckpumpe	Pressure pump		33 l/min max.	
Förderhöhe (Druck)	Delivery pressure (head)		0,7 bar max.	
Saugpumpe	Suction pump		22 l/min max.	
Saughöhe (Sog)	Delivery suction pressure (head)		0,4 bar max.	
Pumpenanschluß	Pump connection		M16x1	
Max. zulässige kin. Viskosität	Max. permissible kin. viscosity		50 mm^2/s	
Eintauchtiefe	Immersion depth		180 mm	
Abmessungen B x T x H	Overall dimension w x d x h		345 x 190 x 395 mm	
Gewicht, netto	Net weight		13 kg	
Netzanschluss umrüstbar	Power supply requirement converible		230V 1~ 50/60Hz oder/or	
			400V 3~N 50Hz	
Max. Stromaufnahme	max. current		16 A oder/or 10 A	
Absicherung	Fuse		20 A oder/or 10 A	
Bestell-Nr.	**Order-No.**	**684.0011BX**	**684.0012BX**	**684.0013BX**
				V2.4/05

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Accessories and periphery: Adapter nom. dia 8/12 mm*, dummy plugs*, sleeve nuts thread M16x1*, connection tubes, external sensor (not with CC1), connection cable, stopcock, displacement rack for reducement of the bath volume, calibration insert * standard equipment

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Bath and circulation-thermostat. Housing, bath and all moistened parts in stainless steel. With cooling coil for water-cooling, pressure- and suction pump.
Plug & Play Technology – new generation of microprocessor controlled Compatible Control. Easy handling and digital display. All controllers (CC1, CC2, CC3) are backward compatible, interchangeable and calibrateable. 3 year warrantee (registration necessary).
Simple user guidance using a 4-line graphical LCD display and a high performance processor. Administrator level and user levels are available.
CC1: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2, similar to **CC1**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function.
CC3, similar **CC2**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko).

Technische Daten	**Technical Data**	**CC1-304B vpc**	**CC2-304B vpc**	**CC3-304B vpc**
Arbeitstemperaturbereich	Operating temperature range		28...300°C	
Mit Kühlgerät	With refrigeration chiller		-20..300°C	
Temperaturkonstanz bei 70°C	Temperature stability at 70°C		0,02 K	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Absolutgenauigkeit	Absolute accuracy		kalibrierbar / setup for calibration	
Temperaturfühler	Temperature sensor		Pt100	
Externer Temperaturfühler	External Temperature sensor	---	Pt100	Pt100
Analoges Interface In/Out	Analogue interface	---	4-20mA alternativ 1-5V	
Sicherheitsklasse	Safety classification		FL	
Heizleistung	Heating capacity		2,0 kW	
Druckpumpe bei Anschluss 12mm	Pressure pump with 12 mm connection		33 l/min max.	
Förderhöhe (Druck)	Pressure		0,7 bar max.	
Saugpumpe bei Anschluss 12mm	Suction pump with 12 mm connection		22 l/min max.	
Saughöhe (Sog)	Suction		0,4 bar max.	
Pumpenanschluss	Pump connection		M16x1	
Max. zulässige kin. Viskosität	Max. permissible kin. viscosity		$50mm^2/s$	
Füllvolumen max.	Bath capacity max.		4 lit.	
Badöffnung B x T / Tiefe	Bath opening w x l / depth		100 x 130 / 145 mm	
Abmessungen B x T x H	Overall dimension w x l x h		210 x 310 x 390 mm	
Arbeitshöhe Bad	Heith of bath opening		190 mm	
Gewicht, netto	Net weight		13 kg	
Netzanschluss	Power supply requirement		230 V 1~ 50/60 Hz	
Max. Stromaufnahme	Max. current (1 Phase)		10 A	
Absicherung	Fuse		10 A	
Schutzklasse	Protection class		IP20	
Min./max. Umgebungsbedingungen	Min./max. ambient temperature		5°C / 32°C	
Bestell-Nr.	**Order-No.**	**679.0001**	**679.0002**	**679.0003**

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Accessories and periphery: Adapter nom. dia 8/12 mm*, bath cover 304 vpc complete*, dummy plugs*, sleeve nuts thread M16x1*, connection tube silicon*, power supply cable*, connection tubes, external sensor (not with CC1), connection cable, stopcock.
* standard equipment

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Refrigerated-Bath Thermostat with air-cooled refrigeration unit. Powerful pressure pump, moistened parts in stainless steel or high-resistant plastics.
Plug & Play Technology – with a new generation of microprocessor compatible controller. Easy handling and digital display.
All controllers are backwars compatible, interchangeable and calibrateable. 3 years guarantee (registration necessary).
Much improved user guidance and expanded with a LCD graphics display and a powerfull microprocessor. Administrator level and several user levels included.
CC1-K6: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2-K6, similar to **CC1-K6**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function.
CC3-K6, similar to **CC2-K6**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko).

Technische Daten	**Technical Data**	**CC1-K6**	**CC2-K6**	**CC3-K6**
Arbeitstemperaturbereich	Operating temperature range		-25...200 °C	
Temperaturkonstanz bei 70°C	Temperature stability at 70°C		0,02 K (DIN 58966)	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Absolutgenauigkeit	absolute accuracy	kalibrierbar / setup for calibration		
Temperaturfühler	Temperature sensor		Pt 100	
Analoges Interface In/Out	Analogue interface in/out	---	---	4-20mA alternativ 1-5V
Sicherheitsklasse	Safety classification		FL	
Heizleistung	Heating capacity		2,0 kW	
Kälteleistung bei 0°C	Cooling capacity at 0°C		0,15 kW	
bei -10°C	at -10°C		0,10 kW	
bei -20°C	at -20°C		0,05 kW	
Kältemittel	Refrigerent		R134a	
Füllvolumen max.	Bath capacity max		4,5 l	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		150 x 140 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		200 x 370 x 535 mm	
Druckpumpe Anschluß 12mm	Pressure pump (adapter nom 12 mm)		10 l/min max.	
Anschluß 8 mm	(adapter nom 8 mm)		7 l/min max.	
Förderhöhe (Druck)	Delivery pressure (head)		0,2 bar max.	
Pumpenanschluss	Pump connection		M16x1	
Pumpenanschluss für Schlauch	Pump connection for hose		NW8/12	
Gewicht, netto	Net weight		20 kg	
Netzanschluss	Power supply requirement		230 V 1~ 50 Hz	
max. Stromaufnahme	max. current		10 Amp.	
Absicherung	fuse		16 A	
Bestell-Nr.	**Order-No.**	**666.0021**	**666.0022**	**666.0023**
Gültig ab Fert.Nr.	**From serial no.**		**60825**	**V1.3/05**

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Accessoires and periphery: adapter nom. dia 8 mm*/ 12 mm*, dummy plugs* and sleeve nuts thread M16x1*, external sensor (not with cc1), connection tubes, connection cable * standard equipment

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Heating- bath thermostat consisting of immersion circulators CC-E, transparent baths and bath bridge.
Baths made of polycarbonate with temperature range up to 100°C. Bath bridge for 111A and 118A with bore hole for cooling probe (e.g. for immersion cooler TC45-TC100E).
name:
example: CC2-105A tansparent bath with 5 l bath capacity, in combination with Compatible Control CC2
CC3-111A tansparent bath with 11 l bath capacity, in combination with Compatible Control CC3

Technische Daten	**Technical Data**	**CC1-105A**	**CC2-105A**	**CC3-105A**
Füllvolumen.	Bath capacity		5 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		120 x 110 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		132 x 280 x 320 mm	
Gewicht, netto	Net weight		5 kg	
Testglaseinsätze Typ A bis H	tube racks typ A to H		1 St.	
Bestell-Nr.	**Order-No.**	**405.0001**	**405.0002**	**405.0003**

Technische Daten	**Technical Data**	**CC1-106A**	**CC2-106A**	**CC3-106A**
Füllvolumen.	Bath capacity		6 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		120 x 210 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		132 x 380 x 320 mm	
Gewicht, netto	Net weight		5 kg	
Testglaseinsätze Typ A bis H	tube racks typ A to H		2 St.	
Bestell-Nr.	**Order-No.**	**405.0011**	**405.0012**	**405.0013**

Technische Daten	**Technical Data**	**CC1-107A**	**CC2-107A**	**CC3-107A**
Füllvolumen	Bath capacity		7 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		120 x 310 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		132 x 480 x 320 mm	
Gewicht, netto	Net weight		6 kg	
Testglaseinsätze Typ A bis H	tube racks typ A to H		3 St.	
Bestell-Nr.	**Order-No.**	**405.0021**	**405.0022**	**405.0023**

Technische Daten	**Technical Data**	**CC1-111A**	**CC2-111A**	**CC3-111A**
Füllvolumen.	Bath capacity		11 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		302 x 165 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		333 x 313 x 325 mm	
Gewicht, netto	Net weight		6,5 kg	
Testglaseinsätze Typ 1 oder 2	tube racks typ 1 or 2		2 St.	
Bestell-Nr.	**Order-No.**	**405.0031**	**405.0032**	**405.0033**

Technische Daten	**Technical Data**	**CC1-118A**	**CC2-118A**	**CC3-118A**
Füllvolumen	Bath capacity		18 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		302 x 325 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		333 x 473 x 325 mm	
Gewicht, netto	Net weight		8 kg	
Testglaseinsätze Typ 1 oder 2	tube racks typ 1 or 2		4 St.	
Bestell-Nr.	**Order-No.**	**405.0041**	**405.0042**	**405.0043**

Zubehör	**Accessoires:**	**105A**	**106A**	**107A**	**111A**	**118A**
Entleerungshahn	drain cock	6026	6026	6026	6026	6026
Badbrücke	bath bridge	6309	6309	6309	6310	6310

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Bath Thermostat consisting of immersion circulators CC-E, insulated stainless steel bath with stainless steel housing and bath bridge. Bath bridge with borehole for cooling probe (e.g. for immersion cooler TC45-TC100E).
Name:
Example: CC2-208B stainless steel bath with 8 l bath capacity, in combination with CC2-E
CC3-225B stainless steel bath with 25 l bath capacity, in combination with CC3-E

Technische Daten	**Technical Data**	**CC1-208B**	**CC2-208B**	**CC3-208B**
Arbeitstemperaturbereich	Operating temperature range		-30...200°C	
Füllvolumen.	Bath capacity		8,5 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		235 x 165 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		290 x 350 x 375 mm	
Gewicht, netto	Net weight		10 kg	
Bestell-Nr.	**Order-No.**	**407.0001**	**407.0002**	**407.0003**

Technische Daten	**Technical Data**	**CC1-212B**	**CC2-212B**	**CC3-212B**
Arbeitstemperaturbereich	Operating temperature range		-30...200°C	
Füllvolumen	Bath capacity		12 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		295 x 193 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		350 x 375 x 375 mm	
Gewicht, netto	Net weight		8,5 kg	
Testglaseinsätze Typ 1 bis 4	tube racks typ 1 to 4		2 St.	
Bestell-Nr.	**Order-No.**	**407.0011**	**407.0012**	**407.0013**

Technische Daten	**Technical Data**	**CC1-215B**	**CC2-215B**	**CC3-215B**
Arbeitstemperaturbereich max	Operating temperature range max		-30...200°C	
Füllvolumen.	Bath capacity		15 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		295 x 193 / 200 mm	
Abmessungen B x T x H	Overall dimension w x l x h		350 x 375 x 420 mm	
Gewicht, netto	Net weight		12 kg	
Testglaseinsätze Typ 1 bis 4	tube racks typ 1 to 4		2 St.	
Bestell-Nr.	**Order-No.**	**407.0021**	**407.0022**	**407.0023**

Technische Daten	**Technical Data**	**CC1-220B**	**CC2-220B**	**CC3-220B**
Arbeitstemperaturbereich max	Operating temperature range max		-30...200°C	
Füllvolumen.	Bath capacity		20 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		295 x 370 / 150 mm	
Abmessungen B x T x H	Overall dimension w x l x h		350 x 555 x 375 mm	
Gewicht, netto	Net weight		11 kg	
Testglaseinsätze Typ 1 bis 4	tube racks typ 1 to 4		4 St.	
Bestell-Nr.	**Order-No.**	**407.0031**	**407.0032**	**407.0033**

Technische Daten	**Technical Data**	**CC1-225B**	**CC2-225B**	**CC3-225B**
Arbeitstemperaturbereich max	Operating temperature range max		-30...200°C	
Füllvolumen.	Bath capacity		25 lit.	
Badöffnung B x T / Badtiefe	Bath opening w x l / bath depth		295 x 370 / 200 mm	
Abmessungen B x T x H	Overall dimension w x l x h		350 x 555 x 420 mm	
Gewicht, netto	Net weight		14,5 kg	
Testglaseinsätze Typ 1 bis 4	tube racks typ 1 to 4		4 St.	
Bestell-Nr.	**Order-No.**	**407.0041**	**407.0042**	**407.0043**

Zubehör	**Accessoires:**	**208B**	**212B**	**215B**	**220B**	**225B**
Entleerungshahn	drain cock	6026	6026	6026	6026	6026
Badbrücke	bath bridge	6303	6308	6308	6308	6308
variabl. Stellboden	var. Bottom		6297	6297	6298	6298
Bad-Deckel vorn	bath cover front	6214	6023	6023	6023	6023
Bad-Deckel hinten	bath cover back				6024	6024
Bad-Deckel einteilig	bath cover 1-piece				6025	6025

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Circulation cooler with air-cooled refrigerating unit and circulation pump.Evaporator (cooler), tank and housing of stainless steel. Pressure- suction pump of plastic.
Plug and Play technology –new generation of microprocessor controlled Compatible Control.
Simple operation with a rotary knob and digital display. Alarms and warnings are given both visually and acoustically. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and can be calibrated. 3 year warrantee. Simple user guidance using a 4 line graphical LCD display and high performance processor. Administrator level and multiple user levels are available.
CC1: Maximum and minimum setpoint for additional safety.
CC2, like **CC1:** With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with a maximum of 5 segments).
CC3; like **CC2:** Including however the complete functions: Programmer (50 segments, may be split into 10 programs), interactive, contains a digital RS232 or RS485 interface as well as a 4 to 20mA analogue interface for bi-directional communication and a potential free alarm contact (Poko).

Technische Daten	**Technical Data**	**minichiller cc1**	**minichiller cc2**	**minichiller cc3**
Arbeitstemperaturbereich	Operating temperature range		-20...40°C	
Temperaturkonstanz bei -10°C	Temperature stability at -10°C		1 K	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Auflösung der Anzeige	Resolution		0.1 K	
Temperaturfühler	Temperature sensor		Pt100	Pt100
Anschluss externer Fühler	Connection external sensor	---	Pt100	Pt100
Analoges Interface In/Out	Analogue interface in/out	---		4-20 mA alternativ 1-5V
Sicherheitsklasse	Safety classification		NFL	
max. Förderleistung der Umwälzpumpe	max. capacity circulation pump		12 l/min	
Förderhöhe (Druck) max.	Pressure of pump max.		0,2 bar	
Förderleistung der Saugpumpe max.	Suction pump max.		12 l/min	
Förderhöhe (Saug) max.	Suction max.		0,1 bar	
Pumpenanschluss	Pump connection		M16x1	
Pumpenanschluss für Schlauch	Pump connections for hose		NW8/12	
max. zulässige kin. Viskosität	max. permissible kin. viscosity		50 mm^2/s	
min. Füllvolumen	min. filling capacity		1,4 l	
Ausdehnungsvolumen	Volume of expansion		2,6 l	
Abmessungen B x T x H	Overall dimension w x l x h		225 x 360 x 380 mm	
Netzanschluss	Power supply requirement		230 V 1~ 50/60 Hz	
max. Stromaufnahme	max. current (1 phase)		2,8 A	
Absicherung	Fuse		6 A	
max./min. Umgebungstemperatur	max./min. ambient temperature		32°C / 5°C	
Gewicht	Weight		23kg	
Bestell-Nr.	**Order-No.**	**466.0011**	**466.0012**	**466.0013**
				V1.1/04

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Accessories and periphery: Adapter nom. dia 8 mm*/ 12 mm*, sleeve nuts thread M16x1*, connection tubes *standard equipment

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Circulation cooler with air-cooled refrigerating unit and circulation pump. Equipped with an 2Kw heating, extended temperature range to 100°C. Evaporator (cooler), tank and housing of stainless steel. Pressure- suction pump of plastic. Plug and Play technology –new generation of microprocessor controlled Compatible Control.
Simple operation with a rotary knob and digital display. Alarms and warnings are given both visually and acoustically. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and can be calibrated. 3 year warrantee. Simple user guidance using a 4 line graphical LCD display and high performance processor. Administrator level and multiple user levels are available.
CC1: Level protection and adjustable overtemperature protection. Maximum and minimum setpoint for additional safety.
CC2, like **CC1:** With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with a maximum of 5 segments).
CC3; like **CC2:** Including however the complete functions: Programmer (50 segments, may be split into10 programs), interactive, contains a digital RS232 or RS485 interface as well as a 4 to 20mA analogue interface for bi-directional communication and a potential free alarm contact (Poko).

Technische Daten	**Technical Data**	**minichiller-1-H1**	**minichiller-2-H1**	**minichiller-3-H1**
Arbeitstemperaturbereich	Operating temperature range		-20...100°C	
Temperaturkonstanz bei -10°C	Temperature stability at -10°C		1 K	
Temperatureinstellung	Temperature adjustment		digital	
Temperaturanzeige	Temperature indication		digital	
Auflösung der Anzeige	Resolution		0.1 K	
Temperaturfühler	Temperature sensor		Pt100	Pt100
Anschluss externer Fühler	Connection external sensor	---	Pt100	Pt100
Analoges Interface In/Out	Analogue interface in/out	---		4-20 mA alternativ 1-5V
Sicherheitsklasse	Safety classification		FL	
Heizleistung	Heating capacity		2kW	
max. Förderleistung der Umwälzpumpe	max. capacity circulation pump	12 l/min		
Förderhöhe (Druck) max.	Pressure of pump max.		0,2 bar	
Förderleistung der Saugpumpe max.	Suction pump max.		12 l/min	
Förderhöhe (Saug) max.	Suction max.		0,1 bar	
Pumpenanschluss	Pump connection		M16x1	
Pumpenanschluss für Schlauch	Pump connections for hose		NW8/12	
max. zulässige kin. Viskosität	max. permissible kin. viscosity		50 mm^2/s	
min. Füllvolumen	min. filling capacity		1,4 l	
Ausdehnungsvolumen	Volume of expansion		2,6 l	
Abmessungen B x T x H	Overall dimension w x l x h		225 x 360 x 380 mm	
Netzanschluss	Power supply requirement		230 V 1~ 50/60 Hz	
max. Stromaufnahme	max. current (1 phase)		6,4 A	
Absicherung	Fuse		10 A	
max./min. Umgebungstemperatur	max./min. ambient temperature		32°C / 5°C	
Gewicht	Weight		23kg	
Bestell-Nr.	**Order-No.**	**466.0011H1**	**466.0012H1**	**466.0013H1**

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Circulation cooler with air-cooled refrigerating unit and circulation pump.Evaporator (cooler), tank and housing of stainless steel. Pressure-suction pump of plastic.
Plug and Play technology –new generation of microprocessor controlled Compatible Control. Simple operation with a rotary knob and digital display. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and can be calibrated. 3 year warrantee (registration necessary). Simple user guidance using a 4 line graphical LCD display and high performance processor. Administrator level and multiple user levels are available.
CC1: Maximum and minimum setpoint for additional safety.
CC2, like **CC1:** With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with a maximum of 5 segments).
CC3; like **CC2:** Including however the complete functions: Programmer (50 segments, may be split into 10 programs), interactive, contains a digital RS232 or RS485 interface as well as a 4 to 20mA analogue interface for bi-directional communication and a potential free alarm contact (Poko).

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Circulation cooler with air-cooled refrigerating unit and circulation pump.Evaporator (cooler), tank and housing of stainless steel. Pressure-suction pump of plastic.
Plug and Play technology –new generation of microprocessor controlled Compatible Control. Simple operation with a rotary knob and digital display. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and can be calibrated. 3 year warrantee (registration necessary). Simple user guidance using a 4 line graphical LCD display and high performance processor. Administrator level and multiple user levels are available.
CC1: Maximum and minimum setpoint for additional safety. **CC2,** like **CC1:** With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with a maximum of 5 segments). **CC3;** like **CC2:** Including however the complete functions: Programmer (50 segments, may be split into 10 programs), interactive, contains a digital RS232 or RS485 interface as well as a 4 to 20mA analogue interface for bi-directional communication and a potential free alarm contact (Poko).

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Chiller with air-cooled refrigerating unit and circulation pump. Evaporator (cooler) and housing made of stainless steel. Pump made of brass. Atmospheric open expansion tank made of stainless steel.
Plug & Play Technology – new generation of microprocessor controlled Compatible Control. Easy handling and digital display. Alarms and warnings are given both visually and acoustically. All controllers (CC1, CC2, CC3) are backwards compatible, interchangeable and calibrateable. 3 year warrantee (registration necessary). Simple user guidance using a 4-line graphical LCD display and a high performance processor. Administrator level and user levels are available. **CC1**: Level protection and maximum and minimum setpoint for additional safety. **CC2,** similar to **CC1**: With however an external temperature sensor connection, external temperature control and a temperature programmer (1 program with max. of 5 segments) and easy ramping function. **CC3,** similar to **CC2**: With however all functions: Programmer (50 segments, divisible into 10 programmes), interactive, contains digital RS232/RS485 interfaces as well as a (4...20mA) analogue interface for bidirectional communication and a potential free alarm contact (Poko). This allows it to be easily integrated into a laboratory automation system.

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Accessories and periphery: Adapter nom. dia 12*, sleeve nuts thread M16x1*, cover expansion tank*, stopcock, connection tubes, external sensor, connecting cable (not with CC1)

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Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP