ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions


Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Mettler Toledo AutoChem LC-RS9* at our website: *Click HERE*

Operating instructions
Bedienungsanleitung
Mode d'emploi
Instrucciones de manejo
Istruzioni d'uso

METTLER TOLEDO
LC-RS25/LC-RS9/LC-RSopen cable
Kabel LC-RS25/LC-RS9/LC-RSopen
Câble LC-RS25/LC-RS9/LC-RSopen
Cable LC-RS25/LC-RS9/LC-RSopen
Cavo LC-RS25/LC-RS9/LC-RSopen

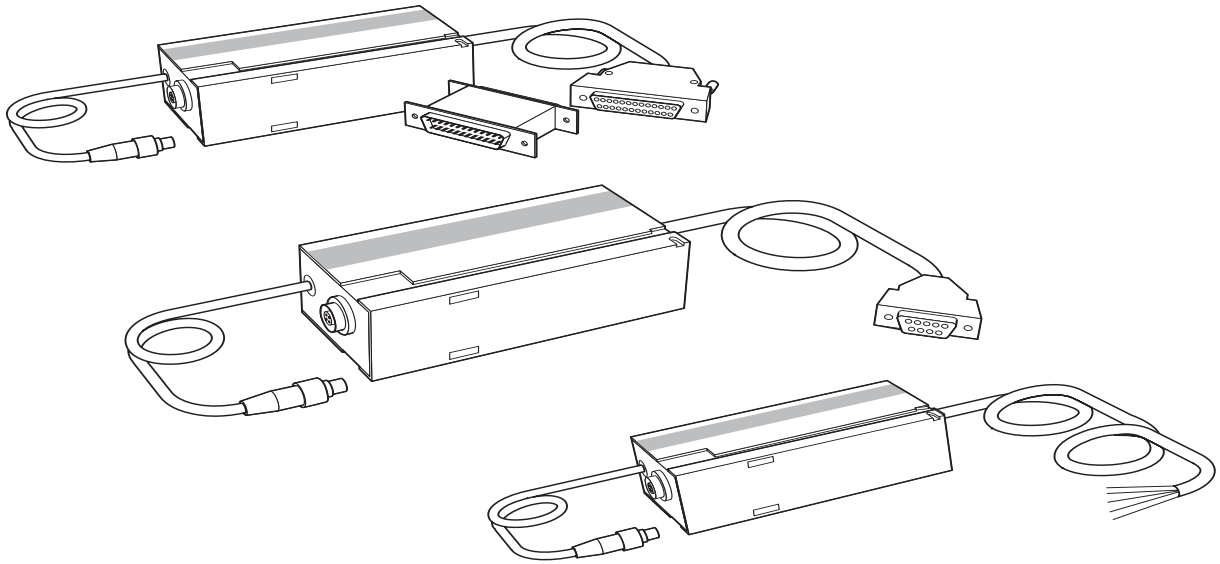
© Mettler-Toledo GmbH 1998
705157A Printed in Switzerland 9806/2.45

Mettler-Toledo GmbH, Laboratory & Weighing Technologies, CH-8606 Greifensee, Switzerland
Phone +41-1-944 22 11, Fax +41-1-944 30 60, Internet: http://www.mt.com

P705157

Matching the cable to peripherals

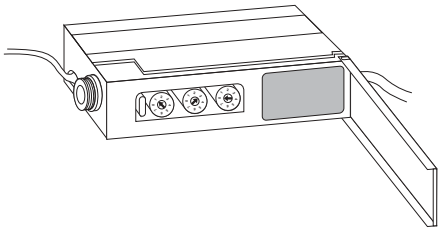The cable is configured in the factory for attachment of a computer (printed in bold).
The following protocol must be set at the peripheral device: 2400 baud, 7 bits, even parity, CTS/DTR.
If other devices need to be attached or work is performed with a different communications protocol, the cable must be appropriately configured using the 3 switches.

Switch left			**Switch middle**		**Switch right**		
Pos.	Device		Pos.	Baud rate	Pos.	Protocol	
0	**Host (PC)**	Computer, database	0	300	**0**	**7 Bit, Even, CTS/DTR**	
1	Printer (24)	Strip printer, 24 characters/line	1	600	1	7 Bit, Odd, CTS/DTR	
2	Printer (80)	Page printer A4, 80 characters/line	2	1200	2	7 Bit, Mark, CTS/DTR	
3	2nd balance	Second balance, reference balance	**3**	**2400**	3	8 Bit, NoPar, CTS/DTR	
4	ID unit	Barcode reader, terminal	4	4800	4	7 Bit, Even, Xon/Xoff	
5	Gen. I/O	I/O module	5	9600	5	7 Bit, Odd, Xon/Xoff	
6	Spec.	I/O special applications	6	19200	6	7 Bit, Mark, Xon/Xoff	
			7	unassigned	7	8 Bit, NoPar, Xon/Xoff	
7	Alt.	other devices, older METTLER TOLEDO devices; in this configuration, special settings apply to the right switch; for operation with 8 bits, only baud rate 9600 is possible.	0…6	300… 19200	0	7 Bit, Even, Syn-Ack	CL handshake
					1	7 Bit, Even, CTS/DTR	CL free mode
					2	7 Bit, Even, Pause 1 s	GA44 protocol
					3	7 Bit, Even, CTS/DTR	T, S, SI, SIR
					4	7 Bit, Even, CTS/DTR	Continuous (T, S, SI, SIR)
					5	7 Bit, Even, CTS/DTR	Option 03
				9600 (fixed)	6	8 Bit, No, CTS/DTR	Auxiliary display
				9600 (fixed)	7	8 Bit, No, Xon/Xoff	Reserved

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
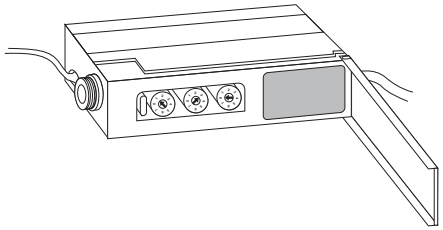
Pin assignment LC-RS25
25-pin connector, female
GENDER CHANGER 25-pin, male

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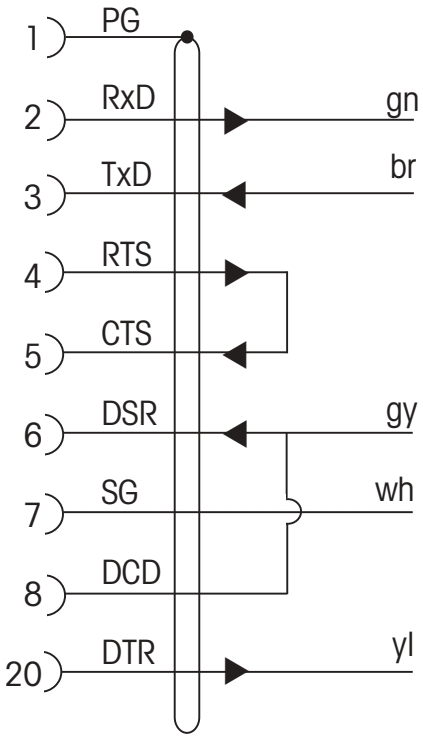
Pin assignment LC-RS9
9-pin connector, female

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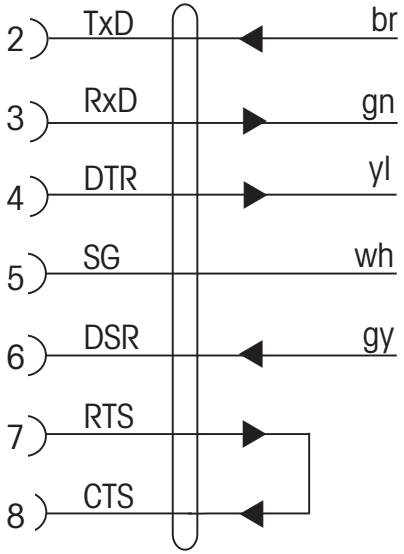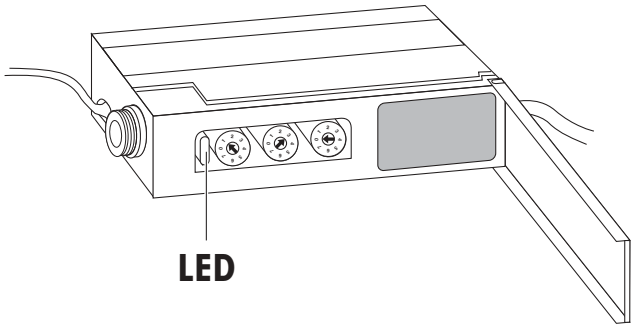
Status displays
A two-color LED is integrated in the cable and shows the status of the cable.

LED lit up green	Data are transferred
LED flashes green	Wrong interface parameters
LED lit up red briefly	Transmission error
LED flashes red	Cable faulty
LED lit up red	Self-test ended with fault, cable faulty
LED not lit up	No power supply or cable faulty

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Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP