ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *EMD Millipore XX6700P10* at our website: *Click HERE*

Technical Assistance

For more information, contact the Millipore office nearest you. In the U.S., call **1-800-MILLIPORE** (1-800-645-5476). Outside the U.S., see your Millipore laboratory catalogue for the phone number of the office nearest you. You can reach us by e-mail at tech_service@millipore.com or visit our web site (www.millipore.com).

Millipore Corporation is pleased to provide internet access to Material Safety Data Sheets (MSDS) for its products that contain hazardous materials. To obtain any MSDS documents that may be associated with this product, go to the MSDS page of our website (www.millipore.com/msds.nsf/home).

Standard Warranty

Millipore Corporation ("Millipore") warrants its products will meet their applicable published specifications when used in accordance with their applicable instructions for a period of one year from shipment of the products. **MILLIPORE MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** The warranty provided herein and the data, specifications and descriptions of Millipore products appearing in Millipore's published catalogues and product literature may not be altered except by express written agreement signed by an officer of Millipore. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and if given, should not be relied upon.

In the event of a breach of the foregoing warranty, Millipore's sole obligation shall be to repair or replace, at its option, the applicable product or part thereof, provided the customer notifies Millipore promptly of any such breach. If after exercising reasonable efforts, Millipore is unable to repair or replace the product or part, then Millipore shall refund to the customer all monies paid for such applicable product or part. **MILLIPORE SHALL NOT BE LIABLE FOR CONSEQUENTIAL, INCIDENTAL, SPECIAL OR ANY OTHER INDIRECT DAMAGES RESULTING FROM ECONOMIC LOSS OR PROPERTY DAMAGE SUSTAINED BY ANY CUSTOMER FROM THE USE OF ITS PRODUCTS.**

Millipore is a registered trademark of Millipore Corporation or an affiliated company.
Viton and Teflon are trademarks of E.I. DuPont de Nemours & Co., Inc.
P36367 Rev. -, 2/01

MILLIPORE

Dispensing Pressure Vessel Fitting Kit
User Guide

XX67 000 PK

Introduction

The Millipore pressure vessel fitting kit is designed for use with Millipore's dispensing pressure vessels:

- 5 liter, catalogue number XX67 00P 05
- 10 liter, catalogue number XX67 00P 10
- 20 liter, catalogue number XX67 00P 20
- 1 gallon, catalogue number XX67 00P 01

The fitting kit is used to control fluid delivery to the filter holder and regulate the gas supply from the pressure source. The quick disconnect fittings facilitate attachment and detachment to and from the vessel.


This fitting kit contains:

- Vent /relief, qty 2
- Pressure guage, qty 1
- Hex nipple, qty 2
- Ball valve, qty 2
- Quick disconnect nipple and coupling, qty 2
- Hose connector, qty 2
- Teflon™ tape, qty 1


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Using the Pressure Vessel Fitting Kit

1. Attach a ball valve to one end of a hex nipple and the other end to the vessel's outlet elbow fitting.
2. Attach the quick-disconnect fittings to the ball valves.

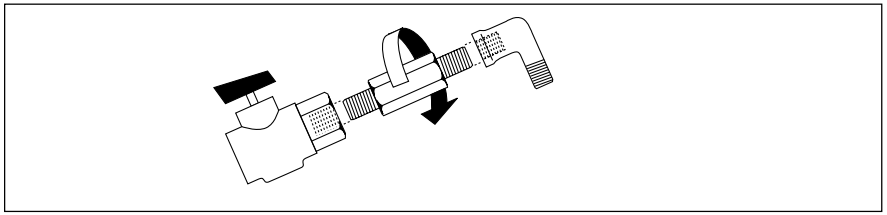


3. Attach the barbed hose connectors to each quick-disconnect fitting.

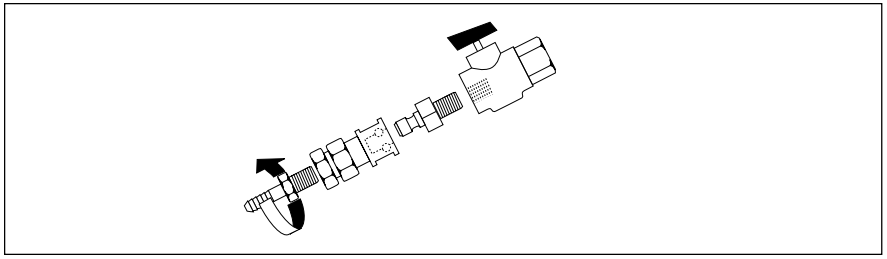

Product Ordering Information

This section lists the catalogue numbers replacement parts for the dispensing pressure vessel fitting kit. See the Technical Assistance section for information about contacting Millipore. You can also buy Millipore products on-line at www.millipore.com/purecommerce.

Product Ordering Information, continued

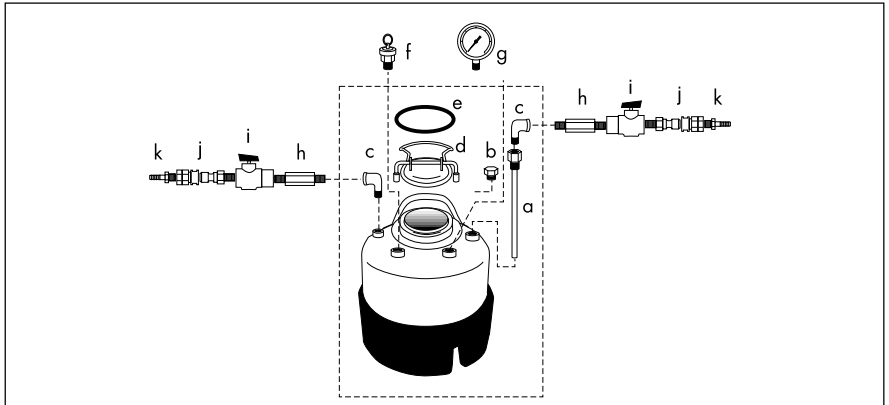


Description	Key	Catalogue Number
Replacement Parts		
Vent/relief valve, 7 bars (100 psig) ASME	f	XX67 00L 24
Pressure gauge, 0-7 bars (160 psig) ASME, 316 ss	g	XX67 00L 15
Pressure gauge, 0-7 bars (0-100 psig)	g	YY13 010 15
Hex nipple, 1/4" NPTM, 51 mm (2 in) length	h	XX67 001 25
Ball valve, 1/4" NPTF	i	YY20 293 48
Quick-disconnect nipple and coupling, 1/4" NPT	j	XX67 000 30
Hose connector, 1/4" NPTM, 316 ss, 2/pk	k	XX67 00L 11
Teflon tape		TP00 013 26
Accessories		
Outlet dip-tube	a	
5 liter		6977
10 liter		6978
20 liter		6979
Pipe plug, 1/4" NPT, 316 ss 2/pk	b	YY13 010 09
Elbow, 1/4" NPTM/Female	c	XX67 001 04
Vessel cover, cam-lock	d	6976
O-ring, Viton™-A	e	XX67 000 59


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Using the Pressure Vessel Fitting Kit

1. Attach a ball valve to one end of a hex nipple and the other end to the vessel's outlet elbow fitting.
2. Attach the quick-disconnect fittings to the ball valves.

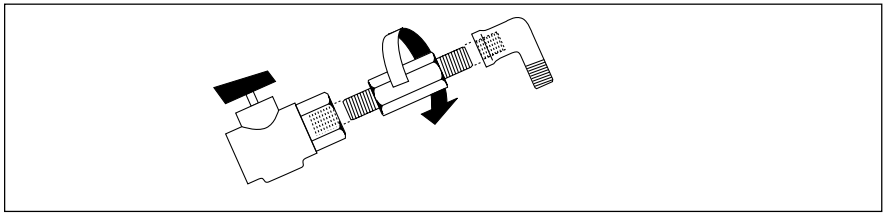


3. Attach the barbed hose connectors to each quick-disconnect fitting.

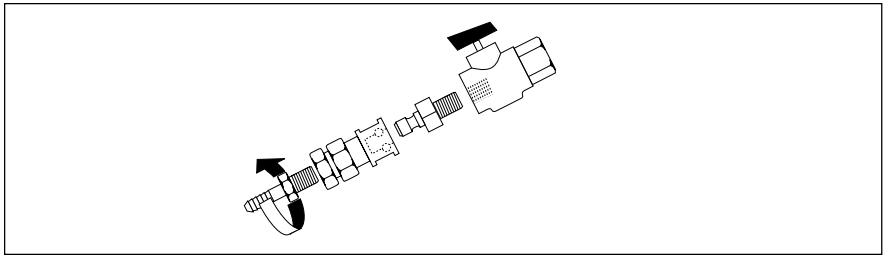

Product Ordering Information

This section lists the catalogue numbers replacement parts for the dispensing pressure vessel fitting kit. See the Technical Assistance section for information about contacting Millipore. You can also buy Millipore products on-line at www.millipore.com/purecommerce.

Product Ordering Information, continued

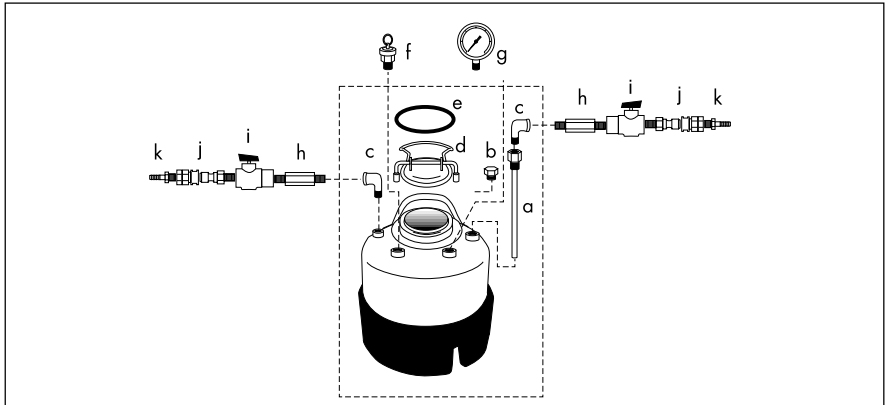

Description	Key	Catalogue Number
Replacement Parts		
Vent/relief valve, 7 bars (100 psig) ASME	f	XX67 00L 24
Pressure gauge, 0-7 bars (160 psig) ASME, 316 ss	g	XX67 00L 15
Pressure gauge, 0-7 bars (0-100 psig)	g	YY13 010 15
Hex nipple, 1/4" NPTM, 51 mm (2 in) length	h	XX67 001 25
Ball valve, 1/4" NPTF	i	YY20 293 48
Quick-disconnect nipple and coupling, 1/4" NPT	j	XX67 000 30
Hose connector, 1/4" NPTM, 316 ss, 2/pk	k	XX67 00L 11
Teflon tape		TP00 013 26
Accessories		
Outlet dip-tube	a	
5 liter		6977
10 liter		6978
20 liter		6979
Pipe plug, 1/4" NPT, 316 ss 2/pk	b	YY13 010 09
Elbow, 1/4" NPTM/Female	c	XX67 001 04
Vessel cover, cam-lock	d	6976
O-ring, Viton™-A	e	XX67 000 59

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Technical Assistance

For more information, contact the Millipore office nearest you. In the U.S., call **1-800-MILLIPORE** (1-800-645-5476). Outside the U.S., see your Millipore laboratory catalogue for the phone number of the office nearest you. You can reach us by e-mail at tech_service@millipore.com or visit our web site (www.millipore.com).

Millipore Corporation is pleased to provide internet access to Material Safety Data Sheets (MSDS) for its products that contain hazardous materials. To obtain any MSDS documents that may be associated with this product, go to the MSDS page of our website (www.millipore.com/msds.nsf/home).

Standard Warranty

Millipore Corporation ("Millipore") warrants its products will meet their applicable published specifications when used in accordance with their applicable instructions for a period of one year from shipment of the products. **MILLIPORE MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** The warranty provided herein and the data, specifications and descriptions of Millipore products appearing in Millipore's published catalogues and product literature may not be altered except by express written agreement signed by an officer of Millipore. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and if given, should not be relied upon.

In the event of a breach of the foregoing warranty, Millipore's sole obligation shall be to repair or replace, at its option, the applicable product or part thereof, provided the customer notifies Millipore promptly of any such breach. If after exercising reasonable efforts, Millipore is unable to repair or replace the product or part, then Millipore shall refund to the customer all monies paid for such applicable product or part. **MILLIPORE SHALL NOT BE LIABLE FOR CONSEQUENTIAL, INCIDENTAL, SPECIAL OR ANY OTHER INDIRECT DAMAGES RESULTING FROM ECONOMIC LOSS OR PROPERTY DAMAGE SUSTAINED BY ANY CUSTOMER FROM THE USE OF ITS PRODUCTS.**

Millipore is a registered trademark of Millipore Corporation or an affiliated company.
Viton and Teflon are trademarks of E.I. DuPont de Nemours & Co., Inc.
P36367 Rev. -, 2/01

MILLIPORE

Dispensing Pressure Vessel Fitting Kit
User Guide

XX67 000 PK

Introduction

The Millipore pressure vessel fitting kit is designed for use with Millipore's dispensing pressure vessels:

- 5 liter, catalogue number XX67 00P 05
- 10 liter, catalogue number XX67 00P 10
- 20 liter, catalogue number XX67 00P 20
- 1 gallon, catalogue number XX67 00P 01

The fitting kit is used to control fluid delivery to the filter holder and regulate the gas supply from the pressure source. The quick disconnect fittings facilitate attachment and detachment to and from the vessel.

This fitting kit contains:

- Vent /relief, qty 2
- Pressure guage, qty 1
- Hex nipple, qty 2
- Ball valve, qty 2
- Quick disconnect nipple and coupling, qty 2
- Hose connector, qty 2
- Teflon™ tape, qty 1

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP