

ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

**Artisan Technology Group**
(217) 352-9330 | sales@artisantg.com | artisantg.com

**Full-service, independent repair center**
with experienced engineers and technicians on staff.

**We buy your excess, underutilized, and idle equipment**
along with credit for buybacks and trade-ins.

**Custom engineering**
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

**Expert team | Trust guarantee | 100% satisfaction**

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Newport / New Focus 0921* at our website: *Click HERE*

205x 10MHz Adj Rcvr revA.fm Page 1 Monday, January 10, 2005 9:29 AM

USER'S GUIDE

# 10-MHz Adjustable Photoreceivers

***Models 2051 & 2053***




2584 Junction Avenue • San Jose, CA 95134-1902 • USA
phone: (408) 919–1500 • e-mail: contact@newfocus.com • www.newfocus.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

## Warranty

New Focus, Inc. guarantees its products to be free of defects for one year from the date of shipment. This is in lieu of all other guarantees, expressed or implied, and does not cover incidental or consequential loss.

Information in this document is subject to change without notice.
Copyright 2004, New Focus, Inc., a division of Bookham Technology plc. All rights reserved.

The New Focus A Division of Bookham logo and NEW FOCUS, Inc. are trademarks or registered trademarks of Bookham Technology plc in the U.S.A or other countries. Products described in this document may be covered by one or more patents in the U.S.A. and abroad.

Document Number 200329 Rev. A

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

205x 10MHz Adj Rcvr revA.fm Page 3 Monday, January 10, 2005 9:29 AM

# Contents


**Operation** 5
- Introduction . . . . . 5
- Using the Photoreceiver . . . . . 7
- Checking the Batteries . . . . . 7

**General Features & Principles** 9
- Photoreceiver Circuitry . . . . . 9
- Optical Power and Output Voltage . . . . . 11

**Frequency Response and Noise** 13
- Measuring Bandwidth . . . . . 13
- Measuring Noise . . . . . 13
- Performance Data for Frequency Response . . . . . 16
- Performance Data for Noise . . . . . 18

**Characteristics** 21
- Physical Specifications . . . . . 21
- Model 2051 Specifications . . . . . 22
- Model 2053 Specifications . . . . . 23

**Customer Service** 24
- Technical Support . . . . . 24
- Service . . . . . 24




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# Operation

## Introduction

The Model 205X is a general-purpose photoreceiver with adjustable gain and bandwidth. These receivers can be powered by batteries or by an external ±15-V power supply. There are two models available, each based on a different photodetector. Free-space (FS) and fiber-coupled (FC) versions are available for each model:

| Model | Wavelength | Diode Type | Active Area |
|---|---|---|---|
| 2051-FC | 300–1070 nm | silicon | 0.8 $mm^2$ |
| 2051-FS | 300–1070 nm | silicon | 0.8 $mm^2$ |
| 2053-FC | 900–1700 nm | InGaAs | 0.0078 $mm^2$ |
| 2053-FS | 900–1700 nm | InGaAs | 0.08 $mm^2$ |


**Note:** *Complete specifications begin on page 21.*

The 10-MHz three-stage transimpedance amplifier includes selectable gain and selectable low- and high-pass filters for easy signal optimization.



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

**Figure 1:** Typical responsivities of the Model 2051 & 2053 photodiodes

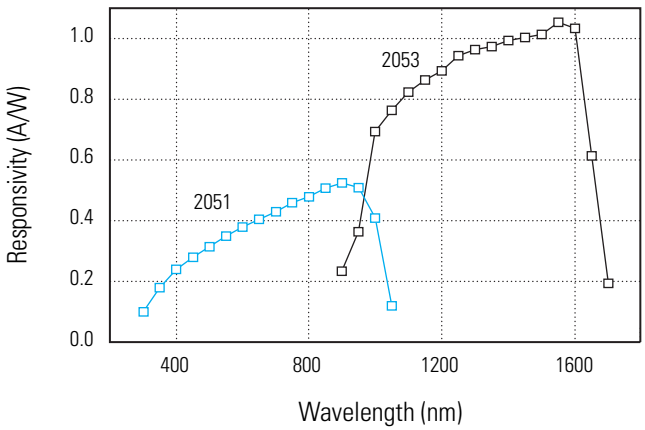


**Note:** *To obtain the value of the "response factor" in V/mW, divide the photodiode responsivity by 1.5. For more information on frequency response and noise, see page 13.*



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

## Using the Photoreceiver

1. **Mount the photoreceiver.** Use the 8-32 thread (M4 for metric versions) on the bottom of the casing to mount the photoreceiver to a post or pedestal.
2. **Supply power.** Power the Model 205X using either two 9-volt alkaline batteries or a ±15-V low-noise linear power supply (such as the New Focus Model 0901).
3. **Connect the receiver output.** Connect your voltmeter, oscilloscope, or other instrument to the Output SMA connector on the receiver.

Note: *If you wish to connect to a BNC cable, you can purchase a BNC-to-SMA adapter such as the New Focus Model 1225.*

4. **Turn on the photoreceiver power.** For external power, use ±15 VDC ON; for battery, use Batt Mode ON.
5. **Align the optical beam onto the detector.** The photodiode is not very large, so take care when aligning the beam.
6. **Adjust the gain.** Use the knob and rocker switch on the receiver to set the gain. The bandwidths vary with the gain setting (see table on page 10).
7. **Adjust the filters.** Select low-pass and high-pass corner frequencies using the knobs on the receiver.
8. **Turn off the photoreceiver power.** When you are finished with the receiver, place the power switch in the ±15 VDC ON position and switch off or unplug the external power supply.

## Checking the Batteries

The Model 205X can be powered by two standard 9-volt alkaline batteries. Under normal operating



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

conditions with low light levels and a high impedance load attached to the BNC connector, the photoreceiver draws about 20 mA from the batteries, and the battery lifetime is approximately 24 hours.

To check the condition of the battery:

1. Turn on the photoreceiver using the power switch.
2. Set the **Low Frequency** adjustment to **DC**.
3. Set the **Gain** to **3x10$^4$**.
4. Focus at least 1 µW of optical power on the detector (or place the detector in front of a desk lamp).

The output should be greater than 7 V. If it is not, replace the batteries with fresh ones.

## Replacing the Batteries

The Model 205X is shipped with two fresh 9-V batteries installed. To avoid confusion due to low batteries, replace the batteries on a monthly basis when the receiver is in frequent use, or use an external linear power supply such as the New Focus Model 0901.

1. Turn off the receiver using the power switch.
2. Use a Phillips-head screwdriver to remove the two screws on the back panel of the photo-receiver.
3. Remove the back panel.
4. Replace the used 9-V batteries with fresh ones.
5. Replace the back panel and the two screws.
6. Recheck the battery level as described above.




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# General Features & Principles

## Photoreceiver Circuitry

The circuitry inside the Model 205X consists of a photodiode followed by a three-stage transimpedance amplifier. The gain can be adjusted from 626 V/A to $18.8 \times 10^6$ V/A in 5-dB steps. The low-noise amplifier design is optimized to maximize bandwidth at each gain setting. At the higher gain settings, the bandwidth is limited by amplifier gain-bandwidth product. The plots of Figure 3 show the typical frequency responses for the different gain settings.

**Figure 2:** Functional schematic of the Model 205X circuitry

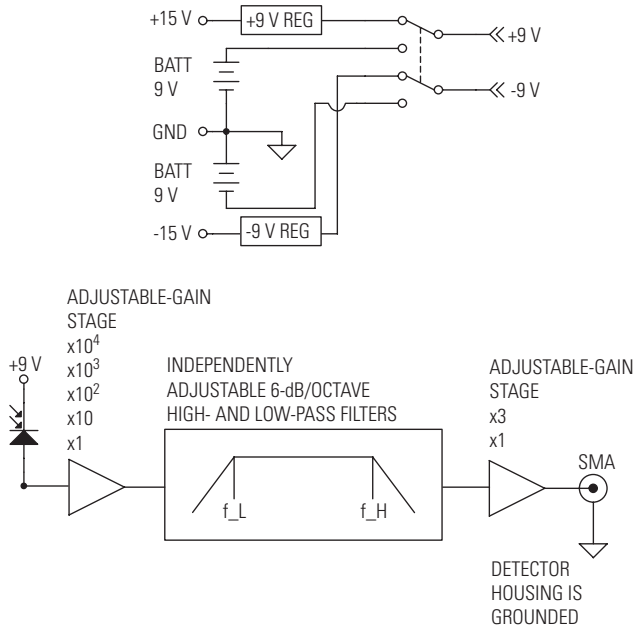




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

The following table summarizes the bandwidth at each gain setting. The bandwidth on the 3x settings is somewhat lower than the 1x settings, and significantly decreases at the highest gain settings. There is little difference in frequency response between the visible (Model 2051) and IR (Model 2053) models. The plots of Figure 3 show the frequency-response details for each gain setting.

| Gain Setting | Specification | Typical Performance |
|---|---|---|
| 1x1 | 10 MHz | 12 MHz |
| 3x1 | NA | 6 MHz |
| 1x10 | NA | 12 MHz |
| 3x10 | NA | 6 MHz |
| $1x10^2$ | NA | 8 MHz |
| $3x10^2$ | NA | 6 MHz |
| $1x10^3$ | NA | 700 kHz |
| $3x10^3$ | NA | 700 kHz |
| $1x10^4$ | NA | 250 kHz |
| $3x10^4$ | 150 kHz | 250 kHz |

## Optical Power and Output Voltage

The typical operating range for these receivers is from a few nanowatts up to 2 to 5 mW (depending on the model and gain setting). Be careful to keep the optical power below the maximum optical power of 10 mW to avoid damaging the photoreceiver.

To compute the approximate output voltage for a given input optical power use the relationship

$$V_{out} = P \cdot R \cdot G,$$



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

where $P$ is the input optical power in Watts, $R$ is the photodetector's response factor in V/mW, and $G$ is the amplifier's gain setting.

**Note:** *Estimate the value of the response factor by dividing the responsivity shown in Figure 1 by 1.5.*

For example, the Model 2051 on the $1x10^3$ gain setting and with 10 µW of optical power at 900 nm on the photodiode will have an output voltage of approximately

$$(0.01\ \text{mW})\cdot(0.35\ \text{V/mW})\cdot(1x10^3) = 3.5\ \text{V}.$$

The maximum differential optical power that can be detected by the photoreceiver is determined by the input optical power at which either stage of the transimpedance gain saturates. We can calculate the saturation power at 900 nm for the Model 2051 at its maximum output voltage of ±7 V with fresh batteries or operating from an external ±15 VDC power supply.

Using the expression $7\ \text{V} = P_{sat}\cdot R\cdot G$, the Model 2051 has a differential saturation power of 20 mW for the lowest gain setting up to 0.7 µW for the highest gain setting. At other wavelengths where the responsivity is lower, the saturation power increases inversely with response factor.



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

205x 10MHz Adj Rcvr revA.fm Page 13 Monday, January 10, 2005 9:29 AM

# Frequency Response and Noise

## Measuring Bandwidth

The frequency response and noise characteristics of the adjustable photoreceiver depend on the selected gain. The figures beginning on page 16 give the typical frequency response and noise behavior for the photoreceivers at each of the gain settings. The frequency response of the transimpedance gain is plotted using the expression

$$20 \cdot \log[\text{Gain}(f)/\text{Gain}(0)],$$

where $f$ is the frequency and Gain(0) is the gain at DC. The photoreceiver's bandwidth is defined as the frequency where the gain has decreased by 3 dB, or a factor of $\sqrt{2}$.

## Measuring Noise

The photoreceiver noise is characterized using the noise equivalent power (NEP), which is a measure of the weakest optical signal that the photoreceiver can detect. The NEP is the optical power which will produce a signal-to-noise ratio of 1 in a 1-Hz bandwidth. The minimum detectable optical power can be found using the relationship

$$\text{Minimum Optical Power} = NEP \cdot \sqrt{\text{BW}},$$

where *BW* is the bandwidth. Note that NEP is a wavelength-dependent quantity that changes with the photodetector's responsivity.



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Another way to characterize the noise is with the photocurrent noise ($I_n$), which is related to NEP by

$$I_n = R \cdot NEP,$$

where $R$ is the photodetector's responsivity (in A/W). The photocurrent noise is independent of wavelength because it gives the noise of the photoreceiver with the photodetector's responsivity factored out.

To characterize the noise of the photoreceiver, the output electrical noise spectrum is measured with a spectrum analyzer. This voltage noise spectrum is converted to an equivalent optical photocurrent noise by dividing the voltage noise by the transimpedance gain (V/A). The photocurrent noise, $I_n(f)$, has units of pA/$\sqrt{\text{Hz}}$ and is plotted in Figure 3 and Figure 4 using the expression $20 \cdot \log[\text{In}(f)/1\ \text{A}]$.

## Calculating NEP

The noise equivalent power (NEP) can be calculated by dividing the photocurrent noise by $R$, the detector's responsivity (see page 6).

From DC to 150 kHz the average photocurrent noise for the Model 2051 on the high gain setting is about 0.34 pA/$\sqrt{\text{Hz}}$, corresponding to an average NEP at 900 nm of 0.68 pW/$\sqrt{\text{Hz}}$. The integrated noise equivalent power from DC to 150 kHz is then obtained by multiplying the average NEP by $\sqrt{\text{BW}}$, the square root of the bandwidth.

The expression $BW = 2\pi f_{3\text{-dB}}/4$ for a one-pole low-pass filter is useful for calculating the equivalent noise bandwidth. Using the high-pass filter set 1 decade below the low-pass cutoff reduces noise-equivalent bandwidth by approximately 10 %. For the Model 2051 with a 3-dB bandwidth of 150 kHz, the equivalent noise bandwidth is 235 kHz. This gives an optical noise equivalent power of about 330 pW, so the minimum detectable optical signal at 900 nm (with a signal-to-noise ratio of 1) for the Model 2051 on the highest gain



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

setting is 330 pW when operating at full detector bandwidth.

You can further improve your signal-to-noise ratio by using optical modulators or choppers with lock-in amplifiers to limit the detection bandwidth. Using such techniques you can reduce equivalent bandwidth to 1 Hz or less.

## Calculating Output-Voltage Noise

The output-voltage noise can be calculated from

$$G \cdot R \cdot NEP \cdot \sqrt{BW},$$

where *G* is the gain (V/V), *R* is the photodiode response factor (V/mW), *NEP* is the average noise equivalent power, and *BW* is the bandwidth. This gives an output noise voltage for the Model 2051 on the high gain setting of

$$(3\text{x}10^4\ \text{V/V}) \cdot (0.35\ \text{V/mW}) \cdot (0.68\text{x}10^{-9}\ \text{mW}/\sqrt{\text{Hz}}) \cdot \sqrt{\frac{2\pi}{4} \cdot 150 \times 10^3\ \text{Hz}} = 3\ \text{mV}_{\text{rms}}.$$

CAUTION!

*The Johnson noise at the input of a 100-MHz bandwidth oscilloscope with 1-MΩ input impedance is 1.6 $mV_{rms}$. This is often the limiting factor in broadband measurements.*

## Summary

With the Model 2051 on the highest gain setting the minimum NEP is 0.68 pW/$\sqrt{\text{Hz}}$, and this yields an output noise voltage of 3 $mV_{rms}$. Viewed another way, for operation at the peak responsivity wavelength of 900 nm and for the high gain setting, you will achieve a signal-to-noise ratio of unity if the input power is 330 pW.

For the Model 2053 with an InGaAs photodiode, the NEP at peak response wavelength of 1500 nm is 0.34 pW/$\sqrt{\text{Hz}}$ over the 150-kHz bandwidth. The full



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


bandwidth signal-to-noise ratio of 1 is achieved around 120 pW.

Note that this assumes operation without any post-photoreceiver filtering and with the full photoreceiver bandwidth. By using the built-in electronic band-pass filter or an optical chopper and a lock-in amplifier, the receiver can detect significantly weaker optical signals.

## Performance Data for Frequency Response

The 3-dB frequency bandwidth is defined as the frequency where the photoreceiver's transimpedance gain has decreased by a factor of $\sqrt{2}$. The typical frequency responses for the Model 2051 and Model 2053 are shown in the following figures.


**Figure 3:** Typical frequency response for Model 205X at each gain setting

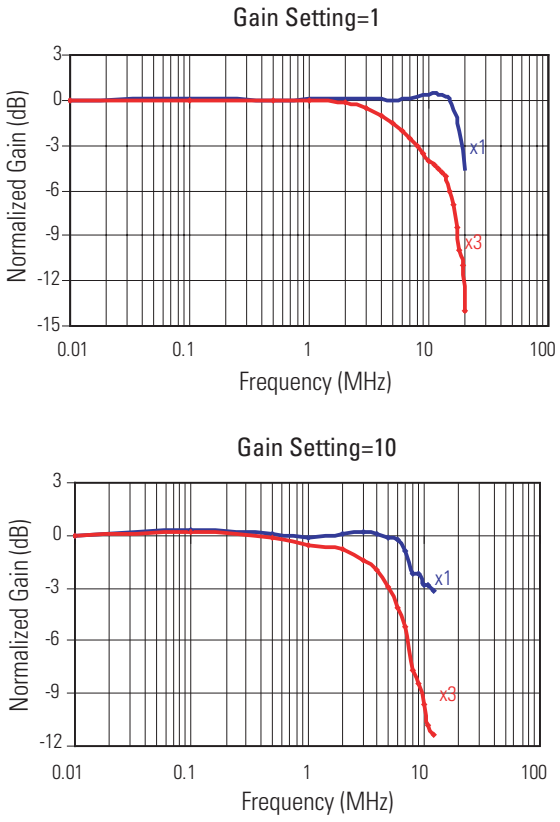





Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


Gain Setting=$10^2$

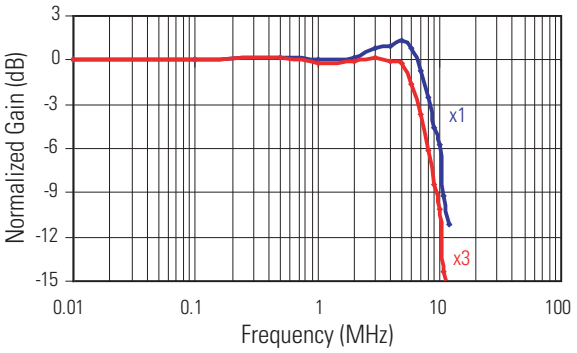


Gain Setting=$10^3$

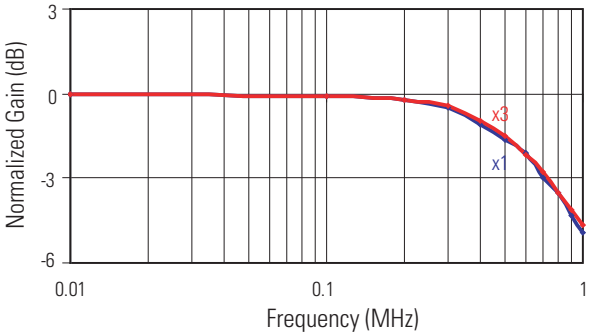


Gain Setting=$10^4$

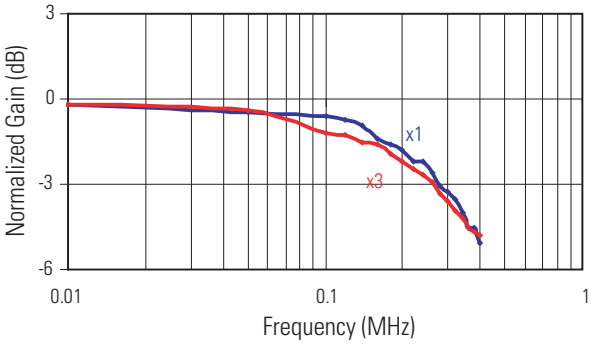




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

## Performance Data for Noise

Figure 4 shows the typical noise spectrum expressed as photocurrent noise for Model 205X photoreceivers on the highest gain setting.

To derive the receiver's Noise Equivalent Power (NEP), divide the photocurrent noise by the photodiode responsivity. To convert to output voltage noise (RMS), multiply the photocurrent noise by the gain setting from the 205X front label, then by 630 V/A (the scaling factor between the gain setting labels and the actual amplifier transimpedance gain).

For example, the output voltage noise (RMS) for Model 2053 in the $3x10^3$ setting is approximately:

$$0.34\ pA/\sqrt{Hz} \times 3 \times 10^3 \times 630\ V/A = 0.65\ \mu V_{rms}/\sqrt{Hz}.$$

For the 700 kHz of amplifier bandwidth in the $3x10^3$ gain setting, the equivalent noise bandwidth is:

$$( 2 \times \pi/4 ) \times 700 \times 10^3\ Hz = 1.1\ MHz,$$

so the predicted output noise voltage is approximately


$$0.65\ \mu V_{rms}/\sqrt{Hz} \times \sqrt{1.1 \times 10^6 Hz} = 0.7\ mV_{rms}.$$

Because the NEP is listed at the highest gain setting, some additional considerations add to the NEP at lower gain settings. First, the noise spectrum (Figure 4) is not flat, rising at frequencies above 100 kHz. This contributes an extra 20% to the output noise voltage in the $3 \times 10^3$ setting compared to $3 \times 10^4$. Also, as the output noise voltage approaches 1 $mV_{rms}$, the Johnson noise limit of your measurement instrument will become important. Note that the Johnson noise for an oscilloscope with 100-MHz bandwidth (assuming perfect roll off) and 1-MΩ input impedance is 1.2 $mV_{rms}$.



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

**Figure 4:** Typical noise spectrum for Model 205X

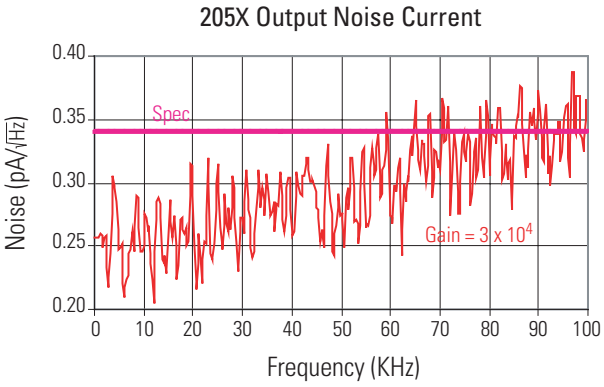




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

205x 10MHz Adj Rcvr revA.fm Page 20 Monday, January 10, 2005 9:29 AM




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# Characteristics

## Physical Specifications

**Figure 5:** Mechanical drawing of the Model 205X casing

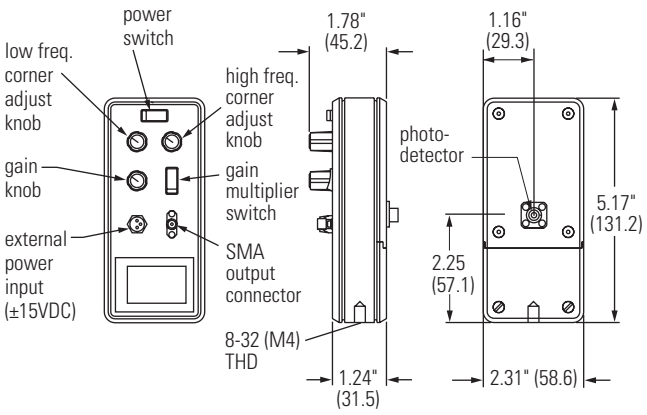




Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# Model 2051 Specifications

| | Model 2051 |
|---|---|
| Wavelength Range | 300–1070 nm |
| 3-dB Bandwidth | 10 MHz, 5 MHz, 150 kHz |
| Rise Time | 80 ns |
| Peak Conversion Gain | $9.4 \times 10^6$ V/W |
| Typical Max. Responsivity | 0.5 A/W |
| Max. Transimpedance Gain | $18.8 \times 10^6$ V/A |
| Output Impedance | 16 Ω |
| Minimum NEP | 0.68 pW/ $\sqrt{Hz}$ |
| CW Saturation Power | 20 mW @ 850 nm |
| Max. Differential Power | 20 mW @ 850 nm |
| Max. Power per Photodiode if balanced (damage threshold) | 20 mW @ 850 nm |
| Detector Material/Type | Si/PIN |
| Detector Active Area | 1.0 mm x 0.8 mm |
| Optical Input | FC or Free Space |
| Electrical Output | SMA |
| Power Requirements | ±15 VDC <150 mA External Power Supply or Two 9-V Batteries |



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# Model 2053 Specifications

| | Model 2053 |
|---|---|
| Wavelength Range | 900–1700 nm |
| 3-dB Bandwidth | 10 MHz, 5 MHz, 150 kHz |
| Rise Time | 80 ns |
| Peak Conversion Gain | $18.8 \times 10^6$ V/W |
| Typical Max. Responsivity | 1 A/W |
| Max. Transimpedance Gain | $18.8 \times 10^6$ V/A |
| Output Impedance | 16 Ω |
| Minimum NEP | 0.34 pW/ $\sqrt{Hz}$ |
| CW Saturation Power | 10 mW @ 1600 nm |
| Max. Differential Power | 10 mW @ 1600 nm |
| Max. Power per Photodiode (damage threshold) | 10 mW @ 1600 nm |
| Detector Material/Type | InGaAs/PIN |
| Detector Active Area | 0.3-mm diam. (FS)<br>0.1-mm diam. (FC) |
| Optical Input | FC or Free Space |
| Electrical Output | SMA |
| Power Requirements | ±15 VDC <150 mA<br>External Power Supply<br>or Two 9-V Batteries |



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

# Customer Service

## Technical Support

Information and advice about the operation of any New Focus product is available from our applications engineers. For quickest response, ask for "Technical Support" and know the model and serial numbers for your product.

**Hours:** 8:00–5:00 PST, Monday through Friday (excluding holidays).

**Toll Free:** 1-866-NUFOCUS (1-866-683-6287)
(from the USA & Canada only)
**Phone:** (408) 919-1500

Support is also available by fax and email:

**Fax:** (408) 980-8883
**Email:** techsupport@newfocus.com

We typically respond to faxes and email within one business day.

## Service

In the event that your photoreceiver malfunctions or becomes damaged, please contact New Focus for a return authorization number and instructions on shipping the unit back for evaluation and repair.



Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

**Artisan Technology Group** is an independent supplier of quality pre-owned equipment

## Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

## We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

## Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP