ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Optronics 60670* at our website: ***Click HERE***

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DEI-750 CE
Model S60671/S60670
PAL/NTSC
CE: Medical Device Directive, low voltage and EMC directives, UL-544

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Important information for users of OPTRONICS instruments

Thank you for your expression of confidence in the OPTRONICS brand name. Like all of our other products, this product is the result of years of experience and great care in manufacture. You and your organization have decided in favor of a modern, high quality, item of equipment from OPTRONICS.

This instruction manual is intended to serve as an aid in the proper setup, installation, and operation of the DEI-750 CE. All essential details of the equipment and all actions required on your part are clearly presented and explained. We ask that you read this manual carefully before proceeding to work with the equipment. Keep it available for ready reference in a convenient and conspicuous location near the equipment.

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. Sale restrictions apply to medical device applications only.

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Warnings and cautions

Please read this manual and follow its instructions carefully. The words **WARNING**, **CAUTION**, and **NOTE** convey special meanings. When they are used throughout this manual, they should be carefully reviewed to ensure the safe and effective operation of this product.

WARNING A WARNING indicates that the personal safety of the user may be involved. Disregarding a WARNING could result in injury to the user.

CAUTION A CAUTION indicates that particular service procedures or precautions must be followed to avoid possible damage to the product.

NOTE A NOTE indicates special information to improve the ease of maintaining the product, or to clarify important information.

The symbol of an exclamation mark within a triangle is intended to alert the user to the presence of important operating and maintenance instructions in the product's accompanying documents.

WARNING Read this instruction manual thoroughly and be familiar with its contents prior to using this equipment. Read the entire instruction manual before assembling or hooking up the camera system.

WARNING The electrical installation of the relevant location must comply with the applicable IEC, CEC, and NEC requirements.

CAUTION To reduce the risk of electrical shock, do not remove cover of equipment. Refer servicing to qualified personnel. Removal of covers by unauthorized personnel will void the equipment's warranty.

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Operating elements, displays, connections and their uses

① Power switch
② Membrane panel with:
- Color bars
- White balance
- Weighting
- Exposure
- Auto/manual
- Down/up arrows

③ Camera cable receptacle
④ RGB-Sync color video output connectors
⑤ BNC video output connectors (comp. video, B/W video)
⑥ Y/C video output connectors
⑦ Keyboard receptacle
⑧ Power fuses
⑨ Power supply cord receptacle
⑩ Potential equalization connector, per DIN 42801
⑪ Back focus adaptor
⑫ DEI-750 camera head
⑬ DEI-750T cooled camera head
⑭ Relay lens (1.0x)
⑮ Reduction lens (0.5x)
⑯ Keyboard (optional)

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Farbbalken

Symbols employed

Read the instructions carefully before operating the equipment

Power on

Power off

Equipotentiality

DANGER: Risk of explosion if used in the presence of flammable anesthetics.

CAUTION: To reduce the risk of electrical shock, do not remove cover. Refer servicing to qualified service personnel.

Protective earth (ground)

Alternating current

Replace fuse with F2.0A, 250V,
5 x 20mm fuse, IEC 127

White balance

Receptacle for camera head cable connector

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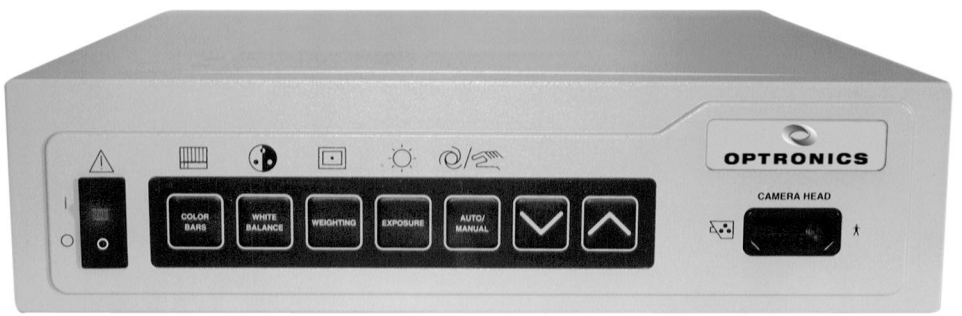
Contents

Warnings and cautions (3)
Pictures of the instrument (5)
Operating elements, displays, connectors and their uses (8)
Symbols employed (9)

Safety instructions

Intended use (12)

Installation and operating instructions

Unpacking the equipment (15)
Basic equipment (15)
Installing the DEI-750 CE (16)
Operation (16)
White balance functions (19)
Exposure control functions (21)
Color shading correction feature (24)
Freeze frame (24)
Enhancement functions (25)
Titling functions (28)
Color bars (30)
Camera reset (30)
Function exit (31)
Front panel controls (31)
Color adjustment (32)
Summary of ALT commands (33)
System configurations (34)

Description of the instrument

Cleaning and storage after use (41)
Service (43)
Warranty (43)

Technical description

Troubleshooting (44)
Video signal termination (48)
Technical data (50)

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Safety Instructions

Please read these safety instructions carefully before using the instrument.

Intended use

The DEI-750 CE is a video camera for attaching to industrial, scientific, or medical/surgical microscopes. Use of the instrument in fields other than that indicated above is not allowed for safety reasons.

The DEI-750 CE may only be used with accessories and items which have been designated by Optronics as suitable for the instrument or the safe use of which is proven.

Unauthorized conversions or modifications to the instrument are not allowed for safety reasons.

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
For an Optronics system to meet the requirements of the Low Voltage, EMC Directives, or Medical Device Directive, both the camera head and CCU must be marked with the CE mark. If any one or more components are not marked with the CE marking, then the combination of components do not meet the Low Voltage, EMC Directives, or Medical Device Directive.

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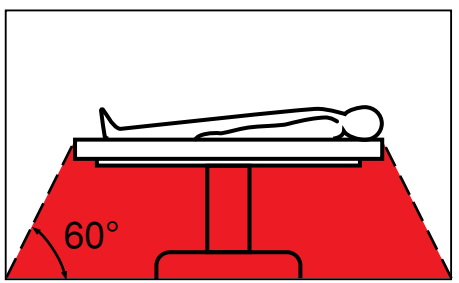Safety precautions at the site of medical installations

The unit may only be used by physicians and medical assistants who have a corresponding specialized qualification and who have been instructed in the use of the equipment.

The unit may be used in medical rooms installed according to national standards.

It is not intended for use in hazardous zones. This means, for example, that when using easily combustible and explosive inhalation anesthetics or mixtures thereof, the instruments must not be operated inside the demarcated hazard zone. Examples of such substances are: anesthetic ether (diethyl ether, cyclopropane) as well as combustible, volatile skin cleansers and skin disinfectants which may create an explosive atmosphere (e.g. detergent ether, petroleum ether).

The instrument is equipped with a connector for attaching a ground line. It should be connected according to national standards.

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Safety precautions when operating the instrument

It is the user's responsibility to make sure the instrument is safe and operates properly before using the instrument.

During treatment with the DEI-750 CE the patient must be treated and kept under observation with the usual medical care. This includes keeping a check on the progress of treatment, as well as monitoring sterile application conditions where required by the type of intervention.

This equipment has been tested and found to comply with the limits for the EMC Directive 89/336/EEC with amendments 92/31/EEC and 93/68/EEC. These limits are designed to provide reasonable protection against harmful interference in a typical medical installation. The equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to other devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference with other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving device
- Increase the separation between the equipment
- Connect the equipment into an outlet on a circuit different from that to which the other device(s) is connected
- Consult the manufacturer or field service technician for help.

WARNING: This is a Class A product. In a domestic environment this product may cause radio interference in which case the user may be required to take adequate measures.

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Unpacking the instrument

Carefully unpack the DEI-750 CE and remove the unit and its accessories from their packaging. Check for missing items and evidence of shipping damage.

File any complaints with the manufacturer or supplier immediately.

If possible, retain the original packing materials for later use; these can come in handy when the unit must be transported.

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Basic equipment

1 Camera head
1 Lens/adaptor
1 Camera control unit (CCU)
1 Set of RGB connecting cables
1 Y/C connecting cable
1 BNC connecting cable
1 Power cord
1 Instruction manual
1 Keyboard

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The equipment is equipped with a connector ⑩ for attaching a ground line.

The ground line should be installed by a qualified electrician.

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Medical installations

Note: The DEI-750 CE may be used only in medical facilities having electrical installations conforming to applicable national, state, and local electrical codes.

This instrument is not intended for use in hazardous zones. Do not operate the instrument within demarcated hazard zones while explosive anesthetic gases are in use.

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Installing the DEI-750 CE

Set the CCU on a flat surface.

Make sure there is sufficient distance on all sides from other equipment and objects.

Before plugging in the equipment, make sure that the voltage on the nameplate corresponds to the voltage of the local power line.

WARNING: Test the video camera before each procedure. Ensure that the proper video image appears on all video monitors before beginning each procedure.

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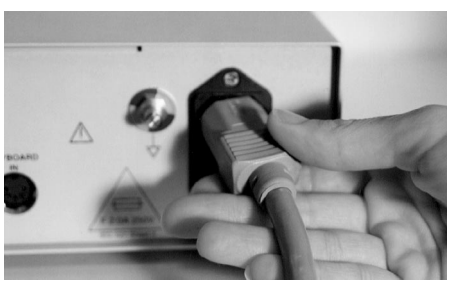
Connect power cord. Insert power cord into power cord receptacle ⑨ as far as it will go.

WARNING: For medical applications always use hospital grade power cord with this camera.

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
Plug camera head connector into receptacle ③.

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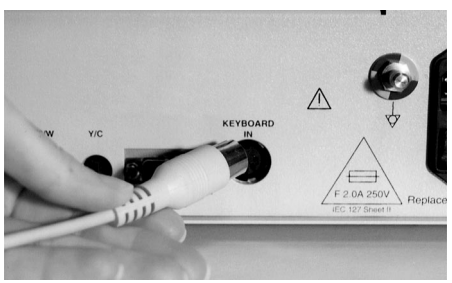Switch on the camera control unit and the video monitor.

Each time the CCU is switched on, a color test-pattern will appear on the video monitor screen. If the camera head is not connected to the CCU upon power up, the color bars will remain on until the head is plugged in.

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All camera functions are controlled from the keyboard. Only the white balance function can be controlled from the camera controller front panel. Most of the keyboard functions can be found on the label above the top row of keys. In general, pressing one of these keys will enter the camera into the corresponding mode. The current function and its corresponding setting will be displayed on the monitor as a line of text along the bottom. Many of the functions have associated parameters or selections which can be changed using the up and down arrow keys in the lower right corner of the keyboard ([↑], [↓]), or the arrow keys on the camera processor front panel. To exit a function, press the same function key again, press [Alt] [SPACE], or press another function key.

Keyboard function keys are indicated in this document within brackets, such as: [EXPOSURE].

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White balance functions

Prior to using the camera, it is important to perform a white balance operation to ensure accurate color reproduction. Each component of a camera/ microscope system can have a significant effect on the white balance. For this reason, it is recommended that the entire optical and illumination system be assembled before white balancing the camera. Be certain that no filters other than neutral density filters are in the optical path of the microscope prior to white balancing.

White balance setting [WHITE BALANCE]

The WHITE BALANCE key displays the current mode for setting white balance. The arrow keys are used to select from 5 modes described below.

AWB
This mode allows the camera to be white balanced to a white object. This is generally the most accurate method for white balancing when using a white light illumination.

Before starting a white balance, set exposure mode to “Auto” (do not use Auto+ or Auto-) and arrange a white, non-shiny object to fill most of the camera’s field of view. This can be a blank place on a slide, a piece of white paper, etc. Next press the enter button. During the cycle the message

“White Balancing...” will be displayed. If the white balance is successful, the message will change to “White Balance OK.”

If the white balance cycle fails, one of these messages will appear:

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White Bal Fail – Color
This means that the color temperature of the light is out of the 2000K to 6000K range of the camera. Try turning up the illumination and check for colored filters in the optical path.

White Bal Fail – Dark
This means that the image is too dark for an accurate white balance. Try turning up the illumination, increasing the camera exposure time, opening an iris, or removing filters in the optical path.

White Bal Fail – Bright
This means that the image is too bright for an accurate white balance. Try turning down the illumination, decreasing the camera exposure time, closing an iris, or adding neutral density filters to the optical path.

The white balance setting is stored in nonvolatile memory in the camera control unit. This means that the camera should not need to be re-white balanced when the camera is turned back on, provided nothing in the system has changed (light source, microscope, etc.).

Please note that most illuminator controls are electrical “dimmers,” and changing their settings has a significant effect on the illumination color temperature. Also, camera color accuracy is best when the illuminator color temperature is high. Therefore, it is recommended that the illuminator be set to maximum at all times. If the illuminator setting is changed, a new white balance operation is recommended.

Note: If the white balance function is activated when the camera is set for a long exposure time (1/4 second or more), a considerable amount of time may elapse before the cycle is complete. If necessary, the cycle can be aborted by turning the camera off, then on again.

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2800K
This mode sets the camera white balance for most normal incandescent lighting such as household light bulbs. This setting is also appropriate for high intensity bulbs, which are being dimmed, such as microscope illuminators.

3200K
This mode sets the camera white balance for high intensity incandescent lighting such as quartz halogen bulbs, and most microscope illuminators at maximum brightness.

5600K/fluorescence
This mode sets the camera white balance for daylight illumination. It is also the setting for use with fluorescent subjects for best color rendition.

Manual white balance
This mode allows the user to manually change the white balance of the camera by using the red gain (Alt-red) and blue gain (Alt-blue) adjustments. See “Set red, blue, gains” section on page 26.

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Exposure control functions

The camera exposure functions determine how much light is recorded by the CCD, as well as the amount of gain that is applied to the signal before it is displayed. There are two basic ways in which the camera exposure is set. In manual mode, the user manually sets the shutter speed. The gain is preset to a fixed level. In automatic mode, the camera analyzes the image and automatically finds the “best” combination of shutter speed and gain. In automatic mode, the camera is constantly monitoring the scene content, and will quickly accommodate changes in scene illumination such as when switching from darkfield to brightfield.

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Set exposure [EXPOSURE]

The EXPOSURE key will enter the camera into the set exposure mode. The monitor will display the current exposure mode (auto or manual) along with the current shutter speed. The arrow keys have two effects, depending on whether the camera is in auto or manual exposure mode. In auto mode, the arrow keys will increase or decrease the exposure by one stop, resulting in a brightening or darkening of the image. The 5 auto settings are indicated on the monitor with “Auto,” “Auto - ,” “Auto +,” “Auto - -,” and “Auto ++.” In manual mode, the arrow keys allow the user to select from any of the available shutter speeds. Press the arrow keys repeatedly until the desired shutter speed is selected.

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Select auto/manual exposure [AUTO/MANUAL]

The AUTO/MANUAL key toggles the exposure mode between auto and manual. After pressing the key, the current exposure mode and shutter speed will be displayed. The arrow keys can then be used to change the exposure.

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Set shutter speed limit [Alt] [L]

The Alt-L key combination enters the camera into the set shutter speed limit mode. This mode allows the user to change the longest allowed shutter speed when the camera is in auto exposure mode. While in this mode, the monitor will display the current limit, which can then be changed with the arrow keys.

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Set auto exposure weighting [Alt] [W]

The Alt-W key combination enters the camera into the set auto exposure weighting mode. This mode allows the user to change the method the camera uses to determine the best exposure in auto exposure mode. While in this mode, the monitor will display the current method, which can then be changed with the arrow keys.

– FULL SCREEN uses the average of the total screen content. This is best for general use.
– WINDOW uses the average of a central area of the screen. This is useful if the outer boundary is either unimportant or obscured by the optical system.
– PEAK uses the brightest feature within the entire scene. This is useful when imaging small bright features such as fluorescent structures or reflective metal surfaces.

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Set gamma [Alt] [G]

The Alt-G key combination enters the camera into the set gamma mode. This mode allows the user to change the gamma curve used in the camera processor. While in this mode, the monitor will display the current gamma, which can then be changed with the arrow keys. A gamma of 1.0 gives the camera a linear response, which is useful for quantitative analysis. A gamma of 0.7 gives the camera the proper nonlinear response to correct for the non-linearity of the monitor. The 0.7 gamma should be used for most visual and qualitative work.

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Color shading correction feature

The color shading feature corrects pink to green background shading on low contrast subjects. This allows the camera to be used with any reduction lens video adaptor without compromising color accuracy.

Color shading correction [Alt] [C]

White balance the camera. Press Alt-C on the keyboard. A color correction message with a value between -100 and +100 will appear on the monitor. To increase or decrease the value, use the arrow keys on the keyboard. Once the desired value is achieved, turn off the monitor message by pressing Alt-C. White balance the camera again.

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Freeze frame

Freeze frame can be used to capture transient events, or to minimize bleaching of fluorescent subjects. Frozen images can be subsequently modified with any of the enhancement functions: contrast, brightness, sharpness, red offset, blue offset.

Note: When using freeze frame, the subject must be stationary in order to obtain stable, flicker-free images.

Freeze frame on/off [FREEZE FRAME]

Pressing this button allows the user to “freeze” the displayed image on the monitor. Pressing the button again will restore the display to “live” video. When freeze frame is on, the letters “FF” will appear along with the current exposure when the Exposure function is turned on.

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Enhancement functions

The enhancement functions are used to provide additional clarity to existing images. They are all digital processing functions, and can be applied to either live or frozen images.

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Set contrast [CONTRAST]

This key enters the camera into the set contrast mode. This mode allows the user to change the contrast slope of the image. While in this mode, the monitor will display the current contrast setting, which can then be changed with the arrow keys. The normal contrast setting is 0. The contrast can be increased in steps up to a setting of +16 or decreased to a setting of -16. Increased contrast is useful for improving visibility of details or contours. Decreased contrast is useful for subjects with excessive contrast, or as a means to accentuate color differences rather than brightness differences. In general, the contrast function will need to be used in conjunction with the brightness function to achieve the optimal image (see below).

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Set brightness [BRIGHTNESS]

This key enters the camera into the set brightness mode. This mode allows the user to change the overall brightness of the image. While in this mode, the monitor will display the current brightness setting, which can then be changed with the arrow keys. The normal contrast setting is 0. The brightness can be increased in steps up to a setting of +32 or decreased to a setting of -32. Increased brightness is useful for improving visibility of features in dark areas. Decreased brightness is useful for suppressing unimportant features in dark areas of the image. In general, the brightness function will need to be used in conjunction with the contrast function to achieve the optimal image (see above).

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Set sharpness [SHARPNESS]

This key enters the camera into the set sharpness mode. This mode allows the user to change the overall detail enhancement of the image. While in this mode, the monitor will display the current sharpness setting, which can then be changed with the arrow keys. The normal sharpness setting is 0. The sharpness can be increased in steps up to a setting of +20. Increased sharpness is useful for improving visibility of fine details and contours in the image.

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Set red, blue offsets [RED], [BLUE]

These keys enter the camera into the set red or blue offset modes. These modes allow the user to change the overall color bias of the image. While in these modes, the monitor will display the current setting, which can then be changed with the arrow keys. The normal setting is 0. The selected color can be increased in steps up to a setting of +32 or decreased to a setting of -32. Positive red settings will make the image appear more red, while negative red settings will make the image appear more green. Positive blue settings will make the image appear more blue, while negative blue settings will make the image appear more yellow.

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Set red, blue gains [Alt] [RED], [Alt] [BLUE]

These keys enter the camera into the set red or blue gain modes. These modes allow the user to change the white balance of the camera manually. This is useful in situations where the camera cannot be automatically white balanced or where the user would like to modify the color rendition of the camera to improve the appearance of the image. (See “White balance setting” section on page 20 for information about white balance). While in these modes the monitor will display the current setting, which can be changed with the arrow keys. The values range from 0 (minimum) to 255 (maximum).

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Set color enhancement [Alt] [WHITE BALANCE]

This key combination enters the camera into the set color enhancement mode. When color enhancement is on, colored objects will appear brighter and more vivid (saturated) while unsaturated objects (whites and grays) will appear darker. This is useful when the subject matter uses color as a means of differentiating features, such as stained cells and multi-dyed fluorescence subjects. Use the arrow keys to turn color enhancement on and off.

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Enhancement reset [RESET]

This key sets the contrast, brightness, sharpness, red offset, and blue offset to their default settings (0). Pressing the key a second time will restore the settings to their prior values. This function enables the user to quickly toggle between enhanced and unenhanced operation in order to better evaluate the effects of the enhancements.

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Titling functions

The camera processor has a built-in titling generator, which allows the user to annotate images with alphanumeric text. When in the title edit mode, the keyboard is used as a typewriter to enter desired text. There are five “pages” of text which can be written to and displayed. Each page can hold 10 lines of text, 24 characters per line.

Title edit/title display/title off [TITLE]

Pressing the TITLE key causes the camera to enter title edit mode, which allows the user to enter text into the current title page. When in title edit mode, the message “Title Edit” will appear at the bottom of the monitor along with the current page number. A flashing cursor will appear in the text area, indicating where the text will be placed. To enter text, simply type the desired message from the keyboard.

Pressing the TITLE key again switches the camera into title display mode, which only displays the current text and does not allow the text to be changed.

Pressing the TITLE key a third time turns the title display off.

Cursor movement [↓], [↑], [←], [→]

In title edit mode, the four arrow keys in the lower right corner of the keyboard can be used to position the cursor anywhere in the text area.

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Page up/page down [Fn] [PgUp], [Fn] [PgDn]

The Fn-PgUp and Fn-PgDn key combinations increment or decrement the current page number when in title edit or title display mode. There are a total of five pages.

Note: The Fn key is in the lower left corner of the keyboard, and the PgUp and PgDn keys are also the up and down arrow keys in the lower right corner of the keyboard.

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Replacement mode/insert mode [INSERT]

The Insert key toggles the text editor between replacement mode and insert mode when the camera is in title edit mode. In replacement mode, the cursor appears as a flashing white square, and typing a character simply replaces the character under the cursor with the new one. In insert mode, the cursor appears as a flashing underscore, and typing a character inserts the new character immediately to the left of the cursor and shifts the characters to the right of the insertion point.

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Delete current character [DELETE]

The Delete key causes the character underneath the cursor to be deleted when in title edit mode.

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Delete previous character [BACKSPACE]

The Backspace key causes the character to the left of the cursor to be deleted when in title edit mode.

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Clear current page [Alt] [DELETE]

The Alt-Delete key combination causes all text in the current page to be deleted when in title edit mode.

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Clear all pages [Ctrl] [DELETE]

The Ctrl-Delete key combination causes all text in all pages to be deleted when in title edit mode.

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Color bars

The camera control unit has a built in color bar generator. This is useful for adjusting the monitor (see “Adjusting Color Monitor” section in this manual), providing a leader for videotape recording, or for troubleshooting a video system. The pattern used is the SMPTE color bar pattern which includes bars of each primary and secondary color and a pluge signal, which is used to set the brightness level of the monitor.

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Color bars on/off [COLOR BARS]

This key toggles the color bar test pattern on and off. When color bars are on, all camera functions are disabled except for titling.

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Camera reset

The camera reset function sets all camera parameters to default states. This is a convenient way to put the camera into a known state if the settings are unknown.

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Standard reset [Alt] [RESET]

The Alt-RESET key combination initiates a standard reset. All enhancements are turned off, the camera is placed in auto exposure mode, the weighting is set to full screen, the white balance is set to 2800K, the color bars are turned off, freeze frame is turned off, and titling is turned off.

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Function exit

A currently active function will be indicated along the bottom line of the monitor (“EXPOSURE,” “AUTO,” etc.). A function can be exited by pressing its function key again, pressing another function key, or by using the general exit function.

General exit [Alt] [SPACE]

The Alt-SPACE key combination allows the user to quickly exit any function, and simultaneously clear the text at the bottom of the screen ([SPACE] is the spacebar at the bottom of the keyboard).

Default exposure mode

If no text appears on the bottom of the monitor image, the camera defaults to exposure set mode. Pressing the up or down keys on the keyboard will change the exposure accordingly.

Front panel controls

The camera processor front panel has buttons for Color Bars, White Balance, Weighting, Exposure, Auto/Manual, Up Arrow, and Down Arrow which are functionally identical to their counterparts on the keyboard. These controls can be used if the keyboard is malfunctioning or missing. If the Exposure or Weighting functions are enabled, use the up and down arrow buttons on the front panel to change the settings. Press the up and down arrow buttons simultaneously to initiate a standard reset of the camera.

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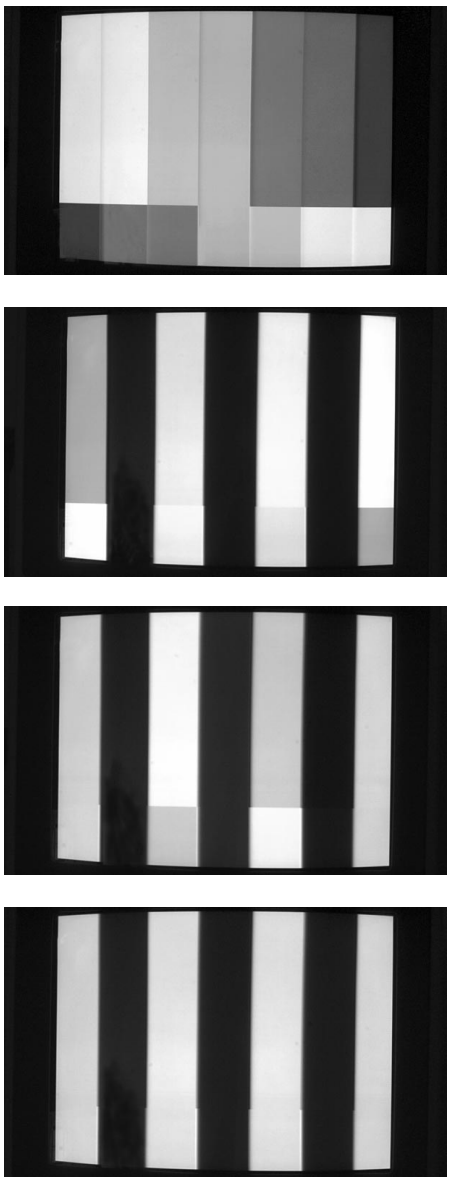
Color adjustment using SMPTE-type color bars

NOTE: The adjustments below are for composite or Y/C video only. Most monitors in RGB mode have preset color by the manufacturer, and color intensity or hue cannot be manually adjusted. Only the BRIGHTNESS and CONTRAST adjustment can be performed in RGB mode.

1) Activate the SMPTE-type color bar generator by pressing the COLOR BARS key.
2) *Contrast adjustment*
 Set contrast for brightest possible image without blooming or distortion.
3) *Brightness adjustment*
 Set the brightness control on the monitor until the right hand “pluge” bar is just visible, as shown.
4) For the next two adjustments (chroma and phase), press the *blue only* button on the monitor, if so equipped. Otherwise, view the monitor through a blue Wratten 47B filter.
5) *Chroma (color) adjustment*
 Adjust the chroma (or color) control on the monitor until the indicated bars match in intensity (*blue only* on, or viewed through blue 47B filter).
6) *Phase (hue, tint) adjustment*
 Adjust the phase (or hue, tint) control on the monitor until the indicated bars match in intensity (*blue only* on, or viewed through blue 47B filter).
7) When the correct results have been achieved, press the *blue only* button again (or set aside the blue filter). The colors should now be accurate.
8) Press the COLOR BARS key again to return to normal camera operation.

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Summary of ALT commands

[Alt] [W]
Set auto exposure weighting

[Alt] [L]
Set auto exposure shutter speed limit

[Alt] [G]
Set gamma

[Alt] [WHITE BALANCE]
Set color enhancement on/off

[Alt] [RED]
Set red gain

[Alt] [BLUE]
Set blue gain

[Alt] [DELETE]
Clear all text in current title page

[Ctrl] [DELETE]
Clear all text in all title pages

[Alt] [RESET]
Initiate standard reset of camera

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
Videosystem
Tastatur
Y/C Recorder
Composite-Videomonitor

System configuration 1

Video system
Keyboard
Y/C recorder
Composite video monitor

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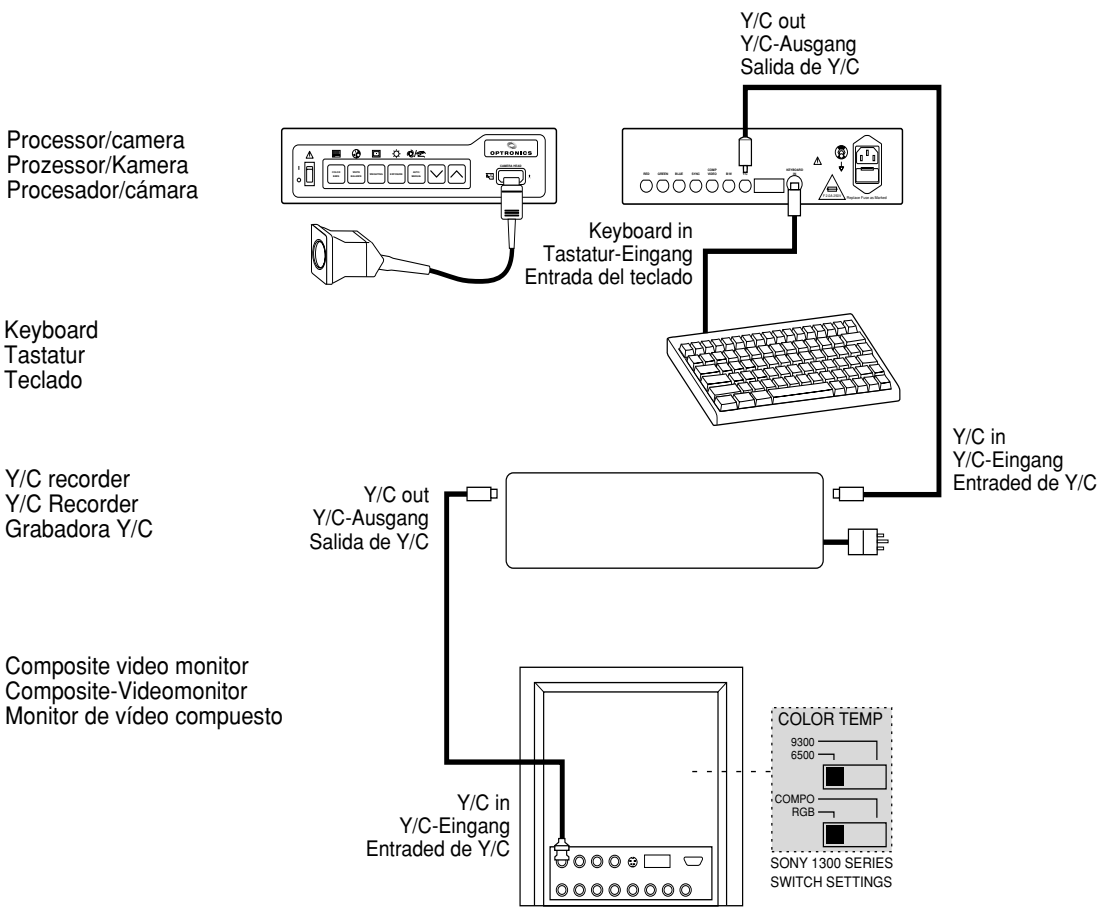
Videosystem
Tastatur
RGB Videomonitor

System configuration 2

Video system
Keyboard
RGB video monitor

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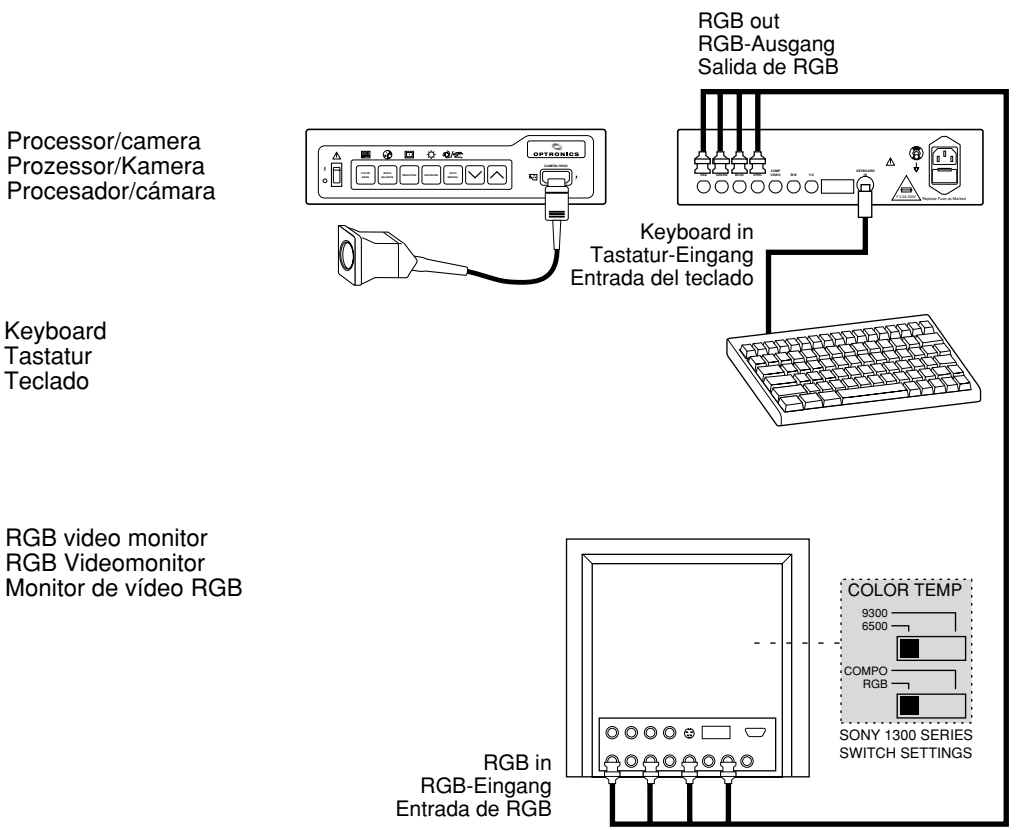
Videosystem
Tastatur
Videorecorder
Composite-Videomonitor

System configuration 3

Video system
Keyboard
Video cassette recorder
Composite video monitor

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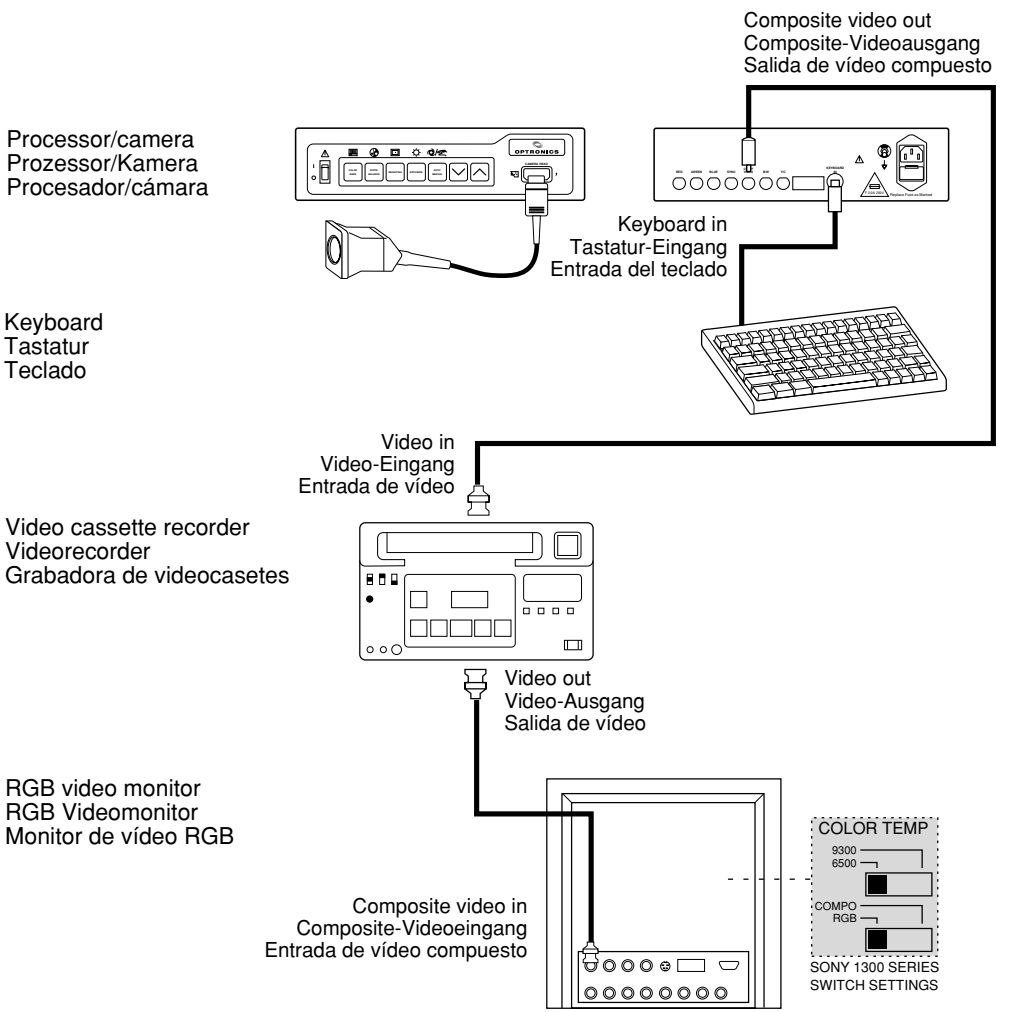
System configuration 4

Video system
Keyboard
Video cassette recorder
Video printer or disk recorder
Composite video monitor

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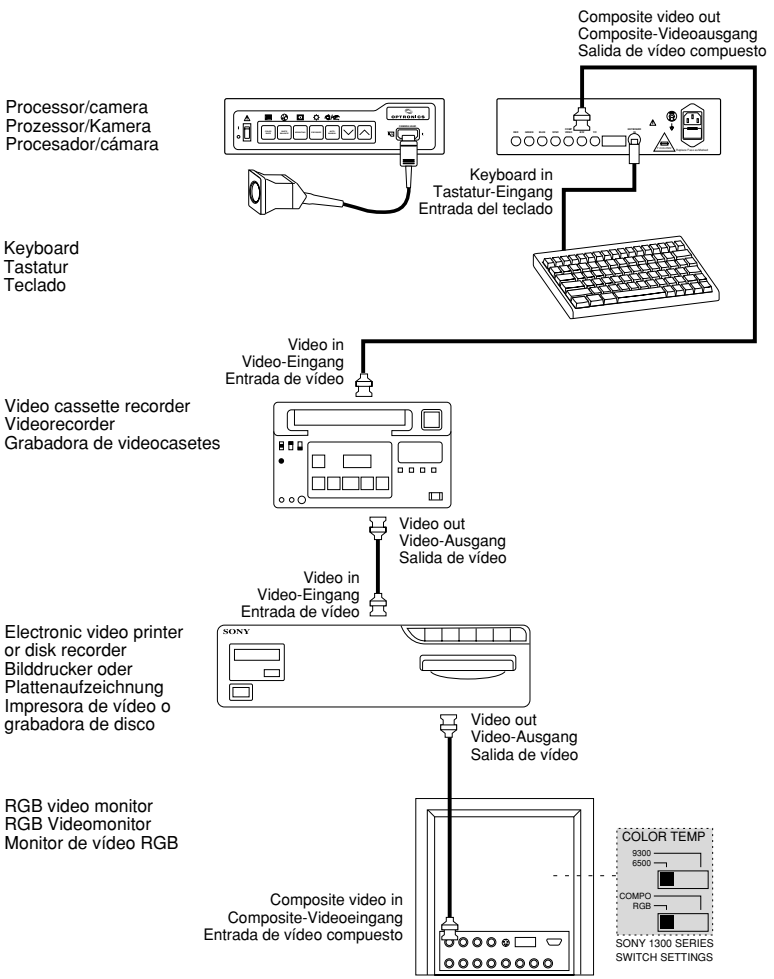
System configuration 5

Video system
Keyboard
Composite video monitor
Additional video monitor

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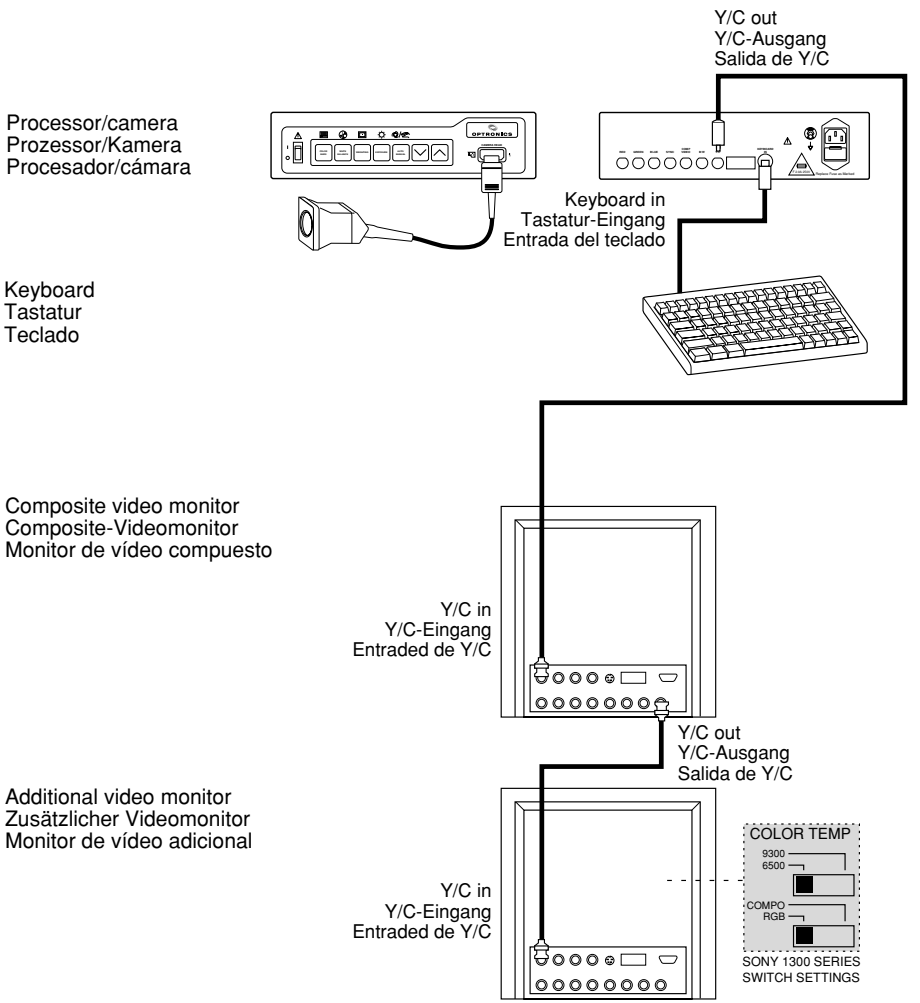
Videosystem
Tastatur
Sony UP5000 Bilddrucker
RGB Videomonitor

System configuration 6


Video system
Keyboard
Sony UP5000 video printer
RGB video monitor

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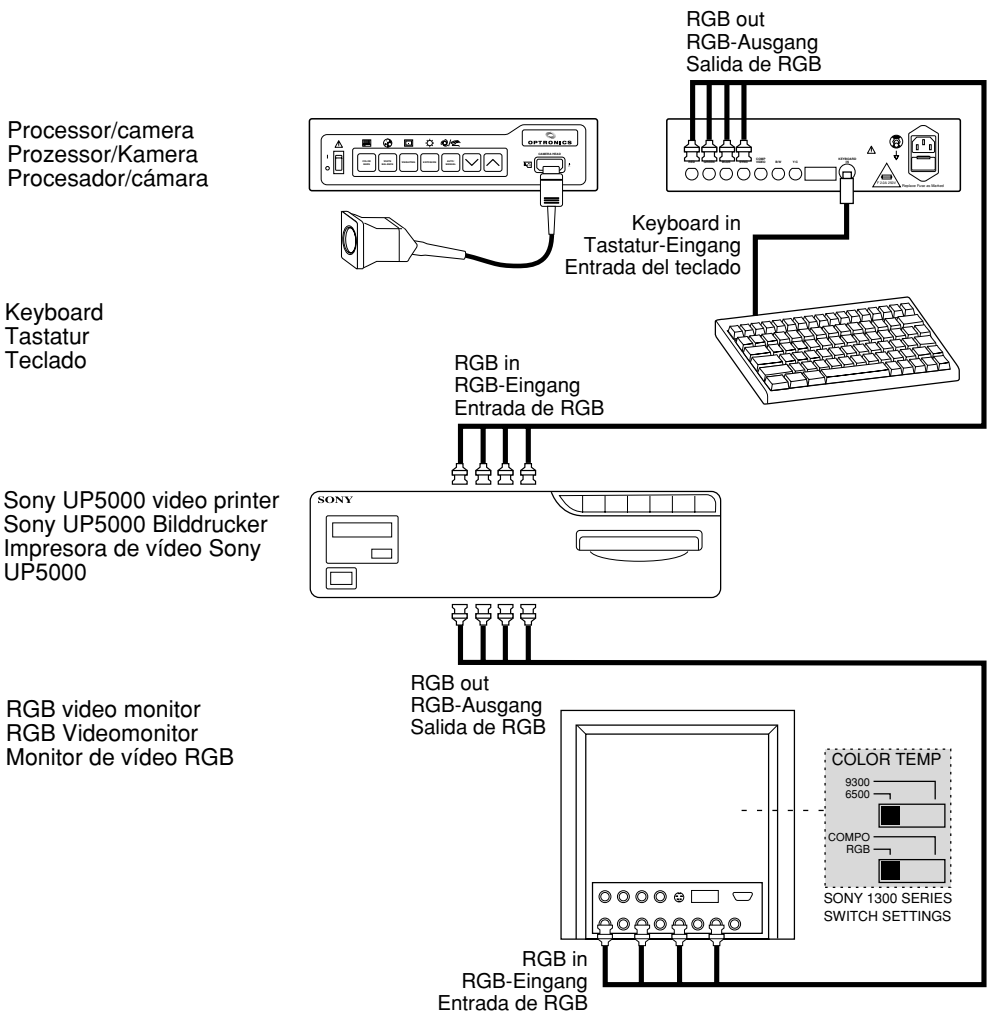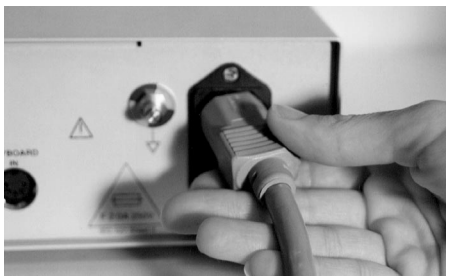Fuse replacement


Switch off the instrument and remove the power plug from the electrical outlet.

CAUTION: Use only fuses approved to IEC 127 Sheet II, Type F for replacement.

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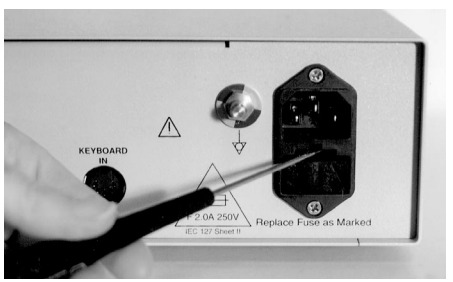Remove the power fuse holder ⑧ with a screwdriver or other tool.

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Insert new fuses of the appropriate rating.

Caution: Only use fuses of the correct rating.

	100–240 VAC
Power fuses	2 x F2.0A

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Replace and tighten the caps of the fuse holder ⑧.
Reconnect the power cord.
Test the instrument for proper operation.

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Cleaning and storage after use

Camera head

1. Clean the camera head with a damp cloth.
2. Dry thoroughly with a soft towel.
3. Replace the camera head into its box. Do not coil the camera cable tightly.

Camera processor

1. Clean the processor cover with a damp cloth.
2. Dry thoroughly with a soft towel.
3. Replace the camera processor into its box.

CAUTION: Do not store the camera in direct sunlight or excessive heat.

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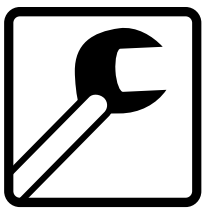
Maintenance operations

Performance of preventive maintenance is not essential. Regular maintenance can, however, contribute to identifying potential problems before they become serious, thus enhancing the instrument's reliability and extending its useful operating life. Maintenance services can be obtained from your local representative or from the manufacturer.

Information on maintenance can be obtained from the manufacturer.

WARNING: To ensure safe operation of the video camera, you should perform the following procedures every 12 months:

1. Check leakage current from chassis < 100µA.
2. Check ground impedance < 0.1 ohms.
3. Check power consumption ≤ rated power.

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Service

Refer all servicing and adjustments to qualified service personnel. DO NOT open camera head or processor chassis as this will result in damage and VOIDING OF WARRANTY.

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Warranty

The cameras are covered by a one year warranty applicable from the date of sale against defects in materials and workmanship or malfunctions under normal use and service. This warranty does not apply to any unit which has been subject to misuse, abuse, neglect, improper installation or operation or that has been altered, adjusted, or tampered with by any person other than authorized service personnel.

Should a malfunction occur, the unit must be returned to **Optronics** or to its authorized representative. Before returning the unit, the customer must first call **Optronics** to obtain a "Return Materials Authorization" number. **Optronics** will repair the unit and return it to the customer.

If, upon examination by authorized service personnel, it is determined that the malfunction is due to misuse or abuse, warranty provisions will not apply. An estimate of the cost of repair work will be given to the customer prior to servicing and repairing the unit.

Optronics is responsible for the effects on safety, reliability and performance of the equipment only if it is serviced by authorized personnel and used in accordance with the Operating Manual Instructions.

If the equipment is used in a manner not specified by **Optronics**, the protection provided by the equipment may be impaired.

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Troubleshooting

Symptom:
No video – power light off.

Possible causes:
Power cord is not properly connected to power cord connector.
Camera processor not plugged in.
Defective power fuse.

Remedy:
Push power cord firmly into power cord connector.
Check that there is electricity to the wall outlet.
Check that processor switch is turned on.
Exchange fuses.

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Symptom:
No video – power light on.

Possible causes/remedy:
Check that monitor is plugged in and turned on.
Plug in monitor and turn on.
Check that video (comp, Y/C, RGB) cable from camera processor to monitor is connected.
Check that camera head is plugged in.
Check that monitor is appropriately switched (correct video in channel or TV/line switch turned to “line”).
Check that all accessories (VCR, video printer, etc.) are turned on and switched to “line?”
Check that VCR record button is on. Some VCRs require the record button to be depressed to operate.
If using a TV receiver, check that the video cable from processor is connected to “video” inputs (not antenna inputs).
Trace video cables through accessories, i.e., video out to video in, video out to video in, etc.

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Symptom:
Intermittent video.

Possible causes/remedy:
Check all connections: camera head cable to processor and video cables to and from all accessory components.
Wiggle camera head cable at processor and camera head – if picture flashes, refer to Service and Warranty for information on how to return a malfunctioning unit.
Check video cables to determine if bad.
See if problem exists with another video source (for example, remove camera from system and play back tape on VCR and monitor).
Try to isolate problem by finding one component that works and testing all other components with it.

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Symptom:
Bad color.

Possible causes/remedy:
Repeat the white balance procedure with correct illumination, with complete system hooked up, pointing the camera at a white object.
Adjust monitor, starting with controls in detent or center position.
Terminate monitor if it is the last component in the system. If it is not the last component, leave the monitor unterminated and terminate the last component in line (refer to section Video Signal Termination).
Check that processor is in a cool, well-ventilated area.
If no color, turn up color on monitor.
If too much color (or color too intense), turn down color on monitor.
The Sony UP5000 video printer can introduce color shifts to the system.
See System Configuration 6 to isolate the camera from printer-related problems.
Check rear panel switch settings on SONY 1300 series monitors. Switches should be set to RGB and 6500.

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Symptom:
Picture too dark or too bright.

Possible causes/Remedy:
Adjust level of light source.
Adjust monitor brightness.
Adjust contrast, picture, or sharpness controls on monitor, starting with detent or center position.
If desired, the automatic exposure system may be overridden by pressing the manual shutter key, and shuttering through the various speeds until desired exposure is achieved.
The Sony UP5000 video printer can introduce interference to the system. See System Configuration 6 to isolate the camera from printer-related problems.

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Symptom:
General.

Possible causes/Remedy:
Try to isolate problem equipment (monitor, camera, VCR, etc.) before returning camera for abnormal image (noise, interference, bad color, etc.). If the camera does seem to be the problem, refer to Service and Warranty for information on how to return a malfunctioning unit.
Video tape the problem and include the tape when returning the equipment for repair. Do not attempt to repair equipment. There are no user-serviceable parts inside the camera.

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Video signal termination

In a video system, a signal is generated and transmitted through various electronic devices. This signal is always carried on a transmission wire (video coaxial cable) having a characteristic impedance of 75 ohms.

The end point of the transmitted video signal must be terminated, otherwise signal reflections may occur, causing incorrect brightness and/or ghosting effects to appear on the monitor image. Termination is accomplished by applying a 75 ohm resistance to the signal at its end point.

The monitor receives and *passively* generates a video signal. It is terminated only when the monitor is the last device hooked up consecutively in a video system receiving the signal. Many devices such as VCRs and video printers receive as well as *actively* regenerate video signals.

Termination is accomplished automatically through the device's internal circuitry or manually by the operator flipping a rear switch on the device to the *Term* or *75 ohm* position.

To determine if a device in your video system requires manual termination, check the rear of the unit. If a termination switch exists, manual termination is indicated and a decision whether or not to terminate must be made. If a switch is not found, the device will automatically apply termination as needed.

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Monitors

Some monitors are automatically terminating, others are not. Check the rear of the monitor for a termination switch to verify. If a switch exists and the monitor is the last device hooked up to the system, set the switch to terminate. If the monitor is not last, the signal is being passively looped through and termination is not required. Set the switch to the unterminated position. Note: on automatically terminating monitors, the termination switch is located in the video input receptacles and activated when a coaxial cable plug is inserted. Alternatively, the video output receptacles deactivate termination when a plug is inserted. If termination is required, ensure that the connecting coaxial cable is plugged into the video in receptacle and not the video out receptacle.

Video cassette recorders/printers

Some video cassette recorders/printers are automatically terminating, others are not. Check the rear of the device for a termination switch to verify. If one exists, set the switch to terminate regardless of whether or not the VCR or printer is the last device hooked up to the system. The only exception is if T-connectors are used to loop the signal through to a successive device. In this case, the signal is being passively generated or looped through and termination is not required. Set the switch to the unterminated position.

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Technical data

DEI-750 CE
Image sensor: 3 x 1/2" CCD-chip, one per color channel
Picture format: 6.4 mm x 4.8 mm
Picture elements (PAL):
752 x 582 pixels per chip, 1,312,992 in total
Picture elements (NTSC):
768 x 494 pixels per chip, 1,138,176 in total
Auto exposure: Microprocessor controlled, Proprietary Continuous Ratchet Algorithm
Auto white balance: 2280° K — 4420° K
Internal lens: f 1.4, primary color (RGB) separation prism
Instrument coupling: "C"-mount
Video output (PAL): PAL 4.43 signal to BNC socket, Y/C signal to S-VHS socket, RGB signal to BNC sockets
Video output (NTSC): NTSC 3.58 signal to BNC socket, Y/C signal to S-VHS socket, RGB signal to BNC sockets
Input: 5-pole DIN connector to keyboard socket
Fuse: 2 x F2.0 A, 250V
Power requirements: 100 – 240V~, 50/60 Hz, 0.5A, 100 VA
Camera head: 37.5 x 47.5 x 60.0 mm (w x h x d)
Weight: 0.37 kg (13.1 oz)
Cooled camera head: 114 x 173 x 106 mm (w x h x d)
Weight: 2.5 kg (5.5 lb) with cable
Camera control unit (CCU):
305 x 89 x 322 mm (w x h x d)
Weight: 3.8 kg (8.4 lb)
Storage conditions:
Ambient temperature: -10° to +60° C (14° to +140° F)
Humidity: 0 to 80% (non-condensing)
Atmospheric pressure: 50 kPa to 106 kPa
Operating conditions:
Indoor use maximum altitude: 2,000 m
Temperature range: +5° to +35° C (41° to 95° F)
Humidity: 0 to 80% (non-condensing)
Mains supply voltage fluctuation ± 10%
Installation overvoltage category II
Pollution degree: 2

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Technical description

General classification

Type of protection against moisture:
Ordinary

Type of protection against electric shocks:
Protection Class 1

Mode of operation:
Continuous

Laboratory classification

Mains supply voltage fluctuation: ± 10%
Installation overvoltage category II
Pollution degree: 2

Equipment test certificates

Instruments are manufactured in accordance with UL-544, EN60601-1 (IEC 601-1).

Equipment test reports

Instruments are designed in accordance with EN61010-1 (not agency certified).

Technical documentation

On request, the manufacturer will provide those circuit diagrams, itemized parts listings, descriptions, sets of adjustment instructions, and other items of available documentation to suitably qualified user personnel duly authorized by the manufacturer for their use in repairing those components of the instrument that have been designated by their respective manufacturers as repairable.

Supply of such technical documentation relating to the instrument shall not be construed as constituting manufacturer's authorization of user's personnel, regardless of their levels of technical training, to open or repair the instrument. Explicitly exempted herefrom are those maintenance and repair operations described in this manual.

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Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

OPTRONICS
Picture This™

DECLARATION OF CONFORMITY

We,

Optronics
175 Cremona Drive, Goleta, California 93117 - USA

declare under our sole responsibility that the product(s) to which this declaration relates is in conformity with the following EU Directives, standard(s) or other normative document(s):

Standards:	**Directives:**
EN 60601-1	MD (93/42/EEC)
EN 60601-1-2	

Performed on units bearing CE marking.

Goleta, CA. USA

30 APRIL 2001
(Place and date of issue)

T.P. Fernandez
Director, Regulatory and Standards Compliance

(Name and signature or equivalent marking of authorized person)

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

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The authorized EU Regulatory Representative for Optronics is:
KARL STORZ GmbH & Co.
Mittelstrasse 8, D-78532 Tuttlingen
Postfach 230, D-78503 Tuttlingen, Germany
Regulatory Department
Fax: +49/7461/708 304

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Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP