ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Fisher Scientific / Hamilton 950P201GY* at our website: *Click HERE*

FORENSICS CABINET
INSTALLATION, OPERATION, MAINTENANCE

INSTALLATION

TWO CABINETS ARE AVAILABLE, FLOOR AND BENCH-TOP MODELS. ON FLOOR UNITS, EFFORTS SHOULD BE MADE TO LEVEL AND PLUMB THE UNIT PRIOR TO MAKING EXHAUST CONNECTION.

1) BOTH UNITS REQUIRE CONNECTION TO A 120 VOLT, 20 AMP CIRCUIT AND TO AN EXTERNAL EXHAUST SYSTEM. SEE ROUGH-IN DETAILS AND WIRING DIAGRAM ON FOLLOWING PAGES.

2) CONNECT EXHAUST STACK OF CABINET TO AN EQUAL OR LARGER DIAMETER EXTERNALLY VENTED SYSTEM. EXHAUST SYSTEM MUST MEET CFM VALUES SPECIFIED ON FOLLOWING PAGES.

OPERATION

FOLLOW OPERATION LABEL AFFIXED TO FRONT PANEL OF UNIT. SEE COPY AT RIGHT.

DAMPER OPERATION
PULL DAMPER TO "OPEN"
PUSH DAMPER TO "CLOSE"

FUMING PROCEDURE

1 INSERT EVIDENCE INTO A CLEAN CHAMBER AND CLOSE THE DOOR AND DAMPER.

2 AFTER EVIDENCE IS PROPERLY EXPOSED, OPEN DAMPER AND EVACUATE THE CHAMBER WITH THE DOOR CLOSED FOR 5 MINUTES.

3 OPEN CHAMBER DOOR AND REMOVE EVIDENCE.

4 KEEP DAMPERS CLOSED ON CHAMBERS WHEN NOT IN USE.

5 CLEAN GLASS AND CHAMBER INTERIORS WEEKLY OR AS NEEDED.

THIS PANEL MAY BE REMOVED TO ACCESS VENT DAMPER PLENUM AND GUIDES FOR CLEANING OF INTERNAL COMPONENTS.

MAINTENANCE

1) PERIODIC CLEANING OF INTERIOR SURFACES IS RECOMMENDED TO MINIMIZE BUILD-UP OF DEVELOPING AGENTS. LAMINATE CLAD INTERIOR SURFACES AND GLASS DOOR PANELS MAY BE CLEANED WITH SOLVENTS AS RECOMMENDED BY DEVELOPING AGENT MFGR. AVOID SOLVENT CONTACT ON WINDOW GLAZING AND VENT PORT GROMMETS.

2) RELAMP BY UNSCREWING INTERIOR GLASS GLOBE FROM LAMP HOUSING. DO NOT USE BULBS IN EXCESS OF 100 WATTS.

3) WIPE DOWN EXTERIOR SPILLS AS SOON AS POSSIBLE.

4) VERIFY PROPER AIRFLOW VOLUMES ANNUALLY. SEE FOLLOWING PAGES FOR SPECIFICATIONS.

5) REPLACE 'V'-SEAL DOOR GASKETS AS NECESSARY. SEAL IS AVAILABLE AT MOST HARDWARE STORES.


6) LUBRICATE DAMPERS AS REQUIRED FOR SMOOTH OPERATION.


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PRODUCT DETAILS

BENCH-TOP UNIT 950P200

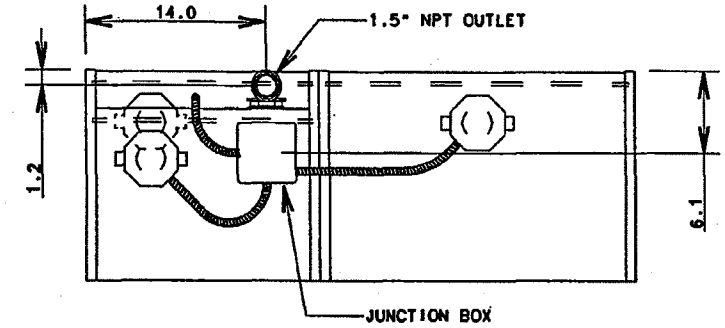

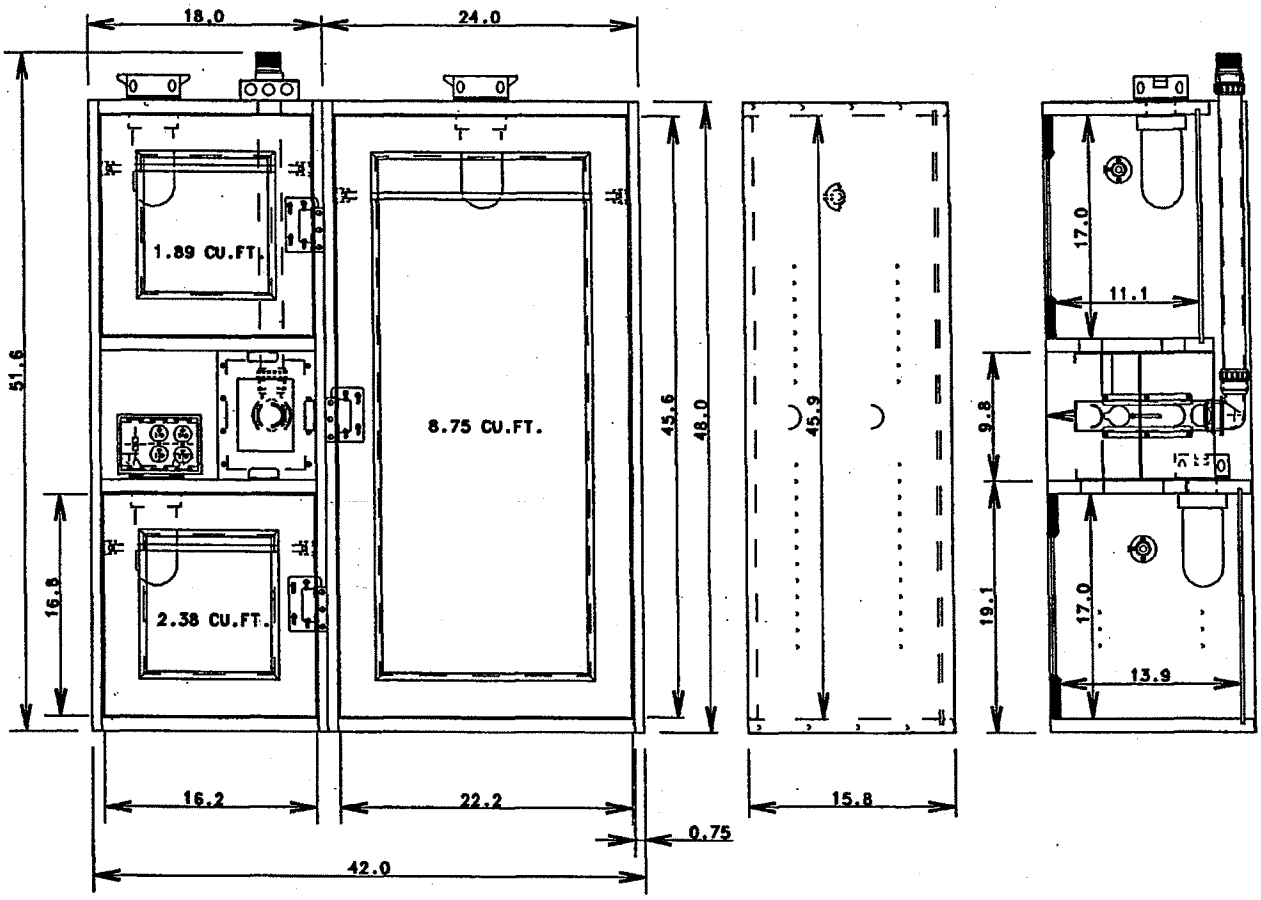

MINIMUM AIR FLOW - 13 CFM


CONNECT TO 1.5" I.D. MINIMUN EXHAUST LINE

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PRODUCT DETAILS

FULL HEIGHT UNIT 950P201

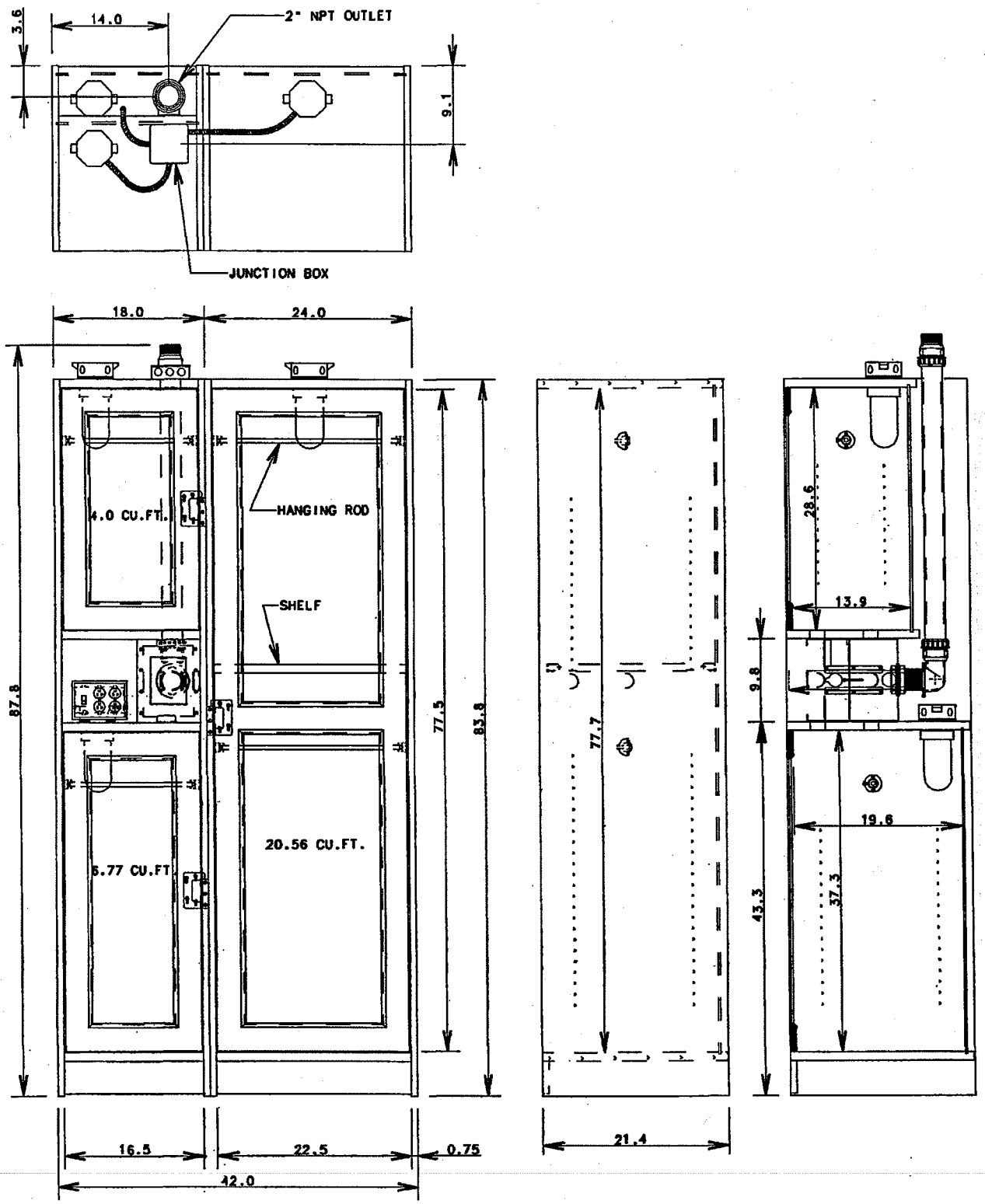

MINIMUM AIR FLOW - 32 CFM


CONNECT TO 2.0" I.D. MINIMUN EXHAUST LINE

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

WIRING DIAGRAM

CONNECT TO 20 AMP CIRCUIT

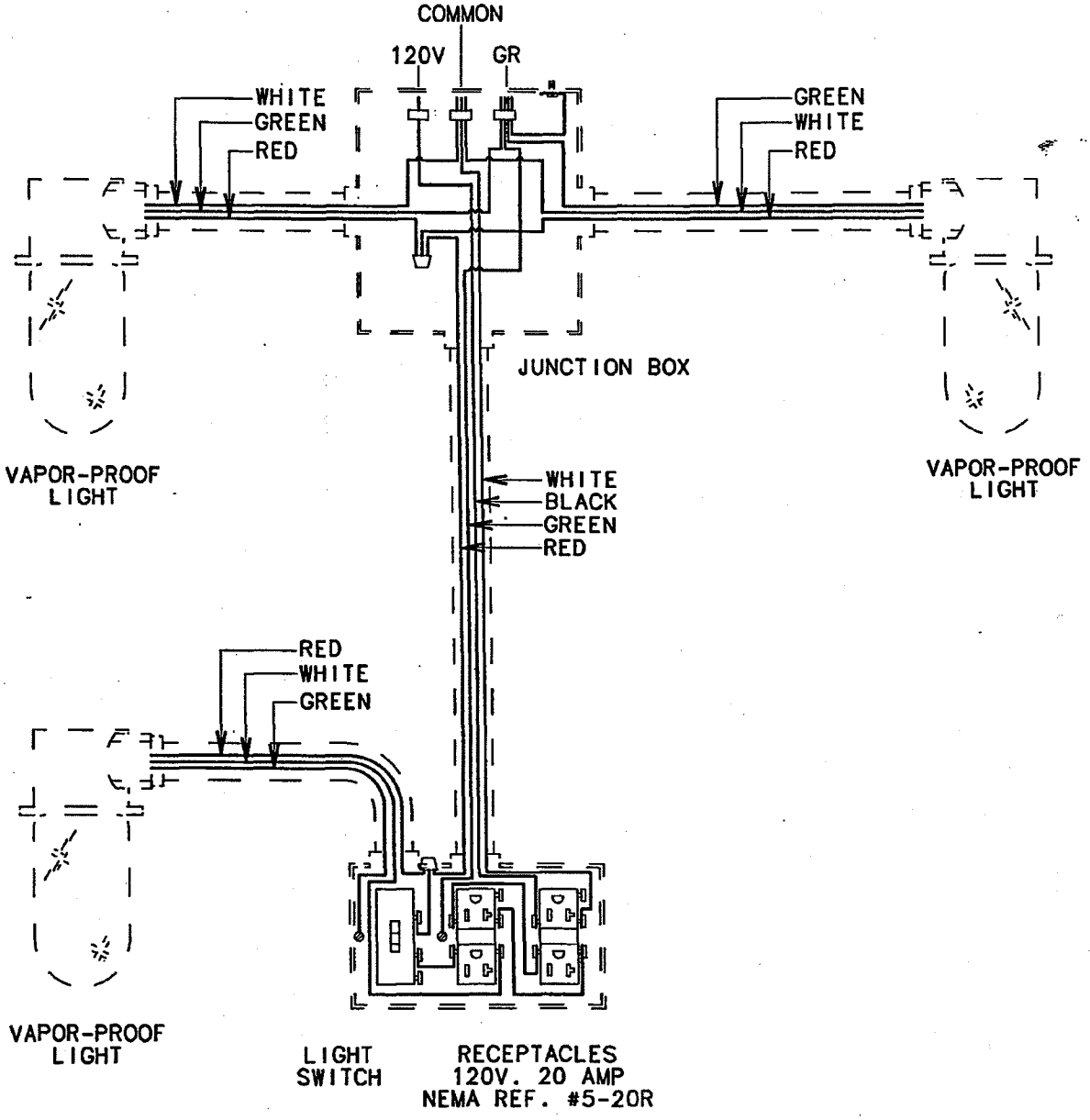

Fisher Hamilton

Fisher Hamilton Inc.
1316 18th Street
Two Rivers, WI 54241
Tel: 920-793-1121
Fax: 920-793-3084

PL-641 Jan98 P/N 70776


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

950P200G-Y

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

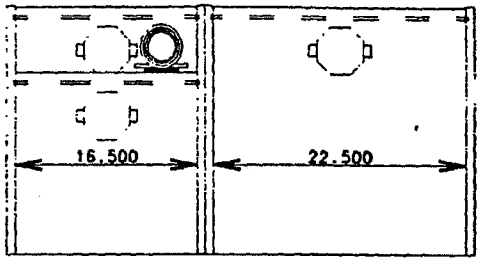

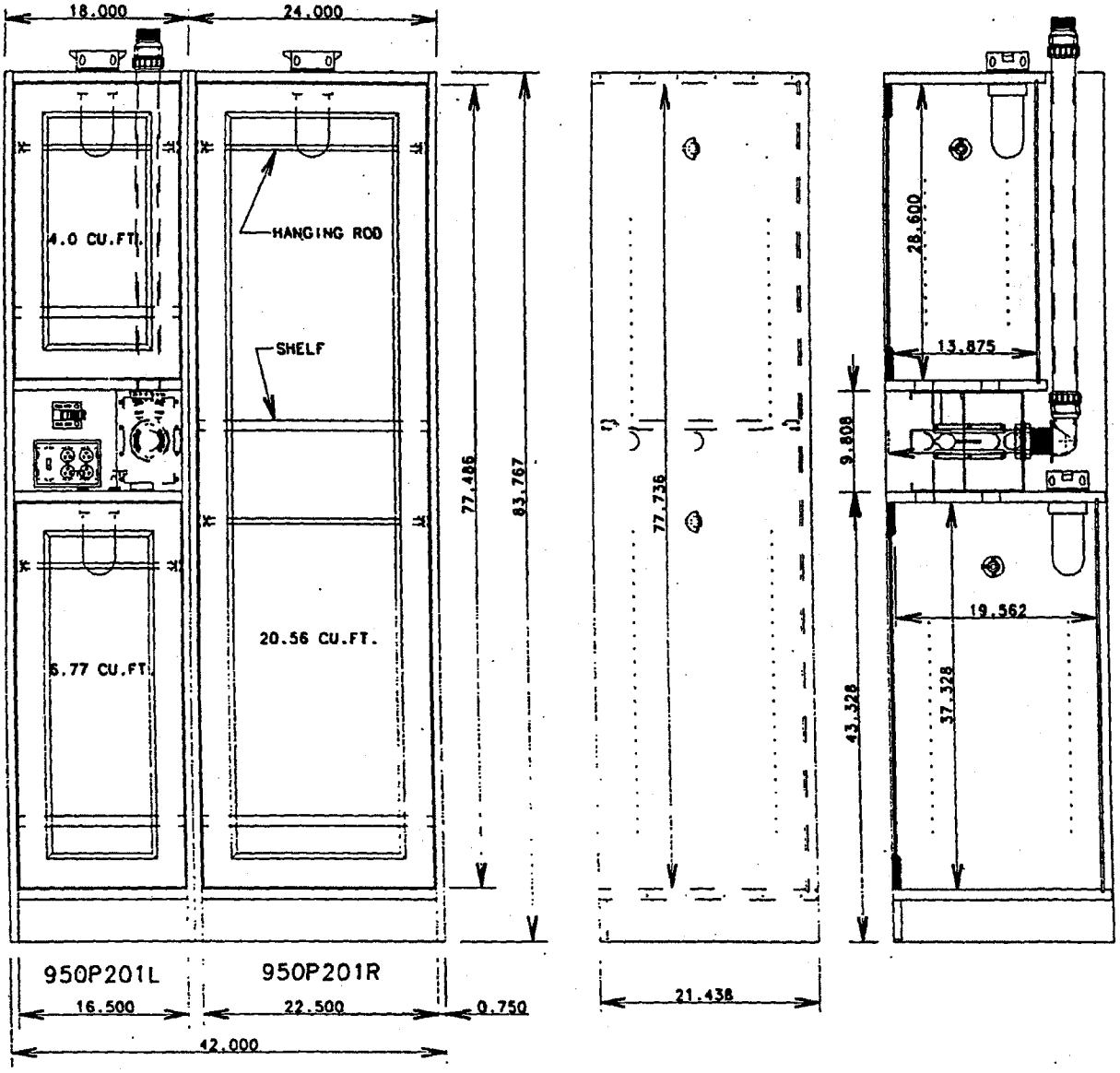

950P201

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP