ARTISAN®
TECHNOLOGY GROUP

Your **definitive** source for quality pre-owned equipment.

Artisan Technology Group
(217) 352-9330 | sales@artisantg.com | artisantg.com

Full-service, independent repair center
with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment
along with credit for buybacks and trade-ins.

Custom engineering
so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

Find the *Gardner Denver Welch Duo-Seal 1400* at our website: ***Click HERE***

Welch
Rietschle Thomas
A Thomas Industries Company

OWNER'S MANUAL
FOR
DUOSEAL® VACUUM PUMPS

MODELS:
1400, 1405, 1402,
1376, 1397, 1374,
1399, 1380, 1392



Contents:

Installation
Operation
Maintenance
Specification
Dimensional Drawings
Parts List
Accessories
MSDS Sheet for DUOSEAL® OIL



WARNING
Never Block the Exhaust Port.
If the exhaust is blocked, pressure will build-up in the pump with the potential of the pump body bursting and causing possible injury to personnel in the area.

Welch Rietschle Thomas
7301 North Central Ave.
Skokie, IL 60077
Phone: (847) 676-8800
Fax: (847) 677-8606 (Technical Support)
Fax: (920) 451-4397 (Ordering)
Web Page: www.welchvacuum.com
E-mail: welchvacuum@thomasind.com

Part No. 67-0777R2.9
Printed in the U.S.A.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

CONTENTS

SECTION		PAGE
1.	Installation	3-5
2.	Operation	6-8
3.	Maintenance	9-11
4.	Trouble Shooting Guide	12
5.	Quick Reference Charts	13-14
6.	Pumping Speed Curves	15
7.	Dimensional Drawings (with rotation)	16-17
8.	Exploded Parts Listings	18-31
9.	Accessory Section	32-33
10.	MSDS Sheet For DUOSEAL® OIL	34

PLEASE READ BEFORE OPERATION

While reading your manual, please pay close attention to areas labeled **WARNING AND CAUTIONS**. The description of each is found below.

WARNING
Warnings are given where failure to observe instruction could result in injury or death to people

CAUTION
Cautions are found where failure to observe the instruction should result in damage to the equipment, associated equipment and process.

These units conform to the SI International system of units of measurement.
The following symbols (with recommendations of IEC1010) of warning will be found on the pump.

Caution- refer to accompanying documents 

Caution- risk of electrical shock 

Caution- hot surface 

WARNING
Motor includes a self resetting thermal cutout and the pump could restart without actuation under fault condition.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 1: INSTALLATION

1.1 Introduction

This manual has been complied not only for the care and maintenance of the DUOSEAL pump now in your possession but as a helpful reference and guide for many problems which are usually associated with mechanical vacuum pumps. Take time to read these instructions carefully and preserve this manual for future reference; we think it will be useful to you.

1.2 Unpacking

Carefully remove the pump from the shipping case and unfasten and remove the wooden skid. Preserve all paper work and inspection tags for future reference. If damage has occurred from shipment a claim must be filed with the carrier immediately; preserve the shipping container for inspection by the carrier. If you are required to communicate with your dealer or with Welch Vacuum be sure to include your order numbers for quick identification. Do not return the pump to the factory with out first calling for a returned goods number.

1.3 Pump Mounting

1.3a Mounted Pumps

Rubber bumpers are supplied with most of our mounted pumps, either loosely or attached. Bumpers are excellent for applications involving a semi-flexible surface such as a bench top; they help to isolate noise and eliminate creeping. For more rigid requirements, the pump base may be bolted directly to a firm foundation with or without the bumpers. All DUOSEAL pumps should be mounted in a horizontal plane.

1.3b Unmounted Pumps

If you have purchased an unmounted pump, refer to parts list for information concerning the motor, motor pulley and belt necessary to drive your particular pump at the recommended speed.

1.4 Pump Location

The pump should be located preferably in a clean and well ventilated area and adequate space should be provided wherever possible for routine maintenance such as changes of oil and belt adjustments and replacements. Above all, the pump should be located as closely as possible to its system in order to utilize it most efficiently. Its location should include such determining factors as the length and size of connections, the number of bends and the type of exhaust connections.

1.5 Exhaust Provisions (See Accessory section in back of this manual.)

Exhaust connections will be determined by the type of system to be exhausted and the desired cleanliness of the atmosphere surrounding the pump. Under normal conditions of mild evacuation nothing more than the dust cap will be necessary to cover the port. Where relatively high gas flows are involved or where the presence of oil vapor is objectionable an exhaust filter may be fastened to the exhaust port in place of the dust cap. The exhaust filters used on our pumps are capable of absorbing and restricting any vapor particle larger than 0.1 micron.



WARNING

Never block the exhaust port. If the exhaust port is blocked, pressure will built-up in the pump with the potential of the pump body bursting and causing possible injury to personnel in the area.

Where extreme exhaust conditions are encountered it is best to pipe the exhaust direct out of the building. Welch recommended wire reinforced PVC, metal pipe or thick walled rubber hose be used as exhaust lines to avoid potential of line becoming crimped or collapsing resulting in the exhaust port being clocked. Be sure to call Welch technical service prior to start-up at 847-676-8800 if you have any questions.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

1.6 Electrical Power

1.6a Power source Review

Review the power source and the motor rating to be sure they agree in voltage, phase and frequency. On three-phase applications the direction of rotation of the motor must be considered. Make a momentary check rotation at the time of power installation and wiring. Momentary backward rotation of the pump is not harmful. Check the layout drawings for proper direction of rotation.



CAUTION
Make certain the power settings on the pump match your power source before attempting to operate the pump



1.6b Overload Protection

Motor thermal overload protection is made available by the motor manufacturer as an aid to minimizing motor failure. Overload protection is a standard feature on all single-phase 60Hz motors. Single-phase motors will normally have automatic overload protection. Motors of 1-1/2 horsepower or larger supplied with DUOSEAL pumps contain no overload protection. Installations of such equipment must comply with local electrical codes which dictate appropriate starter and protection devices. It is strongly suggested that you familiarize yourself with the protection supplied with your motor so that you may react accordingly in the event of an emergency. Automatic reset protection is designed to reset itself after a predetermined cooling period. If the fault to the drive remains unaltered, the motor will cycle on and off until the fault is corrected. The motor data plate will indicate the presence of thermal protection.

1.7 Vacuum Connections (See Accessory Section in back of this manual.)

1.7a Choice of Connections

The choice of connections and fitting can have a very marked effect on the pumping speed at the vacuum chamber. Any connection placed between the pump and the chamber creates an impedance to the flow of gas. This is particularly true at low pressures in the millitorr range where the gas flow is substantially molecular in character. The gas flow is then dependent upon the kinetic activity of the molecules to bring it to the intake of the pump.

1.7b The Effects of Condutance

It has been shown that the conductance of a tube is proportional to the cube of its radius and inversely proportional to its length. Therefore it is imperative that the connecting lines be as large in diameter and as short in length as practical. For best results the diameter of the connecting tube should be at least as large as the diameter of the pump intake. To avoid a large reduction in pumping speed at the vacuum chamber, it is clear that the conductance of the line must be considerably greater than the speed of the pump.

1.7c Metal Joints

If metal piping or tubing is used, it is preferable to solder or braze all of the connections. Where threaded joints must be used, coat the threads with LocTite® Thread Sealant with Teflon®, or Leak Lock and screw together tightly. Flanged connections with elastomer gaskets make excellent demountable joints. Modular vacuum piping and fittings are now extensively used.

1.7d Rubber Tubing Joints

Where metal tubing is used between the system and the pump intake, joints can be make by butting the ends of the two sections together in a short section of vacuum hose.Worm-screw band clamps are useful for securing the hose to the tubing.Whatever the joint you choose to use, cleanliness should be of utmost importance.

1.7e Valves and Stopcocks

Metal valves or stopcocks may be used in the connecting line between the system and the pump to provide a means of isolating the pump from the system. To minimize the impedance of flow, the valve openings should be as large as possible. Lubricate the rotating plug of the stopcock with a film of vacuum grease sufficiently thick enough to prevent seizure.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

1.8 Vacuum Gauges (See Accessory section in the back of this manual.)

The type of vacuum gauge to be used is determined largely by the pressure range to be measured. Pressures in the ranges produced by DUOSEAL pumps can be covered by McLeod, Pirani or Thermocouple gauges. The McLeod gauge is used where high accuracy of measurement is required. The Pirani and Thermocouple gauges are electrical and give continuous readings of the total pressure. There are preferred where rapid pressure changes occur. The McLeod gauge does not measure condensable vapors; therefore, if vapors are present it will generally read lower in pressure than electrical gauges. For higher vacuums in systems employing diffusion, turbo-molecular or ion pumps, a hot or cold cathode gauge is used.

1.9 Traps (See Accessory section in the back of this manual.)

1.9a The need for a Trap

Where corrosive vapors or large quantities of condensable vapors are evolved from vacuum processing, a cold trap may be used in the connecting line to the pump. It will help prevent damage to the pump mechanism and reduce oil contamination. The cold trap, immersed in a suitable Dewar flask, is installed to that the vapors may come in contact with the surfaces of the trap and condense. Commonly used refrigerants are liquid nitrogen or dry ice and acetone. The refrigerant to be used depends upon the freezing point of the contaminations. A variety of cold traps are available from Welch Vacuum, Thomas Industries, inc.

1.9b The care of a Trap

When using a cold trap the refrigerant should be maintained at a high level in the flask to keep the trap at a uniformly low temperature. If the trap is re-warmed it may allow re-evaporation of the condensate. The refrigerant add tube on the liquid nitrogen trap should not be obstructed as the refrigerant boil-off can produce dangerously high pressures. If the trap becomes saturated it should be disconnected from the system, drained and cleaned. An increase in pressure in the vacuum system will normally indicate that the trap has become saturated. To clean the trap, remove the trap from the system, allow the trap to warm up and rinse off the condensate with a suitable solvent in a fume hood. Thoroughly clean and dry the trap before reinstalling in the system.

1.10 Types of Lubricants (See Accessory section in the back of this manual.)

All DUOSEAL mechanical vacuum pumps are normally tested with DUOSEAL® oil and shipped with a full charge to prevent unnecessary contamination. An additional supply of oil is furnished with each pump with instruction to drain and discard the oil contaminated in the pump and replace with the fresh oil. DUOSEAL oil has been especially prepared and is ideally suited for use in mechanical vacuum pumps because of its desirable viscosity, low vapor pressure and chemical stability.

The vacuum guarantee on all DUOSEAL pumps applies only when DUOSEAL oil is used. Other lubricants for special applications are available including various lubricants compatibility, lubricants for use with diffusion pumps as well as other special requirements.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 2: OPERATION

2.1 Starting Procedures

2.1a Starting a DUOSEAL Pump

Before attaching the pump to a system it is well to familiarize yourself with the function and action of the pump which you have now required.Remove the intake and exhaust port plugs and temporarily provide a stopper for the intake and a dust cap for the exhaustReview the power requirements as described in Section 1.6.



CAUTION

Do not run this pump with the intake open directly to the atmosphere. The pump will overheat, excessive oil mist will be emitted from the exhaust and the pump will eventually seize.

2.1b Starting of the Model 1392 with Diffusion Pump.

1. Turn on Mechanical pump ONLY until a pressure of 100 millitorr or less has been reached.



WARNING

Do NOT operate diffusion pump until a pressure of 100 millitorr has been reached.Turning on the diffusion pump before a pressure of 100 millitorr has been reached will damage your diffusion pump and may cause serious damage to surroundingequipment and personnel.

2. After a pressure of 100 millitorr had been reached by the mechanical pump, plug in the diffusion pump. Make sure diffusion pump has oil in it at an appropriate level before turning it on. Check the oil indicator sight glass to confirm oil level. Also, make sure that fan is on when operating diffusion pump. Note that it will take the diffusion pump a few minutes to warm up.

Refer to the Varian diffusion pump manual for more information.

2.1c Cleanliness

Take every precaution to prevent foreign from entering the pump. A fine mesh screen is provided for this purpose in the intake passage of all DUOSEAL pumps.

2.1d Oil Level Determination

The amount of oil suitable for efficient and satisfactory performance should be determined after the pump has reached its operating temperature. Initially, however, the pump should be filled with fresh oil while the pump is idle. Fill the pump until the oil level falls half way of the oil level window. If after a short period of operation the level should fall, it is likely the result of oil entering some of the interior pockets of the pump. If the oil level rises, this signifies oil has drained into the pump cavity while pump was idle. Shut off pump, then drain oil down to proper level.

If a gurgling sound occurs, additional oil must be added. Mechanical pumps will gurgle in varying degrees under four conditions of performance: [a] when operating at high pressure as in the beginning cycles of evacuation of a chamber; [b] when the oil level in the pump reservoir is lower than required; [c] when a large leak is present in the system; and [d] when the gas ballast is open. Awareness of the possibilities will save time in setting up a system. Best performance of a mechanical pump is generally obtained after sufficient time has been allowed for the pump to come to operating temperature.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

2.2 Leak Detection

2.2a Large Leaks

The importance of eliminating all leaks in a vacuum system is obvious when it is realized that a leak into the system, at atmospheric pressure, expands in volume by a factor of 750,000 to 10,000,000 or more. The pump must remove this added volume to maintain the desired vacuum. Fortunately a number of effective techniques for leak detection have been developed. Large leaks can be located by pressurizing the system and painting the suspected area with a thick soap solution. Escaping air will produce soap bubbles.

2.2b Small Leaks

Small leaks may also be detected by spraying a suspected area with acetone or gases rich in hydrogen, and observing a sudden change in pressure on an electrical gauge. The difference in calibration of these gauges, for air and other gases, will produce a distinct change in the pressure reading. To use this method of detection, the system must be under vacuum and the gauge sensing tube must be located between the pump and the area to be probed. Use extreme caution, as these materials are highly flammable!

2.2c Fine Leaks

Locating very fine leaks requires a helium-sensitive, mass-spectrometer leak detector. This instrument will locate leaks which cannot be detected by any other method. Numerous fine leaks can have the total effect of a large leak.

2.3 High Pressure Operation

DOUSEAL two-stage pumps are designed to be most efficient when operated at or near their ultimate blankoff pressure. The pump is designed to remove air from a vacuum system at atmospheric pressure to a deep vacuum. When the pump is operated at elevated pressures up to about 10 Torr for long periods of time, the pump will run hotter. At elevated ambient temperatures under these conditions, the thermal protection switch on single phase motors may cute out. Use of an exhaust filter with coalescing element is required.



CAUTION

Do not run this pump for prolonged periods of time at or above 10 Torr. The pump will overheat, excessive oil mist will be emitted from the exhaust and the pump will eventually seize.



Leaky vacuum systems, too large a vacuum chamber for pump size or high gas loads to the vacuum systems are causes for extended runs of the pump at elevated pressure. If you need additional information, please contact Welch Vacuum at 847-676-8800, extension 1, for further information and precautions.

2.4 Shutdown Procedures

2.4a DuoSeal Shutdown

A few simple precautions are all that is necessary when shutdown is in order. If a gauge is connected to the system, first isolate the gauge, then turn off the power and open the system to atmosphere. If the pump is removed from the system, cover the intake port with a rubber stopper or suitable cover to protect the pump against contamination and loose particles. If the pump has been contaminated in service and is going to be shelved for a prolonged period it is best to drain the oil refill with a fresh charge.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

2.5 the Principle of Gas Ballast

2.5a the Effects of Unwanted Vapor

systems which contain undesirable vapors cause difficulty both from the standpoint of attaining desirable ultimate pressures as well as contamination of the lubricating medium. A vapor is defined as the gaseous form of any substance which is usually a liquid or a solid. Water, oil and a mercury vapors are three of the more common vapors encountered in typical vacuum systems. When such vapors exist in a system, the vapors or mixtures of gas and vapor are subject to condensation within the pump; the precipitated liquid may thus ultimately dissolve or become emulsified with the lubricating medium. This emulsion is recirculated to the chambers of the pump where it is again volatilized causing increased pressure within the system.

2.5b the Presence and Removal of Condensate

Condensation takes place particularly in the compression stroke of the backing or second stage of a two-stage pump. The compression stroke is that portion of the cycle during which the gas drawn from the intake port is compresses to the pressure necessary to expel it past the exhaust valve. Condensation takes place when the ratio between the initial pressure and the end pressure of the compression is high, that is, when the mixture of vapor and gas drawn from the intake port is compressed from a low pressure to high pressure. By adding air though the gas ballast valve to the mixture of vapor and gas being compressed, the pressure required for delivery past the exhaust valve is reached with a considerably smaller reduction of the volume of the mixture; this, depending upon the amount of air added, condensation of the vapor is wither entirely avoided or substantial reduced.

2.5c Pump function With Gas Ballast

In a pump functioning on a contaminated system and operating without the gas ballast, compression within the stage takes place in the normal manner until the saturation pressure of the contaminating vapor contained within the mixture of gas and vapor is reached. The saturation pressure of water vapor is that pressure and corresponding temperature at which the dew point of the vapor is reached and condensation occurs. The saturation pressure of water vapor at an ambient temperature of 20°C is 17.5 Torr, while at 60°C, the approximate operating temperature of a pump, the saturation pressure is 149 Torr. The external side of the exhaust valve is subjected to atmospheric pressure. Consequently a compressive force somewhat greater than atmospheric pressure is required to open the valve and permit expulsion of the gas. Sometime during increased compression of the mixture of gas and vapors, the saturation pressure of 149 Torr for the water vapor is reached and the vapor condenses. The condensate is then allowed to emulsify with the oil which is recirculated within the pump stages this providing continued contamination of the system.

2.5d Pump function with Gas Ballast

On the other hand, when ballast air at atmospheric pressure is supplied to the compression stroke by means of the gas ballast, the partial pressure of the unwanted vapor becomes a very small part of the total pressure of the mixture of gas, vapor and newly supplied air. The vapor is thus prevented from reaching its saturation pressure corresponding to the temperature of the pump and is finally expelled form the pump as a vapor.

2.5e Controlled Ballast Flow

Some degree of variation in ballast flow may be obtained by the amount of opening applied tot he gas ballast. Two or more turns of the gas ballast are sufficient to open it wide. With the gas ballast open, the sound of the exhaust is similar to that of a pump operating against a large leak. Because of the increased pressure introduced into the compression stroke, the pump must work a little hard to function, thus resulting in an increased operating temperature of approximately 8°C over a prolonged period of time. Tests have shown that continuous and prolonged operation for several weeks under these conditions is not injurious tot he pump.

2.5f Other Forms of contamination control

The application of the gas ballast is a moderate and very successful method for the removal of condensable vapors. For very heavily laden systems, other means of removal such as oil separators may be required. For mild cases of contamination the simple expedient of a cold trap or a change of oil serve the purpose.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 3: MAINTENANCE

3.1 Vacuum Problems

3.1a Pressure Determinations

Leakage, contamination and unusual outgassing are the genera causes of problems associated with poor vacuum. To operate a maximum efficiency a system must be thoroughly clean. If the system is completely clean and free from leaks, and unwarranted vacuum problems still exist, the pump should be checked. A simple criterion for the condition of a mechanical pump is a determination of its ultimate pressure capability. This can be accomplished by attaching a gauge directly to the pump. The gauge may be any suitable type provided consideration if given tot e limitations of the gauge being used. Refer to Paragraph 1-8 for further suggestions. If the pressure is unusually high, the pump may be badly contaminated, low on oil or malfunctioning. On the other hand, if the pressure is only slightly higher than the guaranteed pressure of the pump, an oil change may be all that is required.

3.1b Oil Contamination

the most common cause of a loss in efficiency in a mechanical pump is contamination of oil. It is caused by condensation of vapors and by foreign particles. The undesirable condensate emulsifies with the oil which is recirculated and subject to re-evaporation during the normal cycle of pump activity this reducing the ultimate vacuum attainable. Some foreign particles and vapors may form a sludge with the oil, impair sealing and lubrication and cause eventual seizure. A gas ballast calve is helpful in removing vapors, especially water, but it is not equally effective on all foreign substances; therefore, periodic oil changes are necessary to maintain efficient operation of the system. The required frequency of changes will vary with the particular system,. Experiences with the process will help you determine the normal period of operation before an oil change required.

3.1c Oil Overheating

This pump is designed to operate continuously below 10 Torr. Continuous operation of this pump above 10 Torr will lead to overheating and eventual pump failure. See Section 2.3. High Pressure Operation.

3.2 Maintenance of the Model 1392 with Diffusion Pump

3.2a Determination of an Oil Change

The Model 1400 pump of the diffusion-mechanical pump combination may be repaired and maintained as outlined in the preceding paragraphs, The pump fluid of the diffusion pump should be inspected periodically for color change and odor. If the fluid is slightly darkened and had no odor, a simple change of fluid is normally all that is needed. If the fluid had developed a noticeable odor, the pump should be cleaned. Follow the instructions in the Owner's Manual for the Diffusion Pump.

3.2b Diffusion Pump Cleaning Procedures

To clean the diffusion pump, first turn off the power. Turn off the cooling fan and disconnect the power plug after the diffusion pump has cooled. Disconnect the inlet and foreline connections. Do not scratch or damage sealing surfaces. Remove the cold cap and jet assembly and all O-rings and gaskets. Drain oil from the diffusion pump. Thoroughly clean the diffusion pump body interior and the jet assembly with acetone followed by an isopropyl alcohol rinse. The inside of the pump can be dried using dry nitrogen. Install jet assembly with he cold cap inside of the pump body. Fill the pump with 30ml of fresh diffusion pump fluid. Reinstall the pump in the system. For more detailed information and how to replace a burned heater red instructions on the AX-65 Diffusion Pump Instruction Manual supplied with this pump.



WARNING

Solvents, such as acetone, used to clean different pumps are extremely hazard ous. Inhalation may be fatal or harmful. Acetone is also extremely flammable. Never use acetone near heat or open flame. Use only in a well ventilated area.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

3.3 Oil Changes and Oil Level

3.3a Developing a Maintenance Schedule

After studying many examples of pump failure, Welch has found the most common reason is poor condition of the oil. This is why a regular maintenance schedule for the oil is critical to obtain the longest service life out of your DUOSEAL pump. Welch recommends that you examine the condition of the oil on a daily basis in the early days of a new process or experiment. You want to look for discoloration of the oil and whether the oil level is rising. The discoloration can indicate deterioration of the oil and a rising oil level can indicate condensation of vapors is occurring in the pump. When changes occur, the oil needs to be change. If no changes in the oil level or color are observed, extended pump service life is obtained if the oil is changed every three to four months.

3.3b Forced Oil Flushing

When you drain oil through the drain valve, you are not removing the oil and contaminants that are inside the pumping mechanism. You are removing oil only from the oil case. Welch recommends a forced oil flush of DUOSEAL pumps be performed at the regular maintenance oil change. The procedure for the forced oil flush is given below.

1. Check the oil level.

If the oil level is well above the fill mark, (This may indicate either the pump has been overfilled with oil or has ingested a liquid or a large amount of vapor water or organic solvents.) Please go to step 2.
If the oil level is even with the fill mark and you do **NOT** suspect corrosive gases or particulates (hence forth called **contaminants**) ingested have damaged the mechanism, run the pump for 15 minutes to allow the pump oil to warm up before going to step 2.

2. Turn off motor for the vacuum pump. Drain the oil into a clear plastic container (may need pliers to open drain valve). Look for contaminations settling to the bottom of container. If you see contaminants, you will need to repeat step 3 through 5 several times until the oil comes out clear. The oil you drained from the pump came from the oil case only. There may be contaminants in the pumping mechanism. To be sure all contaminants have been removed, the pump mechanism needs to be flushed.

3. Make sure the belt guard is installed before proceeding further**.** Attach a short hose to the drain valve which runs into a clear plastic container.

4. Flushing the pump is carried out by adding a cup of DUOSEAL vacuum pump oil through the intake port (IN) while the pump is turned on for 15 to 20 seconds. While adding the pump oil, the palm of your hand is placed lightly over the exhaust port (OUT). Look for water coming out of the drain. Turn off the pump.

5. Repeat step 4 until clean oil comes out of the drain hose.

6. Fill the pump with the amount of DUOSEAL vacuum pump oil your pump needs.

7. Plug the intake (IN) port with a rubber stopper. Turn the pump on and run the pump for 10 minutes. Close the gas ballast.

8. Check the vacuum reading of the pump by connecting a thermocouple gauge tube to the pump's intake. If the pump is running nearly as good as when it was new, the total pressure reading you will read on wither of these two gauges will be at least 10 micron.

A simple way to connect the gauge tube to the pump is to run the threaded tip of tube through a hole in a rubber stopper. Use pump oil as a lubricant for inserting the tube. The stopper chosen should be bigger than the outer diameter of the intake flange.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

3.3c Refilling The Pump

After you are satisfied that the pump has been thoroughly flushed, refill the pump by pouring new DUOSEAL oil into the exhaust port. Fill to the indicated level and start the pump with the intake closed. A gurgling noise is characteristic when high pressure air is drawn through the pump. It should disappear quickly as the pressure within the pump is reduced. If gurgling continues, add sufficient additional oil through the exhaust port until gurgling ceases.

3.4 Shaft Seal Replacement

To replace the shaft seal of a pump, drain the oil and remove the pump pulley and key. Remove the screws securing the old seal and pry it loose with a screwdriver or similar wedge, being careful not to mar the surface of the pump body against which the seal fits. Discard the seal and its gasket, inspect all surfaces and repair any damages with a fine abrasive stone. Wipe all sealing areas clean and place a film of DUOSEAL oil on both the shaft and the inside bore of the new shaft seal. Using a new gasket, carefully slide the new seal into position and center it on the shaft. It is not necessary to apply any sealant to the gasket. Tighten the mounting screws uniformly and refill the pump with DUOSEAL oil. Follow instructions included in repair kit.

3.5 Repairing Vacuum Leaks

3.5a Location, Cause and Effect

Oil leaks may develop wherever two mating faces are sealed with a gasket. Such seams may fail as the result of deterioration of the gasket material, loosening of the screws caused by temperature variations, or improper care as the result of previous reassembly. Typical gaskets seams in a mechanical pump are located at the oil level window, the shaft seal, the oil drain, and the mating faces of such mechanical surfaces as the intake chamber cover, the oil case, and the exhaust chamber cover. The importance of a gasketed seam is determined principally by its function. If it is a vacuum seal, the ultimate performance of the pump is dependent upon it. If it is an oil seal, the pump may be operated satisfactorily for some time without loss of function. Eventually, of course, a great loss of oil may cause harmful damage.

3.5b Repairing Techniques

An oil seam may be sealed by any of several methods. When an O-ring is employed, the surfaces of the O-ring and its groove should be wiped clean. If the O-ring is not badly deformed or scratched it may be reused by sealing with a slight film of vacuum oil or vacuum grease. Thin composition gaskets are generally used for large irregularly shaped areas. A replacement joint of this type should be thoroughly cleaned of all previous gasket material and the mating surfaces cleaned of any nicks.

3.6 Repairing Vacuum Leaks

3.6a Surface Preparation

Good Vacuum seals are an essential and important attribute of a good mechanical pump. A good seal is dependent upon the quality of the mating surfaces as well as the sealant and its preparation. The mating faces should be carefully inspected for any projections or foreign particles which might interfere with proper mating. Slight projections such as nicks and burrs are most easily removed by rubbing with a fine abrasive stone. The surface of the mating parts may be washed with a solvent or alcohol after which they must be thoroughly dried.

3.6b Temporary Repair

Temporary vacuum repairs are often made by covering the known leak with an industrial sealant such as Loctite with Teflon. Such a practice, however, is not recommended for seals of a permanent nature.

3.7 Drive Problems

If for any reason the pump will not operate, turn off the power and check the fuse and electrical connections. Then try the power to the motor only by removing the belt. If the motor operates properly try hand-rotating the pump in the proper direction with the pump intake port open. If both turn freely then replace the belt and check the belt tension. The tension should be sufficient to drive the pump without visible slippage. Any greater tension will cause noise and possible damage to the bearings of both the motor and pump. Make certain that both pulley grooves are clean and free from oil. The pulleys must be fastened securely on their respective shafts, and in parallel alignment.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 4: TROUBLE-SHOOTING GUIDE

4.1 Poor Vacuum Reading

1. Disconnect vacuum pump from your vacuum system.

2. Obtain 1 gallon of DUOSEAL oil, 1 gallon or larger clear plastic container, clear drain hose (1/4 inch ID by at least 2.5 feet long), rubber gloves, goggles, apron and pliers.

3. Carry out a forced oil flush (See Section 3.3b) if motor can turn pumping mechanism without excessive noise. If pumping mechanism cannot be turned by motor, have pump sent into factory repair center.

4. Check ultimate vacuum with thermocouple gauge. If reading is still unsatisfactory, consult with Welch Technical Service after checking for leaks in your connections.

4.2 Noisy Pump

1. Locate noise source on your vacuum pump (motor, belt, pump).

2. If motor makes noise:
   . Check for loose bolts on motor mount.
   . Motor may need to be replaced

3. If belt makes noise:
   . Pulley may be bent. If so, replace.
   . Loose belt, tighten by loosening bolts on motor and moving motor away from pump until belt exhibits a 1/2 inch deflection.

4. If pump makes noise: vanes may be sticking, carry out a forced oil flush.

5. If noise level is still unsatisfactory, consult with WELCH Technical Service at (847)676-8800, Extension 1.

4.3 Excessive Oil Mist

1. Check oil level. Maintain the level between the two lines on the oil window during operation. Too high a level will cause excess oil mist.

2. Determine the vacuum level in your system with a thermocouple gauge or manometer. **Be sure** the pump is operating at a pressure of 10 Torr or lower for continuous operation.

CAUTION: if not, discuss application with Welch Technical service before proceeding to prevent oil from overheating and breaking down resulting in the pump failing.

3. A system leak(s) or degassing can cause excessive oil mist. Locate and seal the vacuum leaks. High vapor loads may require a larger pump to compensate.

4. An oil eliminator can be installed for coalescing the oil mist from the exhaust gases. If an oil mist eliminator is already installed, the element may need to be replaced if saturated with oil. Open up oil mist eliminator and examine the element; replace if needed.

5. Check the oil condition for any visible change in color or consistency. contaminants in the oil can add to the amount of oil exhausted. Carry out a forced oil flush if contamination is discovered. (See Section 3.3b.)

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 5: QUICK REFERENCE CHART

Motor Options	1400	1405	1402	1376	1397	1374	1399	1380
Wired for 115V, 60Hz, 1Ph	1400B-01	1405B-01	1402B-01	1376B-01	1397B-01	1374B-01	1399B-01	1380B-01
Wired for 230V, 60Hz, 1Ph	1400C-01	1405C-01	1402C-01	N/A	N/A	N/A	N/A	N/A
Wired for 230V, 50Hz, 1Ph	1400C-02*	1405C-02*	1402C-02*	1376C-03*	1397C-03	1374C-02*	1399C-02	N/A
Wired for 230V, 60Hz, 3Ph	N/A	N/A	1402M-01	1376M-01	1397M-01	1374M-01	N/A	N/A
Wired for 230V, 50Hz, 3Ph	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Replacement Parts								
Shaft Seal: Lip Seal	1401E -01	1401E -01	1401D -01	1401D -01	1401F -01	1401F -01	1401E -01	1401D -01
Shaft Seal: Mechanical Seal	1401E	1401E	1401D	1401D	1401F	1401F	1401E	1401D
Dust Cap	41-1345	41-1345	41-0612	41-0612	41-1501	41-1501	41-1345	41-0612
Air Filter Screen	41-0890	41-0891	41-0660	41-1692	41-0937	41-0937	41-0306	41-0660
Intake Chamber Cover Gasket	41-0383	41-1308	41-0234	41-0234	41-1495	41-1495	41-2375	41-0234
Oil Case Gasket	41-1052	N/A	41-0403	41-0403	41-1496	41-1496	41-1052	41-0403
Drain Plug	61-2278	61-2278	61-2278	61-2278	61-2278	61-2278	61-2278	61-2278
Oil Window	41-3753	41-3753	41-3753	41-3753	41-3753	41-3753	41-3753	41-3753
Large Vane	41-1059	41-1294	41-0959	41-1684	41-1477	41-1477	41-1059	41-0959
Small Vane	41-1059	41-1294	41-1487	41-1487	41-1487	41-1487	N/A	N/A
Large Vane Spring	41-1304	41-1304	41-1304	41-1304	41-1488	41-1488	41-1304	41-1304
Small Vane Spring	41-1304	41-1304	41-1304	41-1304	41-1304	41-1304	N/A	N/A
Large Vane Spring Holder	41-1063	41-1303	41-1685	41-0696	41-1489	41-1489	41-1063	41-1685
Small Vane Spring Holder	41-1063	41-1303	41-0696	41-0696	41-0696	41-0696	N/A	N/A
Shaft	41-1060	41-1293	41-2093	41-1683	41-1484	41-1484	41-2374	41-1773
Upper Exhaust Valve	41-1239	41-1239	41-2153	41-2153	41-2161	41-2161	41-1239	41-0695
Lower Exhaust Valve	N/A	N/A	41-2154	41-2154	N/A	N/A	N/A	41-1775
Pressure Release Valve	N/A	N/A	41-0992	41-0992	41-0992	41-0992	N/A	N/A
Repair Kits								
Minor Repair Kit	1400K-03	1405K-05	1402K-05	1376K-05	1397K-07	1397K-07	1399K-03	1380K-04
Major Repair Kit	1400K-04	1405K-06	1402K-06	1376K-06	1397K-08	1397K-08	1399K-04	1380K-03
Base and Belt Guard Kit	1371E	1371A-01	1371A-01	1371B-01	1371F	1371F	1371E	1371A-01

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

	1400B-01	1405B-01	1402B-01	1376B-01	1397B-01	1374B-01	1399B-01	1380B-01	1392B-01 MECH.	1392B-01 DIFFUSION
Free Air Displacement, L/M (CFM)	25 (0.9)	91 (3.2)	160 (5.6)	300 (10.6)	500 (17.7)	650 (23)	35 (1.2)	160 (5.6)	25 (0.9)	30L/Sec.
Gas Ballast	Yes	Yes	Yes	Yes	Yes	Yes	No	Yes	No	-
Ultimate Pressure, Torr (mbar)	1 (1.3)x10-4	1 (1.3)x10-4	1 (1.3)x10-4	1 (1.3)x10-4	1 (1.3)x10-4	1 (1.3)x10-4	15 (20)x10-4	5 (20)x10-4	1 (1.3)x10-4	<5 (6.5)x10-8
Rotational Pump (RPM)	580	525	525	525	400	510	750	525	580	-
Oil Capacity, quarts (Liters)	0.62 (0.59)	2.25 (2.1)	2.25 (2.1)	2.5 (2.4)	1.25 (1.2)	1.25 (1.2)	0.5 (0.47)	3.7 (3.5)	0.62 (0.59)	30ml
Net Wt. Unmounted Pump, lb (Kg)	33 (15)	82 (37)	82 (37)	107 (52)	138 (63)	138 (63)	28 (13)	57 (26)	33 (15)	8 (3.6)
Net Wt. Mounted Pump, lb (Kg)	58 (26)	112 (51)	112 (51)	156 (71)	206 (93)	218 (100)	49 (22)	91 (41)	71 (32)	-
Optional Exhaust Filter	1417	1417P-10	1417P-10	1417P-10	1417P-20	1417P-20	1417	1417P-10	1417	-
Optional Exhaust Filter Element	1417L	1417G	1417G	1417G	1417H-01	1417H-01	1417L	1417G	1417L	-
Pump Pulley	41-2191	41-2194	41-2074	41-2074	41-1492	41-1492	41-2191	41-2074	41-2191	-
Standard Belt	1400A	1405A	1405A	1405A	1397A	1374A	1399A	1405A	1400A	-
Standard Hose Barb Fitting	1393J	1393J	1393K	1393K	1393L	1393L	1393J	1393K	N/A	-
Tubing Needed (I.D. Inches)	7/16	7/16	13/16	13/16	1 5/8	1 5/8	7/16	13/16	N/A	-
Standard Flange Adapter Fitting	1393F	1393F	1393G	1393G	1393H	1393H	1393F	1393G	N/A	-
Standard Flange Size	NW16	NW16	NW25	NW25	NW40	NW40	NW16	NW25	NW16	NW50
Intake Thread Size	3/4-20	1-20	1-20	1-20	1 3/4 -20	1 3/4 -20	3/4-20	1-20	3/4-20	-
Exhaust Thread Size	3/4-20	1-20	1-20	1-20	1 3/4 -20	1 3/4 -20	3/4-20	1-20	3/4-20	NW16
Electrical Requirements 1	115V	115V/230V	115V/230V	115V/230V	115V/230V	115V/230V	115V	115V	115V	115V
Electrical Requirements 2	60Hz	60Hz	60Hz	60Hz	60Hz	60Hz	60Hz	60Hz	60Hz	60Hz
Electrical Requirements 3	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase	1 Phase

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


Section 6: PUMPING SPEED CURVES


HIGH EFFICIENCY PUMPING CHARACTERISTIC OF DUOSEAL® PUMPS


The low RPM design of DuoSeal vacuum pumps makes it possible to maintain high pumping efficiencies into the low micron range. High RPM direct drive vacuum pumps, on the other hand, experience a rapid fall off of pumping speed below 100 microns.

A comparison of the DuoSeal pumping speed curves shown below, versus high RPM pump speed curves conclusively illustrates the superior performance of the DuoSeal pump.

What does that mean for the vacuum pump user? DuoSeal pumps provide considerable faster pumpdown in the low micron range.

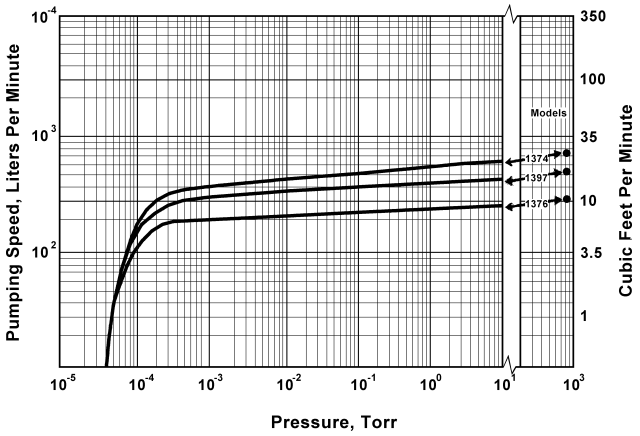

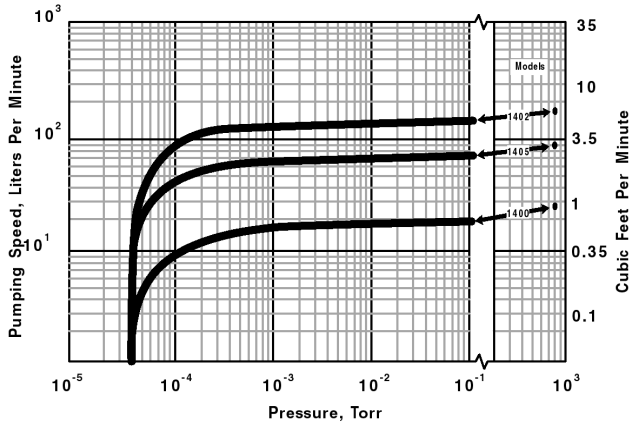

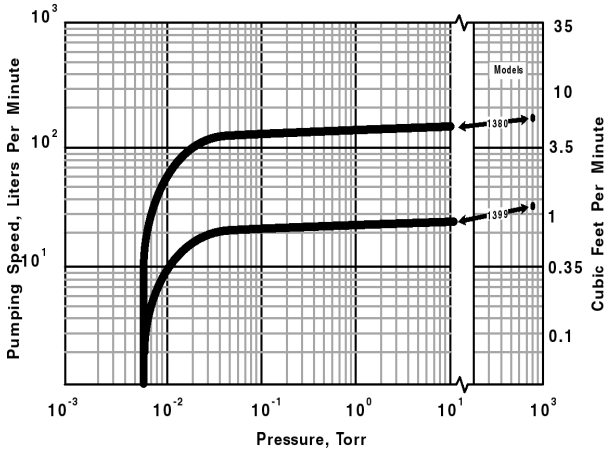


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


Section 7: DIMENSIONAL DRAWINGS

Mounted Pumps


1374B-01

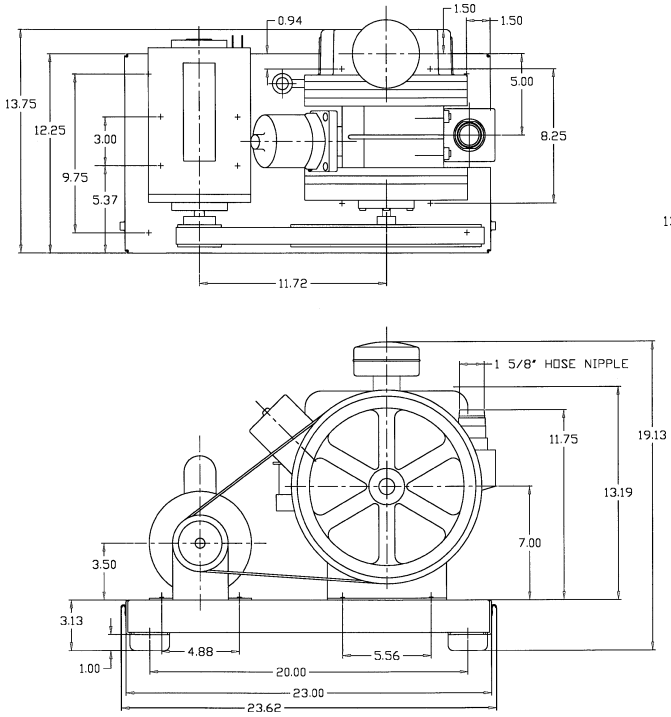


1376B-01

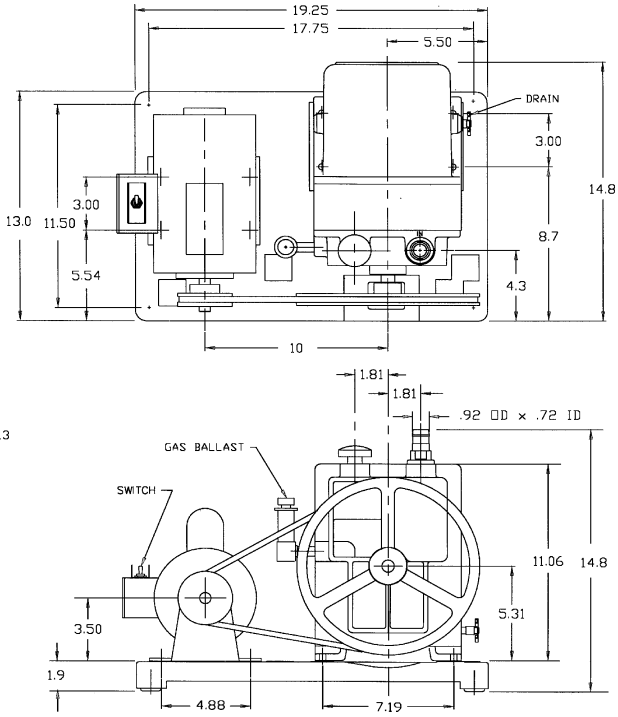

1380B-01

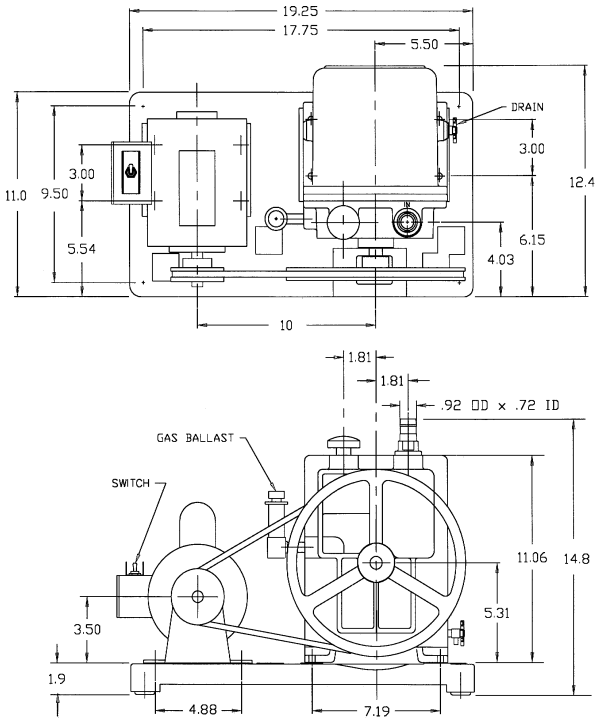


1392B-01

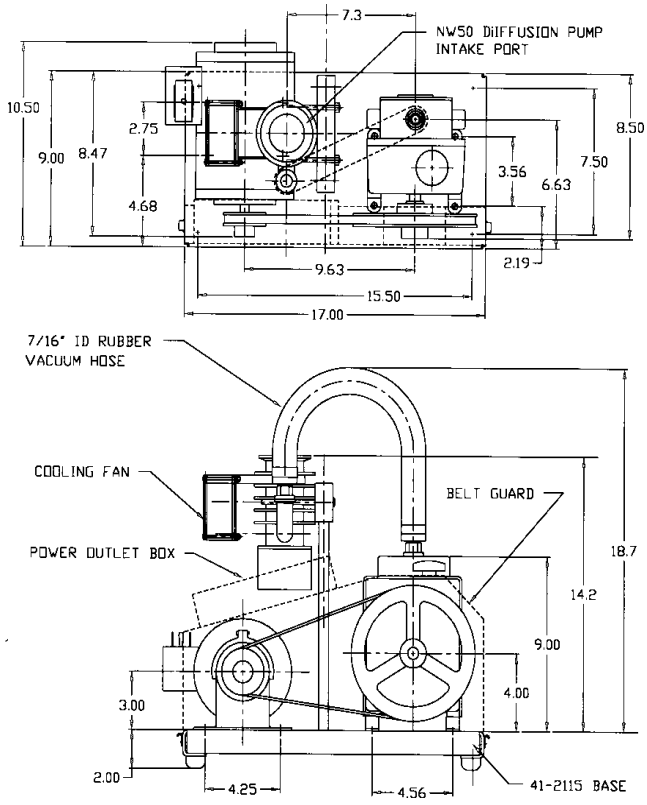


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


1397B-01

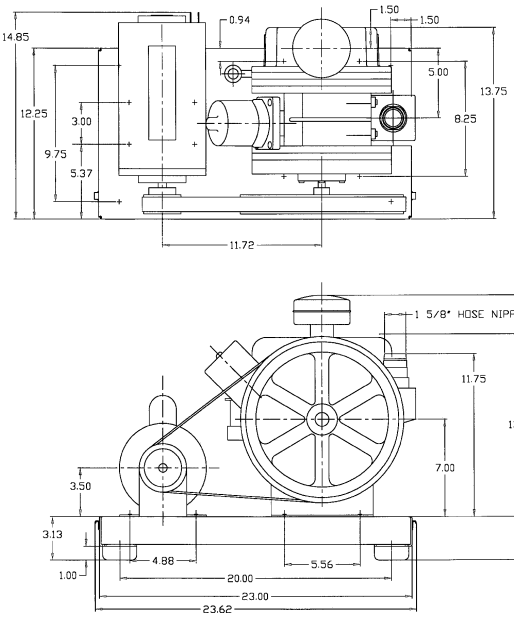


1399B-01

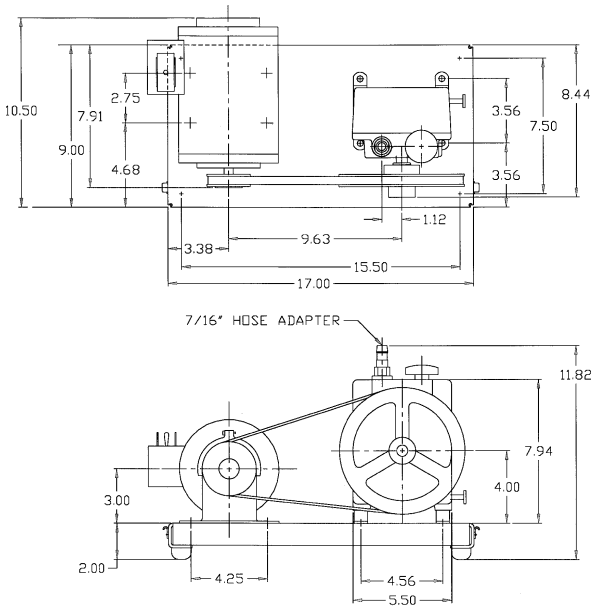


1400B-01

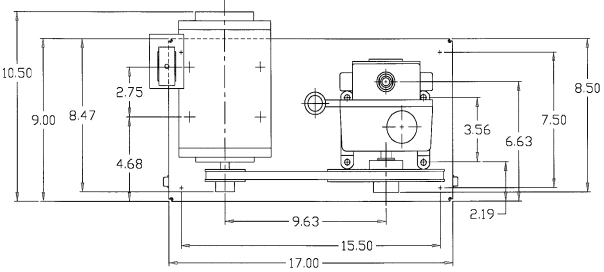

1402B-01

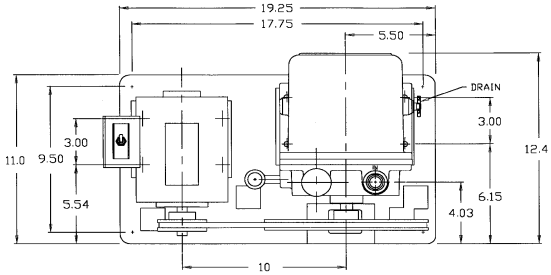

1405B-01

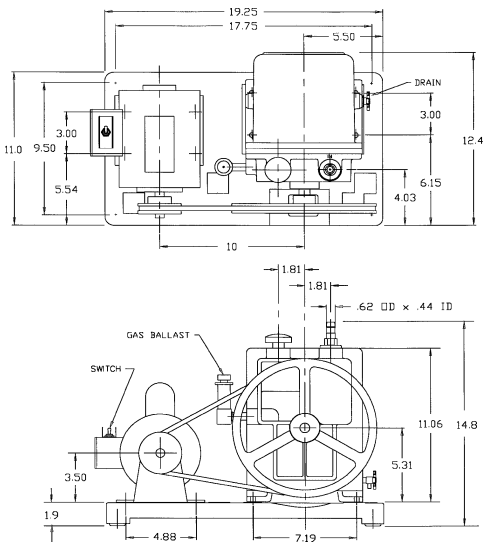

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


PARTS LIST
MODEL 1376 VACUUM PUMP

ITEM	QUAN	PC. NO.	A*	B*	DESCRIPTION
01	1	41-2074			Pulley incl. item 2
02	1	2-01-9306			Socket Head Set Screw 5/16-18x3/8
03	1	61-8378A			Shaft Lip Seal Assy
04	1	41-0643			Seal Gasket
05	3	2-00-2706			Fil. Head Screw 10-32x3/8
06	1	61-8685			Large End Plate
07	1	61-8671			Intake Cover
08	1	41-0234	1	1	Intake Cover Gasket
09	4	2-01-6112			Socket Head Screw 1/4-20x3/4
10	1	41-1692	1	1	Intake Screen
11	1	61-8329A			Intake Screen Bracket
12	6	41-2363			Washer
13	3	2-01-6316			Socket Head Screw 5/16-18x1
14	1	41-0612			Dust Cap
15	1	61-8468A			Adapter Hose
16	1	61-9284A	1	1	O-Ring, Fluoroelast. #214
17	2	41-0403	2	2	Oil Case Gasket
18	10	62-0015			Socket Head Screw 1/4-20x3-1/8
19	1	41-0508	1	1	Thrust Washer
20	1	41-1683			Shaft
21	1	41-0624	1	1	Pulley Key
22	1	41-1688			Large Ring
23	1	41-1682			Large Rotor
24	2	41-1684	2		Large Vane
25	3	41-0696	3	3	Spring Holder
26	6	41-1304	6	6	Vane Spring
27	1	61-8250C			Center Plate
28	2	2-36-9900			Lock Nut 1/4-20
29	2	2-01-6320			Socket Head Screw 5/16-18x1-1/4
30	1	41-1518			Exhaust Rotor
31	2	41-1487	2		Small Vane
32	2	41-0613	2	2	Woodruff Key
33	1	61-8108D			Exhaust Ring
34	1	41-2045			Small End Plate
35	1	41-0672			End Cap
36	1	61-8310A			Pressure Release Tube Part 2
37	1	4-40-1200	1	1	Steel Ball
38	1	41-0992	1	1	Pressure Release Valve
39	1	61-8309B			Pressure Release Tube Part 1
40	3	62-0259			Self Tap Screw #6x1/4 Type 'A'
41	6	2-01-0356			Hex Head Screw 5/16-18x3-1/2
42	1	41-2154	1	1	Lower Exhaust Valve
43	7	2-61-0000			Steel Washer 1/4
44	2	2-01-8128			Set Screw 1/4-20x1-3/4
45	6	2-31-0112			Hex Nut 1/4-20x3/16x7/16
46	4	61-8314A	4	4	Block Valve
47	4	41-2158	4	4	Valve Spring
48	2	41-2156			Valve Block
49	2	41-2159	2	2	Valve Block Spacer
50	1	61-8102A			Flange for Gas Ballast Valve
51	3	2-00-6712			Socket Head Screw 10-32x3/4
52	1	41-1686			Baffle Plate
53	2	41-2175			Set Screw (Special)
54	1	61-8687			Oil Case Spacer
55	1	61-8101			Oil Case Assy incl.items 63,67,72,74,75
56	2	61-8362			Socket Head Screw 8-32x1
57	1	61-8100A			Gas Ballast Valve Adapter
58	1	61-8103	1	1	O-Ring, Fluoroelast. #905-75
59	1	41-1736			Gas Ballast Valve incl. item 64
60	6	61-8456A			Washer Copper 5/16
61	2	41-1056			Washer
62	10	2-63-0193			Split Lock Washer 1/4
63	2	61-2229A			Glass Window
64	1	41-1742	1	1	O-Ring, Fluoroelast. #011
65	3	41-1766			Thrust Disk
66	1	61-8101A			Gas Ballast Connector
67	1	61-2278			Drain Valve
68	1	41-2757			Lock Nut
70	1	2-01-5106			Round Head Screw 1/4-20x3/8
71	1	41-2153	1	1	Upper Exhaust Valve
72	2	61-2370			Cover, Window
73	3	2-01-0316			Hex Head Screw 5/16-18x1
74	2	61-2231	2	2	O-Ring, Fluoroelast. #127
75	4	62-1037			Flat Hd Screw 8-32x7/16
		1-99-5282	1	1	Sealant, 6 cc
		1401D	1	1	Mechanical Seal option**

A*: MAJOR REPAIR KIT CAT # 1376K-06
B*: MINOR REPAIR KIT CAT # 1376K-05
** 1401D is mechanical seal kit includes mechanical seal, gasket and 3 screws
*** MANOR REPAIR KIT CAT # 1376K-11 FOR REFRIGEATION SERVICE PUMPS USES NEOPRENE SEALS

67-0894
Sheet 1 of 2
Rev. 6 01/04

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

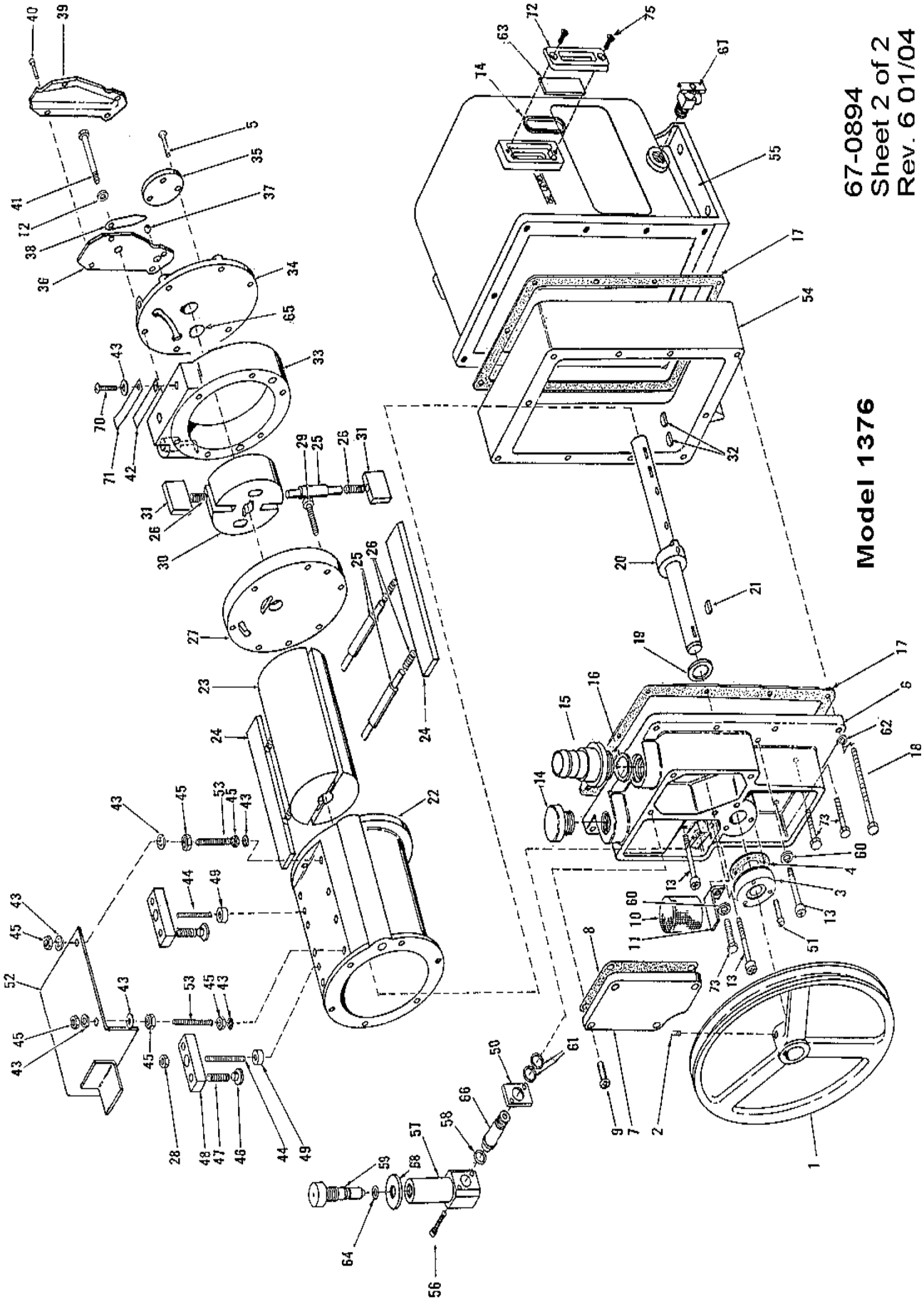

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


PARTS LIST
MODEL 1380 VACUUM PUMP

ITEM	QUAN	PC. NO.	A*	B*	DESCRIPTION
1	1	41-2074			Pulley incl. item 2
2	1	2-01-9306			Soc, Hd. Setscrew 5/16-18x3/8
3	1	61-8671			Intake Cover Chamber
4	4	2-01-6112			Soc Hd Screw 1/4-20x3/4
5	1	41-0234	1	1	Intake Cover Gasket
6	1	41-0660	1	1	Air Filter
7	6	2-01-0316			Hex Hd Screw 5/16-18x1
8	12	41-2363			Steel Washer
9	1	41-0663			Air Filter Bracket
10	3	2-00-2706			Fil Hd Screw 10-32x3/8
11	1	61-8378A			Shaft Lip Seal Assy.
12	1	41-0643			Seal Gasket
13	10	2-01-0114			Hex Hd. Screw 1/4-20x7/8
14	10	2-63-0193			Split Lock Washer 1/4
15	1	61-8505D			Large End Plate
16	1	41-0403	1	1	Oil Case Gasket
18	1	61-8468A			Intake Nipple
19	1	61-9284A	1	1	O-ring #214
20	1	41-0508	1	1	Thrust Washer
21	1	4-06-0754			Retaining Ring
22	2	41-0624	2	2	Woodruff Key
23	1	41-1773			Shaft
24	1	41-1766	1	1	Thrust Disk
25	1	61-8101			Oil Case incl. items 56, 59, 60, 61, 62
26	1	41-1736	1	1	Gas Ballast Valve incl. item 27
27	1	41-1742			O-ring #011
28	2	61-8362			Soc Hd Screw 8-32x1
29	2	41-1056			Aluminum Washer
30	3	2-00-6712			Soc Hd Screw 10-32x3/4
31	1	61-8102A			Flange For Gas Ballast
34	1	41-1772			Ring
35	3	41-1775	3	3	Lower Exhaust Valve
36	3	41-0695	3	3	Upper Exhaust Valve
37	4	2-61-0000			Steel Washer 3/16
38	2	2-01-5106			Rd Hd Screw 1/4-20x3/8
39	1	61-8507B			Buffle
40	3	2-31-0112			Hex Nut 1/4-20
41	1	41-2175			Special Set Screw 1/4-20x2-1/4
42	2	41-0959	2		Vane
43	4	41-1304	4	4	Vane Spring
44	1	41-1771			Rotor
45	2	41-1685	2	2	Vane Spring Holder
46	1	41-2045			Small End Plate
47	1	41-1005			Pressure Release Tube No.2
48	1	4-40-1200	1	1	Steel Ball 3/8
49	1	41-0992	1	1	Pressure Release Valve
50	6	2-01-0318			Hex Hd Screw 5/16-18x1-1/8
51	3	62-0259			Bind Hd Sheet Met Screw #6x1/4
52	1	41-0986			Pressure Release Tube No.1
53	1	41-0672			End Cap
54	1	41-2757			Lock Nut
55	1	61-8103	1	1	O-ring #905
56	2	41-3753			Oil Level Sight
57	1	61-8100A			Gas Ballast Valve Adapter
58	1	61-8101A			Gas Ballast Connector
59	1	61-2278			Drain Valve
60	2	61-2231	1	1	O-Ring Tetr. #127
61	2	61-2228A			Cover, Oil Window
62	4	61-2232			Flat Hd Screw #8-32x1/2
		1-99-5282	1	1	Sealant, 6 cc
		1401D	1	1	Mechanical Seal Option**

A*: MAJOR REPAIR KIT CAT # 1380K-04
B*: MINOR REPAIR KIT CAT # 1380K-03 .

** 1401D is mechanical seal kit includes mechanical seal, gasket and 3 screws

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

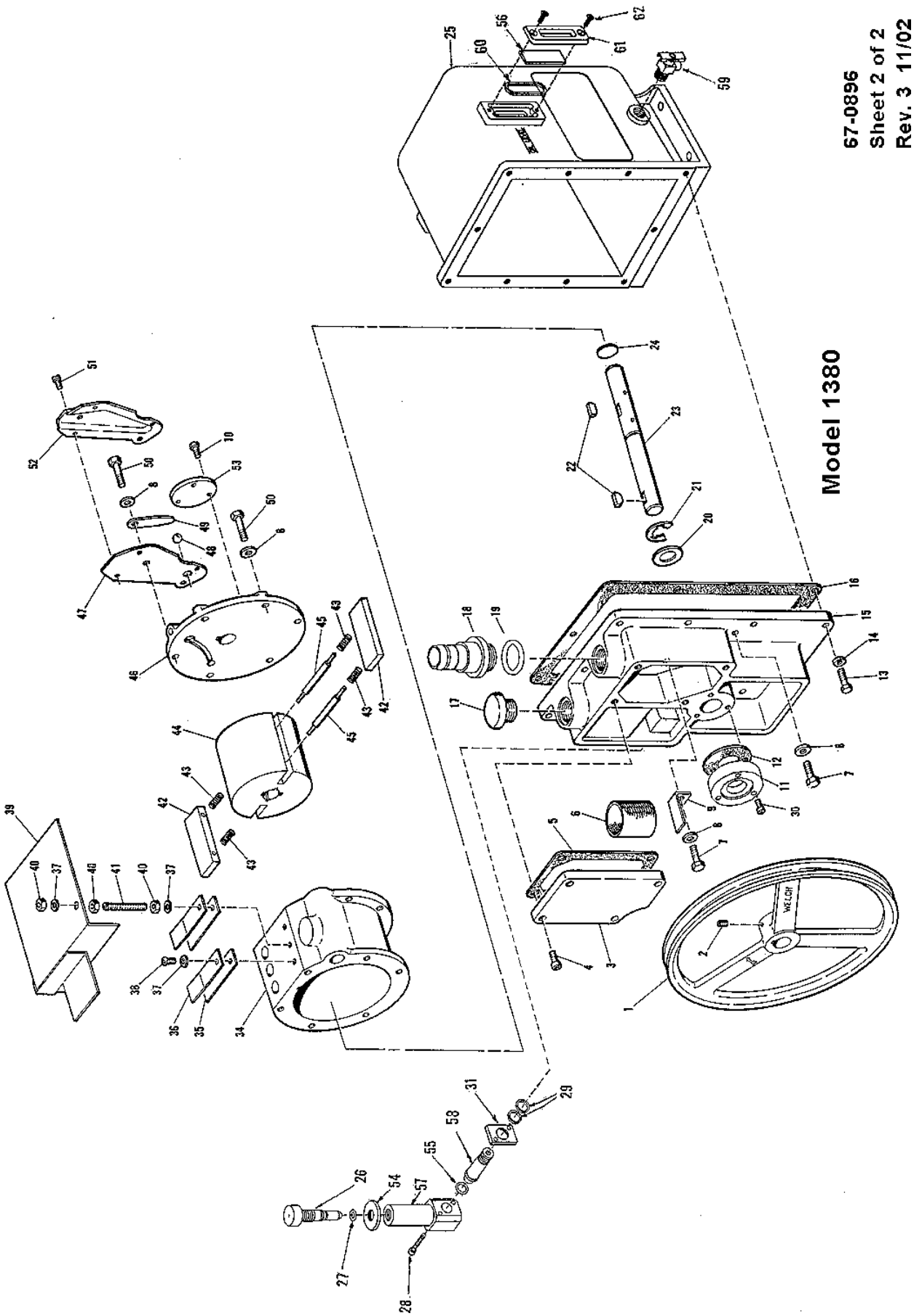

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PARTS LIST MODEL 1397 / 1374 VACUUM PUMP


ITEM	QUAN	PC. NO.‘	A*	B*	DESCRIPTION
1	1	2-01-9306			Set Screw 5/16-18x3/8
2	1	41-1492			Pulley incl. Item 1
3	4	2-01-6116			Soc Hd Screw 1/4-20x1
4	1	61-5074			Shaft Lip Seal Assy.
5	1	41-1494			Seal Gasket
6	2	41-1490	2	2	Woodruff Key
7	1	41-1484			Shaft & Collar
8	2	41-1486			Leg
9	7	2-01-6324			Soc Hd. Screw 5/16-18x1-1/2
10	38	41-2363			Steel Washer
11	3	41-1137			Plug Washer
12	2	41-1136			Plug
13	1	61-8770			Large End Plate
14	1	41-0055			Spacer
15	2	2-00-2704			Fil Hd Screw 10-32x1/4
20	4	2-01-6316			Soc Hd Screw5/16-18x1
21	2	2-01-6320			Soc Hd Screw5/16-18x1-1/4
23	1	61-8880			Cover Plate for Conversion Valve
24	1	41-1508	1	1	Valve Gasket
25	1	41-1736			Gas Balast Valve incl. Item 45
26	1	41-2258			Hardened Sleeve 1/4 Size
27	1	61-8538A	1	1	Valve Backer
28	8	2-61-0000			Steel Washer ¼-20 Bolt Size
29	4	2-31-0112			Hex Steel Nut 1/4-20x7/16
30	1	41-2253			Discharge Valve Cover
31	1	41-1465			Center Plate
32	1	41-2161	1	1	Exhaust Valve
33	2	2-01-0120			Hex Hd Screw 1/4-20x1-1/4
34	1	41-2164			Cover
35	1	2-01-8132			Headless Screw 1/4-20x2
36	1	61-8310A			Pressure Release Tube Part 2
37	1	41-1461			Deflector
38	1	61-8309B			Pressure Release Tube Part 1
39	3	62-0259			Bd Hd Screw #6x1/4
40	1	41-1496	1	1	Oil Case Gasket
41	1	41-1501			Dust Cap
42	1	61-8539A	1	1	Valve Flapper
43	18	2-01-6316			Soc Hd Screw 5/16-18x1
45	1	41-1742	1	1	O-ring #011 Fluoroelastomer
46	2	41-3753			Oil Sight
47	2	61-8362			Soc Hd Screw 8-32x1

A*: MAJOR REPAIR KIT CAT # 1397K-08
B*: MINOR REPAIR KIT CAT # 1397K-07

ITEM	QUAN	PC. NO.	A*	B*	DESCRIPTION
48	1	61-8101A			Connector for GB Valve
49	5	62-0222			Soc Hd Screw 5/16-18x2-3/4
50	1	41-1482			Oil Case Assy. incl Item 54 & (2) Items 46
51	4	62-0213			Soc Hd Screw 5/16-18x1-3/4
52	1	61-8102A			Flange For GB Valve
53	1	61-8103	1	1	O-Ring #905 Fluoroelastomer
54	1	61-2278			Drain Valve
55	1	2-01-0344			Hex Hd Screw 5/16-18x2-3/4
56	3	2-00-2706			Flat Hd Screw 10-32x3/8
57	1	41-0672			Shaft End Cap
58	1	41-0992	1	1	Pressure Relief Valve
59	1	4-40-1200	1	1	Steel Ball 3/8 Dia
60	1	41-1500			Small End Plate
61	1	41-1497			Exhaust Ring
62	1	41-1499			Exhaust Rotor
63	1	41-0696	1	1	Small Spring Holder
64	2	41-1304	2	2	Small Vane Spring
65	2	41-1487	2		Small Vane Metal
66	2	2-01-0316			Hex Hd Screw 5/16-18x1
67	1	61-8469A			Intake Nipple
68	1	61-8475	1	1	O-Ring #224 Fluoroelastomer
69	4	2-01-6314			Soc Hd Screw 5/16-18x7/8
70	1	61-8693			Intake Chamber
71	1	41-1495	1	1	Intake Chamber Gasket
72	1	41-0937	1	1	Screen Filter
73	1	61-8691			Large Ring
74	1	41-1478			Large Rotor
75	2	41-1489	2	2	Large Spring Holder
76	4	41-1488	4	4	Large Vane Spring
77	2	41-1477	2		Large Vane Metal
81	1	41-1516	1	1	Key for Small Rotor
82	4	41-1056	4	4	Aluminum Washer
83	1	41-1469	1	1	Spacer Disk
84	1	2-63-0193			Split Lock Washer 1/4
85	1	41-2256			Coil Spring
86	1	41-2257			Hardened Sleeve 5/16 Size
87	1	41-2553	1	1	Spacer
89	1	61-8100A			GB Valve Adapter
90	1	2-69-7007			SS Washer 11/32x11/16x1/16
91	1	41-2259			Coil Spring
92	1	41-2757			Lock Nut
93	1	41-2165			Guard Stud
94	2	61-8712	1	1	O-Ring #372 Fluoroelastomer
		1-99-5282	1	1	Sealant, 6 cc
		1401F	1	1	Mechanical Seal option**

** 1401F is mechanical seal includes mechanical seal, gasket and 4 screws

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

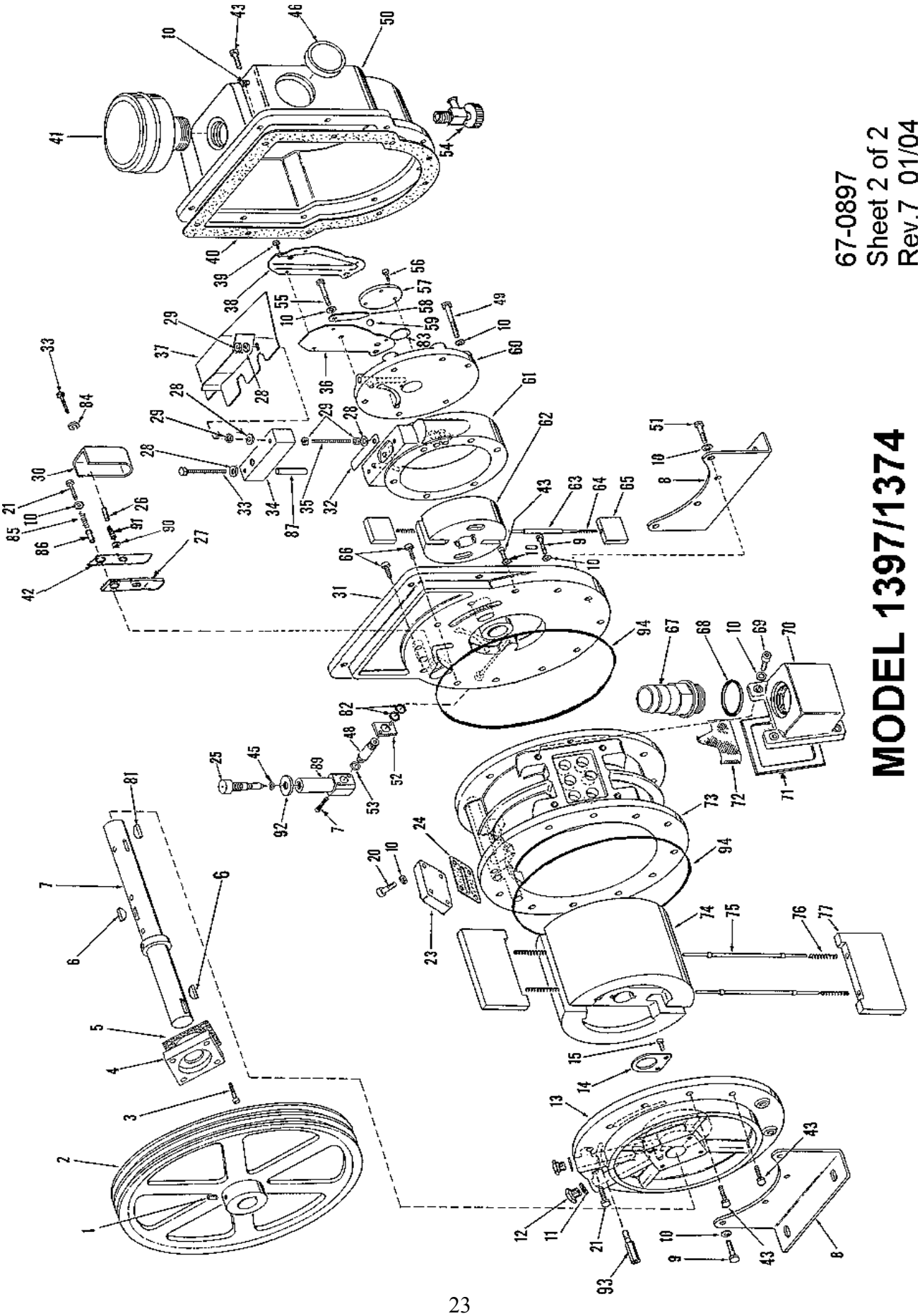

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PARTS LIST
MODEL 1399 VACUUM PUMP

ITEM	QUAN	A*	B*	PC. NO.	DESCRIPTION
1	4			2-01-0308	Hex Hd Screw 5/16 - 18x1/2
2	2			41-1049	Leg
9	1			61-2278	Drain Valve
11	1			41-3753	Oil Sight
12	1	1	1	41-1766	Thrust Disk
13	1			41-2378	Oil Case Assy. incl items 1,2,9,11
14	1	1	1	41-1052	Oil Case Gasket
15	4	1	1	2-00-2706	Fil Hd Screw 10-32x3/8
16	1			41-0672	Shaft End Cap
17	1	1	1	41-1285	Steel Washer
18	4			2-01-0332	Hex Hd Screw 5/16 - 18x2
19	4			41-2363	Washer
20	1			41-2373	End Plate
21	1			41-1042	Inner Ring
22	1	1	1	41-1239	Exhaust Valve
23	2	1	1	2-71-0040	Brass Washer #10
24	1	1	1	41-2370	Valve Cover Stud
25	1	1	1	2-65-0000	Lock Washer #10 Ext.
26	1			41-2368	Valve Cover
27	2	2		41-1059	Vane
28	2	2	2	41-1304	Spring
29	1			41-1040	Rotor
30	1	1	1	41-1063	Spring Holder
31	1	1	1	41-0974	Retaining Ring
32	1	1	1	41-0613	Woodruff Key
33	1			41-2374	Shaft
34	1	1	1	41-0624	Pulley Key
35	1			41-1345	Dust Cap
37	1			61-8467A	Intake Nipple
38	7			2-01-6112	Soc Hd Screw 1/4-20x3/4
39	1			41-2367	Intake Chamber
40	1	1	1	41-2375	Cover Gasket
41	1			41-2369	Intake Chamber Cover
42	6			2-00-0716	Hex Hd Screw 10-32x1
43	1			41-0643	Seal Gasket
44	3			2-00-6712	Soc Hd Screw 10-32x3/4
45	1			41-0306	Air Filter
46	1			41-2371	Air Filter Bracket
47	1	1	1	1401E	Shaft Seal incl. item 43 & (3) 44 **
48	1			2-01-9306	Soc Hd Set Screw 5/16-18x3/8
49	1			41-2192	Pulley incl. item 48
50	1			61-2159A	O-Ring, #116, Fluoroseal
		1	1	1-99-5282	Sealant, 6 cc

A*: MAJOR REPAIR KIT #1399K-04
B*: MINOR REPAIR KIT #1399K-03

** 1401E is mechanical seal which fits all model 1399 serial numbers. There is a lipseal seal on S/N 25773 and above. It is available as a kit. part no. 1401E-01

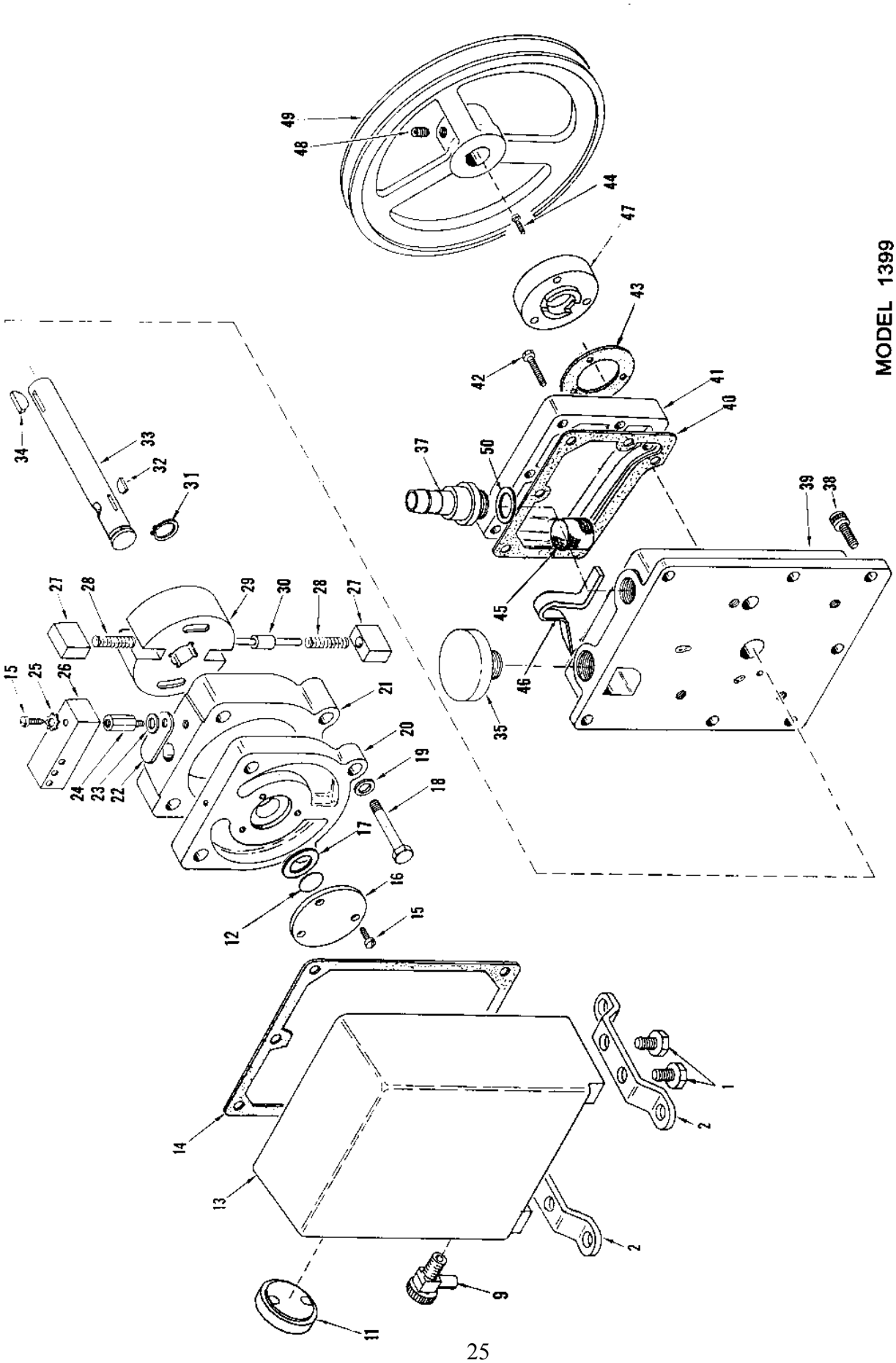

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PARTS LIST
MODEL 1400 VACUUM PUMP

ITEM	QUAN	PC. NO.	A*	B*	DESCRIPTION
1	1	2-01-9306			Soc Hd Set Screw 5/16 - 18x3/8
2	1	41-2191			Pulley include. Set Screw , Item 1
3	1	41-3696			Shaft Seal
4	1	41-0643			Shaft Seal Gasket
5	1	41-0624	1	1	Pulley Key
6	1	41-1060			Shaft
7	1	41-1150	1	1	Retaining Ring
8	2	41-0613	2	2	Woodruff Key
10	2	41-1049			Leg
11	6	41-1138			Washer, Copper, 7/16 OD
12	1	41-3753			Oil Sight
13	1	41-1742	1	1	O-Ring #011
14	1	61-8540C			Oil Case*
15	1	41-1345			Dust Cap
16	1	41-1052	1	1	Oil Case Gasket
17	9	2-01-0332			Hex Hd Screw 5/16 - 18x2
18	9	41-2363			Washer
19	1	41-1045			Exhaust End Plate
20	1	41-1040			Exhaust Rotor
21	1	41-1239	1	1	Exhaust Valve
22	2	2-71-0040	2	2	Brass Washer #10
23	8	41-1056			Aluminum Washer
24	1	41-1039			Vented Exhaust Plug
25	1	61-8679			Center Plate
26	1	41-1285			Steel Washer
27	1	41-0890	1	1	Air Filter
28	1	41-0383	1	1	Gasket
29	1	41-1048			Intake Chamber**
30	1	61-8467A			Nipple
31	2	2-01-0524			Hex Hd Screw 3/8 - 16x1-1/2
32	1	61-8702			Intake End Plate
33	1	61-8673			End Cap
34	5	2-00-2706			Hex Hd Screw 10 - 32x3/8
35	2	2-31-2521			Hex Nut 3/8x16
36	1	61-8675			Intake Ring
37	1	41-1041			Intake Rotor
38	2	41-1063	2	2	Spring Holder
39	4	41-1304	4	4	Vane Spring
40	4	41-1059	4		Vane (Metal)
41	5	2-01-0112			Hex Hd Screw 1/4 - 20x3/4
42	1	41-1042			Exhaust Ring
43	3	2-00-6712			Soc Hd Screw 10-32x3/4
44	1	41-0578	1	1	Seal
45	1	61-8102A			Flange for G.B. Valve
46	1	61-8103	1	1	O-Ring, Fluoroseal
47	2	61-8362			Soc Hd Screw 8 - 32x1
49	1	61-2278			Drain Valve
50	1	2-01-0120			Hex Hd Screw 1/4 - 20x1-1/4
51	1	41-1736			Vented Exhaust Valve
52	4	2-01-0308			Hex Hd Screw 5/16 - 18x1/2
53	2	2-61-0571			Steel Washer .071 Thick x 3/8 ID
54	1	61-2159A			O-Ring, #116, Fluoroseal
55	1	41-1069	1	1	Baffle Stud
56	1	61-8347B			Exhaust Baffle
57	1	61-8100A			Valve Adapter
58	1	61-8101A			Connector
59	1	41-2757			Locknut
60	1	61-8665	1	1	O-Ring, M2.5x86, Fluoroseal
61	1	2-00-2704			Fil Hd Screw 10-32x1/4
62	1	61-9378	1	1	O-Ring, #152, Fluoroseal
63	1	66-0183	1	1	O-Ring, #120, Fluoroseal
		1-99-5282	1	1	Sealant 6 cc
		1401E	1	1	Mechanical Seal Option***

* Oil Case Assembly PT #61-8540C includes Legs item 10, Oil Sight, item 12, oil case, Item 14, Screws, Item 52 & Drain Valve, item 49

** Intake Chamber Assembly PT #41-1048 includes intake chamber, item 29, nipple, Item 30 & O- ring, item 54.


A*: Major Repair Kit #1400K-04

B*: Minor Repair Kit #1400K-03

*** 1401E is mechanical seal kit, it includes mechanical seal, gasket, and 3 screws

Rev. 13 8/01

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

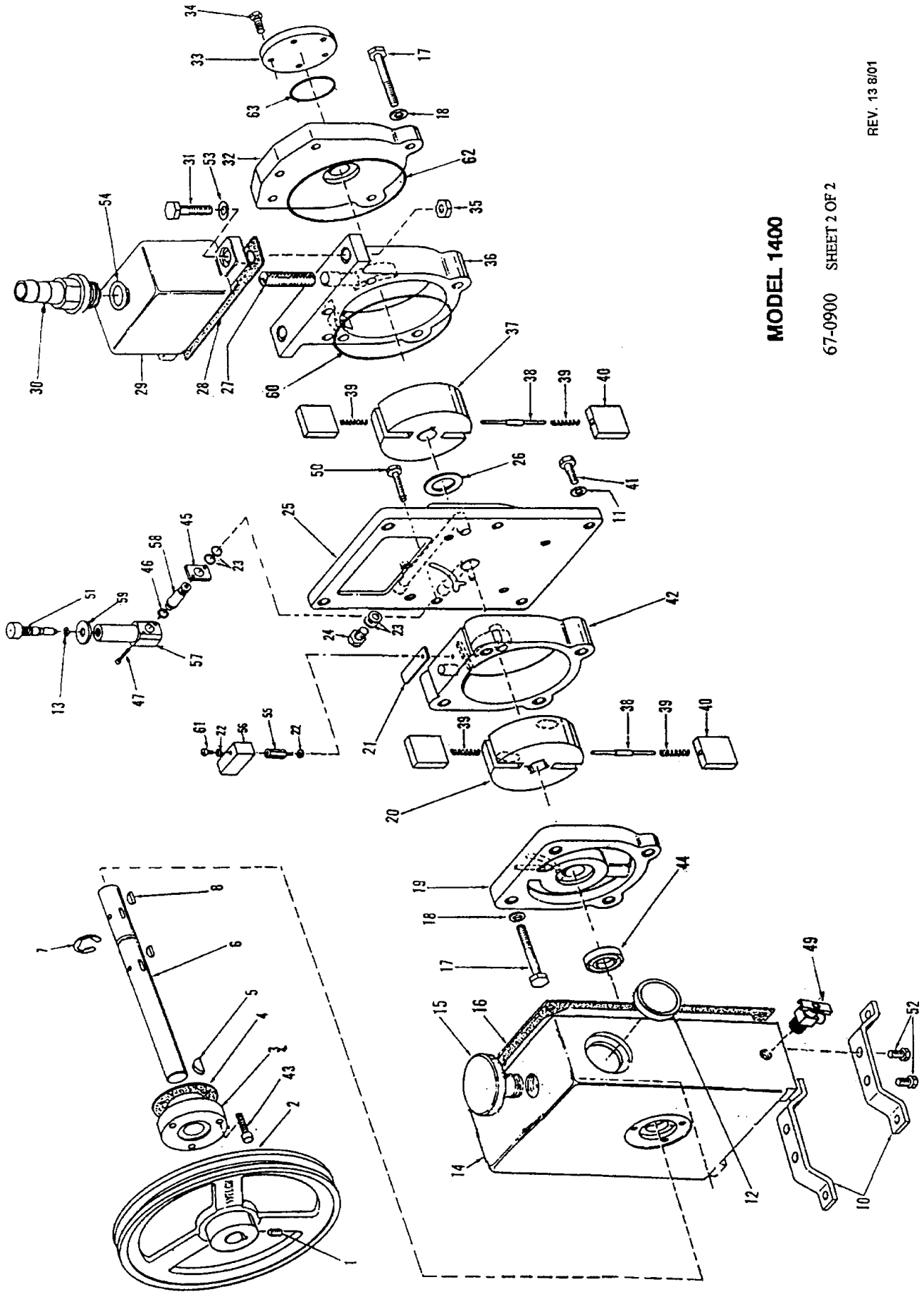

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com


PARTS LIST
MODEL 1402 VACUUM PUMP

ITEM	QTY	PC. NO.	A	B	DESCRIPTION
01	1	61-8113C			Intake Ring
02	3	2-61-0000			Steel Washer 3/16 size
03	3	2-31-0112			Hex Steel Nut
04	1	41-2159	1	1	Valve Block Spacer
05	2	61-8314A	2	2	Valve
06	2	41-2158	2	2	Valve Spring
07	1	2-01-8128			Headless Steel Setscrew
08	1	41-2156			Valve Block
09	1	41-1742	1	1	O-Ring, Fluoroelast. #011
10	1	41-2160			Baffle Plate
11	1	2-36-9900			Nut, Lock 1/4-20
12	3	2-00-2706			Fil. Head Screw 10-32x3/8
14	1	41-2175			Headless Steel Setscrew
15	1	41-2092			Intake Rotor
17	1	61-8250C			Center Plate
19	2	2-01-6320			Socket Head Screw 5/16-18x1-1/4
20	1	41-2154	1	1	Lower Exhaust Valve
21	1	61-8310A			Pressure Release Tube, Part 2
22	1	41-0992			Pressure Release Valve
23	6	41-2363			Washer
24	6	2-01-0356			Hex Head Screw 5/16-18x3-1/2
25	3	62-0259			Self Tap Screw #6x1/4 Type "A"
26	1	61-8309B			Pressure Release Tube Part 1
27	1	41-0672			End Cap
28	1	4-40-1200	1	1	Steel Ball 3/8
29	1	41-2045			Small End Plate
30	1	61-8108D			Exhaust Ring
31	1	41-1518			Exhaust Rotor
32	6	41-1304	6	6	Vane Spring
33	2	41-1487	2		Small Vane
35	2	41-0959	2		Large Vane
36	1	41-2074			Pulley inc. setscrew, Item 37
37	1	2-01-9306			Socket Head Set Screw
38	1	61-8671			Large End Plate Cover
39	1	41-2153	1	1	Upper Exhaust Valve
40	1	2-01-5106			Rd. Head Steel Screw 1/4-20x 3/8
41	4	2-01-6112			Socket Screw 1/4-20x3/4
42	1	41-0234	1	1	Intake Chamber Cover Gasket
43	1	61-8329A			Air Filter Bracket
44	6	2-01-0316			Hex Head Cap Screw 5/16-18x1
45	1	41-0612			Dust Cap
46	1	61-8468A			Intake Nipple
47	1	61-9284A	1	1	O-Ring, Fluoroelast. #214
48	1	41-0403	1	1	Oil Case Gasket
49	2	61-2231	2	2	O-Ring, Fluoroelast. #127
50	1	61-8102A			Flange For G.B. Valve
51	1	41-1736			Gas Ballast Valve
52	2	61-8362			Socket Head Screw 8-32x1
53	1	61-8101			Oil Case Assy, incl. 49,56,59,76 ,77
54	1	41-2757			Lock Nut
55	1	61-8103	1	1	O-Ring, Fluoroelast. #905-75
56	2	61-2229A			Window Glass
57	1	61-8100A			Gas Ballast Valve Adapter
58	1	61-8101A			Gas Ballast Valve Connector
59	1	61-2278			Drain Valve
60	6	61-8456A	6	6	Washer Copper 5/16
61	2	41-1056			Washer
62	2	41-0613	2	2	Woodruff Key
63	1	41-2093			Shaft
64	2	41-0624	2	2	Woodruff Key
65	1	41-0508	1	1	Thrust washer
66	1	61-8668			Large End Plate
67	10	2-63-0193			Split Lock Washer
68	10	2-01-0114			Hex Head Cap Screw 1/4-20x3/4
69	1	41-0660	1	1	Air Filter
70	1	41-0643			Seal Gasket
71	1	61-8378A			Shaft Lip Seal Assy.
72	3	41-1685	3	3	Spring Holder - Intake
73	3	41-1766	3	3	Thrust Disc
74	1	4-06-0754	1	1	Truarc Ring
75	3	2-00-6712			Socket Head Screw 10-32x3/4
76	2	61-2370			Cover, Oil Window
77	4	62-1037			Fl. Head Screw 8-32x7/16
		1-99-5282	1	1	Sealant, 6 cc
		1401D	1	1	Mechanical Seal option**

A*: MAJOR REPAIR KIT CAT # 1402K-06
B*: MINOR REPAIR KIT CAT # 1402K-05
** 1401D is mechanical seal kit includes mechanical seal, gasket and 3 screws

67-0905
Sheet 1 of 2
Rev. 6 11/02

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

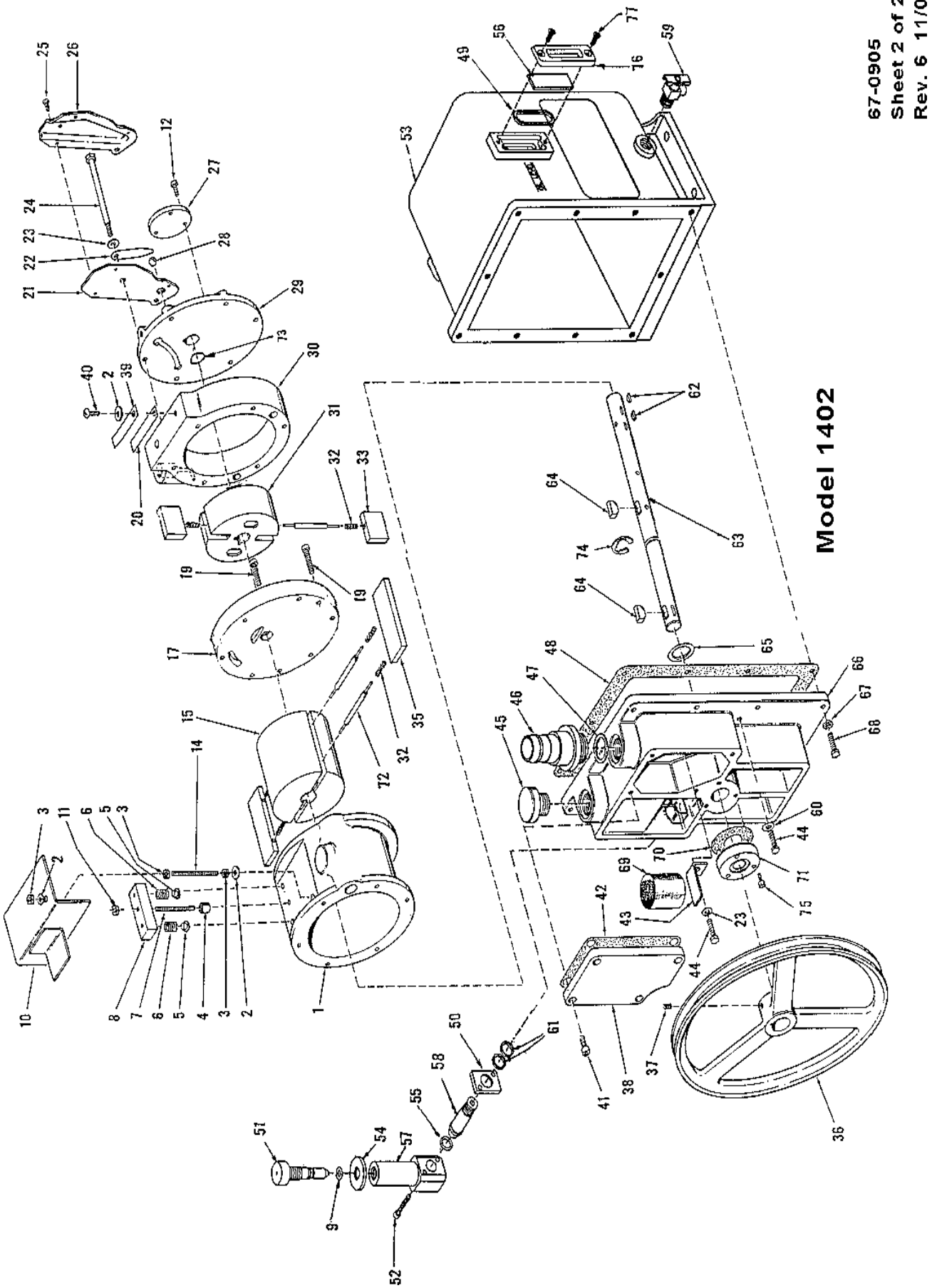

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

PARTS LIST
MODEL 1405 VACUUM PUMP


ITEM	QUAN	PC. NO.	A*	B*	DESCRIPTION
1	1	61-8651			Intake Ring
2	3	2-61-0000			Steel Washer 3/16 size
3	3	2-31-0112			Hex. Steel Nut
4	1	41-2159	1	1	Valve Block Spacer
5	2	41-2157	2	2	Valve
6	2	41-2158	2	2	Valve Spring
7	1	2-01-8128			Headless Steel Set Screw
8	1	41-2156			Valve Block
9	1	41-1742	1	1	O-Ring #011
10	1	41-2160			Baffle Plate
11	1	2-36-9900			Nut, Lock 1/4-20
12	3	2-00-2706			Fil. Hd. Screw 10-32x3/8
14	1	41-2175			Headless Steel Set Screw
15	1	41-2092			Intake Rotor
17	1	61-8250C			Center Plate
19	2	2-01-6320			Socket Hd. Screw 5/16-18x1-1/4
20	1	41-2154	1	1	Lower Exhaust Valve
21	1	41-1005			Pressure Release Tube, Part 2
22	1	41-0992	1	1	Pressure Release Valve
23	6	41-2363			Washer
24	6	2-01-0356			Hex. Hd Screw 5/16-18x3-1/2
25	3	62-0259			Self Tap Screw #6x1/4 Type "A"
26	1	41-0986			Pressure Release Tube Part 1
27	1	41-0672			End Cap
28	1	4-40-1200		1	1Steel Ball 3/8
29	1	41-2045			Small End Plate
30	1	61-8108D			Exhaust Ring
31	1	41-1518			Exhaust Rotor
32	6	41-1304	2	2	Vane Spring
33	2	41-1487	2		Small Vane
34	1	41-0696	1	1	Vane Spring Holder-Exhaust
35	2	41-0959	2		Large Vane
36	1	41-2074			Pulley inc. Set Screw, Item 37
37	1	2-01-9306			Socket Head Set Screw
38	1	61-8671			Large End Plate Cover
39	1	41-2153	1	1	Upper Exhaust Valve
40	1	2-01-5106			Rd. Hd. Steel Screw 1/4-20x 3/8
41	4	2-01-6112			Socket Screw 1/4-20x3/4
42	1	41-0234	1	1	Intake Chamber Cover Gasket
43	1	41-0663			Air Filter Bracket
44	6	2-01-0316			Hex Head Cap Screw 5/16-18x1
45	1	41-0612			Dust Cap
46	1	61-8641			Intake Nipple, (7/16 Hose I.D.)
47	1	61-9284A	1	1	O-ring, #214
48	1	41-0403	1	1	Oil Case Gasket
49	2	61-2231	1	1	O-ring Tetr. #127
50	1	61-8102A			Flange For G.B. Valve
51	1	41-1736			Gas Ballast Valve
52	2	61-8362			Soc. Hd Screw 8-32x1
53	1	61-8101			Oil Case Assy, incl. Items 49,56,59,76 ,77
54	1	41-2757			Lock Nut
55	1	61-8103	1	1	O-ring
56	2	61-2229A			Window Glass
57	1	61-8100A			Gas Ballast Valve Adapter
58	1	61-8101A			Gas Ballast Valve Connector
59	1	61-2278			Drain Valve
60	6	61-8456A	6	6	Washer Copper 5/16
61	2	41-1056			Washer
62	2	41-0613	2	2	Woodruff Key
63	1	41-2093			Shaft
64	2	41-0624	1	1	Woodruff Key
65	1	41-0508			Thrust washer
66	1	61-8668			Large End Plate
67	10	2-63-0193			Split Lock Washer
68	10	2-01-0114			Hex Head Cap Screw 1/4-20x3/4
69	1	41-0660	1	1	Air Filter
70	1	41-0643			Seal Gasket
71	1	61-8378A			Shaft Lip Seal Assy..
72	2	41-1685	2	2	Spring Holder - Intake
73	3	41-1766	3	3	Thrust Disc
74	1	4-06-0754	1	1	Truarc Ring
75	3	2-00-6712			Soc Hd Screw 10-32x3/4
76	2	61-2228A			Cover, Oil Window
77	4	61-2232			Fl Hd Screw 8-32x1/2
		1-99-5282	1	1	Sealant, 6 cc
		1401D	1	1	Mechanical Seal Kit option **

A*: MAJOR REPAIR KIT CAT # 1402K-06

B*: MINOR REPAIR KIT CAT # 1402K-05

** 1401D is mechanical seal kit includes mechanical seal, gasket and 3 screws.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

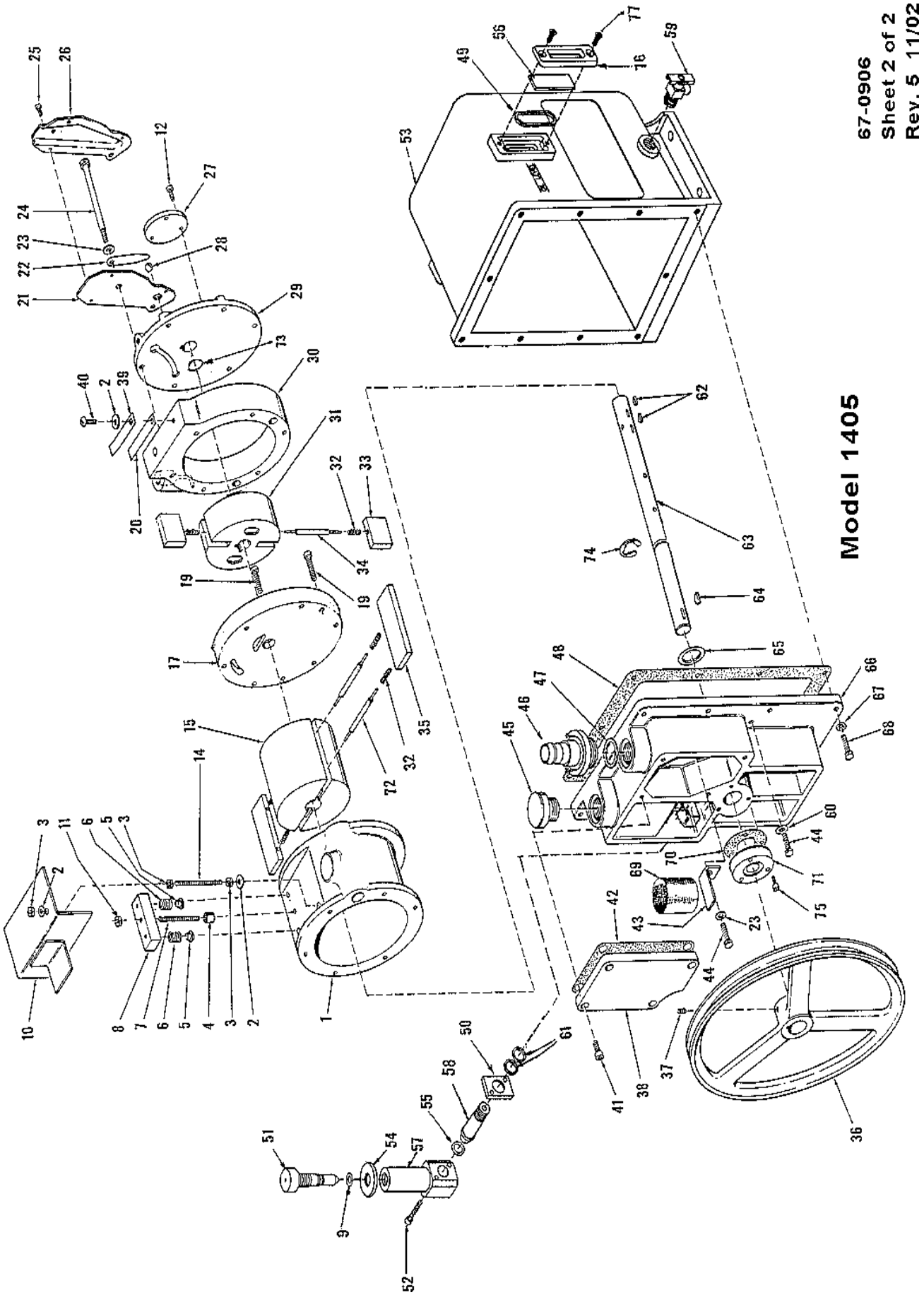

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Section 9: ACCESSORY SECTION

Exhaust Filters

- Filters oil mist from pump exhaust
- Filters to 0.3 micron particle size
- Easy replacement of element
- Reduces pump noise

Standard Exhaust Filters				
	Dimensions		Catalog Number	
Use with Pump Model	Diameter	Height	Filter	Replacement Element
1392, 1399, 1400	2.5	2.5	1417	1417L
1376,1380,1402,1405	5.0	6.25	1417P-10	1417G
1397, 1374	5.0	9.0	1417P-20	1417H-01

Directional Exhaust Filters

- Large, high capacity
- 360° Swivel outlet with tubing connector
- Screw-in type

Directional Exhaust Filters					
	Dimensions			Catalog Number	
Use with Pump Model	Diameter	Height	Outlet Tube O.D.	Filter	Element
1376, 1380, 1402, 1405	5.0	7.5	1.0	1417A	1417G
1374, 1397	6.5	10.0	2.0	1417B	1417H
1373	5.0	7.5	1.0	1417T	1417G

Exhaust Mist Eliminator

Exhaust oil recyclers capture oil mist expelled from the pump, separate air from oil, and return collected oil to the pump via a line connected to the gas ballast inlet. This process effectively reduces oil consumption and minimizes maintenance requirements by extending the time between oil changes.

Catalog No.	Duoseal Model No.
1416H-01	1400
1416C-01	1405, 1376, 1380, 1402

- Requires gas ballast connection kit: Cat. No. 141601 - for models 1400, 1405, 1376, 1380, 1402

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Coaxial Traps

- Reduces oil back streaming
- Stainless steel construction
- Easy replacement of elements

ISO Size	Capacity	Diameter	Height	A[2]	Cat. No. [1]
NW16	To 4 CFM (140 L/min)	3.75	5.25	8.25	541001
NW 25	To 12 CFM (340 L/min)	5.50	5.25	8.25	541002
Copper Element for ISO NW 16 Trap [1.]					541911
Stainless Steel Element for ISO NW 16 Trap [1.]					541921
Copper Element for ISO NW 25 and NW 40 Traps [1.]					541912
Stainless Steel Element for ISO NW 25 and NW 40 Traps [1.]					541922

1. Element must be ordered with trap.
2. Distance between ISO NW inlet and outlet connections.

Modular Sieve Traps

- Blocks backstreaming by absorbing hydrocarbons.
- Regenerates charge using built-in heater.
- Removes trace amounts of water vapor.

ISO Size	Capacity	Diameter	Height	A[2]	Cat. No. [1]
NW16	To 8 CFM (140 L/min)	4.50	6.38	9.00	543001
NW 25	To 8 CFM (226 L/min)	4.50	6.38	9.00	543002
NW40	To 17 CFM (481 L/min)	4.50	9.00	9.00	543003
Synthetic Zeolite Charge [1.]					543950
Activated Alumina Charge [1.]					543960

1. Element must be ordered with trap.
2. Distance between ISO NW inlet and outlet connections.

Dry Ice Cold Traps

- Removable 3 quarter center well for dry ice/alcohol slurry (do not use acetone)
- Visible trapping surface
- Cold temperature up to 12 hours depending on vapor load
- Available with straight tube connections

ISO Size	Diameter	Height	A [1.]	I/O Offset	Cat. No.	Repl. Seal Cat.
NW 16	10	8.25	11.06	3.5	545001	545200
NW 25	10	8.25	11.06	3.5	545002	545200
Straight Tube						
½ O.D.	10	8.25	10.75	3.5	1420H-14	545200

1. Distance between ISO NW inlet and outlet connections.

Liquid Nitrogen Cold Traps

- A highly efficient method of trapping condensible vapors
- Produces a lower mechanical vacuum pump base pressure
- Allows quick clean-up of the trapping surface with the two-piece clamped body design

ISO Size	Capacity	Diameter	Height	A [1.]	Cat. No.
NW 16	0.75 Liter	5.5	11.00	7.0	546001
NW 25	0.75 Liter	5.5	11.00	7.0	546002

1. Distance between ISO NW inlet and outlet connections.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Welch Rietschle Thomas
A Thomas Industries Company

MATERIAL SAFETY DATA SHEET

HMIS Rating: Health: 0 Flammability: 1 Reactivity: 0 Special: X

SECTION 1 - IDENTIFICATION

Chemical Name:	Highly-Refined Petroleum Lubricant Oils
Synonyms:	1407K DUOSEAL Vacuum Pump Oil
Formula:	N/A
Manufacturer's Name:	Welch Rietschle Thomas 7301 North Central Ave. Skokie, IL 60077
Date Prepared:	November 11, 2003
Telephone Number for General Information:	(847) 676-8800

SECTION 2 - COMPOSITION

Exposure Limits in Air

Component	CAS No.	OSHA PEL	ACGIH TVL
Highly-Refined Petroleum Lubricant Oils	64741-88-4	NO	NO

To the best of our knowledge, the above listed component is not hazardous according to OSHA (1910.1200) or one or more state right-to-know lists.

SECTION 3 - HAZARDS IDENTIFICATION

Major Route of Entry: Skin contact

Signs and Symptons of Acute Exposure

Inhalation: No significant adverse health effects are expected to occur upon short-term exposure to this product. Asperation of liquid into the lungs can cause severe lung damage or death.

Eye Contact: Minimal eye irritation may result from short term contact with liquid, mist, and/or vapor.

Skin Contact: This material can cause mild skin irritation from prolonged or repeated skin contact. Injection under the skin, in muscle, or into the blood stream can cause irritation, inflammation, swelling, fever, and systemic effects and mild central nervous system depression. Injection of pressurized hydrocarbons can cause severe, permanent tissue damage. Initial symptons may be minor. Injection of petroleum hydrocarbons requires immediate medical attention.

Ingestion: If swallowed, large volumes of material can cause generalized depression, headache, drowsiness, nausea, vomiting and diarrhea. Smaller doses can cause a laxative effect. If aspirated into the lungs, liquid can cause lung damage.

Chronic Health Effects Summary: Contains a petroleum-based mineral oil. Prolonged or repeated skin contact can cause mild irritation and inflammation characterized by drying, cracking, (dermatitis) or oil acne. Inhalation of petroleum-based mineral oils can cause respiratory irritation or other pulminary effects after repeted or prolonged inhalation of oil mists at concentrations above applicable workplace exposure levels.

THE INFORMATION HEREIN IS GIVEN IN GOOD FAITH, BUT NO WARRANTY, EXPRESS OR IMPLIES, IS MADE.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Conditions Aggravated by Exposure: Personnel with pre-existing disorders should avoid repeated or prolonged contact with this product.

Target Organs: Skin

Carcinogenic Potential: This product does not contain any components at concentations above 0.1% which are considered carcinogenic by OSHA, IARC, or NTP.

SECTION 4 - FIRST AID MEASURES

Take proper precautions to ensure your own health and safety before attempting rescue or providing first aid. For more specific information, refer to Exposure Controls and Personal Protection in Section 8 of this MSDS.

Inhalation: Vaporization is not expected at ambient temperatures. This material is not expected to cause inhalation-related disorders underanticipated conditions of use. In case of overexposure, move the person to fresh air.

Eye Contact: Check for and remove contact lenses. Flush eyes with cool, clean, low-pressure water while occasionally lifting and lowering eyelids. Seek medical attention if excessive tearing, readness, or pain persists.

Skin Contact: Remove contaminated shoes and clothing. Wipe off excess material. Wash exposed skin with soap and water. Seek medical attantion if tissue appears damaged or if irritation persists. Thoroughly clean contaminated clothing before reuse. Discard contaminated leather goods. If material is injected under the skin, into muscle, or into the blood stream, seek medical attention immediately.

Ingestion: Do not induce vomiting unless directed to by a physician. Do not give anything to drink unless directed to by a physician. Never give anything by mouth to a person who is not fully conscious. If large amounts are swallowed or irritation or discomfort occurs, seek medical attention immediately.

Notes to Physician: The viscosity range of the product(s) represented by this MSDS is 100 to 400 SUS at 100°F. Accordingly, upon ingestion there is a low to moderate risk of aspiration. Careful gastric lavage may be considered to evacuate large quantities of material. Subcutaneous or intramuscular injections requires prompt surgical debridement.

SECTION 5 - FIRE FIGHTING MEASURES

NFPA Flammability Classification: OSHA/NFPA Class-IIIB combustible liquid. Slightly combustible!

Flash Point/Method: CLOSED CUP: 210°C(410°F). (Pensky-Martens (ASTM D-93)) OPEN CUP: 238°(460°F) (Cleveland.).

Lower Flammable Limit: No data

Upper Flammable Limit: No data

Auto-Ignition Temp.: Not available

Hazardous Combustion Products: Carbon Dioxide, carbon monoxide, smoke fumes, unburned hydrocarbons and trace oxides of Sulfur and/or nitrogen.

Special Properties: When heated above its flash point temperature, this material will release vapors which, if exposed to an ignition source, can ignite. In enclosed spaces vapors can ignite with explosive force. Mists or sprays may burn at temperatures below the flash point.

Extinguishing Media: Use dry chemical, foam, Carbon Dioxide or water fog.

Fire Fighting Protective Clothing: Firefighters must use full bunker gear including NIOSH-approved positive pressure self-contained breathing apperatus to protect against potental hazardous combustion or decomposition products and oxygen deficiencies

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

SECTION 6 - ACCIDENTAL RELEASE MEASURES

Take proper precautions to ensure your own health and safety before attempting spill control or clean-up. For more specific information, refer to Emergency Overview on Page 1, Exposure Controls and Personal Protection in Section 8, and Disposal Considerations in Section 13 of this MSDS.

Do not touch damaged containers or spilled material unless wearing appropriate protective equipment. Slipping hazard; do not walk through spilled material. Stop leak if you can do so without risk. For small spills, absorb or cover with dry earth, sand, or other inert non-combustible absorbent material and place into waste containers for later disposal. Contain large spills to maximize product recovery or disposal. Prevent entry into waterways or sewers. In urban area, cleanup spill as soon as possible. In natural environments, seek cleanup advice from specialists to minimize physical habitat damage. This material will float on water. Absorbent pads and similar materials can be used. Comply with all laws and regulations.

SECTION 7 - HANDLING AND STORAGE

Handling: Avoid water contamination and extreme temperatures to minimize product degradation. Empty containers may contain product residues that can ignite with explosive force. Do not pressurize, cut, weld, brase solder, drill, grind or expose containers to flames, sparks, heat or other potential ignitionsources. Consult approprite federal, state and local authorities before reusing, reconditioning, reclaming, recycling or disposing of empty containers and/or waste residues of this product.

Storage: Keep container closed. Do not store with strong oxidizing agents. Do not store at temperatures above 120°F or in direct sunlight for extened periods of time. Consult appropriate federal, state, and local authorities before reusing, reconditioning, reclaiming, recycling or disposing of empty containers or waste residues of this product.

SECTION 8 - EXPOSURE CONTROLS AND PERSONAL PROTECTION

Engineering Controls: Provide exhaust ventilation or other engineering controls to keep the airborne concentrations of mists and/or vapors below the recommended exposure limits (see below). An eye wash station and safety shower should be located near the work-station.

Personal Protective Equipment: Personal protective equipment should be selected based upon the conditions under which this material is used. A hazard assessment of the work area for PPE requirements should be conducted by a qualified professional pursuant to OSHA regulations. For certain operations, additional PPE may be required.

Eye Protection: Safety glasses equipped with side shields should be adequate protection under most conditions of use. Wear goggles and/or face shield if splashing or spraying is likely, especially if material is heated above 125°F(or 51°C). Have suitable eye wash water available.

Hand Protection: Use gloves constructed of chemical resistant materials such as neoprene or heavy nitrate rubber if frequent or prolonged contact is expected. Use heat-protective gloves when handling product at elevated temperatures.

Body Protection: Use clean and impervious protective clothing (e.g., neoprene or Tyvek®) if splashing or spraying conditions are present. Protective clothing may include long-sleeve outer garment, apron, or lab coat. If significant contact occurs, remove oil-contaminated clothing as soon as possible and promptly shower. Launder contaminated before reuse or discard. Wear heat protective boots and protective clothing when handling material at elevated temperatures.

Respiratory Protection: Vaporization or misting is not expected at ambient temperatures. Therefore, the need for Respiratory protection is not anticipated under normal use conditions and with adequate ventilation. If elevated airborne concentrations above applicable workplace exposure levels are anticipated, a NIOSH-approved organic vapor respirator equipped with a dust/mist

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

prefilter should be used. Protection factors vary depending upon the type of respirator used. Respirators should be usedin accordance with OSHA requirements (29 CFR 1910.134).

General Comments: Use good personal hygiene practices. Wash hands and other exposed skin areas with plenty of mild soap and water before eating, drinking, smoking, use of toilet facilities, or leaving work. DO NOT use gasoline, kerosene, solvents, or harsh abrasive skin cleaners. Since specific exposure standards/control limits have not been established for this product, the "Oil Mist, Mineral" exposure limits shown below are suggested as minimum control guidelines.

Occupational Exposure GuidelinesSubstance: Oil mist, mineral

Applicable Workplace Exposure Levels:
TWA: 5 (mg/M^3) from ACGIH (United States)
STEL: 10 (mg/M^3) ACGIH (United States)
TWA: 5 (mg/M^3) from OSHA (United States)

SECTION 9 – PHYSICAL AND CHEMICAL PROPERTIES

Physical State: Liquid
Color: Light Amber to Amber
Odor: Mild Petroleum Odor
Specific Gravity: 0.88(Water = 1)
pH: N/A
Vapor Density: GT 1 (Air = 1)
Boiling Point/Range: N/A
Melting/Freezing Point: N/A
Vapor Pressure: 0.01mm of Hg (20°C)
Viscosity (cSt @ 40°C): 65
Density: 7.30 Lbs/gal.
Solubility in Water: Insoluble in cold water
Viscosity (ASTM D2161) = 339 SUS @ 100°F
Volatile Characteristics: Negligible volatility
Additional Properties: Gravity, API (ASTM D287) = 30.0 @ 60°F

SECTION 10 – STABILITY AND REACTIVITY

Chemical Stability: Stable

Hazardous Polymerization: Not expected to occur.

Conditions to Avoid: Keep away from extreme heat, sparks, open flame and strong oxidizing conditions.

Materials Incompatibility: Strong oxidizers

Hazardous Decomposition Products: No additional hazardous decomposition products were identified other than the combustion products identified in section 5 of this MSDS.

SECTION 11 – TOXICOLOGICAL INFORMATION

For other health-related information, refer to the Emergency Overview on Page 1, and the Hazards identification in Section 3 of this MSDS.

Toxicity Data:
ORAL (LD50):
DERMAL (LD50):

Distillates, petroleum, solvent-refined heavy paraffinic:
Acute: >5000 mg/kg [Rat].
Acute: >2000 mg/kg [Rabbit].

Distillates, petroleum, solvent-refined heavy paraffinic:
Mineral oil mists derived from highly refined oils are reported to have low acute and sub-acute toxicities in animals. Effects from single and short-term repeated exposures to high concentrations of mineral oil mists well above applicable workplace exposure levels include lung inflammatory reaction, lipoid granuloma formation and lipoid pneumonia. In acute and sub-acute studies involving exposures to lower concentrations of mineral oil mists at or near current work place exposure levels produced no significant toxicological effects. In long term studies (up to two years) no carcinogenic effects have been reported in any animal species tested. Analyses conducted by method IP346 indicate that the polycyclic aromatic concentration of this mineral oil is below 3.0 weight percent.

THE INFORMATION HEREIN IS GIVEN IN GOOD FAITH, BUT NO WARRANTY, EXPRESS OR IMPLIES, IS MADE.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

SECTION 12 – ECOLOGICAL INFORMATION

Eco-toxicity: Analysis for ecological effects has not been conducted on this product. However, if spilled, this product and any contaminated soil or water may be harmful to human, animal, and aquatic life. Also, the coating action associated with petroleum and petroleum products can be harmful or fatal to aquatic life and waterfowl.

Environmental Fate: An environmental fate analysis has not been conducted on this specific product. However, plants and animals may experience harmful or fatal effects when coated with petroleum-based products. Petroleum-based (mineral) lube oils will normally float on water. In stagnant or slow-flowing waterways, an oil layer can cover a large surface area. As a result, this oil layer might limit or eliminate natural atmospheric oxygen transport into the water. With time, if not removed, oxygen depletion in the waterway might be enough to cause a fish kill or create an anaerobic environment.

SECTION 13 – DISPOSAL CONSIDERATIONS

Hazard characteristic and regulatory waste stream classification can charge product use. Accordingly, it is the responsibility of the user to determine the proper storage, transportation, treatment and/or disposal methodologies for spent materials and residues at the time of disposition.

Conditions of use may cause this material to become a hazardous waste, as defined by Federal or State regulations. It is the responsibility of the user to determine if the material is a hazardous waste at the time of disposal. Transportation, treatment, storage and disposal of waste material must be conducted in accordance with RCRA regulations (see 40 CFR 260 through 40 CFR 271). State and/or local regulations may be more restrictive. Contact the RCRA/Superfund Hotline at (800) 424-9346 or your regional US EPA office for guidance concerning case specific disposal issues.

SECTION 14 – TRANSPORT INFORMATION

DOT Status:	Not regulated by the U.S. Department of Transportation as a hazardous material.
Proper Shipping Name:	Not regulated
Hazard Class:	Not regulated
Packing Group(s):	Not applicable
UN/NA ID:	Not regulated
Reportable Quantity:	A Reportable Quantity (RQ) has not been established for any components of this material.
Placards:	Not applicable.
Emergency Response Guide No.:	Not applicable.
HAZMAT STCC No.:	2911990
MARPOL III Status:	Not a DOT "Marine Pollutant" per 49 CFR 171.8

SECTION 15 – REGULATORY INFORMATION

TSCA Inventory: This product and/or its components are listed on the Toxic Substance Control Act (TSCA) Inventory.

SARA 302/304: The Superfund Amendments and Reauthorization Act of 1986 (SARA) Title III requires Facilities subject to Subparts 302 and 304 to submit emergency planning and notification information based on Threshold Planing Quantities (TPQs) and Reportable Quantities (RQs) for "Extremely Hazardous Substances" listed in 40 CFR 302.4 and 40 CFR 355. No components were identified.

SARA 311/312: The Superfund Amendments and Reauthorization Act of 1989 (SARA) Title III Requires

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

facilities subject to this subpart to submit aggregate information on chemicals by "Hazard Category" as defined in 40 CFR 370.2. This material would be classified under the following hazard categories:No SARA 311/312 Hazard categories identified.

SARA 313: This product contains the following components in concentrations above de minimis levels that are listed as toxic chemicals in 40 CFR Part 372 pursuant to the requirements of Section 313 of SARA: No components were identified.

CERCLA: The Comprehensive Environmental Response, Compensation, and Liability Act of 1980 (CERCLA) requires notification of the National Response Center concerning release of quantities of "hazardous substances" equal to or greater than the reportable quantities (RQs) listed in 40 CFR 302.4. As defined by CERCLA, the term "hazardous substance" does not include petroleum, including crude oil or any fraction thereof which is not otherwise specifically designated in 40 CFR 302.4. This product or refinery stream is not known to contain chemical substances subject to this statute. However, it is recommended that you contact state and local authorities to determine if there are any other reporting requirements in the event of a spill.

CWA: This material is classified as an oil under Section 311 of the Clean Air Act (CWA) and the Oil Pollution Act of 1990 (OPA). Discharges or spills which produce a visible sheen on waters of the United States, their adjoining shorelines, or into conduits leading to surface waters must be reported to the EPA's National Response Center at (800) 424-8802.

California Proposition 65: This material may contain the following components which are known to the State of California to cause cancer, birth defects or other reproductive harm, and may be subject to the requirements of California Proposition 65 (CA Health & Safety Code Section 25249.5): Toluene: 0.0008%

New Jersey Right-To-Know Label: Petroleum Oil

Additional Regulatory Remarks: No additional regulatory remarks.

SECTION 16 – OTHER INFORMATION

Revision: 4 **Revision Date:** 11/11/03

ABBREVIATIONS :

AP= Approximately EQ= Equal
GT= Grater Than LT= Less Than
NA= Not Applicable ND= No Data
NE= Not Established ACGIH= American Conference of Governmental Industrial Hygienists
AIHA= American Ind. Hygiene Association IARC= Internal Agency for Research on Cancer
NTP= National Toxicology ProgramNIOSH= National Institute of Occupational Safety and Health
OSHA= Occupational Safety and Health Administration
NPCA= National Paint and Coating Manufactures Association
HMIS= hazardous Materials Information System
NFPA= National Fire Protection Association EPA= Environmental Protection Agency

All statements, information, and data provided in this Material Safety Data Sheet are believed to be accurate and reliable, but are presented without guarantee, warranty, or responsibility of any kind, expressed or implied on our part. Users should make their own investigations to determine the suitability of the information of products for their particular purpose.

Welch and DUOSEAL is a registered trademark of Welch Rietschle Thomas.

THE INFORMATION HEREIN IS GIVEN IN GOOD FAITH, BUT NO WARRANTY, EXPRESS OR IMPLIES, IS MADE.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

NOTES:

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

WARRANTY

UNPACKING

Inspect the pump carefully. If any damage has occurred, file claim with the carrier immediately. Save the shipping container for carrier to inspect.

OPERATING PUMP

Refer to the enclosed Instruction/Operation Manual for all information to properly operate and maintain the pump.

WARRANTY

This Welch Vacuum product is warranted to be free from defects in material and workmanship. The liability of Welch Rietschle Thomas under this warranty is limited to servicing, adjusting, repairing or replacing any unit or component part which in the judgment of Welch Rietschle Thomas has not been misused, abused or altered in any way causing impaired performance or rendering it inoperative. No other warranties are expressed or implied. The method of executing this warranty: servicing, adjusting, repairing or replacing shall be at the discretion of Welch Rietschle Thomas. Vacuum pumps that have been used for any period, however short, will be repaired under this warranty rather than replaced.

The warranty is effective for one year from the date of original purchase when:

1. The warranty card has been completed and returned.
2. The product is returned to the factory or other designated service centers, freight prepaid.
3. The product in our judgment is defective through no action or fault of the user.

If the product has become defective through misuse, abuse, or alteration, repairs will be billed regardless of the age of the product. In this event, an estimate of the repair costs will be submitted and authorization of these charges will be required before the product is repaired and returned.

To reduce additional charges and delays either within or outside of the warranty period, contact Welch Rietschle Thomas @847-676-8800 for a return authorization number. Products without a return authorization number will be refused by our receiving department. Before shipping, properly pack the pump, insure it against loss or damage, and on the outside of the pump packaging and the packing slip write in the return authorization number. Pumps damaged due to improper packaging are the customer's responsibil-

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Sales and Service Offices
For Welch Rietschle Thomas

USA, Canada, Latin America
Welch Rietschle Thomas
P.O. Box 183
Skokie, IL 60076-0183 USA
Phone: 847-676-8800
Fax (Technical): 847-677-8606
Fax (Orders): 920-451-4397
Email: welchvacuum@thomasind.com

Germany
Rietschle Thomas Puchheim GmbH
Siemensstraße 4
Gewerbegebiet Nord
82178 Puchheim
Germany
Phone: +49 89 80 900 136
Fax: +49 89 80 900 137
Email: mgraham@rtpumps.com

Switzerland/Liechtenstein
Rietschle Thomas Switzerland
Frauenfelder Str. 49
CH-8370 Siranch TG Switzerland
Phone: 41 71 969 3500
Fax: 41 71 969 3501
E-mail: info@rtpumps.ch

Australia
Rietschle Thomas Australia, Pty. Ltd.
30 Bearing Road, Seven Hills
New South Wales Australia 2147
Phone: 61 (2) 96 20 70 00
Fax: 61 (2) 96 20 79 55
Email: tiap@rtpumps.com

New Zealand
Rietschle Thomas New Zealand
P.O. Box 45
40 Anvil Road, Silverdale
Auckland, New Zealand 1330
Tel: 64 9 426 0370
Fax: 64 9 426 0371
Email: tiap@rtpumps.com

United Kingdom
Rietschle Thomas UK
Unit 2, Alton Business Centre
Omega Park, Alton Hants
Hampshire GU34 2YU United Kingdom
Phone: 01420 544 184
fax: 01420 544 183
E-mail: ukinfo@rtpumps.com

Denmark
Rietschle Thomas Denmark, A/S
TåstruphØj 11
P O BOX 185
4300 Holbæk, Denmark
Phone: 45 59444050
Fax: 45 59444006
Email: rtpumpsdk@rtpumps.com

France
Rietschle Thomas France, S.A.S.
Zone Industrielle Liesbach
8, rue des Champs
68220 Hesingue, France
Phone: +33 3 89 70 26 76
Fax: +33 3 89 70 91 20
Email: service.commercial@rietschle.fr

Italy
Rietschle Thomas Italia, S.p.A
Via Brodolini, 17
20032 Cormano (Milano), Italy
Phone: +39 02 614 512 1
Fax: +39 02 66 50 33 99
Email: info.it@rtpumps.com

Netherlands
Rietschle Thomas Netherlands B.V.
Bloemendalerweg 52
1382 KC WEESP, Netherlands
Phone: +31 294 418686
Fax: +31 294 411706
Email: verkoop@rtpumps.com

Sweden
Rietschle Thomas Sweden AB
Götlundagatan 2
SE-12471 Bandhagen, Sweden
Phone: 46 8 447 1830
Fax: 46 8 447 1839
Email: infosweden@rtpumps.com

Mexico
Rietschle Thomas Mexico
la Privada Jesus Maria #110
San Jose del Arenal 20130
Aguascalientes, Mexico
Phone: 49 960 581
Fax: 49 960 380
E-mail: adiaz@rtpumps.com

Japan
Rietschle Thomas Japan
1794 Nippa-cho Kohoku-ku
Yokohama-shi, Kanagawa
223-0057 Japan
Phone: 81 45 533-0390
Fax: 81 45 533-0391
Email: tiap@rtpumps.com

Hong Kong/PRC
Rietschle Thomas Hong Kong
Units 1-5, 25th Floor, Metropole Square
No. 2 On Yiu St., Siu Lek Yuen
Shatin, New Territories, Hong Kong, P.R.C.
Phone: 852 2690 3502
Fax: 852 2792 4598
Email: tiap@rtpumps.com

Taiwan
Rietschle Thomas Taiwan
2F, No. 9, Rong 11, Lane 327
Chung Shan Road, Section 2
Chung Ho City, Taiwan
Tel: 886 2 2246 4236
Fax: 886 22 2245 0216
Email: tiap@rtpumps.com

Korea
Rietschle Thomas Korea
Room 501, Joong Ang Induspia Building
517-13, Sang Dae Won Dong, Joong Won Ku
Sung Nam City, Kyung Ki Do
Korea 462-713
Tel: 82 31 740 9533
Fax: 82 31 740 9538
Email: tiap@rtpumps.com

Brazil
Rietschle Thomas Brasilien Comérico Ltda.
Rua Shigeru Hayashi, 49 CIC
81170-640 CURITIBA - PR
Brazil
Tel: 55 41 347 0492
Fax: 55 41 347 0290
Email: mschmidt@rtpumps.com

Vacuum Pump Repair Facility

7301 North Central Ave.
Skokie, Illinois 60077
Phone: (847) 676-8800 Ext. 1
Fax: (847) 677-8606

OWNER'S MANUAL
For
DUOSEAL VACUUM PUMPS
MODELS 1374, 1376, 1392, 1380, 1399, 1400, 1402, &1405
Part No. 67-0777

Copyright®1993-2004 Welch Rietschle Thomas
DUOSEAL and Welch are registered trademarks of Welch Rietschle Thomas.

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at **artisantg.com** for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

ARTISAN® TECHNOLOGY GROUP